MAGNET-RELAX™ SMART EYE MASSAGER

USER MANUAL

BETTER THAN A MASSAGE,
A UNIQUE EXPERIENCE

Magnetinnov
MEDICAL INNOVATION

TABLE OF CONTENTS

EN

DE

FR

IT


ES

Product characteristics	02
Product configuration	03
Product functions	04
Using the product	05
Bluetooth connection and music control	06
Warning	06-07
Maintenance and storage	08
Troubleshooting	09

Product characteristics

EN Automatic mode :

DE
FR
IT
ES

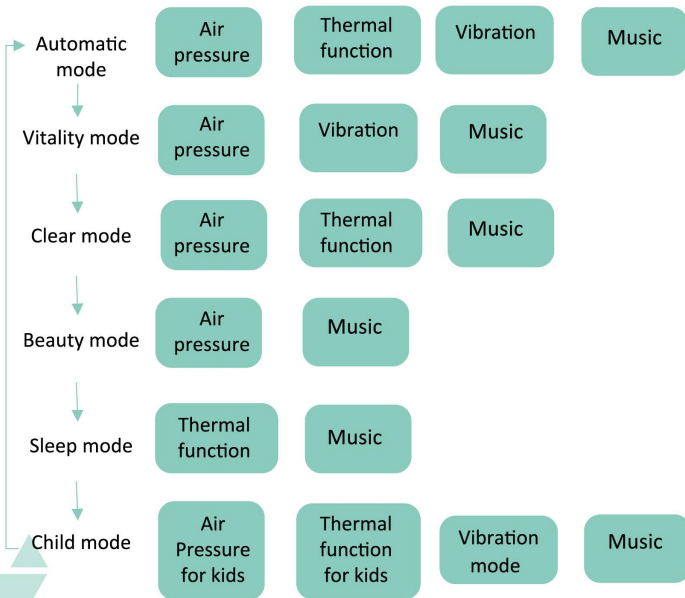Product configuration


EN Your massage device pack includes the following:

DE Eye mask, USB cable, user manual.


FR

IT

ES

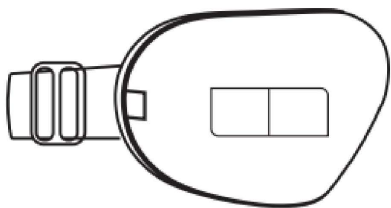
Eye mask

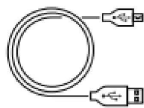USB cable

User manual

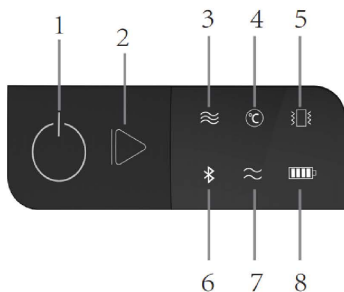1. On/off/mode switch button (long holding for two seconds to turn mask on/off, short holding for mode switch)
2. Volume control, 3 volume levels available with a possibility to fully turn the music off (Default volume level selection at start, long pressing to change the volume level)
3. Air pressure massage function
4. Thermal function
5. Vibration massage function
6. Bluetooth indicator
7. Child mode activation
8. Charging mode

Using the product

EN
DE
FR
IT
ES

Prior to usage

Use the 5V/1A adaptor (not provided with the product) for charging the product. You can also plug the mask with the use of the USB cable provided, into another device, such as a laptop or an external battery.

Turning the device on

Hold the "on/off" button for 2 seconds in order to turn the product ON. After being switched on, the light indicator remains on and the device begins operating. The duration of a session is of 10 minutes. Once the session is over, the device is automatically switched off.

Switching the modes and their sequence

Hold the "on/off" button for less than one second, in order to change the mode. The modes will be displayed in the following sequence: Automatic mode, Vitality mode, Basic mode, Beauty mode, Sleep mode and Child mode.

Charging

When the product indicates that the battery is low, you will need to plug your device in. During the charging period, the red indicator remains on, it switches off once the device is fully charged. The Magnetinnov would require around 2 to 4 hours of charging.

Warning: It is strongly recommended not to use the product while it is charging.

Bluetooth connection and music control

EN
DE
FR
IT
ES

1. Music is not activated on the product by default, you will need to establish the Bluetooth connection after switching it on. Activate Bluetooth on the device and search for the device name EM-1901 in order to finalize the device pairing procedure.
2. Press on the music control button briefly in order to adjust the volume. You can also adjust the volume of the music from your phone, after it has been paired with the device by Bluetooth.

Warning

This device can be used by children starting from the age of 8 years old and by people with reduced physical, sensorial or mental abilities, solely under supervision and by following the instructions on the usage of the device, under the condition of a clear awareness of the following warnings. Read the user manual very carefully prior to using the product. The company is not responsible for the consequences of an improper usage of the device, or from a non-respect of the instructions outlined in this user manual.

1. Do not use while driving.
2. Do not use in environments where the temperature or humidity levels are excessively high.
3. To be used carefully by people who manifest disease symptoms.
4. Do not use the device on the areas of the body that have been recently operated or have had surgery performed on them on previous occasions.
5. Place the device on your eyes at first and then switch it on, in order to avoid damage to the inflatable bag due to overpressure.

- 6.** The usage of the mask is not recommended for individuals that have suffered from a stroke, heart disease, hypertension or hypoglycemia.
- 7.** Do not consume any alcohol prior to using the mask.
- 8.** Remove contact lenses, hair ties and any other accessories from the head and the face in order to avoid any risk of injury.
- 9.** Use the electrical adaptor which conforms to the norms of your country and pay attention to electrical safety (5V/1A, not provided).
- 10.** The product is a home health device. It cannot act as a replacement to a medical treatment by any means.
- 11.** You must close your eyes during the massage.
- 12.** If you are feeling unwell during the usage, stop the session immediately and consult a healthcare professional.
- 13.** Do not use the device if you had surgery performed on your eyes or if you are suffering from any eye diseases. Do not use if your eyes are swollen.
- 14.** The device is equipped with a hot compress function during which the temperature can get reach 42°C. People who are sensible to temperature variations need to be particularly cautious when using the device.
- 15.** If the massage area has got any wounds, redness, inflammation of other physical abnormalities, you must stop the usage of the device immediately.
- 16.** We recommend a daily 10-minute session. The mask would switch off automatically after 10 minutes.
- 17.** Minors, handicapped people and people of third age must use this product under constant supervision.
- 18.** For cleaning, dry the mask with a soft cloth with alcohol and a cleaning product. Do not wet the device.

Warning:

This symbol on the device or its packaging indicates that this electronic device is not to be disposed of in the regular unsorted wastebins. The usage of sorted disposal bins is required. Contact your local authorities in order to obtain information on the available waste collection systems. In the event where electronic devices are disposed of in regular wastebins, dangerous substances may mix with groundwaters and return to the food chain, which could be compromising to your health and the well-being of those around you.

Maintenance and storage

Storage

If the device is being stored for more than 3 months, it needs to be unwrapped and charged prior to using it. The device does not have to be too close to the floor, the wall or the roof. The air ventilation has to be good in order to avoid direct sunlight exposure and an erosion by corrosive gases.

Transportation

The product is equipped with an anti-vibration system, it is well adapted for transportation by air, train and sea. Rain, snow splashes and collisions are to be avoided.

Power verification

Check the electric power. If the electric tension level surpasses the standard limit, do not use the device.

Maintenance

Take the necessary precautions to ensure that the device is protected from humidity, rats and dust.

During cleaning, the device must be switched off and must not be charging. After usage, the device is to be stored in a safe place. Do not place any heavy items on it. In the event of malfunctioning, contact a professional, do not try to disassemble the device on your own.

Troubleshooting

EN If a product develops a defect during its functioning, proceed to a verification of it. If the issue cannot be resolved or if the product is still not functioning, proceed to unplugging the device from electrical power and contact the company's customer service team immediately in order to request a diagnostic procedure.

DE
FR
IT
ES Do not disassemble or try to repair the device on your own.

Issue:	Possible reasons:	Actions and solutions:
The device issues mechanical or other noises.	The normal sound of the inflation of the air pump of the mechanism that produces vibrations.	<ul style="list-style-type: none">• No actions are required.• If the noise of the mask is disturbing you, we suggest you to use the sleep or child mode.
The device switches off automatically. The device does not switch on, the buttons are not responsive. The strength of the massage is weak.	<ol style="list-style-type: none">1- The device has got low or no battery.2- The 10-minute timer has run out.	<ol style="list-style-type: none">1- Charge the device.2- Switch the device off and back on

Technical parameters:

Product name	Eye massage device
Product number	EM-1901
Power	5W
Power adaptor	DC5V 1A
Size	Approx. 220mmX80mm(L*H)
Net weight	Approx. 305g

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Magnet-Innov©