
RENPHO

User Manual

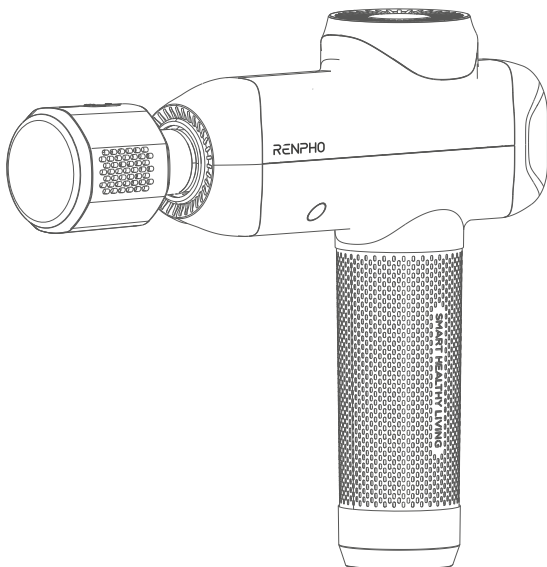Model: R-C001SH

Massage Gun

Thank you for choosing the RENPHO Massage Gun! Please read this user manual carefully and use the device only as directed. RENPHO accepts no responsibility for any cause of injuries due to misuse or negligence to follow the instructions outlined. Please keep this manual for future reference.

Important Safety Instructions

READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.

When using this device, the following basic precautions should be always followed:

1. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
2. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
3. Children should be supervised to ensure that they do not play with the device.

DANGER - To reduce the risk of electric shock:

1. Do not reach for the device that has fallen into water.
2. Do not use while bathing or in a shower.
3. Do not place or store the device in the areas where it may fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING - To reduce the risk of fire, electric shock, or injury to persons:

1. Please consult your doctor before using if you are a user of the following persons:
 - 1) Children under 8
 - 2) Pregnant women
 - 3) Persons with blood stasis
 - 4) Persons with fractures

- 5) Persons with heart disease, with a pacemaker, defibrillator, or metal contents in the body
- 6) If you have experienced joint surgery within 90 days, do not use the device on your body within 3 inches of the surgical position

If you feel unwell or have any unusual discomfort, please stop using the device immediately.

2. Do not use the massage gun on

- 1) Sensitive skin, swollen or inflamed areas, in areas where there is poor circulation, where skin eruptions are present, or in the area of unexplained abdominal pain
 - 2) Head, face, cervical vertebrae, body bones
 - 3) Varicose veins
 - 4) Tumor site
 - 5) Throat area
 - 6) Genital area (do not use this device for sexual purposes. It is a powerful massaging device and may cause severe bodily harm.)
 - 7) An unconscious or sleeping person
 - 8) Any possible fracture site
 - 9) Positions close to the bony areas: the humerus, the back of the foot, the back of the hand, etc
3. The device should never be left unattended when plugged in. Unplug from the outlet when not in use, and before cleaning, putting on or taking off parts.
 4. Do not operate under a blanket or pillow. Excessive heating may occur and cause fire, electric shock, or injury to persons.
 5. Never operate the device with the air openings blocked. Keep the air openings free of lint, hair, etc.
 6. Close supervision is necessary when this device is used by, on or near children, invalids, and disabled persons.
 7. Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 8. Never operate this device if it is not working properly, has been damaged, or dropped into water. Send the device to the service center for examination and repair.
 9. Never drop or insert any object into any opening.

10. Do not operate in the areas where aerosol (spray) products are being used or where oxygen is being administered.
11. Do not keep the device near a magnetic field or in the areas where it can be affected by a magnetic field as it may cause batteries to discharge.
12. Do not use the device before bed as it has a stimulating effect and may delay sleep.
13. Avoid touching the gap between the massage head and body with your finger or hair to avoid pinching.
14. It is recommended to start the massage with speed level 1. Do not use for more than 30 minutes at a time. Excessive use could lead to slight bruises or device overheating that shortens its lifespan.
15. Keep the device away from heat and avoid using it in places where the surrounding temperature is relatively hot.
16. Do not use the device with wet hands.
17. Only use the original massage head officially provided in the package. The use of a third-party massage head may cause damage to the device and void the warranty.
18. Do not drop, apply excessive force or heavy weight over the device to avoid malfunction and/or damage.
19. If a high temperature or foul odor occurs while using the device, stop using the device immediately and unplug it.
20. It is prohibited for non-professional personnel to repair the device or replace its parts.
21. Use only the massage head that produces the best results.
22. Do not press the massage head vigorously or cover the device during a massage.
23. Do not use the device in close proximity to long hair, loose clothing or jewelry.
24. Do not massage the same body part for a prolonged period of time. If you feel unwell, please stop using the device immediately.
25. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
26. The device has a heated surface. Persons insensitive to heat must be careful when using the appliance.

About the Battery

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short-circuited by metal objects.
3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

About the Product

1. What's in the Box

1 × Massage Gun	1 × Type-C Charging Cable
1 × User Manual	2 × Rubber Sleeve
1 × Travel Case	6 × Massage Head

2. Product Parameter

Model: R-C001SH

Input: 5V $\equiv$ 2A

Rated Power: 25W

Battery Voltage: DC 11.1V

Battery Capacity: 2500mAh, 27.75Wh

Amplitude: 12mm

Bluetooth: Maximum RF power: 2.3dBm

Operation frequency: 2402MHz to 2480MHz

Speed: 1-6

RPM: 1850-2600r/min

Charging Port: Type-C


Noise: ≤ 55 dB

Auto-off Timer: 10 minutes

Dimension: 6.7×8.4×2.4inch

3. Product Overview

Replaceable Massage Head

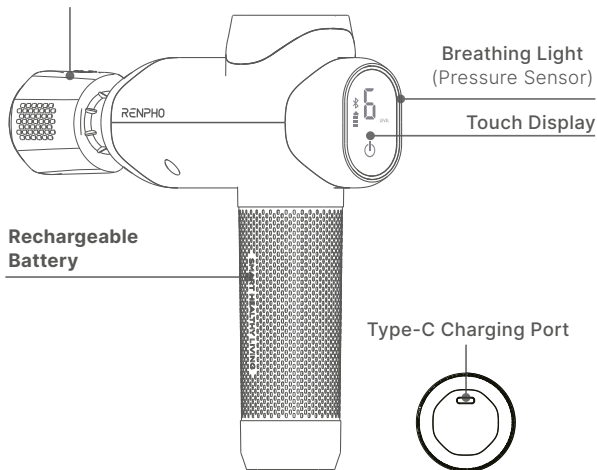**Ball head**

For large muscle groups like quads, glutes, extreme muscles, and joint soreness.

**Air cushioning head**

For sensitive muscles.

**U-shaped head**

For spine, neck and Achilles tendon.

**Shovel-shaped head**

For abdominal muscles and lower back.

**Bullet head**

For joints, deep tissues, trigger points and small muscle areas like feet and wrists.

**Heating & Cooling Head**

For all body parts.

**Rubber Sleeve**

Each massage head is pre-installed with a rubber sleeve that ensures the massage head is securely attached to the massage gun. If it is damaged, replace it with a spare one.


Breathing Light Instructions

Status	Breathing Light (Pressure Sensor)
Ideal pressure applied	Blue
Stall warning: Excess pressure applied	Pink
Stalled: Too much pressure applied	Red
Very low battery	Flashes red 3 times
Charging	Flashing red (Battery level < 50%)
	Flashing blue (50% ≤ Battery level < 100%)
Fully charged	Off

Before Using

Charging the device

***Fully charge the device before the first use.**

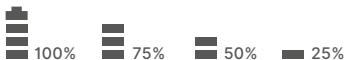1. When the lowest battery level indicator flashes or the breathing light flashes red 3 times, it means that the device is in (very) low battery. Please charge the device using the provided Type-C charging cable and a DC 5V 2A/3A charger (not included).
2. It takes about 3.5 hours to fully charge the device. When charging the device, the battery indicators will flash in sequence. The indicators will turn solid once the charging is complete.

Note: The device cannot be turned on while charging.

Using the Product

Step 1: Choose a preferred massage head and attach it. Ensure the massage head is fully inserted into the slot.

Note:


- a. Ensure the device is off before you swap the massage head.
- b. To remove the air cushioning head, firmly grasp the head attachment joint and pull away from the device.

Bluetooth Indicator
Once BT is connected, its
indicator will be on.

Battery Level Indicator

Power/Speed Button

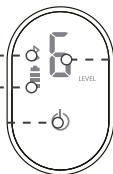Speed Level Display

Step 2: Press and hold  for 2s to turn on/off the device. Once the device is turned on, Bluetooth will be enabled by default with its indicator off.

Step 3: Press  again to cycle through the speed levels among level 1/2/3/4/5/6.

Step 4: Target your muscle groups for massage with the device. It is recommended to rest the massage head perpendicularly over your muscle and begin massaging on the lowest setting and work your way up to higher speeds as you get used to the intensity.


Heating & Cooling Massage Head

- Long press to enable the heating function.
Red indicator will turn on.
- Short press to cycle through heating levels: 1-2-3-4
- Long press again to disable the heating function.

Heating Temperature: 100~117°F/38-47°C.

Timer: 10mins

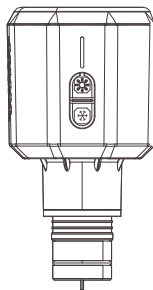- Long press to enable the cooling function. Blue indicator will turn on.
- Short press to cycle through cooling levels: 1-2-3-4
- Long press again to disable the cooling function.

Cooling Temperature: 46~57°F/8-14°C. Timer: 10mins

Note:

- Auto shut-off function: The device and the heating & cooling massage head will shut off automatically after 10 minutes of continuous use.
- Locked-motor protection: If you apply excessive pressure on the device during your session, the device will stall to prevent injury or bruising.
- It is recommended to start every massage session with minimal intensity before you get used to it. Do not use for more than 30 minutes at a time.

App Guide

Step 1: Download the app

Search and download the “RENPHO Health” App from Apple App Store / Google Play or scan the QR code below to download the RENPHO Health App.

RENPHO Health | 🔍

* Google Play is a trademark of Google LLC
Apple and the Apple Logo are trademarks of Apple Inc, registered in the U.S. and other countries

Note: Due to continuous updates and improvements, the “RENPHO Health” App may appear slightly different.

Step 2: Log in/Sign up

- Log in directly with an existing RENPHO App account.
 - If you don't have any RENPHO App account, tap “**Sign up**”.
- Register a new account of RENPHO Health App using your email address and complete your profile information.

Step 3: Add Device

***Enable Bluetooth on your smartphone before adding the device.**


-Method 1:

- Press and hold the Power button for 2s to turn on the device.
- A window will pop up on the phone screen and tap “**Connect**”.

-Method 2:

- Tap "+" on the top right corner of the "Home" page to pair manually.
 - Select "Message" and tap "Message gun".
 - A window will pop up on the phone screen and tap "+".
- Note:** If the device is connected to the App and left idle, it will automatically shut down in 2 minutes.

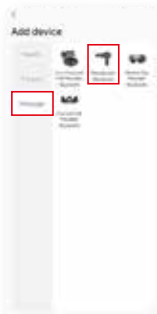Step 4: Message

- Method 1: Message in free mode

- After a successful connection, you will enter the "Free mode" page.
- If you choose "Manual" mode (Default), you can choose the speed levels among level 1-6 manually to get your personalized message.
- If you choose "Auto" mode, the device will automatically increase the speed by 1 gear every 30 seconds. Tap  to reset the time.

- Method 2: Choose the video course

- Tap "Course" and choose "Recommendation" for specific session or choose "Body" for specific muscle group massage.
- Choose your preferred course. Each session consists of corresponding videos course(s) with tips on which massage head to use and duration.
- Tap "Start" to follow it to massage.

*** When the video is playing, you can tap it for additional settings.**

- 1) Tap  to pause the video or choose to exit/continue the video.
- 2) Tap  to freely choose another course.
- 3) Tap Auto  to switch between Auto mode (Default) and Manual mode.
- 4) Tap  to turn on/off the background music or adjust the music volume.

Cleaning and Maintenance

1. Place the device in a cool, dry place when not in use.
2. Always turn off the device before cleaning. Wipe with a soft cloth slightly dampened with non-abrasive cleaners. Keep away from all solvents and harsh detergents.
3. Do not attempt to disassemble or repair the device. The device does not contain user-serviceable parts.

Troubleshooting Guide

Problem	Cause	Solution
The device cannot be turned on.	Short press the Power button.	Press and hold the Power button for 2s to turn the device on.
	The battery is dead.	Charge the device.
The device shuts down suddenly.	Auto-off timer: The device will automatically shut down after 10 minutes of continuous use.	Turn on the device again. It is not recommended to massage for more than 30 minutes at a time.
	Built-in locked-motor protection: When the massage head is pressed with excessive force, the protection will be triggered, and the breathing light will flash red 3 times before the device automatically shuts down.	
	Built-in overheat protection: It will automatically shut off the massage gun when overheating occurs.	Wait for the device to cool down before turning it on again.
No Bluetooth connection.	Bluetooth on the smartphone is not enabled.	Enable Bluetooth on the smartphone.
	The massage gun is connected by another mobile device.	Unbind the massage gun on another mobile device.
The device does not respond while operated via the App.	Bluetooth is disconnected.	Keep your smartphone close to the device and reconnect.

Problem	Cause	Solution
Slow App loading.	Poor Internet.	Operate in a good Internet environment.
The video courses cannot be downloaded.	No/poor Internet.	Operate in a good Internet environment.
If you encounter any other problems, please feel free to contact our customer service.		

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

Email: support-eu@renpho.com (UK&EU)

* For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

3. Personal data

All end-users of WEEE are reminded that you are responsible for deleting personal data from the WEEE to be disposed of.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/53/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

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