
RENPHO

Smart Health Simplified

User Manual

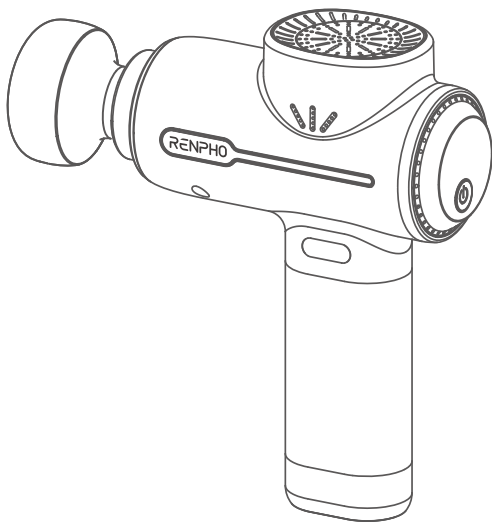Model: R-C002SH

Massage Gun

Contents

English..... EN - 01 ~ EN - 17

Deutsch..... DE - 18 ~ DE - 33

Français..... FR - 34 ~ FR - 49

Italiano..... IT - 50 ~ IT - 65

Español..... ES - 66 ~ ES - 81

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.

When using this device, the following basic precautions should be always followed:

1. The device is intended for household use only and not intended for medical use.
2. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
3. Children should be supervised to ensure that they do not play with the appliance.

DANGER - To reduce the risk of electric shock:

1. Do not attempt to retrieve the device if it falls into water.
2. Do not use while bathing or in a shower.
3. Do not use the device with wet hands.
4. Do not place or store the device in areas where it may fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING - To reduce the risk of fire, electric shock, or injury to persons:

1. Please consult your doctor before using if you are a user of the following persons:

- 1) Children under 8
- 2) Pregnant women
- 3) Persons with blood stasis
- 4) Persons with fractures
- 5) Persons with heart disease, a pacemaker, defibrillator, or metal contents in the body
- 6) If you have experienced joint surgery within 90 days, do not use the device on your body within 3 inches of the surgical site.

If you feel unwell or have any unusual discomfort, please stop using the device immediately.

2. Do not use the massage gun on

- 1) Sensitive skin, swollen or inflamed areas, in areas where there is poor circulation, where skin eruptions are present, or in the area of unexplained abdominal pain
 - 2) Head, face, cervical vertebrae, body bones, varicose veins, tumor site, throat area, any possible fracture site, positions close to the bony areas: the humerus, the back of the foot, the back of the hand, etc
 - 3) Genital area (Do not use this device for sexual purposes. It is a powerful massaging device and may cause severe bodily harm.)
 - 4) Unconscious or sleeping persons
-
3. Do not operate under a blanket or pillow. Excessive heating may occur and cause fire, electric shock, or injury to persons.
 4. Never operate the device with the air openings blocked. Keep the air openings free of lint, hair, etc.
 5. Close supervision is necessary when this device is used by, on or near children, invalids, and disabled persons.
 6. Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 7. Never operate this device if it is not working properly, has been damaged, or dropped into water.
 8. Never drop or insert any object into any opening.
 9. Do not operate in the areas where aerosol (spray) products are being used or where oxygen is being administered.
 10. Do not keep the device near a magnetic field or in the areas where it can be affected by a magnetic field as it may cause batteries to discharge.
 11. Avoid touching the gap between the massage head and body with your finger or hair to avoid pinching.
 12. Keep the device away from heat and avoid using it in places where the surrounding temperature is relatively hot.
 13. Do not drop, apply excessive force or heavy weight over the device to avoid malfunction and/or damage.
 14. If a high temperature or foul odor occurs while using the device, stop using the device immediately.
 15. Do not use the device in close proximity to long hair, loose clothing or jewelry.

16. Avoid applying excessive pressure on the massage head and refrain from massaging the same area continuously for an extended period of time.
17. The device has a heated surface. Persons insensitive to heat must be careful when using the device.

About the Battery

- 1. Do not expose the device or battery to excessive temperatures.**
- 2. Be aware of the risk of terminals of the battery-operated device or battery being short- circuited by metal objects.**
- 3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.**
- 4. Do not dismantle, open or shred battery.**
- 5. Do not subject cells or batteries to mechanical shock.**
- 6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.**
- 7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.**
- 8. Do not leave a battery on prolonged charge when not in use.**
- 9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.**
- 10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).**

SAVE THESE INSTRUCTIONS

WHAT'S IN THE BOX

1 x Massage Gun
5 x Massage Head
1 x Type-C Charging Cable
2 x Rubber Sleeve
1 x Travel Case
1 x User Manual

SPECIFICATIONS

Product Name: Massage Gun

Model: R-C002SH

Input: 5V $\equiv$ 2A

Amplitude: 9mm

RPM: 1800/2050/2300/2500/2800r/min

Auto-off Timer: 10mins

Battery Voltage: DC 7.4V

Battery Capacity: 2000mAh 14.8Wh

Charging Time: $\leq$ 1.5H

Working Time: about 3H

Noise: $\leq$ 55dB

Product Dimensions: 7.3x5.7x2.3 inch/186x147x60mm

Product Weight: 1.5lb/676g

Environmental Conditions:

Temperature: 5°C - 35°C


Relative Humidity: $\leq$ 90%;

Atmospheric Pressure: 70kPa-106kPa

PRODUCT OVERVIEW

Interchangeable Massage Head

Breathing Light

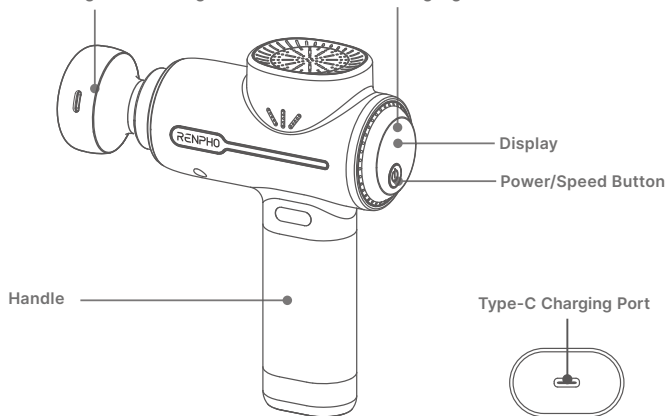MASSAGE HEAD ATTACHMENTS

Ball head

For large muscle groups like quads, glutes, extreme muscles, and joint soreness.

Crescent head

For superficial muscles like brachioradialis and extensor digitorum.

U-shaped head

For spine, neck and Achilles tendon.

Heating head

For all parts of the body.

Bullet head

For joints, deep tissues, trigger points and small muscle areas like feet and wrist.


Rubber Sleeve

Each massage head is pre-installed with a rubber sleeve. If it is damaged, replace with a spare one.

DO NOT use on your head.

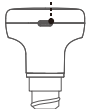Heating Head

- Type-C charging port (Use 5V 1A/2A charger, not included)
When the massage gun reaches a certain amplitude, the massage head will automatically turn on the heating function.


- Built-in 600mAh lithium battery
- Charging time: about 90 minutes
- Battery life: about 90 minutes

* When the red light inside the heating head flashes, please charge it.

CHARGING

Fully charge the device before using it for the first time.

To charge, plug the provided Type-C charging cable to the unit's charging interface and into the mobile power source or mobile charger (not included). Ensure the input should be DC 5V, 1A/2A.

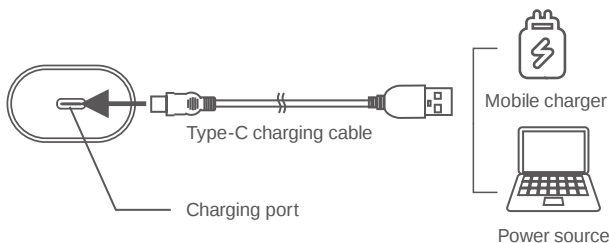
Status	Display
Low battery	Flashes red "⚡".
Normal charging	Shows current battery level
Fast charging	Shows "⚡" and current battery level

Note: The device cannot be turned on while charging.

ATTACHING THE MESSAGE HEAD

Select a preferred message head and screw it fully into the slot of the device until it is securely attached. Ensure the device is turned off when installing the message head.

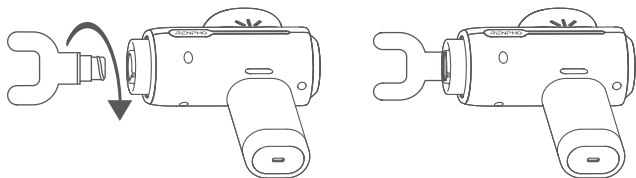USING THE PRODUCT

- Press and hold "  " for 2s to turn on the device. The indicator light will turn on and the display will show the current battery level and then "L0".
- After power on, short press "  " to cycle through speed levels: L1/L2/L3/L4/L5 (RPM: 1800/2050/2300/2500/2800r/min).
- Press and hold "  " for 2s again to turn off the device. The indicator light and the display will turn off.

Note:

- The device will enter standby status if left idle for 120s after power on.
- Auto-off function: The device shuts off automatically after 10 minutes of continuous use.
- Locked-motor protection: If you apply excessive pressure on the device during your session, the device will stall to prevent injury or bruising.
- It is recommended to start every massage session with minimal intensity before you get used to it. Do not use for more than 30 minutes at a time.


Breathing Light Instructions

Status	Breathing Light
Normal operation	Blue
A little much pressure applied	Orange
Stall warning: Excessive pressure applied	Red
Charging	Flashes green slowly when battery level $\geq 20\%$
	Flashes red slowly when battery level $< 20\%$
Fully charged	Off

APP GUIDE

Step 1: Download the App

Search and download the "Renpho Health" App from Apple App Store / Google Play or scan the QR code to download the Renpho Health App.

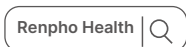Note: Due to continuous updates and improvements, the "Renpho Health" App may appear slightly different.

*Apple and the Apple Logo are trademarks of Apple Inc., registered in the U.S. and other countries.
Google Play is a trademark of Google LLC.

Step 2: Log in/Sign up

- a. Log in directly with an existing RENPHO App account.
 - b. If you don't have any RENPHO App account, tap "Sign up".
- Register a new account of Renpho Health App using your email address and complete your profile information.

Step 3: Add Device

* Enable Bluetooth on your smartphone before adding the device.

- Method 1:




- Press and hold the Power button for 2s to turn on the device.
- A window will pop up on the phone screen and tap "Connect".

- Method 2:

- Tap "+" on the top right corner of the "Home" page to pair manually.
- Select "**Massage**" and tap "**Massage gun**".
- A window will pop up on the phone screen and tap "+".


Note: If the device is connected to the App and left idle, it will automatically shut down in 2 minutes.

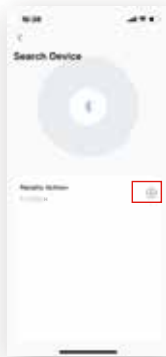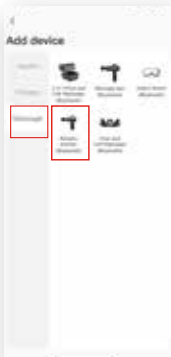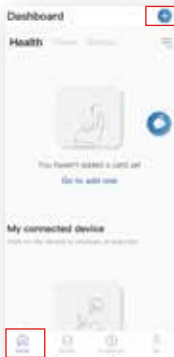Step 4: Massage



- Method 1: Massage in free mode

- After a successful connection, you will enter the "Free mode" page.
- If you choose "Manual" mode (Default), you can choose the speed levels among level 1-5 manually to get your personalized massage.
- If you choose "Auto" mode, the device will automatically increase the speed by 1 gear every 30 seconds. Tap  to reset the time.

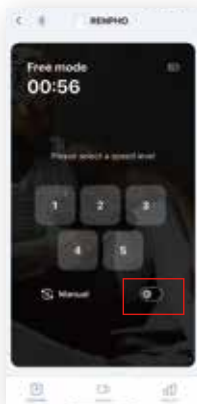
- Method 2: Choose the video course

- Tap "Course" and choose "Recommendation" for specific session or choose "Body" for specific muscle group massage.
- Choose your preferred course. Each session consists of corresponding videos course(s) with tips on which massage head to use and duration.
- Tap "Start" to follow it to massage.

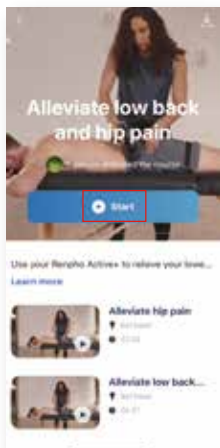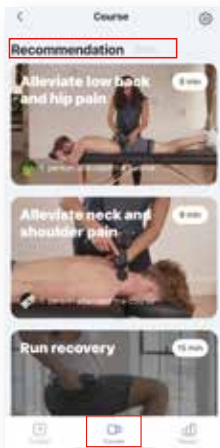*** When the video is playing, you can tap it for additional settings.**

- 1) Tap  to pause the video or choose to exit/continue the video.
- 2) Tap  to freely choose another course.
- 3) Tap  to switch between Auto mode (Default) and Manual mode. 
- 4) Tap  to turn on/off the background music or adjust the music volume 

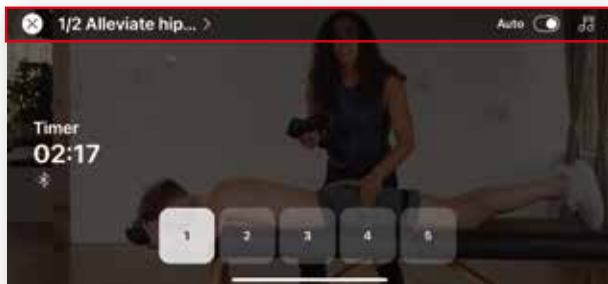CLEANING AND MAINTENANCE

1. Turn off the device prior to cleaning. Use a soft towel or cloth to wipe the outer surfaces clean.
2. Do not attempt to disassemble or repair the device. The device does not contain user-serviceable parts.
3. Place the device in a cool, dry place when not in use.

TROUBLESHOOTING GUIDE

Problem	Cause	Solution
The device cannot be turned on.	Short press the Power button.	Press and hold the Power button for 2s to turn the device on.
	Battery runs out.	Charge the device.
The device stops working suddenly.	Auto-off timer: The device automatically shuts off after 10 minutes of continuous use.	Restart the device.
	Battery runs out.	Charge the device.
E1	The motor is locked.	Restart the device.
E2	The device is overheating.	Wait for the device to cool down before restarting it.
E3	The Power button is pressed while charging.	Do not turn on the device while charging.
Bluetooth connection fails.	Bluetooth on the smartphone is not enabled.	Enable Bluetooth on the smartphone.
	The massage gun is connected by another mobile device.	Unbind the massage gun on another mobile device.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one year limited manufacturer warranty from the date of delivery. For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

-  **TEL: +44(0)7434 666088 (UK ONLY)**
Monday-Friday 9:00AM-5:00PM (GMT)
+1(844) 417 0149 (US&CA)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+49 40 23969285 (DE ONLY)
Monday-Friday 9:00AM-6:00PM (CET)
-  **Email: support-eu@renpho.com (UK&EU)**
support@renpho.com (US&CA)
support-au@renpho.com (AU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitters that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s).

Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

3. Personal data

All end-users of WEEE are reminded that you are responsible for deleting personal data from the WEEE to be disposed of.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/53/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

Maximum RF power: -2.0dBm
Operation frequency: 2402MHz to 2480MHz

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