

## QUICK START

Before use: please check to make sure massager is charged. Massager will not work while in charging mode. The default setting time is 30 minutes. The massager will automatically power down after 30 minutes. To power on/off, press and hold the power button.



1

Separate velcro and unfold product.

Fold product around lower leg until comfortable. Make sure to orientate so power button is on top.

2




3

Press and hold power button to turn on. A light will appear.

Press the power button again to switch between high and low intensities.

4

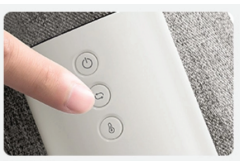


5

Press the second button to switch between three different massage modes.

Press the last button to switch between high, low or zero heat modes.

6

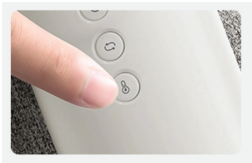

HAND WASH ONLY

## PRODUCT OVERVIEW

Thank you for choosing the OneCompress Air Compression Leg Massager. For ease of use, please read the quick start / online manual thoroughly and understand all before use. The pictures in this manual are for illustration only. Please be sure to read and understand all safety precautions.

## PRODUCT CARE

### Main Unit

The display screen of the main unit should be properly protected. Please use a soft cloth to clean the main unit.

### Fabric Material

Do not put product in washing machine. When cleaning and maintenance are needed, please use alcohol to sterilize and remove stains or wipe it with a clean wet towel.

### Storage Requirements

Please keep product out of reach for children, away from humidity, high temperature, direct sunlight and water splashes. When not in use for a long time, turn off power and keep charging once every two months.

### Troubleshooting

If the product does not work normally, please check whether the following conditions fix the issue. If the suggestions provided cannot solve the problem, please return the product to Onecompress for an exchange.

Issue / Malfunction	Possible Cause	General Solution
The working intensity of massager is weak or cannot be felt	Product cannot start	Hold power button
	Improper tightness	Readjust fit
	Low power	Please charge
The massager cannot be turned on	The battery has ran out of charge or is exhausted.	Please charge. Restart after charge is restored.

Please review safety precautions here:  
[www.onecompress.com/pages/leg-massager-safety](http://www.onecompress.com/pages/leg-massager-safety)

**Questions? Email us at: [support@onecompress.com](mailto:support@onecompress.com)**

