


Massage chair
User operation manual

APEX

APEX

AP-Bonita

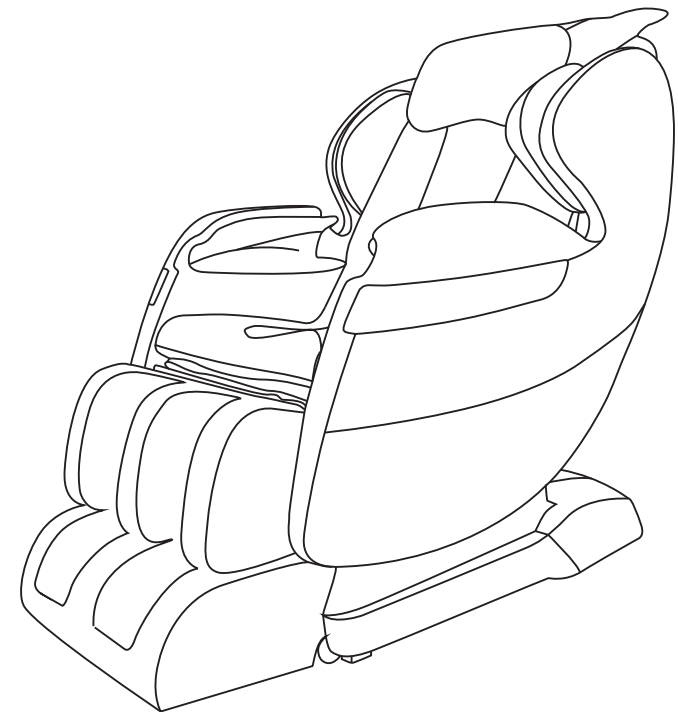

Trouble shoot

In order to prevent accidents, during abnormal conditions, please discontinue use of the chair and disconnect the power. If the following methods are not effective, please contact the service and warranty center.



Possible causes of an abnormal condition	Correct operation
During function, the chair has excessive noise.	The noise is normally caused by the design of the mechanism, not by the function.
When connect the power, can't operate massage chair.	1. Check the power plug to ensure it is well connected from the massage chair to the electrical outlet. 2. Check the power switch on the massage chair to see whether it is 'on' position. The indicator light should also be on.
During use, the controller can be operated, and the massage chair stops functioning.	Long duration operations, the massage chair temperature is over heat. If overheated, the massage chair will enter a protection status automatically. Please turn power off and wait 30 to 50 minutes in order to let the heat dissipate.
During use, controller can't be operated, and the massage chair stops functioning.	After the automatic time out period starts, repress the auto keystroke and the massage chair will begin functioning again.

Specification




Product Name	Apex Massage chair
Model	Apex-Bonita
G.W.	115kgs
Rated Power	180W
Rated Voltage	110V

Safety cautions

- Before use, please read the 'safety cautions', in order to operate correctly.
- This signal calls users attention, help users to use the machine correctly, and avoid users and other people getting hurt or damage.
- To distinguish damage degree, the cautions are divided into 'warnings' and 'attentions'. All are relevant safety conditions, please do obey.

 Warnings	The chair maybe cause minor injury if the user apply incorrect operation.
 Attentions	The chair maybe cause personal injury if the user apply incorrect operation.


< Examples >

	△ This signal stands for warning and attention. (The left signal means normal warning and attention)
	⊘ This signal means that it is not allowed. Picture represents detailed indicative content (Ex: the left picture means disassembly is not allowed)
	● Signal means that must be forced or indicated. Picture shows that detailed indicative content (Ex: the left picture represents taking off the plug from power supply)


After read, please keep the manual for your records and future use.

Maintenance and preservation

Maintenance

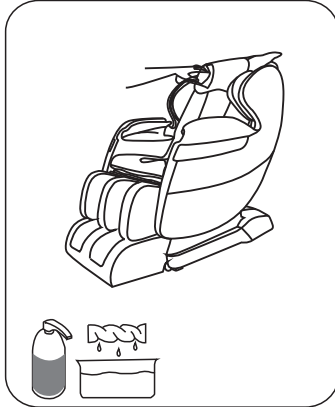


- Use a dry fabric to clean the sweat.
- When cleaning stains, please use a diluted liquid soap solution.
- Use only a damp cloth and softly clean the massage chair surface.

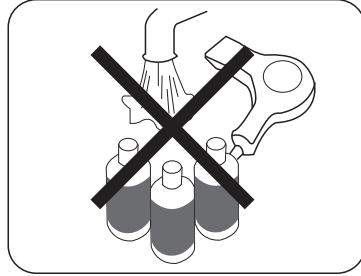

- When cleaning the massage chair, ensure the power is disconnected. Be careful to not allow wet hands to touch the plug.
- Don't clean the massage chair with water, we recommend a dry cloth to clean.
- Don't use a blow dryer to dry the sweat from the massage chair.
- Don't use the petrol, diluent, lotion or bleach.

Preservation

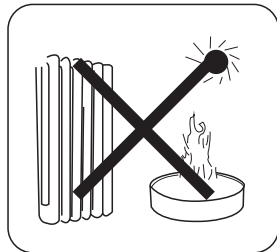

- Don't place the chair in direct sunlight for long periods, to avoid damage or change in color.
- Don't place in high temperature locations, such as heaters, stoves, to avoid damage or change in color.

Warnings


⊘	Please don't use other voltage except rated voltage to avoid fire, electric shock, and breakdown.
	Please stop use immediately if the power cord and plug is broken, or the power socket is loosen, to avoid electric shock and short circuit.
	To avoid accident and injury, people who could not express well, operate the chair is not allowed.
	Please don't allow children to use and play on chair, or sit and stand on the backrest, armrest, and legrest.
	Power cord should not be broken, reprocessed, twisted, tied and pressed by heavy objects (especially not pressed by this machine). Please don't curl the cushions in storage, it may cause short circuit, fire or electric shock.
⊘	Please do not place the machine in damp environments, it may cause electric shock, fire, breakdown, and promote rust.

Attentions

!	Please stop using if your body appears unusual or in an abnormal situation.
	Make sure the power switch and reset function is normal in case of causing accident or injury.
⊘	Please check the massage chair function before use if it stops working after a long time. When the main power is slowing down, please turn the power switch to off and unplug the power cord from the power socket to avoid possible accident occurred when the main power recover.
	If no use, unplug the power cord, put it 30cm above ground to avoid insulation aging caused by dust and moisture, which may cause electric leakage and fire.
⊘	When use the massage chair, please take off hard accessories such like watch, jewelry, etc to avoid injury.
	When massaging hip and leg, please take out hard objects from pockets to avoid accidents and injury.
	Please don't fall asleep during massage, which may cause uncomfortable body or injury
	Please do not use the machine around stove such with fire places, and don't smoke during usage, don't put it on electric blanket in case of causing fire.
	Please don't use the machine after meal in case of causing indigestion
	Please don't use the machine intoxicated in case of causing accident or injury.
	Please don't sit on the chair more than one people at the same time in case of causing accident or injury.
	Please don't use the chair without any clothes, it may cause ache by massaging skin in a strong strength.
Don't sit down or stand up when the chair legrest is lifting, and don't stamp the foot part when the footrest is not totally erect in case of causing accident or injury.	

Product installation

Controller installation method

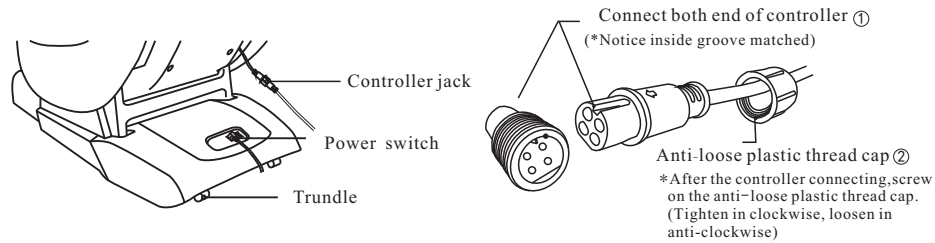


Massage chair setting method

Around the massage chair should has enough space, and lay on smooth ground

Remark: When the backrest is low and calfrest is raised, there needs more than 50cm space in front and more than 10cm space behind the chair to avoid crash.

Remark : Suggest roll out a base rug under the chair to avoid damage of the ground and carpet.


Remark: To avoid leather damage or color change, please don't put the chair to the sun or in high temperature place.

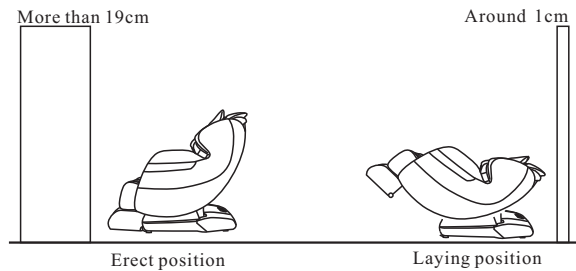

Massage chair moving method

Lift the backrest and push the whole chair by hands

- Attentions**
1. Don't move it with people or object on it, in case of turnover.
 2. Don't move the chair or drag wheels on ground that is easy to scratch.
 3. Please don't take up the seat part, leg part and shoulder part when moving the chair.
 4. Please notice its weight when lifting the calfrest.

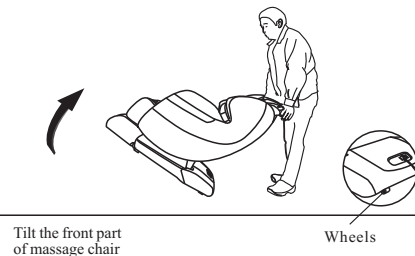

Each time before massage

Connect power


1. Before using, ensure the voltage marked on the power cord whether suit the local using condition, the power cord socket need ground wire.
2. Turn on power switch
The power switch is under the backrest.
After connecting power, the power on/off button flashing.

*Warm Tips:

We recommend to use massage chair at room temperature. If use the massage chair in cold space, please don't increase the temperature suddenly. When removing the chair from cold conditions to warm conditions, please use the chair after one hour, otherwise it will cause malfunction.


Check List

1. Check the location
 - ① Check to ensure there is no persons, pets, or objects around the massage chair
< standby position >
2. Check the massage chair state
 - ① Before using, check the backrest cover and other area's cover if there is any damage or wearout, if (no matter how small it is), please stop use, unplug the power socket and ask for maintenance.
 - ② Check the power cord, controller wire or any object is not in massage chair.
 - ③ Make sure the plug is totally in socket, and don't get metal or dust on plug in case of electric shock, short circuit and fire.

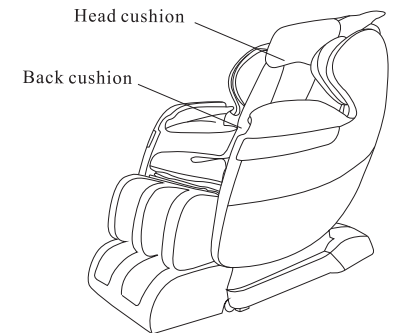

Controller instruction

Controller

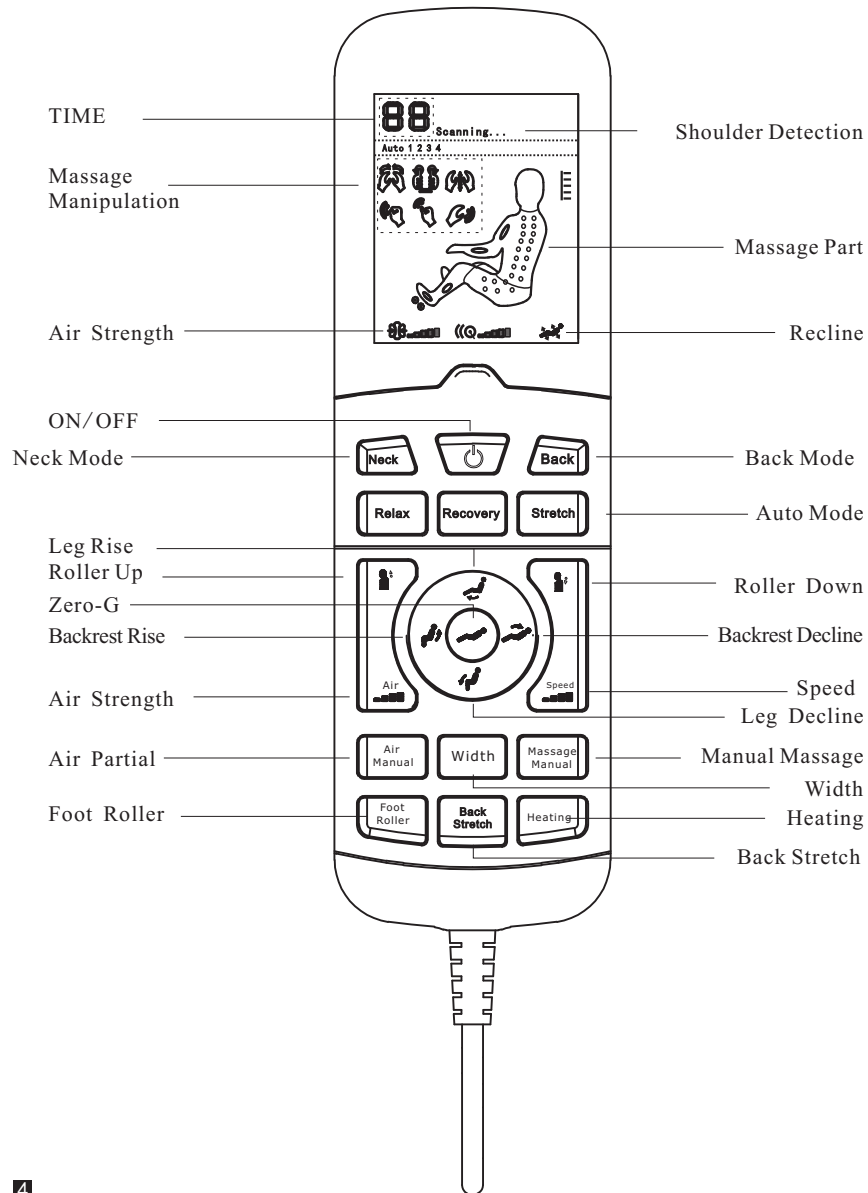

Screen display

