

Osaki
LED Therapy Dome

User's Manual Book



Product Composition



Main device



Special eye protectors
for optical instrument





Charger

Functions on the panel



Product installation

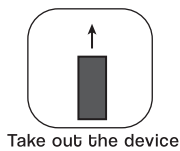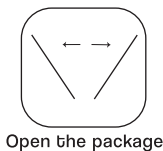

Installation

- Open the box, take the instrument out carefully, hold the bottom with two hands and unfold the instrument to form a semicircle. This product must be used on a flat leveled surface.
- Plug the DC terminal of the adapter into the DC port and make sure the connection is stable.
- Plug the AC terminal of the adapter into the power socket. If the instrument is powered on, there will be a beeping sound and the power key will flash and then proceed to standby mode.

(d) Warning

- When powered off the buttons for the functions are not visible. Once powered on the buttons will light up.
- When you press any of the function buttons you should hear a beeping sound to confirm that the command was received.

Functions



Red Light - Wave length 635nm

Red light targets wrinkles and believed to stimulate skin regeneration.



Orange Light - Wave length 610nm

Orange light helps with revitalization of the skin.



Yellow Light - Wave length 590nm

Yellow light is to address people with facial redness.



Green light - Wave length 525nm

Green light is believed to create a more consistent skin tone.



Blue Light - Wave length 464nm

Blue light is believed to help people with acne bacterial removal, acne inflammation reduction, sebaceous gland control and black & white head care.



Purple Light - Wave length 308nm

Purple light is believed to reduce scarring.

How to use

- (a) The device is designed to use with the person lying down on a flat leveled surface.
 - (b) Turn on the unit and slide your head under the arch. Eye protectors **MUST** be used at all times when using this device.
 - (c) If there is any irritation or any abnormal feeling when using the device immediately stop using and contact a medical doctor.
 - (d) This unit should not be used for more than 30 minutes per day.
- (2) Please be sure to sterilize this product frequently with a not toxic antibacterial wipe that is damp and not wet.

Note: Red light, orange light, yellow light, green light, blue light and purple light can not work at the same time. (When the pre-set time is up, the device will enter into standby mode and all function will be off.)

(d) Timer

The default lasting time for all functions are 30mins and if there two or more functions are working, the longest duration will be shown on the LED screen with counting down. When the screen shows 0, the device will quit working mode and enter into standby mode.

Safety Precautions

1. Please read the “safety precautions” carefully for correct operation before using.
2. Please comply with the relevant content of safety precautions.
3. Please keep the instruction manual well.

***Medical Attention* The following persons can NOT use this product:**

- (a) Pregnant women
- (b) People with terminal cancer or acute disease
- (c) People who are cardiopulmonary dysfunction
- (d) People with mental disability or handicap
- (e) People who are with coronary stent after bypass operation, embedded cardiac pacemaker or implanted silicone prosthesis subcutaneously should use with caution.

About the power supply

1. Power supply DC24V / 2.5A
2. Please use the provided power plug and adapter. If not it may lead to electrical shock and injury.
3. Do not damage, bend, pull or tie the power plug or wire. This may cause fire or electric shock due to damage to the plug or wire.
4. Keep the plug clean and do not let the power cable is has no cracks or exposure to the internal wires. If damaged do not use the plug and contact the manufacturer for a replacement.
5. Do not use wet hands to insert or unplug the power cable, which may cause potential electric shock.
6. For safety reasib unplug this unit from the power outlet when not in use and store away in the provided cover to keep clean and from damage.

About usage

1. Do not dismantle, repair or alter this product for any reason. This may cause electric shock or cause danger due to abnormal operation.

(Please contact our After Service department when you need to repair.)

2. In case this product is defective contact the manufacturer immediately.

※ The company reserves the right to modify the contents of this manual without prior notice.


Grounding

1. In order to prevent electric leakage and electric shock accidents, please take grounding measures.

In case of abnormal conditions of the product: in order to prevent electric leakage and electric shock accidents, please be sure to take grounding measures.

In the case of normal products without grounding measures: minor leakage of electricity can make the body feel uncomfortable.

(be sure to use a grounded power outlet and take grounding precautions.)

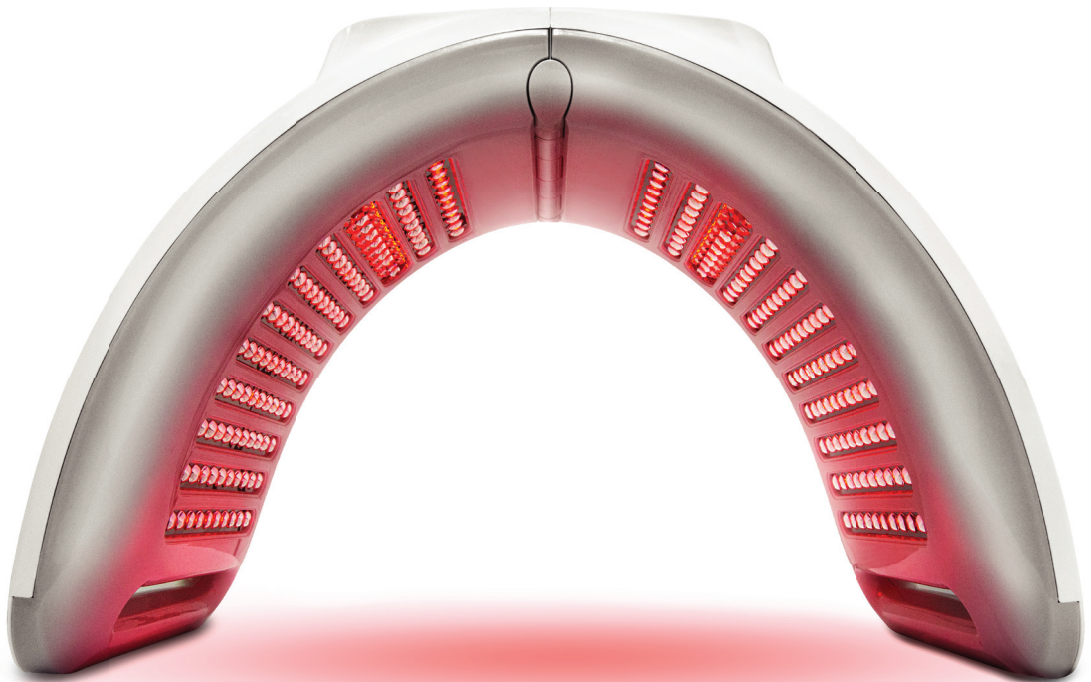

Company Information

Osaki

Contact Number: 1-888-848-2630

Website: www.titanchair.com

Osaki
LED Therapy Dome