

**AmaMedic**


# **NECK TENS MASSAGER**

**WITH HEATING THERAPY**






# Inside the Package

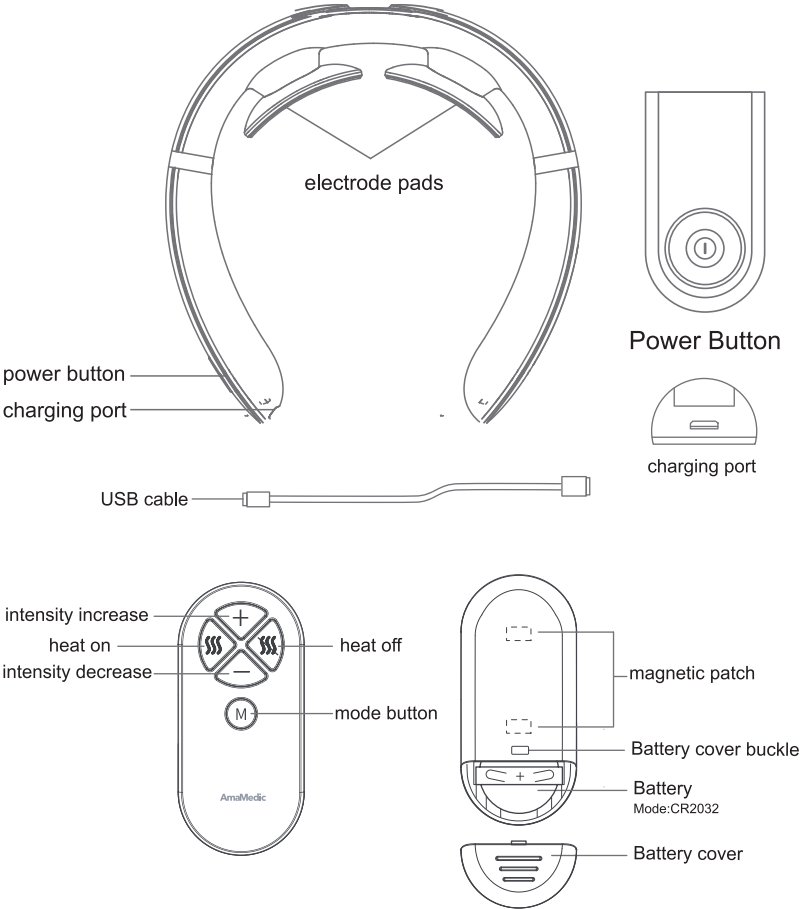

# Operation

---

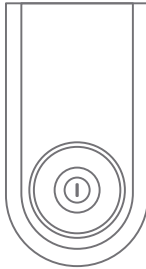

Power Button

## **Power button in the main unit**


Hold down the power button for about 2 seconds until you hear a beep which indicates that the massager has powered ON.

The heating pads will instantly start to warm up. There 15 levels of intensity. Increase the intensity by pressing the button.

To turn off the unit press and hold down the button for about 2 seconds until you hear 2 beeps.

# Buttons in the Remote Control

---

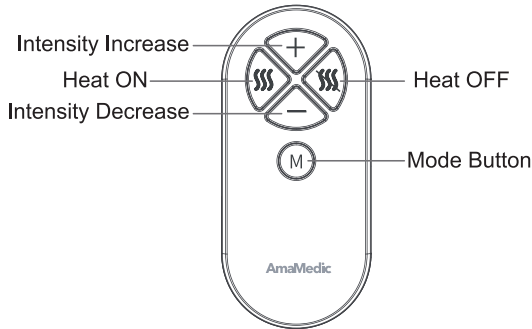

## Message Modes

There are 4 different message modes on this unit.

To change the mode simply press the mode button till reach your desired message mode/style.

## + Button

The + button when pressed will increase the intensity of the massage.

There are a total of 15 levels of intensity. The product will beep twice to indicate it has reached level 15.

## - Button

The - button decreases the massage intensity.

The product will beep twice to indicate that it has reached the lowest setting of level 1.

## Heat on button

Press to turn on the heat function and it will beep once.

The heat function is default to start when the product is turned on.

The temperature of the heat is approximately 42 °C.

## Heat off button

Press to turn off the heat function and it will beep twice to indicate it is off.

---

## Indicator light

The blue light will stay lit when the product is turn on.

The blue light will flash slowly when the product needs charging. Use Micro-USB cable to connect the product with adaptor(5V 1A) for charging.

The red light will flash slowly when the product is charging.

The red light will stop flashing and stay lit when the product is charged.

It needs about 3 hours for full charging.

**Remark:** *Remark: The product does not work when charging.*

## Timer

The product will turn off automatically after 15minutes of use.

\*If the product is on but not making any contact with the body or solid object it will automatically turn off after 2 minutes.

## Safety advise

Read all safety advice and instructions.

Non-observance of safety advice and instructions may cause electric shock, fire and/or serious injuries

Company Information:

## AMAMEDIC

1001 W. Crosby Rd.  
Carrollton TX 75006  
United States



**AmaMedic**