

Titan

Pro-Alpha
Massage Chair
OPERATION INSTRUCTIONS

SAVE THESE INSTRUCTIONS



Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay close attention to the safety warnings and keep this manual book for future use.

Note: Our Company reserves the right to revise the design and description of this product without any further notice. Please refer to the actual color of this product.

CONTENTS

Safety and maintenance	02
Name and function of components	04
Function	05
Method of usage	06
Product specification	13

Safety and maintenance

1. Important Safety Warning.



- Don't allow children to touch the movable parts of this product.
- Please use a well-grounded power supply that is suited for this product.
- Please unplug after use or before cleaning to avoid injury or damage to the product.
- Please only use the product as instructed by the manual.
- Don't use accessories that are not recommended.
- Don't use this product outdoors.
- Please read the manual instruction carefully before operation.
- Any other usage that is not listed in this manual instruction is forbidden.
- 20 minutes of usage each time is recommended.
- Do not use it if the cover or leather is damaged.
- Do not use it if the radiator hole is covered or blocked.
- Do not drop anything into this product.
- Do not fall asleep while using this product.
- Do not use if intoxicated or feel unwell.
- Do not use within one hour after having meal.
- Do not make the massage function too strong to avoid Injury.

2. Environment for usage.

- Please don't use under high temperature and moist environments such as bathrooms.
- Please don't use when the environmental temperature changes sharply.
- Please don't use under heavy-dusty or caustic environment.
- Please don't use where there is not enough space or good ventilation.

3. People who are not suitable to use this product.

- People who are suffering from osteoporosis
- People who are suffering from heart disease or carrying electronic medical equipment such as heart pacemaker.
- People who are suffering from fever.
- Pregnant women or women that are in menses period are not suitable to use this product.
- People who are injured or have skin disease.
- Children under 14 or those who have mental illnesses who require supervision
- People who are told by doctor to have rest or people who feel unwell
- People whose bodies are wet
- The appliance has a heated surface. Persons insensitive to Heat must be careful when using the appliance.

4. Safety.

- Check the voltage to see if it is suitable to the specification of this product.
- Don't use wet hands to unplug.
- Don't allow water go inside this product to avoid electric shock or cause damage to this product.

Safety and maintenance

- Do not damage the wires or change the circuit of this product.
- Do not use a wet cloth to clean the electric parts such as switch and plug.
- Be away from this product under power cut state to avoid injury in case the power come backs suddenly.
- Stop using this product while it works abnormal, and please consult local agent immediately.
- Stop using this product if you feel unwell, and please consult health care doctor.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised when near or using the appliance.
- If the supply cord is damaged, it must be replaced by the manufacturer, its Service agent or similarly qualified persons in order to avoid a hazard.


5. Maintenance.

- This product should be maintained by specialized agents, users are prohibited to disassemble or maintain by themselves.
- Please don't forget to cut the power after usage.
- Don't use this product if the socket is loose.
- If this product will be left unused for long, please curl up the wires and place the product in dry and dustless environment.
- Don't place this product under high temperature or near fire and avoid longtime exposure to direct sunlight.
- If the detachable cord / supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a Hazard.
- Please clean this product with a dry cloth. Don't use thinner, benzene or alcohol.
- The mechanical components of this product are specially designed and made, it is not necessary to be specially maintained.
- Don't use edged res to thrust this product.
- Don't roll or pull this product on uneven ground, it should be lifted before moving.
- Please use it intermittently; don't use this product for long periods of time.

6. Solution of ordinary malfunctions.

- It is normal for the motor to make sounds while in usage.
- If the controller isn't functioning properly, please check the plug and socket and see if it is properly connected and make sure the switch is turned on.
- If the rated working time is over, the product will automatically activate the power-off button; If this product works for an extended period of time, the temperature protector will make the product automatically power-off, it should be used again after half an hour's rest.

Name and function of components

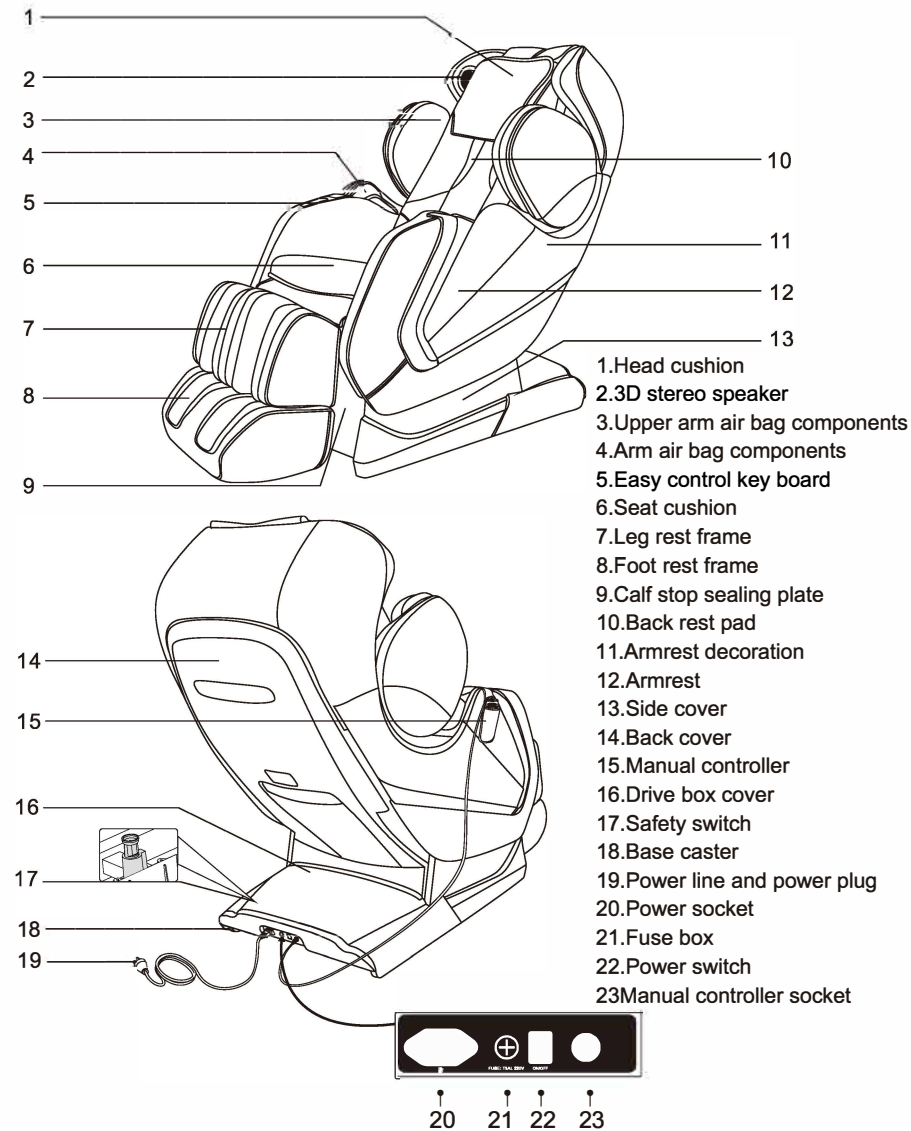

Function

Functions:

- Designed with a set of vertically movable, four wheel driven, muted intelligent massage hands.
- Designed with emulation massage methods (shiatsu, kneading, tapping, knocking, simultaneous kneading and flapping).
- Set with shoulder automatic detection and location.
- Intelligent automatic massage functions: express, wholesome, Oriental, Swedish, classic, Chiro.
- Set with memory set on function M1 and M2
- Manual massage for the upper body with three options: fixed, partial and overall massage. Five levels of speed: flapping, shiatsu and knocking. The width between the two kneading balls is adjustable with three Levels: wide, medium and narrow. On partial and fixed point, the massage hands can move up and down to adjust the massage point accurately.
- Air pressure massage function with three modes: full body air pressure, arms (upper arm set with 4 airbags, arm set with 8 airbags), legs. Set with 16 airbags with three levels of intensity.
- Leg designed with roller massage function at the bottom of foot, two pairs of papillae on toe root rotating massage; on the arch and heel mastoid skin scraping massage.
- The upper arm area has two width adjustment functions (When you adjust to the last grade, push into the arm about a gear distance then lose your hand and arm device will automatically reset, and so on) in order to adapt to the different shoulder breadth of users.
- Back heating function: carbon fiber far infrared heat source.
- The footrest use stepless extended to suit for different body heights.
- Backrest frame down slides used technology to deviate from the wall by keeping zero space, reach the maximum space saving.
- 3D stereo sound speaker installed into armrest.
- Fast-pair blue tooth to play music.
- Easy control key board on right armrest for on/off, decline/recline.
- Leg lifting adjustments, backrest frame lifting adjustments and automatic lay down mode adjustment function.
- L track that fits the curves of the body, can realize the zero gravity when back down.


Method of usage



Note: The specific displays refers to the actual object of this product

Method of usage


1. Switch power supply on

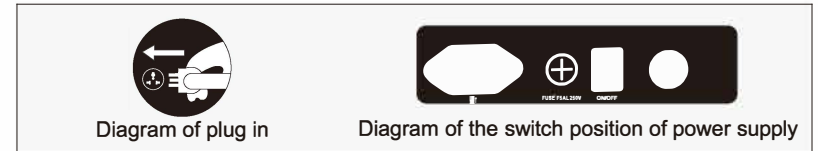

2. Start to massage

- 1). Press power button, the massage chair will automatically recline and start the massage.
- 2). Wait for the automatic shoulder position and massage acupressure point detection.
- 3). If the shoulder massage balls do not align with the shoulders, you can adjust the up and down button to a suitable position. The shoulder position has 10 adjustable grades. It will run the automatic massage if there is no adjustment on the shoulders within 10 seconds.

Button	Description	Display
	Move the massage balls upward to a suitable position	
	Move the massage balls downward to a suitable position	

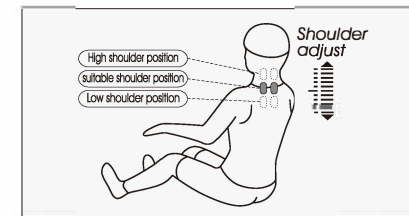

4) Choose massage function to massage

- Memory function massage
Press the "ORIENTAL/ M1 " or "CHIRO/ M2 " button for 3 seconds, then enter into memory massage function "M1" or "M2" ,will revert to the previous memory massage mode.








Function	Button	Display
Memory massage		Menory
		Menory

Method of usage


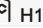
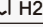
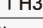


- Characterized auto massage

Function	Button	Description	Display
Characterized auto massage	EXPRESS	Quick fix massage for those who are in a rush	Auto F0
	WHOLESOME	A full body massage with rolling massage on the neck area.	Auto F1
	ORIENTAL/M1	Rolling massage on the neck with fast moving massage on the back.	Auto F2
	SWEDISH	Kneading massage with scrubbing.	Auto F3
	CLASSIC	Decline and recline the back rest and foot rest to enforce lying automatically. Is an auto-cycle mode, to stretch out the waist, meanwhile, the foot airbag, back mechanism massage waist.	Auto F4
	CHIROM2	Emphasizes on massage along the spinal area plus stretching action.	Auto F5




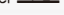
- Upper body manual massage, lower body manual massage
- 1) Upper body manual massage

Function	Button	Adjust Massage Ball	Display
Mode: One button cycle the five manual massage modes, kneading, flap, flap & kneading, shiatsu, tap, and off.	MODE	Knead & flap: Speed	
		Flap: Speed And Width	
		Knead: Speed	
		Shiatsu: Speed And Width	
		Knock: Speed And Width	
Speed	SPEED	Three Grades Speed Adjustable	
Width	WIDTH	Three Grades Widths Adjustable	



Method of usage

Function	Button	Description	Display
Back Stretch	MASSAGE AREA	Fixed position (fixed point) massage back and forth	 H0
	MASSAGE AREA	Small area (partial area) massage back and forth	 H1
	MASSAGE AREA	Overall massage back and forth	 H2
	MASSAGE AREA	Upper body massage back and forth	 H3
	MASSAGE AREA	Lower body massage back and forth	 H4
		Massage point adjust (in fixed point, partial modes you can adjust the massage position)	Massage hand move up
		Massage hand move down	





- Air pressure function massage

Function	Button	Description	Display
Air pressure mode	AIR MODE	Three air pressure modes(full body, arm and leg) for Lower body	Arm  Feet 
Intensity	AIR POWER	With three grades adjustable	Arm  Feet 

- Other manual function






Function	Button	Description	Display
Roller	FOOT ROLLER	Sole roller: sole roller switch, with 3 speed levels available	
Heating	HEAT	Far infrared back heating, after 3 minutes should start heat therapy	

- Other manual function

Function	Button	Description
Calf rise and drop		Calf rest lifting button, press this button, calf restlifting slowly, stop when release it
		Calf rest decline button, press this button, calf rest decline slowly, stop when release it
Backrest rise and drop		Backrest lifting button, press this button to lift the backrest, release to stop
		Backrest decline button, press this button to decline the backrest, release to stop
Auto lay down mode	ZERO	with three modes circulation function, press the button to turn into auto lay downmode.


Method of usage

Shortcut keys

Button	Description
	Pause and run button: to stop massage, press on/off button, if it stops more than 20 minutes, the massage chair will shut down all the massage functions
	Calf rest lifting button, press this button to lift calf rest, release to stop
	Calf rest decline button, press this button to decline calf rest, release to stop
	Backrest lifting button, press this button to lift the backrest, release to stop
	Backrest decline button, press this button to decline the backrest, release to stop

3. Turn off power, stop massage


- 1) Press power button during massage progress, all massage function will stop immediately, the backrest and footrest will return to the original position automatically, meanwhile you can set up a memory function. After rated working time, all massage functions will be finished, meanwhile you can set up a memory function, but backrest and footrest won't return to original position.
- 2) Memory store function: After rated working time or in off mode, the hand controller displays "Memory" and glitter, remind the users whether need to store the massage mode before turn off, you can press the memory function buttons "M1" or "M2" as you like, store the massage mode.
- 3) Cut off the power of the whole machine.
Figure (Switch power off)

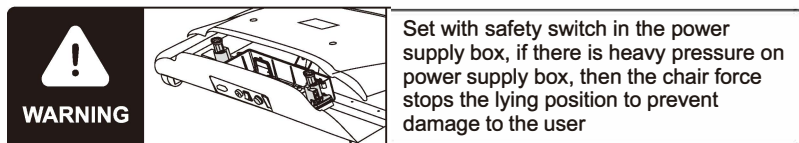


Method of usage

4. Calf installation

- Install the rotating shaft in calf, fixed axis first, make sure the hole is aligned, so that the R shape pin assembly (chart 1)
- After calf is connected with the fixed axis aligned, then insert (chart 2)

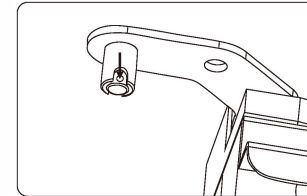


(chart 1)

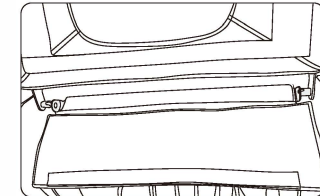


(chart 2)

- Plug R pin, then move to the other side of the calf, and the plug R pin (chart 3)
- Finished (chart 4)

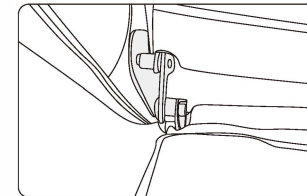


(chart3)

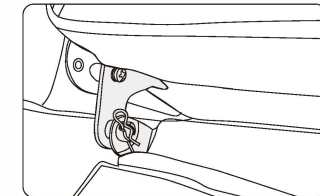


(chart4)

Note

- When the machine under the working state, with lying and calf rest lifting, you could not apply gravity on the calf, to prevent the massage chair does harm to people before dumping(chart1).
- Ensure the mechanism hands not in the seat part ,otherwise will cause uncomfortable in the Buttock when you seating in the chair(chart2).

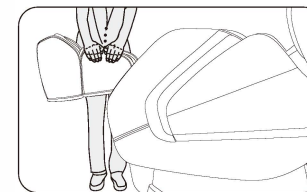

(chart1)

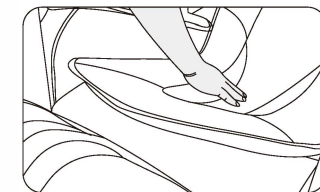


(chart2)

Method of usage


5. Head cushions and backrest pad's dismounting and use manual

- The head cushion has three fold-down layers. The users can choose the number of layers to find the best intensity of back massage.
- Separate the head pillow to reduce the intensity of the back massage.
- Backrest cushion and backrest connected by zipper (1), head cushion and backrest pad connected by thread gluing (2)

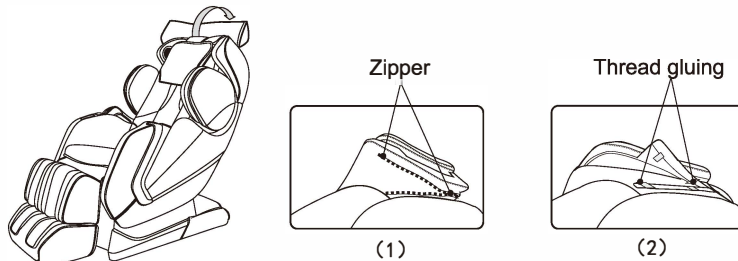

6. Floor Protection

- It is recommended to not place the massage chair directly on wooden floors to prevent damage to the floor, please place a carpet or rug under the chair.

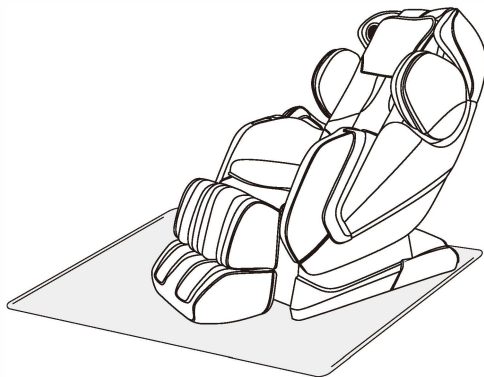


Caution: Do not move the chairs with its wheels on the wood floor, rough ground or in a narrow space. Two people are required to lift the chair. While lifting, take the calf sealing board and backrest cover grooves as the focal point(Please do not lift armrest Directly.)

Method of usage

7. Method of movement

- Make sure that all wires are far above the ground, tilt the backrest backwards (with the gravity center resting on the wheels), push the chair forward or backward with your hands and bring back the chair to the normal position in a slow and gentle manner.

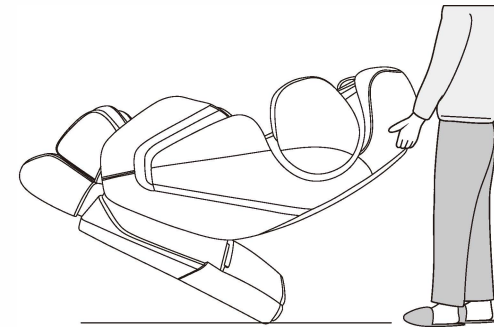

Note: During movement you must turn off the power switch first and then unplug the power line and hand controller cord.

Product specification

Model: Pro-Alpha

Name: Massage Chair

Rated voltage: 110-120V~ 60Hz

220-240V~ 50Hz/60Hz

Rated power input:150W

Rated time: 20Min

Safety structure: Class I

Material introduction: PVC, PA, steel parts and electric &ele-ctron parts

MATERIAL	PU	LEATHER	CLOTH	WOOD