
MOER

HOME, SMART MOES HOME

Instruction Manual

Portable TENS+Infrared Heating Massager

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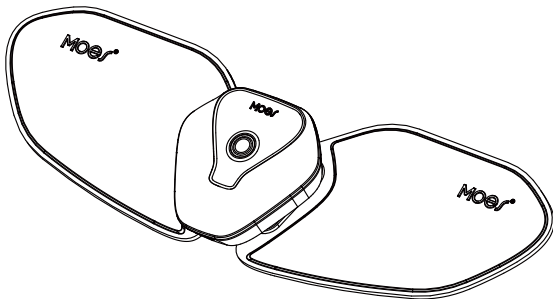
Device Description

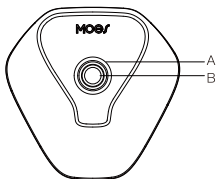A: Host switch button/The aperture indicates the charging status

B: Charge/Low battery indicator

C: Temperature gear button

D: Massage intensity increase button

E: Shoulder massage indicator light

F: Remote control switch button

G: Massage intensity digital display

H: Belly massage indicator light

I: Back massage indicator light


J: Temperature gear button

K: Waist massage indicator light

L: Massage intensity decrease button

M: Leg massage indicator light

N: Charging indicator light

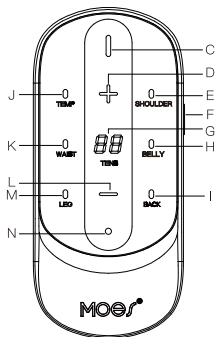Instructions

1. Installation

- ①. Remove the blue protective film of the gel and attach the gel pads to the thin electrical coating on the back of the product.
 - * Please wipe the coating carefully before use, so as not to affect the conductivity due to lack of adhesion caused by human sweat and grease.
 - * Please make the gel pads and electrical coating fit tightly without any bubbles. If there are bubbles please use your finger press the gel pads which are separated by transparent film.
- ②. Correctly identify the left and right directions of wings spread and install the host (Make sure the magnet is aligned).

2.Power on

Remote control pairing: Long press the power button for 6S when the host and remote control are both off respectively. The white light of the host and the switch light of the remote control flashes. When the white light of the host is stay on and the switch light of the remote control switches to white and stay on the pairing successfully.

Normal power on: Long press the power button for 3 seconds when the host and remote control are both off respectively. The white light of the host is stay on, the digital display, the massage part light, the increase and decrease light, the switch light of the remote control is stay on.(After the remote control pairing is successfully connected).

※ Remote control pairing is required when the product is used for the first time. The host can be automatically paired with the remote control after normal power on.

3.Attaching the gel

Remove the white protective film of gel, stick the product onto where the massage is required.

4.Start massage

Press the+ button of remote control, the digital display screen will change from 00 to 01, the host starts massage state.

5.Intensity

Press the "+/-" button to adjust the intensity, choose the most proper intensity you want. The intensity is adjustable from 0 to 15.

※ Long press the "-" button to switch the massage output.

6.Mode

Press the remote control switch button to switch different body parts' massage actions .There are 5 massage modes for different body parts.

7.Heating function

Press the temperature gear button to adjust the heating temperature. The indicator turns green for low temperature, yellow for middle temperature and red for high temperature. The remote control enters the sleep state in one minute when it is under working status, with the light and the digital display off, and the status light displayed normally. It can be awakened by pressing the side key.

※ Preheat 8 seconds after sticking on the body for a better experience.

8. Child lock function for remote control

The remote control without any operation will automatically enter the child lock function after 30s. When you need to use the remote control, you can press the switch button to wake it up, or press the switch button to lock the remote control screen after using.

9. Power off

The device will be automatically closed after 15 minutes, or long press the remote control power button for 3s to power off.

10. Storage

For convenience of next use, please stick the white protective film back to the gel and store it in a cool and dry place.

※ Cover the back of the product completely with a white protective film.

Maintenance

1. In order to prolong the service life of the built-in battery, please store the massager after full charge when not in use for a long time, and charge it once every 2-3 months.
2. If dust or sweat appears on the surface of the massager, please gently wipe with wet wipes (no water should flow out of the wipes).
3. When the gel pads are difficult to attach due to perspiration, please follow the steps below to maintain:
 - ① Please use your finger get a water droplet and put it apply evenly onto the gel or gently wipe off the dust;
 - ② After thoroughly wipe the water, put the gel upward and place it in a good ventilation and shade place drying;
 - ③ Air dry until there is no moisture on the surface, stick the gel pads on the white protective film and store it in a cool and dry place.※ If the stickiness of the gel can not recover, please change a new one.
4. The method of getting the gel off
 - ① Pinch one and of the gel pads with your fingers and tear it off (It is recommended to roll it to tear off);

- ② Please completely remove the residual gel pads on the conductive coating to prevent the bad impact of incomplete clear.(It will be easier to remove by taking use of the torn hydrogel pad to stick the residual gel with a small amount.)
 - ③ Please make sure the machine is power off while tearing the gel pads off, otherwise you may get pained or scalded.
 - ④ After the gel is teared off, please wipe the conductive coating carefully to avoid insufficient adhesion or conductivity due to greasy sweat or body sweat or grease.
5. Please do not use or store this product in the following places:

● Direct sunlight

● Near fireworks

● Hot and humid place

● Dusty area

● A place with water

● Strong electromagnetic site

Suggested massage parts

Shoulder

Back

Belly

Waist


Leg

Do not massage parts

Head

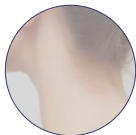Neck

Chest

Sensitive parts

Simple Troubleshooting

Problem	Possible Cause	Solution
The device does not switch on when the ON / OFF button is pressed.	The battery may be flat.	Fully charge the battery . If it is still not possible to switch the device on ,please contact the customer service.
The device does not stick to the body.	The gel pads may be dirty or worn :The treatment area of the skin may not be clean enough.	1.Please clean the adhesive surface of the gel pads using a damp lint-free cloth. The gel pads must be replaced if they still do not stick securely. 2.Clean the skin prior to any application do not use skin care lotions or oil spriror to treatment .
There is no noticeable stimulation.	The gel pads may be dirty or worn;The battery may need charging.	1.Ensure that the gel pads are in firm contact with the treatment area. 2.The battery is almost flat. Charging it fully.
An unpleasant sensation is experienced where the gel pads are applied	The gel pads are worn. This may cause irritated skin , as even distribution of the current across	Replace the gel pads.
Skin in the treatment area turns red	If redness under the gel pads disappears quickly , this is not dangerous , and is caused by the locally stimulated , increased blood flow caused by TENS and heat.	Immediately stop treatment and wait until your skin has returned to its normal condition. However, consult your doctor before you continue treatment if the skin irritation persists and if it is accompanied by an itchiness or inflammation, as this may be caused by an allergic reaction to the adhesive surface.

<p>The skin beneath the surface of the heat turns red.</p>	<p>The reddening of the skin under the heated surface of the massager may be due to intense stimulation of blood flow in the heated parts.</p>	<p>The skin should return to its normal condition a few minutes after treatment. However, consult your doctor before continuing to use the heat function if the redness persists and if you experience any itchness or inflammation.</p>
<p>The massager becomes too hot.</p>		<p>Press the On / Off button to stop the heating function, and just continue using the TENS function.</p>

Warnings And Safety Notes

Prohibited users

1. People has implanted electrical devices(e.g.pacemaker)
2. Patients known with severe arrhythmia, angina pectoris and heart failure.
3. Pregnant women and children under 12 years old.
4. All kinds of pets.

People with caution as advised by doctor

1. Patients who are receiving treatment or taking medication for a serious medical condition.
2. Patients with fever.
3. Patients with acute inflammation, suppurative inflammation and bleeding tendencies.
4. Patients with blood abnormalities, infectious diseases or skin diseases.
5. People with allergic skin.

Notes on use

1. Keep out of reach of children.
2. The elderly and those who cannot take care of themselves should be accompanied when using.
3. Do not put the product or gel pads in the mouth.

4. Please do not use this product during strenuous exercise.
5. Please keep your skin dry when using.
6. Do not stick on the wound.
7. Please avoid metal objects, such as necklaces, bracelets, etc.when using.
8. If you feel uncomfortable during use, please stop using.
9. Do not use continuously for too long (no more than 30 minutes is recommended) to prevent muscle fatigue, skin tingling or redness.

Included in package

Host x1	Remote control x1	Left wing x1	Right wing x1
User manual x1	USB charging cable x1	Gel pads x8	Cleaning wipe x1

※ The package of this product does not contain power adapters. The power adapters applicable to this product are: power adapters of various mobile phones and USB interfaces of computers. Using other types of power adapters may damage the product.

Technical specification

Material: Silicone +ABS+PC	Size:274mm*75mm*22mm	Weight:150g
Rated voltage: 5V	Rated power: 7.5W	Battery capacity: 500mA
Charging voltage: DC5V/3A	Battery charging time: ≈2h	

※ The product does not support external mobile power supply.

SERVICE

Thank you for your trust and support to our products, we will provide you with a two-year worry-free after-sales service (freight is not included), please do not alter this warranty service card, to safeguard your legitimate rights and interests. If you need service or have any questions, please consult the distributor or contact us.

Product quality problems occur within 24 months from the date of receipt, please prepare the product and the packaging, applying for after-sales maintenance in the site or store where you purchase; If the product is damaged due to personal reasons, a certain amount of maintenance fee shall be charged for repair.

We have the right to refuse to provide warranty service if:

1. Products with damaged appearance, missing LOGO or beyond the service term
2. Products that are disassembled, injured, privately repaired, modified or have missing parts
3. The circuit is burned or the data cable or power interface is damaged
4. Products damaged by foreign matter intrusion
(including but not limited to various forms of fluid, sand, dust, soot, etc.)

RECYCLING INFORMATION

All products marked with the symbol for separate collection of waste electrical and electronic equipment (WEEE Directive 2012/19 / EU) must be disposed of separately from unsorted municipal waste. To protect your health and the environment, this equipment must be disposed of at designated collection points for electrical and electronic equipment designated by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences for the environment and human health. To find out where these collection points are and how they work, contact the installer or your local authority.

WARRANTY CARD

Product Information

Product Name_____

Product Type_____

Purchase Date_____

Warranty Period_____

Dealer Information_____

Customer's Name_____

Customer Phone_____

Customer Address_____

Maintenance Records

Failure date	Cause Of Issue	Fault Content	Principal

Thank you for your support and purchase at we Moes, we are always here for your complete satisfaction, just feel free to share your great shopping experience with us.

If you have any other need, please do not hesitate to contact us first, we will try to meet your demand.

FOLLOW US

@moes_smart

MOES.Official

@moes_smart

www.moes.net

@moessmart

@moes_smart

AMZLAB GmbH
Laubenhof 23, 45326 Essen

EVATOST CONSULTING LTD
Address: Suite 11, First Floor, Moy Road
Business Centre, Taffs Well, Cardiff,
Wales, CF15 7QR
Tel:+442921680945
Email:contact@evatmaster.com

WENZHOU NOVA NEW ENERGY CO.,LTD
Address: Power Science and Technology
Innovation Center, NO.238, Wei 11 Road,
Yueqing Economic Development Zone,
Yueqing, Zhejiang, China
Tel:+86-577-57186815
Email:service@moeshouse.com
Made In China

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