# glad skin


Inflammatory Skin Expert

## ACNEFEKT Gel

STAPHEFEKT™  
100  
Sa  
S. aureus killer

**REDUCES ACNE SYMPTOMS**  
ANTI-PUSTULES  
ANTI-COMEDONES  
REDUCES SKIN INFLAMMATION

Salicylic acid-free | Alcohol-free  
**Acne skin**

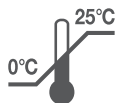Micreos Human Health B.V.  
Bankplein 2  
2585 EV Den Haag  
The Netherlands

Gladskin BV  
Bankplein 2  
2585 EV Den Haag  
The Netherlands

# UK CA

UKRP:  
Advena Ltd  
Pure Offices, Plato Close,  
Tachbrook Park, Warwick,  
CV34 6WE  
United Kingdom

INS-GA-ENNLDEFERSE-V7.2  
REVISION DATE MAY-2023

### ENGLISH

Gladskin Acnefekt Gel is a medical device. Please read this information carefully before using the product, even if you have used it before. That way, you will keep up-to-date on the effects of this product and how to use it properly as we may sometimes need to change the information. Please do not hesitate to contact us if you have any questions after reading this information.

### FOR WHOM IS GLADSKIN ACNEFEKT GEL INTENDED?

Use Gladskin Acnefekt Gel if you suffer from acne. The symptoms of acne are pimples on your face, neck, chest and back, and may also include minor inflammations and red, easily irritated skin. *Staphylococcus aureus* bacteria can cause inflammation of the skin and worsen the acne symptoms.

### TOLERANCE

This formula is fragrance-free and hypoallergenic. It is suitable for face and body.

### HOW DOES THIS PRODUCT WORK?

Gladskin Acnefekt Gel contains the main ingredient Staphefekt™ SA.100. Stemming from Microbiome Science, this unique ingredient has the ability to specifically target the main bacteria responsible for inflammation of the skin, while preserving the other bacteria essential for skin health.

### WHAT DOES THIS PRODUCT DO?

- Reduces acne symptoms
- Anti-pustules
- Anti-comedones
- Reduces skin inflammation

### HOW TO USE THIS PRODUCT

 > For external use. Use the product on intact skin that is free of creams, cleansers, powders and so on.

 > Use regularly two to three times a day.

> Apply a thin layer of the product on and around the areas with acne.

> Let the gel dry for at least 5-10 minutes before you apply any other products on your skin.

When using the Gladskin Acnefekt Gel you may notice a thin film left on your skin. This is remaining gel which can be wiped off after 10 minutes.

### WHEN NOT TO USE THIS PRODUCT

Do not use this product if your skin is not intact or if you are known to be allergic to any of the ingredients.

The product should be a clear gel, free from visible particles. If this is not the case, do not use the product and contact us.

### WHAT CAN YOU EXPECT?

Gladskin Acnefekt Gel can quickly reduce acne symptoms. In some cases it can take a few weeks to see results. Most acne is related to an imbalance in the bacterial composition on the skin (microbiome). If your acne is not due to this imbalance you may not experience an improvement. Some users may experience redness, a burning sensation, swelling or itch after using Gladskin Acnefekt Gel.

If you use Gladskin Acnefekt Gel for the first time your acne symptoms may temporarily increase. This should improve within 10 days.

If your acne does not improve or side-effects are severe or persist, please stop using Gladskin Acnefekt Gel and consult your doctor.

**In case of serious side-effects, please inform us and the national competent authority of the country where you reside.**

### INGREDIENTS

Aqua, Propylene glycol, Hydroxypropyl methylcellulose, Glycerin, Arginine HCl, Sodium chloride, Trometamol, Staphefekt SA.100, Calcium chloride.

Gladskin Acnefekt Gel is hypoallergenic and dermatologically tested. It is a medical device to be used on intact skin. Gladskin Acnefekt Gel does not contain any lanolin, acids, alcohol, perfume or parabens and is non-comedogenic (i.e., will not block your pores).

### SHELF LIFE AND STORAGE

The expiration date of this product is indicated on the dispenser. For an optimal benefit of the main ingredient in Gladskin Acnefekt Gel we recommend you store the product in the refrigerator. Gladskin Acnefekt Gel should not be stored at temperatures below 0°C or above 25°C.

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