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
Paraffinbad ClassicSmooth

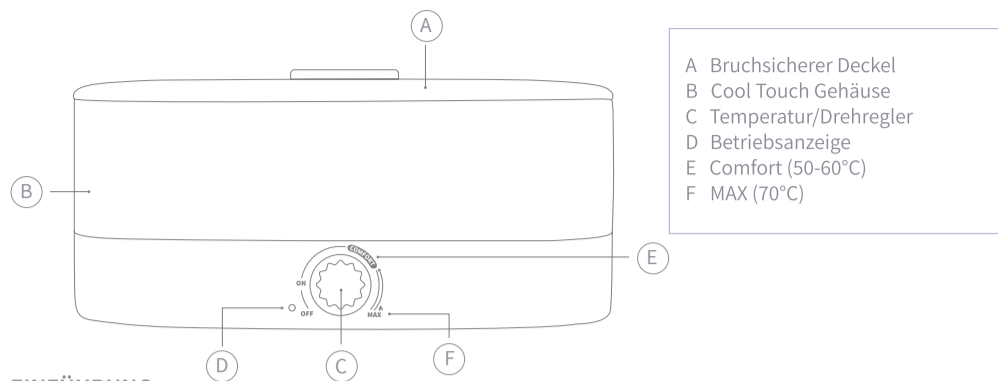EINFÜHRUNG

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
4

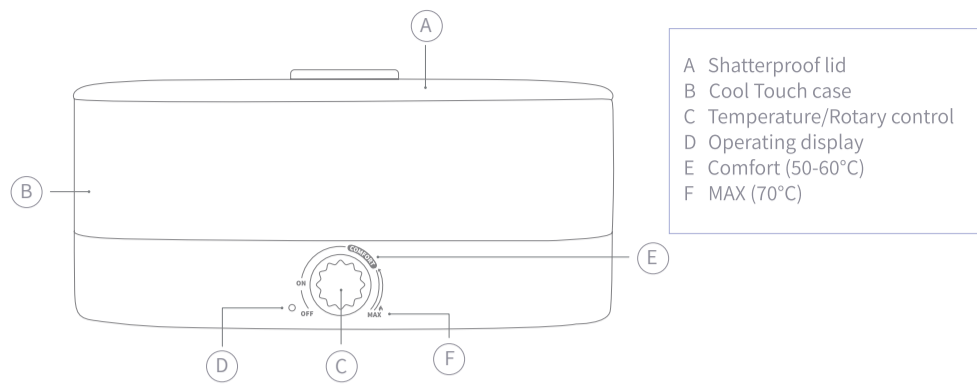INTRODUCTION

We are pleased that you have chosen the WAXSENSATION ClassicSmooth Paraffin Bath, a generous and allergy-friendly device that provides ample space for your hands, feet, and other body parts while pampering and caring for your skin. For skin that feels as good as it looks! This user manual provides you with all the essential information for safe and efficient use of the device. Please take the time to read the instructions thoroughly to get the most out of your new product and to always stay safe.

TECHNICAL DATA

Typ:	Paraffin bath ClassicSmooth
Connection:	220V-240V ~ 50/60Hz
Rated power:	200W
Dimensions:	LxWxH – approx. 37,0 x 22,0 x 16,5cm
Weight:	1.44 kg
Capacity:	3500ml
Temperature range:	30 to 70°C

PACKAGE CONTENTS

- Paraffin bath
- Paraffin wax
- Heat protection grid
- Cotton gloves
- Disposable gloves
- Cleaning brush and spatula

OPERATION

Before starting the device, we recommend you clean the interior carefully with a damp cloth. Also, please remove any protective film from the device. Connect the device to the power supply with the power cord.

- Place the paraffin wax blocks in the designated container. Turn the control knob clockwise to the „MAX“ position. The illuminated operation indicator signifies that the device is now active (For three 450g blocks of paraffin wax, a melting time of about 90 minutes is required. At the MAX temperature, the wax heats up to ~ 70°C).
- When the wax is fully melted, turn the control knob to the „Comfort“ area. Allow the wax to cool down a bit. In the „Comfort“ area, the device has a temperature range of ~ 50°C - 60°C.
- The temperature control can be turned to increase or decrease the temperature.
- When the set temperature is reached, the operation indicator goes out.
- Once the temperature is set, the device operates constantly in this temperature range. If the internal temperature of the container falls below the selected temperature, the device automatically starts the heating process and the operation indicator turns on.
- To switch off the device, turn the knob counterclockwise to the „OFF“ position.

1

SAFETY INSTRUCTIONS

Please follow all instructions carefully. Improper use of the device, which contradicts the safety instructions and operating instructions given in this manual, can lead to injuries or damage to people and property. The manufacturer or seller cannot be held liable for such damages that are due to non-compliance with these instructions.

Always place the heat protection grid in the device before use.

Avoid contact with the bottom of the wax container as it can get very hot during operation.

Avoid using the device at temperatures above 60°C to prevent burns.

It is important that you test the temperature of the wax before dipping your hands or feet into it to ensure that it is not too hot.

Once the device is connected to the power source, it should never be left unattended. Make sure that children do not approach the device without supervision.

Use the device only for the purpose described in this manual.

Only for external use. Never use with sick, injured or inflamed skin.

The heated wax must not be applied near eyes, mucous membranes, rashes, growths, injured or heat-sensitive skin.

Avoid submerging the device in water or other liquids. In case the device falls into the water, do not touch either the water or the device. Unplug the power cord immediately.

Keep the cable away from heated surfaces.

Always unplug the power cord after use.

DISPOSAL

Power tools, accessories and packaging should be recycled in an environmentally sound manner.

FOR EU COUNTRIES ONLY:

Do not dispose of power tools in household waste! According to the European Directive 2002/96/EC on waste electrical and electronic equipment (WEEE) and its implementation in national law, power tools that are no longer usable must be collected separately and recycled in an environmentally sound manner.

Within the EU, this symbol indicates that this product must not be disposed of with household waste. Old devices contain valuable recyclable materials which should be recycled and in order not to harm the environment or human health by uncontrolled waste disposal. Therefore, please dispose of old equipment through appropriate collection systems.

3

APPLICATION

- Remove any jewelry and thoroughly clean the skin in the areas to be treated before application.
- Always check the wax temperature before dipping. If it feels too hot, adjust the temperature down. Once the paraffin has cooled to a suitable temperature, you can start the application.
- Dip your hands, feet, or elbows slowly into the container of the device until they are submerged under the wax surface and covered with wax. Then slowly take them out. Re-dip them after a few minutes. Repeat this process up to five times, until several layers of wax have formed on your hand.
- Now put on disposable gloves and cotton gloves over them to keep the heat and intensify the effect.
- After about 15 minutes of exposure time - the exact time results from the layers of wax and the stored temperature - remove the gloves and dried wax.
- Take a moment now to gently remove any wax from your hands or elbows. You have now successfully completed the treatment! Now you can cream your skin and enjoy the pleasant feeling of softness and suppleness.

IN CASE OF MALFUNCTION

- Check the correct power supply.
- Unplug the power cord from the socket. Plug the power cord back into the socket. The device is now ready for operation again.
- If the problems persist, we recommend that you contact your dealer or service center. There you will receive competent support in troubleshooting and possible repairs.

CLEANING

To guarantee hygiene, it is recommended to disinfect the device externally after use. We recommend disinfecting the device only by wiping.

- We recommend a thorough cleaning after approximately 20 applications.
- Unplug the device and let the wax harden for 12 hours.
- Turn on the device for about 3 minutes before cleaning, then turn it off and unplug it.
- Remove the wax block by gently turning the device upside down so that the wax block cleanly falls out of the device. Remove any residue with the included silicone spatula. Dispose of the wax.
- Never use sharp or pointed objects to remove the wax!
- Wipe the inside of the device with damp paper towels or a cloth.
- Under no circumstances use scouring agents or sharp cleaning agents!

2

LEGAL WARRANTY

This device has been manufactured with all due care and thoroughly tested. The dealer grants the legal warranty of 2 years from the date of delivery according to the following conditions. In principle, there is no liability for age-related defects and wear parts depending on the typical product shelf life.

The legal warranty begins with the delivery of the goods. Warranty claims must be made within the 2 years. Claims occurring after the expiration of the legal warranty cannot be considered. A so-called legal warranty claim occurs if the customer is delivered defective (e.g. manufacturing defects, transport damage), wrong goods or too few goods. In such a warranty case, the buyer is obliged to make the purchased item available to the seller for inspection.

The consumer can choose between the elimination of the defect, e.g. by repair (so-called rectification) or the delivery of defect-free goods (so-called replacement). However, this is only possible if this choice does not cause disproportionate costs.

If rectification (e.g. repair) or replacement delivery is not possible or fails, the customer has the following rights: Withdrawal from the purchase contract, reduce the purchase price, damages or reimbursement of expenses.

Reimbursable costs:

In the event of a legal warranty claim, the seller must bear all costs of subsequent performance, e.g. repair costs, installation and removal costs, transport costs (outward and return), labor costs, material costs, travel costs, postage, telephone costs and storage costs. If the buyer reports a defect although there is no defect at all, he may be liable for damages.

If it comes after a defect examination to the return, the costs for it the buyer carries, if no material defect could be confirmed. Nevertheless, the dealer is obliged to return the goods.

For a smooth complaint processing the following information is essential:

- Original invoice / receipt
- Written description of the defect and photos of the detected defect.

4