

DELUXE SELF-PROPEL COMMODORE



WARNING

- DNR Wheels Pte Ltd shall not be held liable for any damages, injuries or accidents arising from misusing, mishandling and negligence.
- DNR personnel shall demonstrate the basic handling of the commode.
- The safety features for operating the commode must be strictly observed.
- Any defects must be reported immediately to DNR Wheels Pte Ltd.
- When transporting the commode, brakes should be in an unlocked position. The castors are designed for use in an indoor environment and on even surface.
- Castors may be damaged when moving the commode along any rough & uneven surface.

SAFETY INSTRUCTION

- Commode should not be use if found damaged as it may injure the user or caregiver.
- Caregiver should be familiar with the functionality of the commode.
- Commode must not be overloaded (Maximum load capacity: 100kg).
- Ensure all the castors and wheels are locked and leg rest is drawn back before transferring the user in/ out of the equipment.

- If the user is on the commode, castors and wheels must be locked. An unlocked commode poses a risk of falling to the user.
- While using the commode, footrest should be drawn out and user's feet are to be well rested on footrest.
- Bucket can be removed by sliding out from rear end of the commode seat.
- Users with limited physical strength should be supervised or assisted while using commode.
- Do not leave user with poor stability unattended on the commode.

CLEANING

- Use a cloth to wipe the commode dry after use.
- Dry out in the sun if possible.
- Clean out debris from the castors regularly for easy manoeuvring.


Technical Specifications

Seat Width:	43cm (17")	Overall Width:	61.5cm (24.2")
Seat Depth:	43cm (17")	Overall Depth:	102cm (40.2") (without footrests) 109cm (43") (with footrests)
Armrest to Armrest Width:	45cm (17.7")	Overall Height:	94cm (37")
Seat to Armrest Height:	19cm (7.5")	Seat Cushion Opening:	Front Opening: 8cm (3.1") (W) Middle Opening: 18cm (7") (W) Depth: 36cm (14.2")
Seat to Backrest Height:	38.5cm (15.2")	Front Castors:	12.7x3cm (5"x1.18")
Floor to Seat Height:	52cm (20.5")	Rear Wheels:	24"
Footrest to Seat Height:	40 – 50cm (15.7"-19.7")	Material:	Aluminium Grade 6000
Width between Rear Castors:	34cm (13.4")	Max Weight Capacity:	100kg
Floor to Bucket Holder:	42cm (16.5")	Product Weight:	20kg


Flip up the armrest



Flip up the footplate

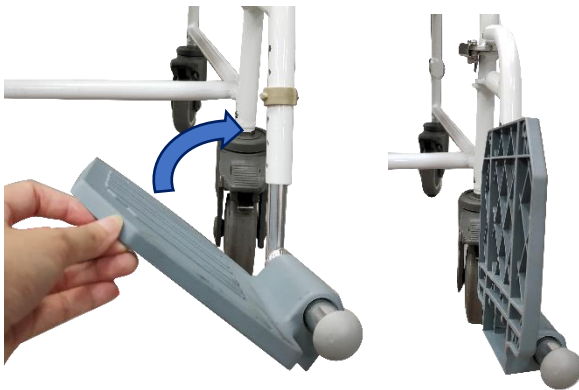

Adjust the height of the footplate

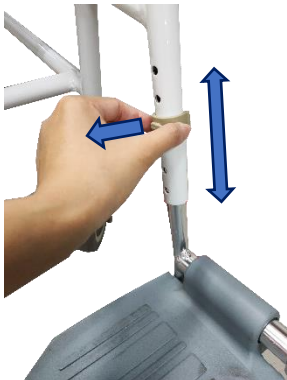

1. Pull out the pin.
2. Adjust the footplate to desired height.

Detaching the footrest



Put on and release the brake



Push forward to put on the brake



Pull backward to release the brake

Lock and unlock the castor



Castor is locked



Castor is unlocked