

Nu Skin Devices

Beautiful skin makes you feel unstoppable. With Nu Skin's at-home skincare devices, you have the freedom to build the routine that helps you get it.



ageLOC LumiSpa & Activating Cleanser



ageLOC Boost & Activating Serum



ageLOC Galvanic Spa & Facial Gels



ageLOC Galvanic Body Spa & Body Shaping Gel

WHAT IS IT FOR?	Thorough cleansing	Brightening, glow	Anti-ageing skincare	Body shaping and toning
USAGE	As a first step in your routine, morning and evening for two minutes.	After cleansing and toning, once per day for two minutes.	After cleansing, 2-3 times per week	3 times per week
BENEFITS	Visibly cleaner, smoother, softer, and radiant skin, with reduced appearance of pores.	Boosts your skin's natural beauty for brighter, plumper, and bouncier skin ¹ .	Firming and lifting your skin while removing impurities and helping to visibly reduce lines and wrinkles ² .	Diminishing the appearance of cellulite and smoothing skin to provide a firmer and more toned appearance.
THE TECH BEHIND IT	Two-Sense Motion Technology	Variable Pulse Microcurrent	Direct/Galvanic microcurrent	Pulsating direct/Galvanic microcurrent

*Source Euromonitor International Limited; Retail Value RSP terms; all channels; 2017 to 2020. Beauty Systems are at-home Skin Care Beauty Devices that are exclusively paired or recommended to be used with a topical consumable of the same brand. Claim verification based on Euromonitor custom research and methodology conducted January-March of 2021. Sales of at-home skin care beauty devices included sales of electrical facial cleansers as defines in Passport database. This category does not include hair care/removal appliances, body shavers, and oral care appliances.

¹ Results obtained from a third-party professional assessment based on a 12-week in vivo, single-blind consumer use study with 32 female subjects aged 25 to 40 who used the ageLOC Boost System once a day for 12 weeks.

² Results obtained from a 2020 third-party professional assessment based on a 4-week in vivo, single-blind, randomised consumer use study with 30 female subjects who used the system three times a week.



ageLOC LumiSpa & Activating Cleanser

1 step. 2 minutes. 7 benefits.




Benefits

- 37% smoother skin after a single use³.
- 37% softer skin after a single use³.
- 2X visibly cleaner skin after 12 weeks⁴.
- 62% more radiance after 12 weeks⁴.
- Reduced appearance of pores.
- Clearer, brighter skin.
- Visibly refreshed complexion.

Activating Cleansers

Formulated for the needs of dry, normal/combo, oily, sensitive, or blemish-prone skin and to enable the precise cushioning and cleansing effect the device needs. Each cleanser contains our ageLOC Blend to target the sources of visible signs of skin ageing.

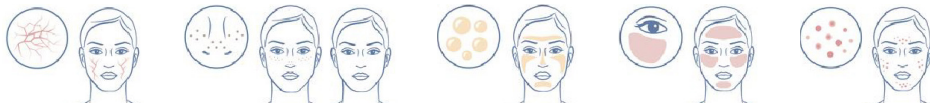

How does ageLOC LumiSpa work?

With its soft silicone head, ageLOC LumiSpa's proprietary Two-Sense Motion Technology enables the device to move at the precise frequency needed to reveal smoother, softer skin after the first use.

How to use

- 1 Moisten your face with warm water and apply the Activating Cleanser that best suits your skin type. Wet the silicone head and turn on your device.



- 2 Glide it slowly across your face in broad strokes. Don't scrub or press too hard. The device will pause if you are misusing it. It will pause every 30 seconds. This is your cue to move to another area of your face. After 2 minutes, your device will switch off.



- 3 Rinse your face with water and gently pat it dry, then apply any further skincare. For optimal results, use the ageLOC LumiSpa system for 2 minutes, twice per day.



- 4 Remove and rinse the Silicone Head. Gently towel dry and store your device in a dry place.

Go further with eye care



Skincare doesn't end with cleansing. Complete your routine by caring for the delicate skin beneath and around your eyes. Use **ageLOC LumiSpa Accent Head & IdealEyes** with your device to exfoliate and invigorate your eye area with the gentle touch your skin needs.

³ Results based on a 4-week in-vivo consumer use study performed on 80 subjects with sensitive skin who used the ageLOC LumiSpa Activating Cleanser for Normal to Combination Skin with the Firm Silicone Head.

⁴ Results based on a 12-week in-vivo consumer use study featuring 150 female subjects with different skin types. Results extracted above correspond to subjects with normal to combination skin. Those subjects were instructed to use the ageLOC LumiSpa with the Normal Silicone Head and ageLOC LumiSpa Activating Cleanser for Normal to Combination Skin in place of their normal facial cleanser, morning and night.



ageLOC Boost & Activating Serum

Enjoy brighter, plumper skin in just two minutes a day.



Benefits

- Boosts your skin's natural beauty for brighter, plumper, and bouncier skin¹.
- Promotes skin cellular turnover.
- Helps ease signs of premature ageing.
- Offers gentle exfoliation.
- Patent-pending Variable Pulse Microcurrent Technology that optimises and activates skin for a more youthful-looking complexion.


Activating Serum

Provides optimal conductivity for the Variable Pulse Microcurrent and skincare benefits from its key ingredients:

- ageLOC Blend
- Brightening Blend: Sea daffodil extract, Vitamine C, Niacinamide
- Caviar Lime Fruit Extract
- Plumping Blend: various fruit extracts, hyaluronic acid


Variable Pulse Microcurrent Technology

Intermittent, variable pulses are sent to the skin. Together with the Activating Serum, it optimises skin for a more youthful-looking complexion.

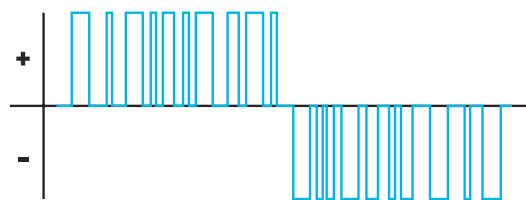

How to use

- 1 After cleansing and toning, turn on device by swiping finger up full length of touch surface.

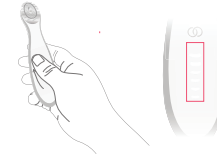

- 2 Dispense three pumps of ageLOC Boost Activating Serum onto the device's conductive surface. Spread over one half of the face.



- 3 Begin application by swiping finger up full length of touch surface. Glide conductive surface over face in an upward/outward motion.



- 4 Repeat Step 2 & Step 3 on other half of the face.



- 5 Massage excess Activating Serum into skin, then follow with your daily skincare regimen.

Before/After picture¹



Before application

After 12 weeks



Before application

After 12 weeks

First visible results in just 7 days¹ and benefits that last throughout the day⁵.

⁵ Results obtained from a third-party professional assessment based on a 6-week in vivo, consumer use study with 20 female subjects, between 25 and 40 years old, who used ageLOC Boost Activating Serum together with the ageLOC Boost Device, 1/day, for 6 weeks.



ageLOC Galvanic Spa & Facial Gels

Looking for the secret to youthful-looking skin? You've found it.



Benefits

- Help to visibly reduce lines and wrinkles and reveal a firmer, more lifted, contoured, appearance⁶.
- Target the sources of ageing to preserve the look of youth and reduce the appearance of ageing.
- Revive and invigorate tired, stressed skin.
- Helps combat the negative effects of oxidative stress caused by environmental stressors.
- Galvanic current supports the supply of beneficial ingredients to the epidermis.


Facial Gels

The secret to those amazing benefits? A combination of your device and these key ingredients:

- Marigold and sea kelp extracts (in the Pre-Treat Gel).
- ageLOC ingredients, Arginine and skin conditioning agents, Green tea extract and Jania rubens extract (CEA) (in the Treatment Gel).

Direct/Galvanic microcurrent technology

Like a set of magnets, the device uses direct currents to support the delivery of ingredients to the skin. Using different settings, a negative current first moves the negatively charged Pre-Treat Gel to the skin, where it binds to impurities. A positive current then withdraws those impurities, while pushing the positively charged ageLOC ingredients in the Treatment Gel to the epidermis.

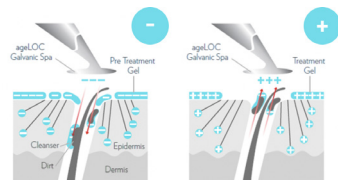


How to use



- 1 After cleansing, evenly apply the contents of one vial of Pre-Treat Gel onto the entire face and neck, avoiding eyes, eyelids, and lips.
- 2 Push the selection button until the number 1 (2 minutes) is lit on the display.
- 3 Thoroughly wet hands with NaPCA Moisture Mist or water. Hold the device with a damp hand while touching the chrome plate on the back and begin to slowly glide the Conductor head over your skin. Be sure to avoid sensitive areas like your eyes, eyelids, and lips.
- 4 When your galvanic application is over, the device will emit a long beep followed by two short beeps and automatically turn off. Remove any remaining gel from your skin and the device with a damp towel.
- 5 Repeat the entire process using a whole vial of ageLOC Treatment Gel and selecting application number 2 (3 minutes) on your device.

Before/After picture⁶

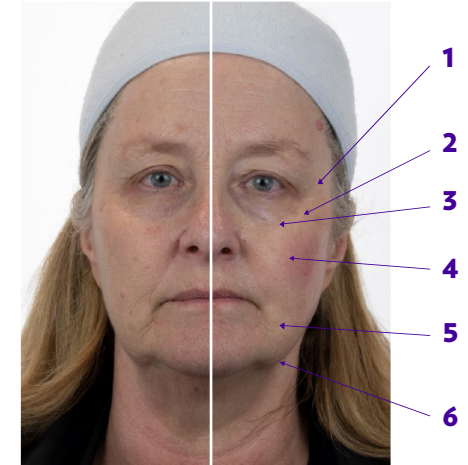

Before application After **one** application

- 1 Reduction in appearance of crow's feet.
- 2 Reduction in appearance of lines and wrinkles in the under-eye area.
- 3 Eye bag lift/improved skin firmness in eye area.
- 4 Reduced appearance of fine lines in cheek area.
- 5 Reduction in visible nasolabial wrinkles ("laugh lines").
- 6 Uplift look in chin area/improved skin firmness.

See 48% improvement to skin elasticity after 4 weeks⁶!

⁶ Results obtained from a 2020 third-party professional assessment based on a 4-week in vivo, single-blind, randomised consumer use study with 30 female subjects who used the system three times a week.



ageLOC Galvanic Body Spa & Body Shaping Gel

Your body optimised.



Benefits

- Target the sources of ageing to preserve the look of youth and reduce the appearance of ageing.
- Helps smooth the appearance of cellulite for a firmer, contoured appearance of your skin.
- Strengthens the extracellular matrix (ECM) of the skin to help build and maintain healthy-looking skin
- Patented pulsating Galvanic current aids in stimulating skin, and supports the supply of beneficial ingredients to the epidermis.


Body Shaping Gel

These great benefits are brought to you by a combination of your device and these negatively charged key ingredients:


- Laminaria digitata extract
- Theobroma cacao extract

Pulsating direct/Galvanic microcurrent technology

Using repelling electrical charges, our Galvanic pulsating microcurrent technology is optimised to support the delivery of key ingredients on the skin of the arms, abdomen, buttocks, and thighs while the microcurrent pulses help stimulate the skin. With ageLOC Body Shaping Gel, this technology helps reduce the visible signs of skin ageing for a smoother, firmer, youthful look.

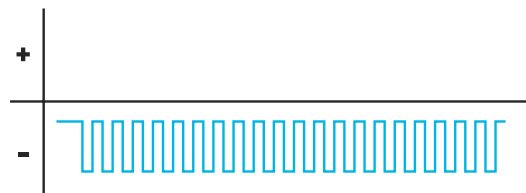

How to use

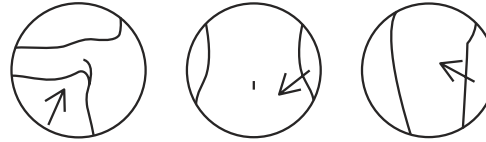

- 1** Apply ageLOC Body Shaping Gel on the area to be immediately targeted.
- 2** Thoroughly moisten your fingers and the palm of your hand with water or NaPCA Moisture Mist.
- 3** Press the on/off button once. The ageLOC rings will appear in the middle of the display.
- 4** Touch the conductive surface to your skin where you have applied ageLOC Body Shaping Gel. You will hear 1, 2 or 3 beeps as the device automatically adjusts the current level for your skin. Begin massaging lightly in a circular motion where the product has been applied.
- 5** At the end of your application, the device will emit a long beep followed by two short beeps indicating that application is complete. Remove excess product from your skin and follow by ageLOC Dermatic Effects for an extension of benefits.

Before/After picture⁷



Before application After 6 weeks

Going further with ageLOC Dermatic Effects



Apply ageLOC Dermatic Effects twice a day for a more contoured, smoother looking you.

82% of users had improvements in the firmness of the skin on their thighs⁸.

⁷ After six weeks of using ageLOC Galvanic Body Spa combined with ageLOC Body Shaping Gels three times a week and ageLOC Dermatic Effects twice a day.

⁸ Results obtained from a third-party professional assessment based on a 12-week in vivo consumer use study with 30 female subjects aged between 18 and 60. ageLOC Body Shaping Gel and ageLOC Galvanic Body Spa were used three times per week for five minutes each on thighs, abdomen, buttocks, and upper arms. ageLOC Dermatic Effects was used twice daily.



Nu Skin Devices Comparison Table



LUMISPA SYSTEM
(Device with ageLOC
Activating Cleansers)



LUMISPA ACCENT
(Device with ageLOC
LumiSpa IdealEyes)



**ageLOC GALVANIC
FACIAL SPA SYSTEM**
(Device with ageLOC
Facial Gels)



**ageLOC
GALVANIC FACIAL
SPA WITH TRU FACE
LINE CORRECTOR**



**ageLOC GALVANIC
FACIAL SPA WITH
ageLOC NUTRIOL
INTENSIVE SCALP &
HAIR SERUM**



**ageLOC GALVANIC
FACIAL SPA WITH
ageLOC BODY
SHAPING GEL**



**ageLOC GALVANIC
BODY SPA WITH
ageLOC BODY
SHAPING GEL**



**ageLOC BOOST
SYSTEM**
(Device with ageLOC Boost
Activating Serum)

MAIN BENEFITS

GALVANIC CURRENTS SUPPORT
THE SUPPLY OF INGREDIENTS
TO THE EPIDERMIS

X

X

X

X

X

WHEN USED WITH THE
SPECIFICALLY FORMULATED
TOPICAL, THE MICROCURRENTS
OPTIMISE AND ACTIVATE SKIN

X

WHEN USED WITH THE
SPECIFICALLY FORMULATED
TOPICAL, THE MICROCURRENTS
OPTIMISE AND STIMULATE SKIN

X

REDUCE VISIBLE
SIGNS OF SKIN AGEING

X

X

X (SCALP)

X

X

TARGET THE PREMATURE
SIGNS OF SKIN AGEING

X

DECREASE APPEARANCE OF
FINE LINES AND WRINKLES

X

X (AROUND EYES/
MOUTH/FOREHEAD)

MORE CONTOURED
APPEARANCE

X*

LIFTED UPPER
EYELIDS

X

FIRMER LOOKING
SKIN

X (AROUND
THE EYES)

X

X (THIGHS & ABDOMEN
& UPPER ARMS)

X (THIGHS & ABDOMEN
& UPPER ARMS)

X

IMPROVEMENT IN
SKIN ELASTICITY

X*

PLUMPER BOUNCIER
LOOKING SKIN

X**

Disclaimer:

We have listed the main benefits, but please note
that many of our devices + cosmetic products have
secondary benefits as well.

*Results obtained from a 2020 third-party professional assessment based on a 4-week in vivo, single-blind, randomised consumer use study with 30 female subjects who used the system three times a week.

**Results obtained from a third-party professional assessment based on a 12-week in vivo, single-blind consumer use study with 32 female subjects aged 25 to 40 who used the ageLOC Boost System once a day for 12 weeks.



Nu Skin Devices Comparison Table



LUMISPA SYSTEM
(Device with ageLOC
Activating Cleansers)



LUMISPA ACCENT
(Device with ageLOC
LumiSpa IdealEyes)



**ageLOC GALVANIC
FACIAL SPA SYSTEM**
(Device with ageLOC
Facial Gels)



**ageLOC
GALVANIC FACIAL
SPA WITH TRU FACE
LINE CORRECTOR**



**ageLOC GALVANIC
FACIAL SPA WITH
ageLOC NUTRIOL
INTENSIVE SCALP &
HAIR SERUM**



**ageLOC GALVANIC
FACIAL SPA WITH
ageLOC BODY
SHAPING GEL**



**ageLOC GALVANIC
BODY SPA WITH
ageLOC BODY
SHAPING GEL**



**ageLOC BOOST
SYSTEM**
(Device with ageLOC Boost
Activating Serum)

MAIN BENEFITS

RADIANCE SKIN /
BRIGHTNESS

X

SKIN TONE
EVENNESS

GENTLE EXFOLIATION TO
PROMOTE SKIN CELLULAR
TURNOVER

SMOOTH & SOFT
LOOKING SKIN

X

X

X

X

THOROUGHLY
CLEANSED SKIN

X

IMPROVE PORE
APPEARANCE

X

REDUCED APPEARANCE OF
EYE PUFFINESS/BAGS
AND DARK CIRCLES

X

APPEARANCE OF VOLUMINOUS,
THICKER, HEALTHIER
LOOKING HAIR

X***

STRENGTHEN
AND FORTIFIES HAIR

X

SMOOTH THE APPEARANCE
OF ORANGE PEEL
LOOKING SKIN

X

X

**Results obtained from a third-party professional assessment based on a 12-week in vivo, single-blind consumer use study with 32 female subjects aged 25 to 40 who used the ageLOC Boost System once a day for 12 weeks.

***Results obtained after the combined use of the ageLOC Nutriol Intensive Scalp & Hair Serum twice daily followed by a two-minute ageLOC Galvanic Spa application together with the ageLOC Nutriol Shampoo and conditioner as the regular hair wash routine.

Disclaimer:
We have listed the main benefits, but please note
that many of our devices + cosmetic products have
secondary benefits as well.

