



Dear Customer,

Many thanks for purchasing our product, please kindly read the manual instructions carefully and operate the system in the correct way.

✧ Working Principle


This system combines functions of Far-infrared ray heating and Airpressure massage therapy.

Airpressure therapy is also called pressotherapy that aims to massage the whole body so that to promote blood circulation、eliminate body muscle tiredness、tighten skin and improve body metabolism.Far infrared heating further helps speeding up body discharges sweat so that to get rid of toxins,quickly burns redundant fatness/cellulite.Both combination, it can achieve the effects of weight loss、body contouring and health-keeping.

✧ Main Effects


- Burning fat、losing weight 、 body contouring and keeping fit.
- Perspiration detoxification,discharge heavy metals along with sweat,promote the metabolism.
- Relieve rheumatism and cold,unclog the meridian,reducing body fatigue and relieving stress.
- Activate,purificate and optimize the blood cells and body tissues,improve body immunity.
- Stimulate body circulation system, promote blood circulation, strengthen and cultivate Viscera.
- Enhancing medication/nutrition lead-in and obsorption.

✧ Panel Introduction

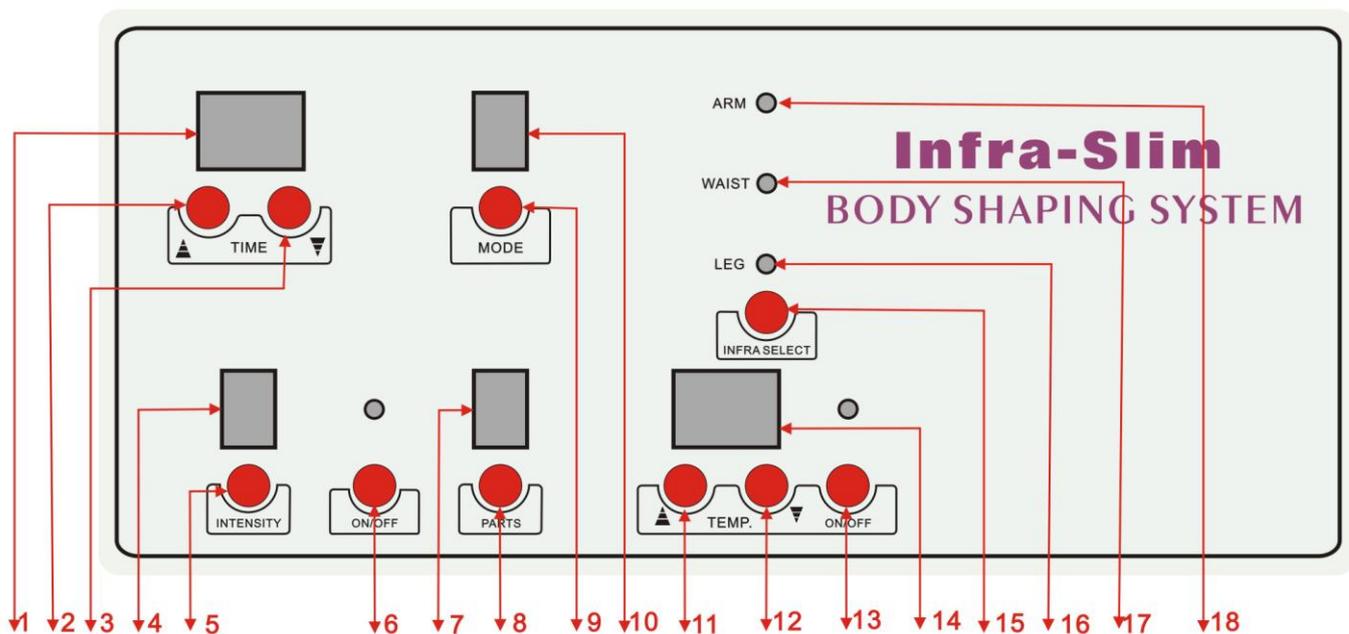

1.TIMER	2,3.TIME UP/DOWN
4.AirPressure INTENSITY	5.INTENSITY select
6.Airpressure ON/OFF	7.Massage PARTS
8.PARTS select	9.MODE select
10.Massage Mode	11,12.Temp. UP/DOWN
13.Infra ON/OFF	14.Temperature Display
15.Infra select	16.LEGS Infrared
17.WAIST Infrared	18.ARMS Infrared

✧ AirPipes Connection

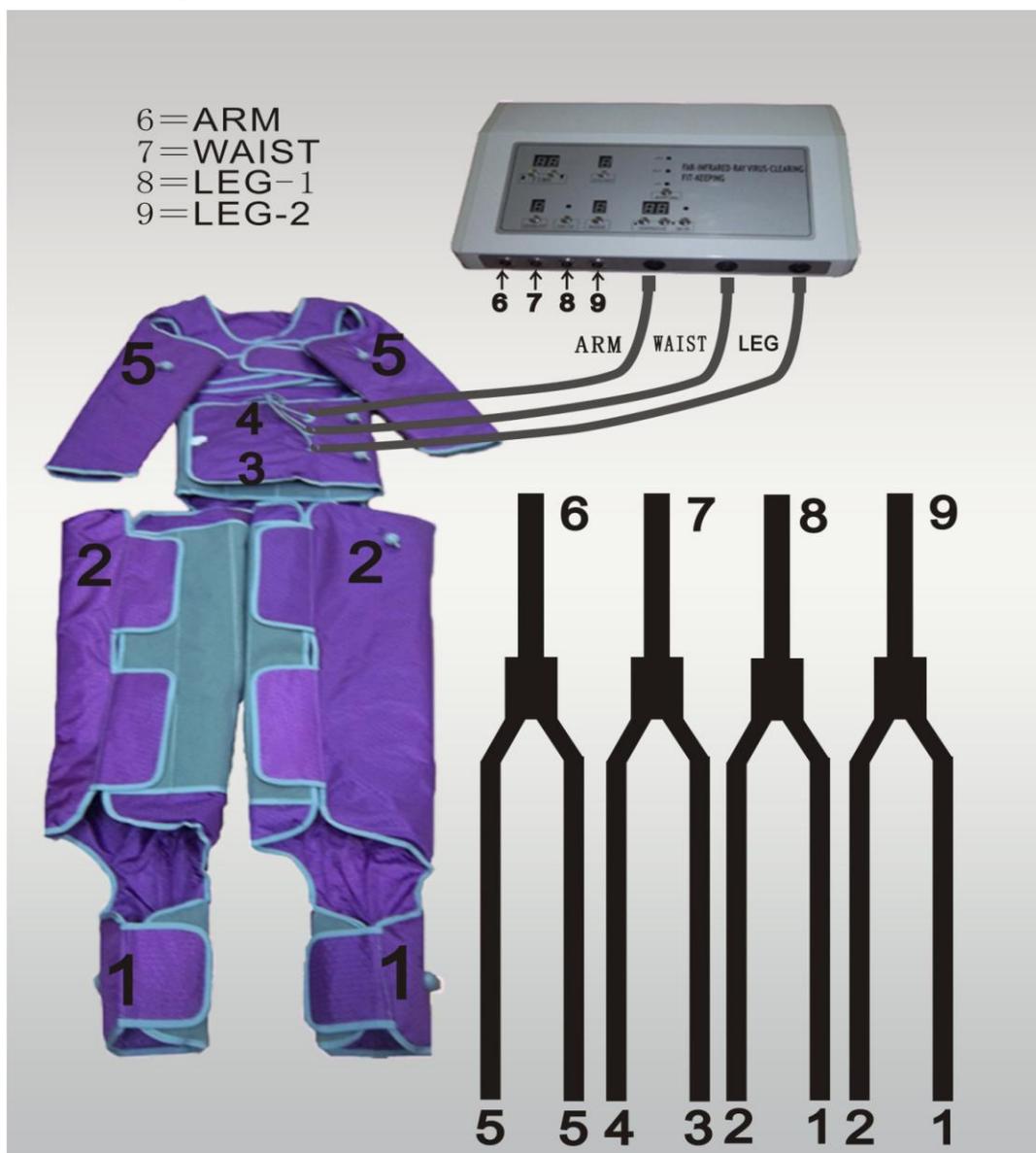


Operation Instructions

◆ Airpressure Steps:

1. Connect airpipes to the according outlet.
2. Connect the power wire, then press power switch to start the machine.
3. Set **TIMER** to the desired time using the ▲ and ▼ arrows.
4. Select massage **PARTS**


“1” : LEGS(FEET)

Airwave massages LEGS and FEET, which tightens legs muscle、 make legs slim
To massage feet, promoting blood circulation of body.

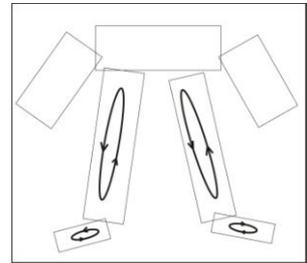


“2” : WAIST

Airwave massages WAIST, which tightens waist muscle 、 make waist slim.

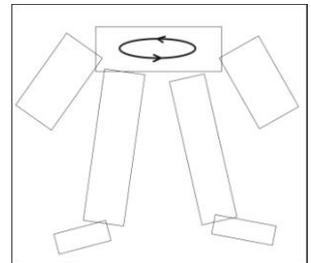


“3” : ARMS

Airwave massages ARMS, which tighten arms muscle and make arms slender.

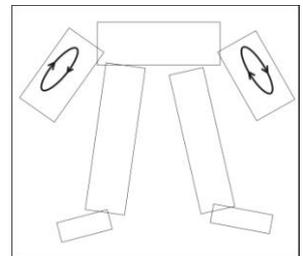

“4” : LEGS(FEET)-WAIST-ARMS

Airwave massages LEGS to WAIST then to ARMS, which makes whole
body relaxed、 tighten body muslce.

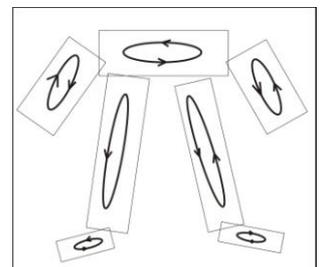

5. Set massage **MODE**(1~6 mode,adjustable).

- Mode " 1" means each inflation lasts 2 seconds,then deflation 1 second.
- Mode " 2" means each inflation lasts 4 seconds,then deflation 1 second.
- Mode " 3" means each inflation lasts 6 seconds,then deflation 1 second.
- Mode " 4" means each inflation lasts 8 seconds,then deflation 1 seconds.
- Mode " 5" means each inflation lasts 10 seconds,then deflation 1 seconds.
- Mode " 6" means each inflation lasts 12 seconds,then deflation 1 seconds.

6. Set the desired **INTENSITY** (1~7 levels)

7.Press **ON/OFF** to begin airpressoure massage

8.Press **ON/OFF** gain to pause massage

◆ Infrared Heat Steps:

1.Connect Infrared Wires to the according outlets.

Notes:

- Have to connect each Infrared Heat wire to the according outlet so that can adjust the temperature.If machine shows "EO" ,it means the wrong connection,then temperature can not be adjusted.
- Each temperature(from 35°C to 80°C) (ARM-1 to LEG-2) can be set freely,and each temperature can be different
- Each part can work independently or can be set to work together.

2. Press **INFRAR SELECT**, **ARM** indicator lamp is on.

Set the desired temperature for ARMS.

Press **ON/OFF** to begin the Infrared heat of ARMS. (Press again **ON/OFF** again to pauze ARMS)

3.Press **INFRAR SELECT** , **WAIST** indicator lamp is on.

Set the desired temperature for WAIST.

Press **ON/OFF** to begin the Infrared heat of WAIST (Press **ON/OFF** again to pauze WAIST)

4. Press **INFRAR SELECT,LEG** indicator lamp is on.

Set the desired temperature for LEGS.

Press **ON/OFF** to begin the Infrared heat of LEGS(Press **ON/OFF** again to pauze LEGS)

5.After doing the above operations,temperature of each part is going up,slowly you can feel the heat.

6.Turn off the power when stop using the system ,pull off the accessories and store the machine well.

✧ Technical Parameter

Voltage Supply	AC220V-240V/50Hz or AC110v/60Hz (Optional)
Power Supply	400W
Max Time	90 mins
Temperature	35°C to 80°C,optional
AirPressure INTENSITY	1-7(Level,optional)
Massage PARTS	LEGS (FEET), WAIST、 ARMS
Massage Mode	6 modes
Machine Dimension	45(L)*30(W)*20(H)
Costume Size	FREE(LARGE)
Net.Weight	12.5KGS
Gross.Weight	14.0KGS(All accessories included)
Packing Dimension	56*37*37(cm)

Caution Items

- The suggested time of using this machine should not be over 30 minutes.
- Never pull out the power plug or other connected pipes as the machine is power on.
- Before doing therapy, customers must take down the sharp objects,which would probably puncture the suit and cause serious consequences;
- Stop using the unit if the suit cracks to avoid the danger of electric shock.

✧ Forbidden Items

- Pregnant women and the women have not completely recovered from childbirth.
- People with high blood pressure,heart disease and serious diabetes.
- People with serious skin diseases.
- People with metal rack inside body.

✧ Packing List

Item Name	Quantity	Remarks
Machine Controller	1 (set)	
Costume	1 (pc)	One-piece
AirPipes	1 (set)	
Power Wire	1 (pc)	
Spare Fuse	2 (pc)	
English Manual Instructions	1 (copy)	