

iSee360

EYE MASSAGER

USER MANUAL

OBM Distribution / Breo-USA
13900 Sycamore Way
Chino, CA 91710
<http://www.breousa.com>
info@breousa.com

breo[®]

**Thank you for purchasing
BREO iSee360.**

Please read this manual carefully before using
the product, and keep it in a safe place for future
reference.

Contents



Introduction	3
Inside the Package	3
Function Deployment	4
Preparation	5
Operation	6
Precautions	7-8
Troubleshooting & Technical Parameters	9

Introduction


iSee360 adopts air pressure and hot compress massaging technology. It rubs, presses, hot compresses, and massages the ocular region in order to relax it. You can enjoy the soothing music while massaging to completely relax your brain.

Inside the Package

The package includes: Main Unit, Controller, AC Adapter, Batteries, User Manual, and Earphones.

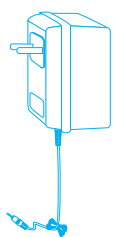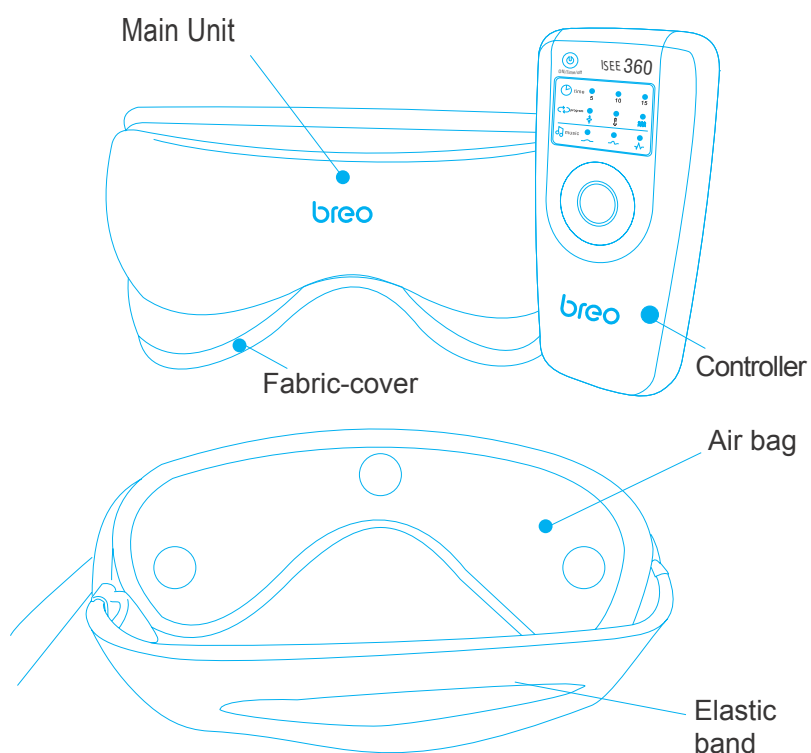


AC Adapter

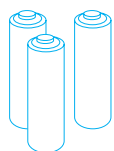

Batteries




User Manual

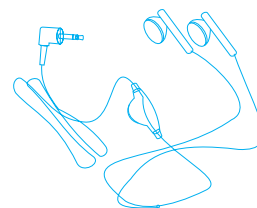

Earphones

Function Deployment

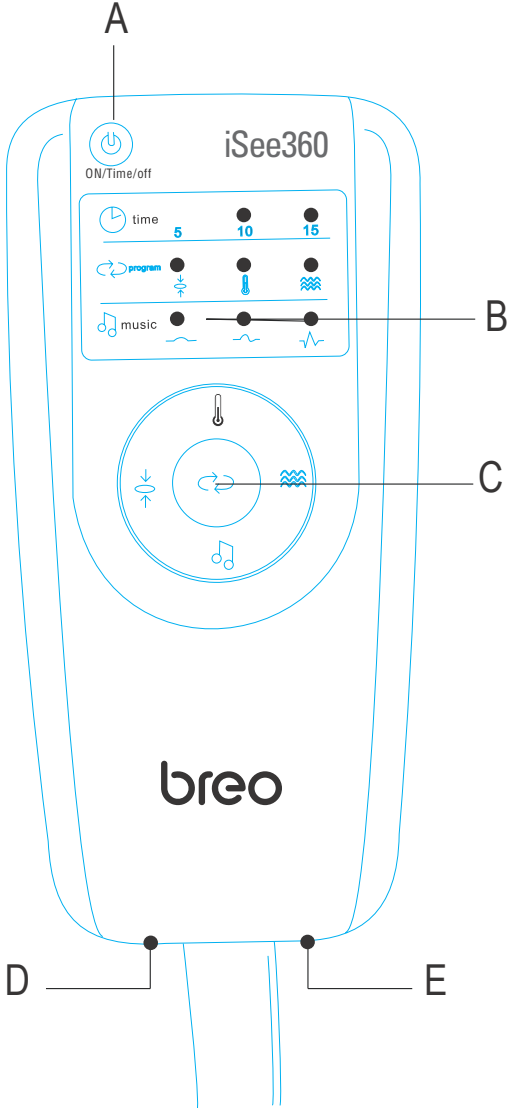

- [A] On/Time/Off
- [B] LED Screen
(From top to bottom: Time / Air pressure, Heat, Vibration / Music)
- [C] Operating Buttons
Clockwise: Heat / Vibration / Music / Air pressure
Center Button: Mode selection
- [D] Earphone Jack
- [E] DC Power Port

Preparation

Prior to use, check the integrity of the machine and accessories:

- a. No damage on connection between controller and main unit.
- b. Battery power is sufficient.
- c. No damage on AC adapter.

Power Options

1. Batteries

- a. Open the battery cover on the back of the controller.
- b. Install 4 AA batteries in accordance with the positive and negative poles.
- c. Close the battery cover.



2. AC Adapter

- a. Insert the plug of AC adapter in the controller's DC port.
- b. Connect the other end to AC power socket.
- c. Insert the earphones' cylindrical plug into the controller's earphone jack.

* Choose power option #1 OR #2.

Wear

This massager is to be used for the eyes only. For ease and safety, adjust buckles until a comfortable fit is achieved.

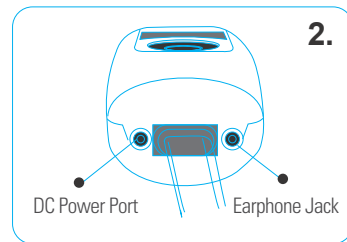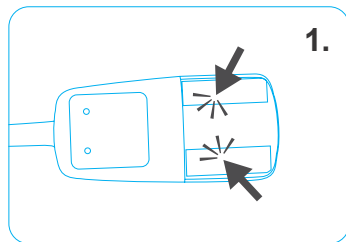

Front view



Side view

Operation

Power on – Press On/Time/Off Button

Preset default is Air pressure + Hot compress / Vibration + Hot compress / Air pressure + Vibration (3 modes circulate), with relaxing music & 5 minutes duration.

Manually operate related function buttons to switch among music, time, vibration, hot compress and air pressure modes.

Power off – Press and hold On/Time/Off Button.

Manual mode

On/Time/Off

Press once – turn ON (initiates default mode)

Press again – TIME (select massage duration: 5 min / 10 min / 15 min)

Press and hold – turn OFF

Mode selection

Press to initiate these modes:

Air pressure + Heat + Vibration

Heat + Vibration

Air pressure + Vibration

Air pressure only

Heat only

Vibration only

Music

Press to select music (Relaxing / Nature) (*Sounds of Nature featured: running water, birds chirping, waves, crickets, dripping water)

Press and hold – STOP the music.

Vibration

Press to initiate vibration modes:

circular motion / strong / medium / weak / light / comfort / strong thump / medium thump / fast thump / vibration left and right alternatingly / fast vibration left and right alternatingly.

Heat

Press to select temperature. (Low-High)

Air pressure

Press to select mode. (Mode 1/Mode 2/Mode3/Mode 4/Loop)

Precautions

1. Please read the manual carefully and operate accordingly.
2. Children under 7 years old cannot use the massager.
3. This massager should not be used by people who suffer from brain disease, eye disease, or another physical or mental disorder.
4. Take off contact lenses / eye glasses before use.
5. Put on / wear the unit first BEFORE turning the power on.
6. Do not use while walking or driving.
7. Keep your eyes closed for a few seconds after use, then open them slowly.
8. Keep the massager away from humidity, heat, and/or corrosive gas.
9. If the massager will not be used for a long time, remove the batteries to avoid battery leakage.
10. This massager is designed as a personal care product. It is not for diagnosis/treatment of any health disorder, nor meant for commercial use.
11. To clean, wipe with aseptic cotton and alcohol. Do not wash or immerse in water.
12. This product is not a substitute for healthcare professionals' advice and treatment. In case of doubt or if necessary, please consult your healthcare professional.

13. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
14. Please note that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
15. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help

Troubleshooting & Technical Parameters

When operational problems occur, please refer to the chart below.
Do not disassemble the machine or try to repair it by yourself.

Problem	Possible Reason	Possible Solution
Power-ON failure	Cable is broken / in bad condition	Check the condition of the cable
	Batteries are low / out	Change batteries
	Batteries are not installed correctly	Check if positive/negative poles of batteries are positioned right Make sure the batteries are snug in place
	AC adapter is not connected properly	Check the connection of the AC adapter to the power socket .
Weak heat or air pressure	Batteries are low / out	Change batteries Make sure the batteries are installed correctly and snug in place

Technical Parameters

Power	DC: 6V
Ideal environment	Temperature: 5°C-40°C
	Relative Humidity: ≤80%
	Atmospheric pressure: 860hpa ~ 1060hpa
Dimension	Main Unit: 230mm × 105mm × 100mm
	Controller: 141mm ¼68.5mm ¼44mm
Net Weight	485g
Power	6W
Heating temperature	30°C-45°C