



# QUICKCHANGE<sup>®</sup>

## INCONTINENCE WRAP



Change  
in 60  
Seconds

### Superabsorbent Core


Locks in moisture and output to keep surrounding area dry including skin and linens

### Unique Shape

Designed to be changed without fully undressing

### Works with Briefs/Diapers

Use inside briefs for quicker changes

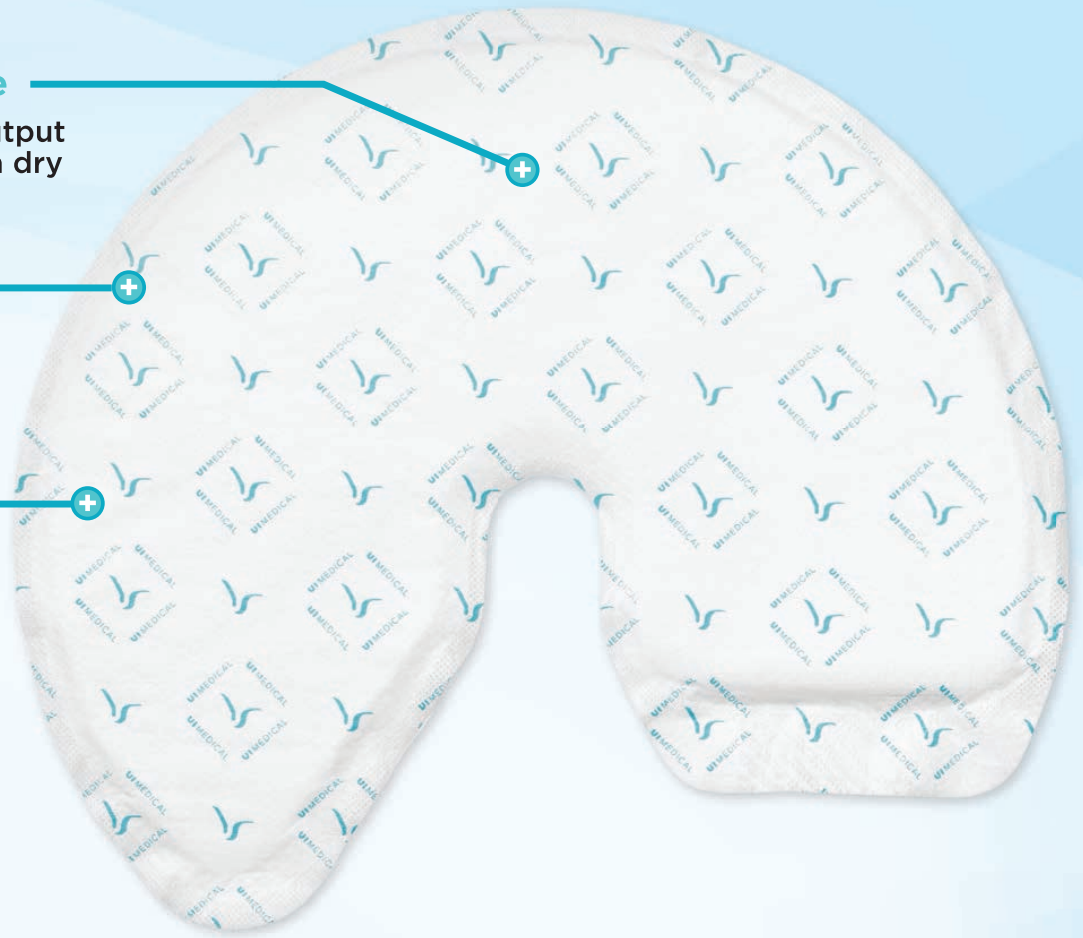

### ACHIEVE MAXIMUM PERFORMANCE

- The base must have a snug seal around the penis. Make sure there are no gaps.
- The wrap is typically folded over twice. It should never be rolled into a cone/funnel.
- The round edge of the folded wrap should point towards the head.
- The hook fastener requires a firm pinch to hold.
- Change the wrap as necessary.
- If using inside a brief/diaper, wraps may be replaced without also replacing an soiled diaper.

### CLIENT PROFILE:

This product is primarily for urinary incontinent, bedbound or chairbound men who do minimal walking or for nighttime use in active men.

### WHY QUICKCHANGE WRAPS?

It's named QuickChange for good reason: it only takes 60 seconds to change a wrap, regardless of a man's size, anatomy or weight. One caregiver, one very quick change, allowing a higher and better level of care...at home.



Easier Changes



Promotes Skin Health



Minimize Sleep Interruptions



Works With Briefs/Diapers



Holds >17 Ounces



No Lifting Required

### SCAN ME

For easy-to-follow pre-recorded training videos and instructions.

Call to schedule a no-obligation, training session at 800.206.2816.

[quickchange.com/how-to-use](http://quickchange.com/how-to-use)





# DIRECTIONS FOR USE

## Standard Penis: Low Urine Pressure (1 wrap)



Place white side up, round edge towards head, and two flaps towards legs.



Fold pointed flap over penis at 45°, then fold straight flap over, slightly pulling upwards.



With two fingers, pinch hook fastener and inner flap together to ensure attachment.

## Standard Penis: Strong Urine Pressure (2 wraps)



Follow steps 1-3 above. Place 2nd wrap underneath and 3-4" higher than 1st wrap.



Pull the top center of the 2nd wrap down and over the 1st wrap, just to the right of center.



Lift the straight flap up and over towards the center, pinching the hook fastener to secure.

## Retracted/Indwelled Penis (1 wrap)



Place white side up, round edge towards head, two flaps pointed towards legs.



Slide pointed flap underneath scrotum, lift and scoop upwards and then fold straight flap over.



Pinch hook fastener and inner flap together to ensure attachment.

## Abdominal Obesity: Hidden/Covered Penis (2 wraps)



Place underneath the scrotum, pull towards head to ensure contact, tuck between the legs to catch runoff.




Place a 2nd wrap, white side up, over belly so that source of the urine is visible through the center opening. Fold hook flap.



Fold the 2nd wrap in half downwards and lengthwise, making sure the center hole is adequately covered.

## Enlarged Scrotum: Medium (1 wrap)

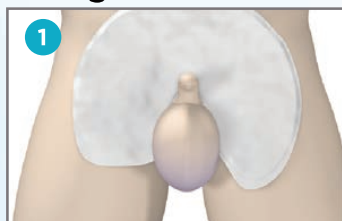

Place white side up, round edge towards head, two flaps towards legs.




Pull pointed flap underneath scrotum creating bowl shape and pull straight flap over.




Pinch hook fastener and inner flap together to ensure attachment.


## Enlarged Scrotum: Large (2 wraps)

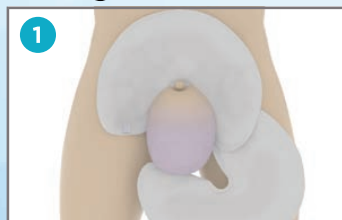

Position two wraps around scrotum & penis; pinch & secure a straight and a pointed flap together as a fulcrum.

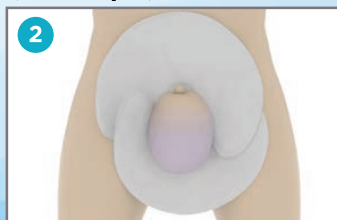


Pull other two flaps towards each other, starting a bowl shape underneath scrotum & pinch two flaps together.

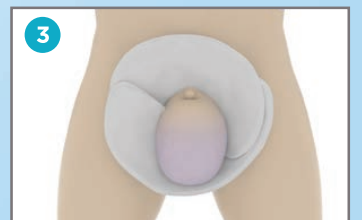


Readjust and tighten the flaps repeatedly, completing the bowl shape. Ensure good contact under scrotum.

DO NOT...

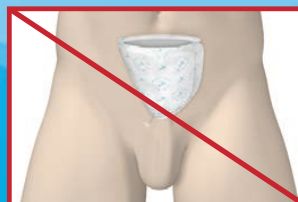

Roll into a tight, cone/funnel shape.

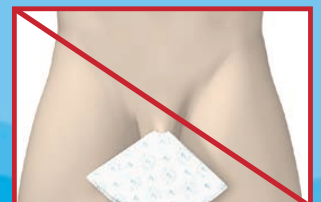

Leave the QuickChange pointed towards the legs.