

- Stand on the scale for a few seconds to get your complete body analysis details on your app.
- Get off the scale when your results are displayed on the screen.
- Pair the scale by following instructions on the app and get off the scale and wait for it to turn off.
- Please do not use this product if you are pregnant or have a pacemaker.

CAUTION

- Please take off shoes and socks and keep your bare feet in touch with electrodes for measuring body composition accurately.
- Do not place the scale on wet and slippery surfaces and should not stand on it with wet feet to avoid slipping.
- The scale is a high precision measuring device. So, never jump or stomp on the scale.
- Please handle it carefully to avoid damage when moving it.
- When weighing, both feet must be in the center of scale to avoid tipping of device.

PRODUCT SPECIFICATIONS

Product Dimensions:	27 X 27 X 2.2 cm
Product Weight:	1.1 kg
Limit of Measurement:	5 kg-180 kg/ or 400 lb/ or 28 st
Units of Weight Measurement:	kg / lb/ st
Connection:	Wireless
Operating temperature:	5 °C - 35 °C
Battery:	Lithium Battery
Display:	LCD digits
Functionality:	Auto On/ Auto Off
Platform:	6mm tempered safety glass with round corners
Division:	100g/0.2lb

MAINTENANCE AND STORAGE

Please use soft tissue or a damp cloth with neutral detergent, alcohol, or glass cleaner to clean the surface of the scale. Avoid using soap or other chemicals as they may lead to damage. The storage environment must be dry. Also, keep it away from water, heat, and extreme coldness.

Disclaimer:

The scale is for family use only. Body Weight, BMI, BFR, Muscle, Body Hydration, Bone mass, BMR, Visceral fat, Protein rate, Body age, Standard weight, Body Fat measurements are for reference only. You should consult a doctor when you undertake any diet or exercise program.

CUSTOMER SUPPORT

CONTACT ADDRESS

USA
Nureca INC,USA
276 5th Avenue, Suite 704-397,
New York (NY) - 10001, USA

INDIA
Corporate Office (Mumbai)
Nureca Limited
Office No. 101, 1st Floor, Udyog Bhavan
Sonawala Lane, Goregaon East
Mumbai, Maharashtra - 400063, India.

Contact us
India: +91-7527013265 / +91-9356658436
Website: www.drtrust.in
Corp Website: www.nureca.com
Email: customer@nureca.com
Connect with us on social networks
Facebook: @drtrust
Instagram: @drtrustsin
Youtube: NurecaUSA

COPYRIGHT © 2022 NURECA LTD ALL RIGHTS RESERVED

ABOUT US

Driven by the passion for innovation, we at Dr Trust endeavour to provide our customers with the latest medical inventions with an objective to promote good health and wellness all around the world. All the medical devices and health monitors provided by Dr Trust are supported by accurate, latest and ground breaking technologies, innovated at our headquarter in NY, USA. All our products adhere to the most stringent CE and FDA guidelines and are strongly recommended by doctors and health practitioners. Our products are designed in the utmost exemplary ways to ensure that their accuracy and convenience are unrivalled. The ease of their use and operation makes them even more suitable for users of all age groups.

Dr Trust strives to enhance the quality of lifestyle by providing with the most trusted and innovative health care and wellness products. Being a renowned global leader in health care products, Dr Trust ensures that our technically efficient team works dynamically and tirelessly to provide the best of the medical devices to our clients. The products that we have to offer are suitably designed for use at homes, laboratories and hospitals.

Our ground breaking solutions allow you to monitor your health in the easiest ways possible. In today's era when all of our lives are too hassled to handle, it becomes a bit difficult to pay attention to our health. But it has now become easier with the coming of the monitoring devices which can be conveniently used at homes and even on the go.

We bring to you a variety of best self medical devices, trusted and used by Doctors, medical professionals and home users all over the world.



Scan to View
Product Demo Video
www.drtrust.in/509

QUICK START GUIDE

Step1
Download the "Dr Trust 360 App".

Search for Dr Trust
360 App in App Store

OR
Scan the QR code to
download the app



drtrust360.com
The Science of Living

Step2
Open the App and complete your registration as a user.

Step3
Step on the scale to connect it to your phone.

Step4
The Dr Trust 360 app will then automatically scan for the scale. Select the scale and touch the "Click to Pair" icon next to it.

Step5
Stand on the scale again and click to begin measurement.

OVERVIEW

The Dr Trust Smart Body Fat Scale 2.0 - 509 is a high-tech smart scale that allows multiple users to use and sync their data to the Dr Trust 360 App for easy tracking. It is a great option for those looking to monitor multiple metrics including body weight, body fat, BMI, and bone and muscle mass.

Built with 4 high-precision sensors, it uses "step on" technology for showing highly accurate results instantly by simply stepping onto the scale. While an LCD backlit display makes it easy to view the results, even in a dim light room.

The smooth surface of this digital scale is made of sturdy strong high-quality tempered glass. It is easy to clean and completely safe with rounded corners. This scale is equipped with 300 mAh internal lithium-ion batteries so that it can conveniently re-energize at any time through various sources.

Key Features

- Easy-to-read display
- Auto calibration
- 4 High-precision Sensors
- Compact, sleek design
- Lightweight
- Convenient for large families
- Easy to connect with the App
- Syncs multiple users' data

MEASURING BODY WEIGHT ONLY

When you stand on the scale, it will measure your body weight.



Instructions

Please follow given steps:

- 1) Place the scale on a hard and flat surface.
- 2) Carefully step on the platform of scale and stand as stationary as possible until the weight reading is shown

and locked on the LCD display.

- 3) The scale will be automatically off without further operation.



Note:

- Do not use it on the carpet or soft surface.

MEASURING BODY COMPOSITION

Operating Instructions

Download App

Apple mobile phones: can download in the App Store;
Android phones: can download in Google Play Store;



drtrust360.com
The Science of Living

Scan the QR code to
download the app



Android 5.0
or above




ISO 9.0 and
above


Pairing the Scale with Your Smartphone

Open Dr Trust 360 App. Follow the on-screen instructions to complete the registration process. After successful activation of your account put the scale on a hard and flat surface and then step on the scale lightly for pairing. Choose the device and then click on "Pair a device" to start pairing your smartphone with scale.

(Following messages will be displayed for the completion of pairing process).

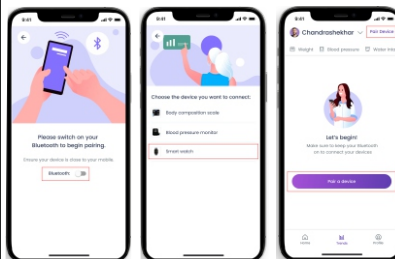

After pairing, stand on the scale again with bare feet (do not move) and click "Go" to begin measurement. The scale shuts off after few seconds. App will show and save all indexes and history record. (Following screens will be displayed for showing the complete data).

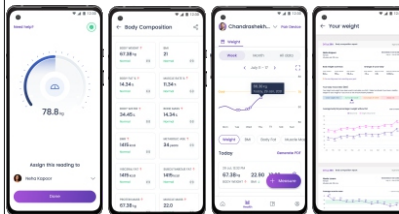

Note: To make a proper measurement of body fat, please keep your bare feet on the electrodes of scale.

OTHER INTELLIGENT FEATURES

- Informative links for each data point
- Link to corresponding article with detailed explanation
- All graphs are scrollable
- Outer segments are also scrollable

- Each data point is clickable to show its data by marker
- Dynamic colour coded marker backgrounds that change according to position
- Graphs of all smart-scale measurement data to bottom-bar
- Instantly compare each data point across time with all other readings
- Progress across all data points are easily trackable by user
- Current selected data pops-up as auto-marked graph

POWER SUPPLY

It must be re-energized when the weight machine starts showing low battery and an error message while taking weight.

ERROR CODES

Overload and Low Power Indication

- 1) LO--When the LCD indicates "Lo", please re-energize the scale.
- 2) Err--Maximum capacity exceeded. When the scale is overloaded (Over 180 kg/400 lb) "Err" will appear on the screen.



COMMON TROUBLESHOOTING

- Please ensure your Smartphones and Dr Trust 360 App both are updated to the latest versions.
- Please ensure that the Dr Trust 360 App has been given necessary permissions to work.
- The scale connects within the App only.
- You need to pair the scale for the first time only.