

## • ABOUT US •

Here at Dr Physio USA, we make the most imaginative, unwinding, and restorative, well-being and individual care Products.

in the hubbub of today's lifestyle, it becomes impossible to nurture yourself at times. In the chaos of maintaining every aspect of life efficiently, we forget that our bodies are not machines. they need absolute care and tenderness. Pain in neck, shoulder, back or foot should never be taken for granted. They can lead to serious ailments in the long run. To help you ease the situations of painfulness, Dr Physio USA brings to you excellent ways of relaxation and pain relief through our premium gadgets. We offer a non-intrusive answer to help renew and reestablish your muscles and relieve throbbing pain with our various products

The variety of items and massagers we offer are intended to target particular areas of the body viz neck, shoulder, back, foot and so on. the adaptable capabilities of our gadgets are meant for utmost body relaxation after a long day at work.

Sit Back, Relax and enjoy the amazing Dr Physio experience by Dr Trust!!

Dr Physio from Dr Trust  
[www.drphys.io](http://www.drphys.io) | [www.drtrust.in](http://www.drtrust.in)

## Dr Physio

### Pillow Massager Pro - 1009

#### • Quick Start Guide •

##### Step 1

Position the cushion on your sofa and bed

##### Step 2

Connect the plug into a power outlet and then to the power point located in the back of massager.

##### Step 3

Turn the speed regulation dial to set the level of massage speed.

##### Step 4

To activate soothing heat with the massage, push the heat button.

##### Step 5

Tapping or kneading massage can be chosen by pushing right side button up and down.

##### Step 6

When finished, turn off massage and heat functions and unplug the pillow massager.




1009

## ● OVERVIEW ●

Dr Physio Pillow Massager is excellent for getting instant relief from neck pain and back pain. It comes with 6 deep-kneading rotating Shiatsu Massage nodes that rotate to deliver a powerful and intense shiatsu massage to relieve muscle pain quickly. This massage pillow can be used simply on the chair in the office, car or kitchen to give yourself an effective massage treatment. This little size massager can be used all over the body hassle free and can be packed easily in a bag for use on the go.

## ● CONTROLLER ●

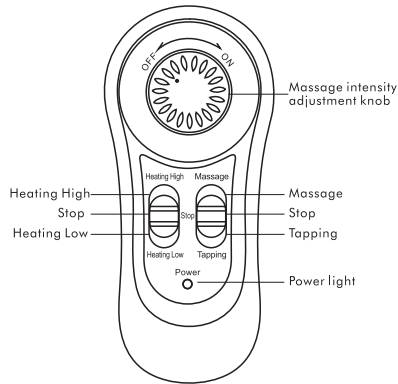

## ● OPERATION ●

Connect the massager to an external power supply. The power indicator light will start flashing on the controller.

### Adjusting Kneading or Tapping massage

1. The button on the right side of the controller helps you to adjust the level of percussion massage. There are 3 gears - "massage", "stop", and "tapping". Just push the button to the corresponding position.
2. The button on the left side of the controller helps you to adjust the heat. There are 3 stalls - "heating high", "stop", and "heating low". Just push the button to the corresponding position to enjoy a soothing warm massage.

### Adjusting the Intensity of Massage

1. Adjust the knob from "OFF" to "ON" in the direction of rotation, open the kneading or tapping massage function, and gradually increase the intensity of the massage.
2. Adjust the knob from "On" to "OFF" in the direction of rotation, gradually reduces the intensity of the massage.



## ● CLEANING AND MAINTENANCE ●

- Firstly, pull out the plug from the socket before to start cleaning the machine.
- Always keep the massager in dry and dust-free area.
- Avoid water or other liquid to contact with the machine.
- Please use only a dry cloth to clean it.
- Avoid water or other liquid to contact with the machine.
- Please do not use volatile oil, Bleach, or other corrosive solvents to wipe the machine.
- For long-term storage, place keep the unit in its box or a safe, dry, and cool place.

## ● TECHNICAL SPECIFICATIONS ●

|               |                |
|---------------|----------------|
| Rated Voltage | AC 220 V-240 V |
| Rated Power   | 25W            |
| Rated Time    | 15 minutes     |

## ● IMPORTANT NOTES ●

1. Do not use product on wrong voltage.
2. This is not recommended for medical use.
3. Do not put the controller in water to avoid controller failure.
4. Do not use the machine in the bathroom, swimming pool, spa, and other humid environments.
5. Do not place the massager too close to heat sources.
6. Do not damage the product casing or jacket, to avoid the risk of injury.
7. Before connecting the power plug to the socket, be sure to check the power switch in the "Off" position.
8. Do not pull the plug wires directly with a wet hand.
9. After use, be sure to push the power switch in the off position and unplug power cord.
10. Do not attempt to repair it, contact to a professional if it is damaged.
11. Pregnant woman, people with limited ability of any health condition and children should not allow to use this product.
12. Allow the unit to cool 30 minutes after use.
13. Do NOT use product on swollen or inflamed areas of the skin.
14. Avoid contact with sharp edges or pointed objects, which might cut or puncture the fabric.



## • KEY MESSAGE AREA •

Dr Physio Pillow massager can be used on the shoulder, back, waist, neck, thigh and calf area etc.

Back



Waist



Neck



Leg(Calf)



## • CUSTOMER SUPPORT •

### CONTACT ADDRESS

USA

Nureca INC.USA  
276 5th Avenue, Suite 704-397,  
New York (NY) - 10001, USA

INDIA

Corporate Office (Mumbai)

Nureca Limited  
128 Gala Number Udyog Bhavan,  
1st Floor Sonawala Lane, Goregaon East  
Mumbai City Maharashtra 400063

Contact us

India: +91-7527013265 / +91-9356658436

Website: [www.drtrust.in](http://www.drtrust.in)

Corp Website: [www.nureca.com](http://www.nureca.com)

Email: [customercare@nureca.com](mailto:customercare@nureca.com)

Connect with us on social networks

Facebook: @drtrust

Instagram: @drtrustisin

Youtube: NurecaUsa

COPYRIGHT©2021 NURECA LTD ALL RIGHTS RESERVED



Scan to View  
Product Demo Video  
[www.drtrustusa.com/1009](http://www.drtrustusa.com/1009)