

● SAFETY INFORMATION ●

- The glass platform is slippery when wet. Make sure the platform is dry before you are stepping on it for measurement.
- This scale is not intended for commercial use.
- To weigh safely and accurately, please stand on the center of the platform and stand still until a result is calculated.
- DO NOT strike or drop the scale.
- DO NOT stand on the edge of the scale or jump on it.
- Always place the scale on a hard, dry, and flat surface before measurement.
- Handle the scale with care, as it is a precision instrument.
- If the scale fails to turn on, check whether the batteries are installed correctly. Replace the batteries if necessary.
- If an error appears on the display or if the scale does not turn off automatically, remove the batteries for 3 seconds, then reinstall them.
- DO NOT overload the scale (max. weight capacity: 400lb / 180 kg).

● CARE AND MAINTAINANCE ●

- Clean the scale with a damp cloth and avoid getting water or any other liquid inside the scale.
- DO NOT use chemical or abrasive cleaning agents to clean the scale.
- Keep the scale in a cool environment with dry air.
- Always keep the scale in a horizontal position; do not store vertically.

● TECHNICAL SPECIFICATIONS ●

Weight Capacity	400 lb (180 kg)
Weight Units	Kg/lb/st
Battery	3 x 1.5V AAA batteries
Platform	Glass Finish with Perimeter Light
Display	LCD with Backlight
Body Type	6mm Tempered Glass
Automatic Shut Off	10 seconds

● CUSTOMER SUPPORT ●

CONTACT ADDRESS

USA
Nureca INC.USA
276 5th Avenue, Suite 704-397,
New York (NY) - 10001, USA

INDIA
Corporate Office (Mumbai)
Nureca Limited
128 Gala Number Udyog Bhavan,
1st Floor Sonawala Lane, Goregaon East
Mumbai City Maharashtra 400063

Contact us
India: +91-7527013265 / +91-9356658436
Website: www.drtrust.in
Corp Website: www.nureca.com
Email: customercare@nureca.com
Connect with us on social networks
Facebook: @drtrust
Instagram: @drtrustisin
Youtube: NurecaUsa



Scan to View
Product Demo Video
www.drtrustusa.com/507

COPYRIGHT©2021 NURECA LTD ALL RIGHTS RESERVED

● ABOUT US ●

Driven by the passion for innovation, we at Dr Trust endeavour to provide our customers with the latest medical inventions with an objective to promote good health and wellness all around the world. All the medical devices and health monitors provided by Dr Trust are supported by accurate, latest and ground breaking technologies, innovated at our headquarter in NY, USA. All our products adhere to the most stringent CE and FDA guidelines and are strongly recommended by doctors and health practitioners. Our products are designed in the utmost exemplary ways to ensure that their accuracy and convenience are unrivalled. The ease of their use and operation makes them even more suitable for users of all age groups.

Dr Trust strives to enhance the quality of lifestyle by providing with the most trusted and innovative health care and wellness products. Being a renowned global leader in health care products, Dr Trust ensures that our technically efficient team works dynamically and tirelessly to provide the best of the medical devices to our clients. The products that we have to offer are suitably designed for use at homes, laboratories and hospitals.


Our ground breaking solutions allow you to monitor your health in the easiest ways possible. In today's era when all of our lives are too hassled to handle, it becomes a bit difficult to pay attention to our health. But it has now become easier with the coming of the monitoring devices which can be conveniently used at homes and even on the go.

We bring to you a variety of best self medical devices, trusted and used by Doctors, medical professionals and home users all over the world.


Dr Trust
Supernova Personal Scale - 507

● QUICK START GUIDE ●


Step 1
Check batteries and place the unit on a hard surface.

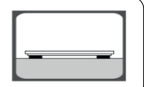

Step 2
Gently step on the scale and it will activate automatically.




Step 3
The scale will display your weight and lock it.

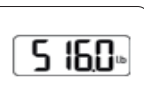

● INTRODUCTION ●

The Dr Trust Supernova Weight Scale- 507 is intended for use in the measurement of bodyweight. It is a completely automatic system that makes body weight measurements easy so that you can manage your weight easily to maintain good health. It has a large display that shows your accurate weight measurements clearly after each use.

● POWER SUPPLY ●

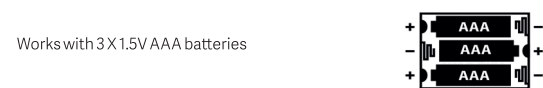


● BATTERY REPLACEMENT ●


- Remove the battery cover
- Insert 3XAAA batteries inside the battery compartment
- Make sure to check the polarity.
- Close the battery cover tightly.

● OPERATION ●

STEP 1. Put the scale on a hard & flat surface. Avoid keeping it on carpet or soft surfaces.

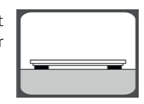


STEP 2. Gently step on the scale and it will activate automatically.




STEP 3. Put your feet again on the scale and remain standing on it without moving or shifting your feet or body. The scale will then display your weight and lock it.

● AUTOMATIC SWITCH OFF ●

When you step off the scale, it will automatically switch Off. Auto shut-off occurs if the display shows "0.0" or shows the same weight reading for approximately 10 seconds.




● TO CONVERT WEIGHT UNIT ●

- Find unit switch on the back of scale.
 - Press the unit conversion button when the scale is power on, LCD shows current weight unit. Press the button again to convert the weight unit between kg-lb-st.
- 

● OVERLOAD INDICATION ●


When the scale is overloaded (over 180kg) the following will appear on the screen:



Note: Please step off to avoid damage.

● LOW BATTERY INDICATION ●

When the battery is low, please open the battery cover, then take batteries carefully out of the battery compartment to replace the old batteries with the new ones.



● RE-TEST ●

When the battery is low, please open the battery cover, then take batteries carefully out of the battery compartment to replace the old batteries with the new ones.

