

• ABOUT US •

Driven by the passion for innovation, we at Dr Trust endeavour to provide our customers with the latest medical inventions with an objective to promote good health and wellness all around the world. All the medical devices and health monitors provided by Dr Trust are supported by accurate, latest and ground breaking technologies, innovated at our headquarter in NY, USA. All our products adhere to the most stringent CE and FDA guidelines and are strongly recommended by doctors and health practitioners. Our products are designed in the utmost exemplary ways to ensure that their accuracy and convenience are unrivalled. The ease of their use and operation makes them even more suitable for users of all age groups.

Dr Trust strives to enhance the quality of lifestyle by providing with the most trusted and innovative health care and wellness products. Being a renowned global leader in health care products, Dr Trust ensures that our technically efficient team works dynamically and tirelessly to provide the best of the medical devices to our clients. The products that we have to offer are suitably designed for use at homes, laboratories and hospitals.

Our ground breaking solutions allow you to monitor your health in the easiest ways possible. In today's era when all of our lives are too hassled to handle, it becomes a bit difficult to pay attention to our health. But it has now become easier with the coming of the monitoring devices which can be conveniently used at homes and even on the go.

We bring to you a variety of best self medical devices, trusted and used by Doctors, medical professionals and home users all over the world.

Dr Trust

Zen Recharge Personal Scale - 503

• QUICK STARTUP GUIDE •

Note:

This scale operates on rechargeable battery so always check charging before proceeding.

Step 1

Place the scale on a hard and flat surface.

Step 2

Gently step on the scale and it will activate automatically.

Step 3

Remain standing evenly on the scale, without moving or shifting your feet or body.

Step 4

The scale will then display your results.




OVERVIEW

Dr Trust Platinum Recharge Personal Scale 501 is a digital rechargeable scale specially designed with patented Accu gauge sensor technology to provide more accurate and precise results. This technology also ensures instant results and lasts long. It has a large display, so the readability of results is simple and easy for all. Also, the platform of the scale is large enough to accommodate large size feet comfortably. It has skid-proof feet to prevent any potential risks. This device is powered by a rechargeable 500 mAH battery and turns off automatically if it is not in use. It has a USB port for easy charging through the tablet, laptop, or any other compatible sources.

OPERATION

- Step 1:** For highly accurate results, place the scale on a hard and flat surface. Avoid putting it on carpet and other uneven flooring.
- Step 2:** When you are ready to measure, gently step on the platform. The scale will activate automatically.
- Step 3:** Remain standing evenly until weight appears. The weight will be locked on the display for some time.
- Step 4:** The scale will automatically turn off after 10 seconds of non-use.

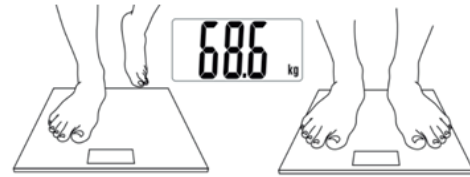

Note:

Please ignore the first-time weight after installing the batteries. As it will consider the weight of scale itself.

CONVERTING WEIGHT UNITS

- Step 1:** Activate the scale by applying pressure onto the platform.
- Step 1:** Find the unit conversion button on the back of the scale. Press the button to convert units from pounds (lb) to kilograms (kg) or stone (st).



Unit button



ERROR CODES

1. **LO**--When battery is low the LCD indicates "Lo".
2. **Err**--Maximum capacity exceeded. When the scale is overloaded (over 180kg/396lbs).

IMPORTANT SAFETY INFORMATION

- The glass platform could be slippery when wet. Make sure the platform is dry before you are stepping on it for measurement.
- To measure your weight safely and accurately, please stand still on the platform until a result is calculated.
- DO NOT strike or drop the scale.
- DO NOT stand on the edge of the scale or jump on it.
- Always place the scale on a hard, dry, and flat surface before measurement.
- DO NOT use chemical or abrasive cleaning agents to clean the scale.
- Keep the scale in a cool environment with dry air.
- Always keep the scale in a horizontal position, do not store vertically.

- Always keep the scale away from strong electromagnetic fields.
- DO NOT overload the scale (max weight 180kg/396lb).
- Before measurement, always remove shoes and socks for best weight result accuracy.
- Do not stack anything on the scale during no use.
- If the scale does not work, please check the charging.
- This is for family use only and should not be used commercially for any reason at all.
- Individuals with the following conditions may not use this scale:
 - Pregnant women
 - Those with symptoms of edema
 - Those on dialysis treatment
 - Those who use heart pacemakers or other implanted medical devices

PRODUCT SPECIFICATION

Max Weight Capacity:	180kg/396lb/28st
Technology:	High precision strain gauge sensor
Graduation:	0.1kg/0.2lb
Weight Units:	kg/ lb/ st
Battery:	450mAH rechargeable battery



Platform:	Tempered Glass
Automatic Functions:	Auto on & Auto off
Display:	LCD
Indications:	Low Power and Overload Indications

Package Contents

- 1 x Digital Bathroom Scale
- 1 x User Manual
- Standard micro USB cable

CUSTOMER SUPPORT

CONTACT ADDRESS

USA

Nureca INC.USA
276 5th Avenue, Suite 704-397,
New York (NY) - 10001, USA

INDIA

Corporate Office (Mumbai)

Nureca Limited
128 Gala Number Udyog Bhavan,
1st Floor Sonawala Lane, Goregaon East
Mumbai City Maharashtra 400063

Contact us

India: +91-7527013265 / +91-9356658436

Website: www.drtrustusa.com

Email: customercare@nureca.com

Connect with us on social networks

Facebook: @drtrust

Instagram: @drtrustisin

Youtube: NurecaUsa

COPYRIGHT©2020 DR TRUST ALL RIGHTS RESERVED



Scan to View
Production Demo Video
www.drtrustusa.com/503