The logo features a red, curved shape above the letter 'i' in 'Urikar'.

Urikar

MASSAGE GUN

Pro 2

User Manual

English	01
Deutsch	08
Français	15
Español	22
Italiano	29
日本語	36


CONTENTS

1. What's in The Box	01
2. Introduction to The Pro 2	01
3. How to Use	02
• Using The Pro 2	
• Choose Angle of Handle	
• Choose Mode	
• Choose Massage Head	
• Charging	
4. Parameter Information	07


What's in The Box

Urikar Massage Gun Pro 2 x1

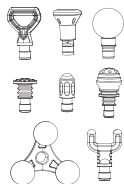Massage Heads x8

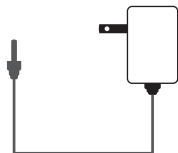Power Charge x1


User Manual x1

Handbag x1

Introduction to The Pro 2


- ① Power Switch
- ② Mode Key | Mode Indicator
- ③ Speed Knob
- ④ Handle-Rotation Button
- ⑤ Speed Level Indicator
- ⑥ Power Indicator & Battery Indicator
- ⑦ DC Input Port
- ⑧ Anti-Slip Handle

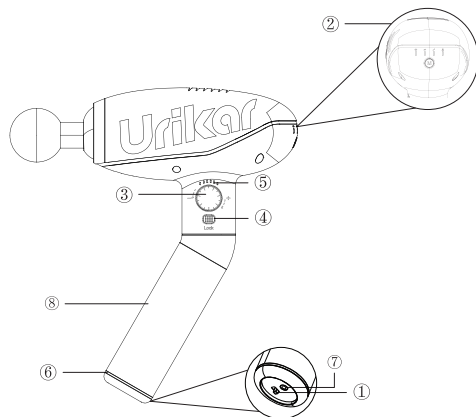
How to Use

● Using The Pro 2

1. Turn on the power switch ① and the message gun is in standby state.

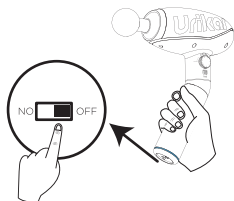
2. Adjust the speed knob ③, and the message gun starts to work (Adjust the mode and speed according to your physique and preferences).

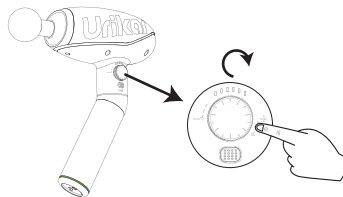● Choose Angle of Handle

Hold handle-rotation button ④ to change the angle of the handle (180°rotatable), which is convenient for massaging the back muscles.

Note: You can hear a click when the rotation is completed.

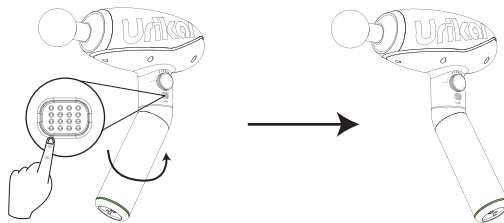
● Chosse Mode

1. Mode Switching


Short press the mode key to switch arbitrarily between constant speed mode, fluctuation mode and stairway mode.

2. Mode Introduction


Constant Speed Mode: Constant speed massage and percussion.

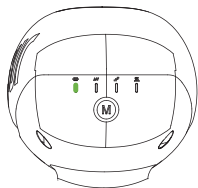
Fluctuation Mode: Switch high and low gear back and forth to massage and percuss.

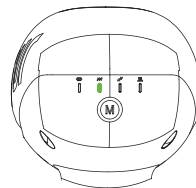Stairway Mode: Advanced massage and percussion from low to high gear.

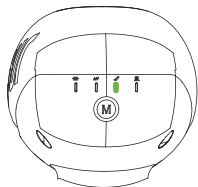Heating Mode: Constant temperature 40-45°C massage and percussion.

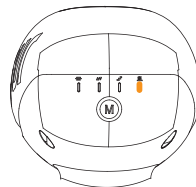
When using the heating mode, please replace the special heating head first, and then long press the mode key to enter the preheating mode (Orange light flashing), after preheating 45s, the heating head enters a constant temperature state (Orange light is always on).

Note: You can't adjust the speed knob in fluctuation and stairway mode.

3. Intelligent Protection Mode


When the machine detects that a continuous opposite external force is applied to the massage head and exceeds the limit set by the machine, the machine will intelligently identify and stop turning, which can effectively protect the motor and prolong the life of the motor.

● Chosse Massage Head

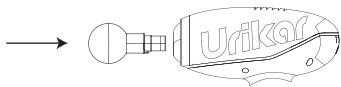
1. Install

Choose the massage head you needed and insert by pressing in, remove by pulling it out.

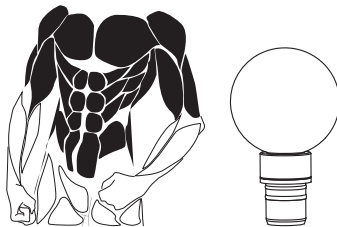
2. Introduction of Massage Head

A. Ball Head: Uniform and soft strength suitable for whole body muscle massage and relaxation.

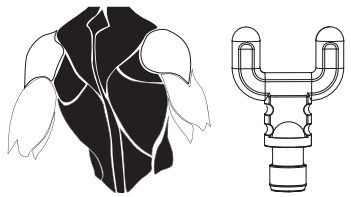
B. U-Shaped Head: Suitable for shoulders, cervical vertebrae and bilateral muscles.


Note: Don't directly hit the cervical spine and spine, it will cause personal injury.

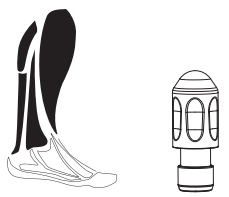C. Bullet Head: Suitable for tendon, acupoint, etc.

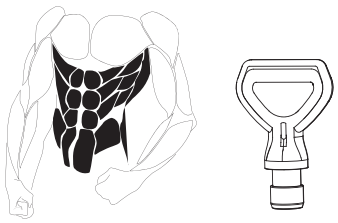
D. D-Shaped Head: Suitable for back, cervical and waist massage.

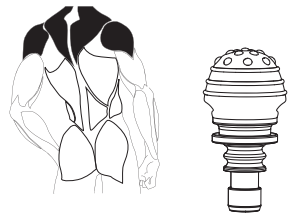
E. Air Plug Head: Suitable for massaging bones and sensitive parts.

F. Flat Head: Suitable for relaxation and plasticity of various muscle parts of the body.

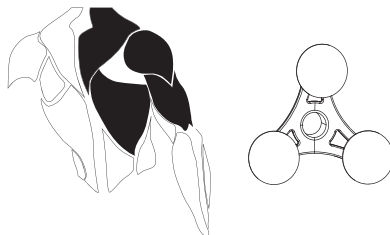G. Triangular Soft Head: Suitable for shoulder massage and relaxation.

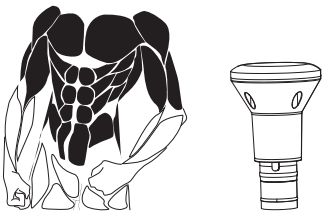
H. Heating Head: It can promote the decomposition of lactic acid after exercise when heated, and is suitable for whole body muscle massage and relaxation.

● Charging

1. How to Charge

Use the included power adapter to connect the charging port and power supply of the massage gun.

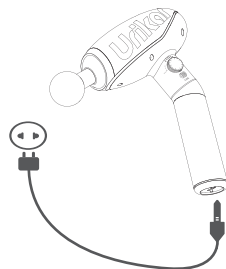2. Description of Battery Indicator

- A. Low Battery Display: Battery <20%, red light is on.
- B. Charging Display: Green light flashes.
- C. Fully Charged Display: The green light is always on.

Note: This product cannot be turned on to use when charging.

Parameter Information

Name: Massage gun

Brand: Urikar

Model: Pro 2

Size: 268*292*93mm

Weight: 1.2kg

Input Voltage: 100-240V, 50/60Hz

Rated Current: About 1000 mA

Rated Power: 30W

Work Mode: Constant speed mode, fluctuation mode, stairway mode, heating mode

Power Supply: 2600mAh rechargeable lithium battery

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