The logo for Urikar, featuring the brand name in a white, sans-serif font. A red curved line is positioned above the letter 'i'.

MASSAGE GUN

Pro 3


User Manual

English	01
Deutsch	09
Français	17
Español	25
Italiano	33
日本語	41


CONTENTS

1. What's in the Box	01
2. Introduction of the massage gun	01
3. How to Use	02
• Using the massage gun	02
• Using the Timing Function	02
• Choose Angle of Handle	02
• Choose Massage Head	03
• Charging	05
4. Parameter Information	06
5. Maintenance	06
6. Precaution for Use	06


What's in the Box

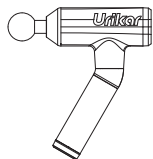Urikar Massage Gun Pro 3 x1

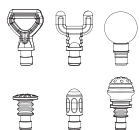
Massage Heads x6

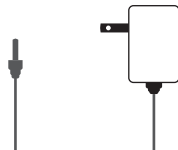Power Charge x1


User Manual x1

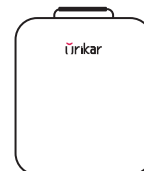Carry Case x1

Introduction of the massage gun

- ① Massage Head
- ② Rotatable handle
- ③ Power Indicator
- ④ Power Switch
- ⑤ DC Input Port
- ⑥ Battery Indicator
- ⑦ Timing Display
- ⑧ Speed Level Display
- ⑨ Speed ^ Button
- ⑩ Speed v Button
- ⑪ Start/Pause Button

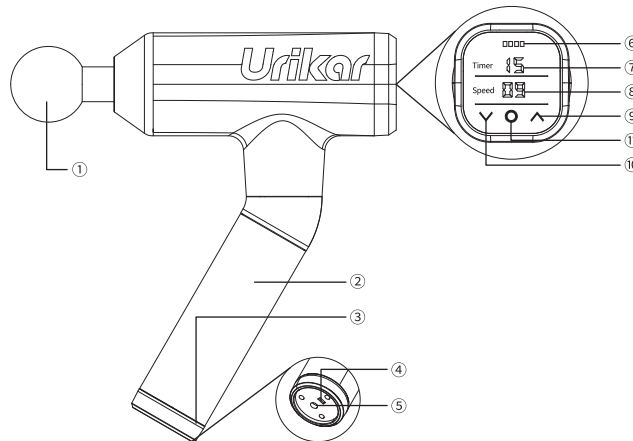How to Use

● Using The Massage Gun


1. Insert a suitable massage head, turn on the power, the display shows two lines of "-.-", the blue power indicator flashes in a circle to indicate that the device is turned on.

2. Short press the "⊙" to start/pause the device.

3. Use the speed button "v" or "A" to adjust the desired speed level(0-30 gear).

Note: Long press "A" to increase the speed, long press "v" to decrease the speed. If you press "v" to the lowest "00", the machine stops working and is in standby mode. The "Timer" on the display screen represents the time from the start to the end of the massage, measured in minutes (it is recommended for general users to use 15min). "Speed" stands for speed gear, the device has 0-30 gears.

4. If it is not used for a long time, please turn off the power switch.

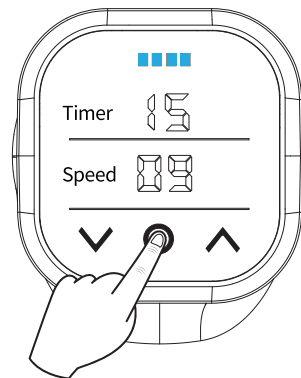● Using Timing Function


1. Massage Gun will be automatically timed (in minutes) after turn on.

2. Turn off the power switch or turn to standby, the time will be reset and displayed as "-.-".

3. The maximum time count is 99 (in minutes) .

● Choose Angle of Handle



As shown in the picture on the right, the handle can be rotated 180° to change the angle, which is convenient for massage of the back muscles.

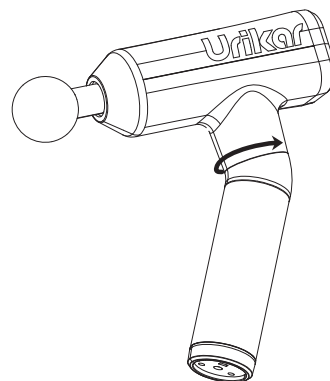● Choose Massage Head


1. Install of the Massage Head

The convex groove of the massage head is aligned with the groove of the gun and insert by pressing in.

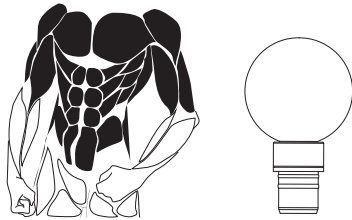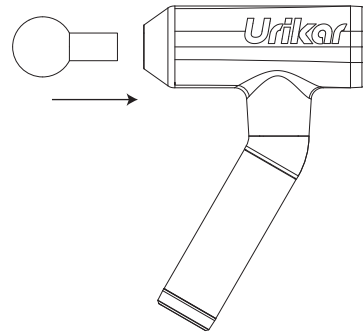
2. Introduction to Massage Head

A. Ball Head: Uniform and soft strength suitable for whole body muscle massage and relaxation.

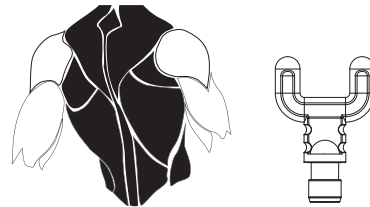
B. U-Shaped Head: Suitable for shoulders, cervical vertebrae and bilateral muscles.


Note: Don't directly hit the cervical spine and spine, it will cause personal injury.

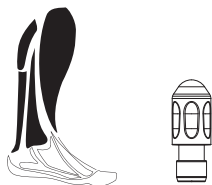C. Bullet Head: Suitable for tendon, acupuncture, etc.

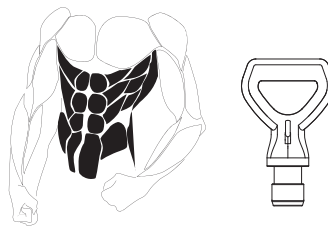D. D-Shaped Head: Suitable for back, cervical and waist massage.

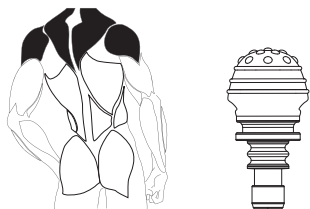E. Air Plug Head: Suitable for massaging bones and sensitive parts.


F. Flat Head: Suitable for relaxation and plasticity of various muscle parts of the body.

● Charging

1. How to Charge

Use the included power adapter to connect the charging port and power supply of the massage gun.

Note: When charging, the blue power indicator light is always on.

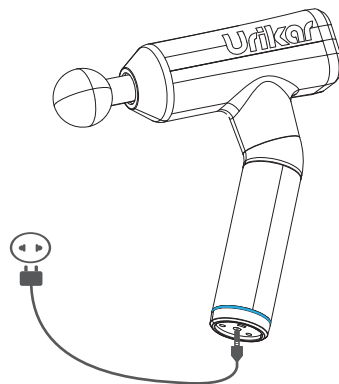2. Description of Battery Indicator


There are four blue grids on the display, and each grid represents 25% power and the corresponding indicator flashes in the charging range.

When the power is less than 10%, the first battery indicator and the power indicator flash red;


When the power is less than 25%, the first battery indicator and the power indicator are always red;

When it is 50% full, the first and second battery indicator are always blue, and the third battery indicator flashes;

When the battery is fully charged, the battery indicator are all blue.
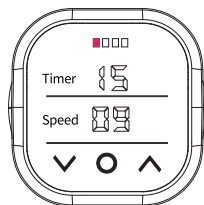
< 25%

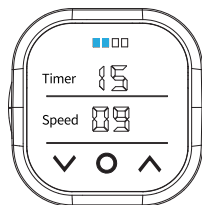50%

100%

Parameter Information

Name: Massage gun

Brand: Urikar

Model: Pro 3

Size: 18.65×26.99×6.8cm / 7.34×10.63×2.68inch

Input Voltage: 100-240V, 50/60Hz

Speed Level: 30 Speed

Power Supply: 2600mAh

Display Method: Touch Sensing

Maintenance

1. Wipe the surface of the massage gun with a slightly damp towel, and then dry it with a soft cloth.
2. When you need to store or travel for a long time, please put the massage gun into the case, and fully charge it next time before you use it.
3. The massage gun should be stored in a dry, well-ventilated, non-corrosive gas place.
4. The massage gun is suitable for transportation by car, train, ship and plane.
5. The product should be fixed during transportation to prevent violent vibration, collision, rolling and falling, and to be upward, light and rainproof.

Precaution for Use

In order to reduce the risk of malfunction, electric shock, injury, fire, death and equipment damage, please read all the instructions and warnings in this manual carefully before using this product.

Forbiddance

1. It is forbidden for minors to use it. If you are injured, do not use this product. Please consult your doctor before using this product when your muscles are damaged.
2. It is forbidden to use on the surface of non-dry, non-clean skin. You can use it across clothes, gently press and move on the skin, each part for about 60 seconds.
3. It is forbidden to use above the neck and parts with less muscle coverage and on bones. Use this product only on the soft tissues of the body without causing pain or discomfort.

4. It is forbidden to use when bruising, pain or discomfort appears.
5. It is forbidden to put any objects into the ventilation holes of the motor.
6. It is forbidden to immerse the product in water.
7. It is forbidden to use other chargers for charging, you must use the charger provided by our company.
8. It is forbidden to disassemble this product in any way.
9. Forbid fingers, hair or other parts of the body near the axis and back of the massage head to prevent pinch and injury.

Note

1. Keep away from liquids or heat sources, keep the vents free of dust.
2. After charging is complete, please unplug the charger.
3. Don't use this product while charging.
4. Please use this product according to the instructions.
5. Before each use, carefully check whether the product can be used normally according to the operation methods of this manual.

Warning

This product cannot be used under any of the following situations without the prior consent of the doctor: Pregnancy, diabetic complications (such as neuropathy or retinal injury), wearing a pacemaker, recent surgery, epilepsy or migraine, herniated disc, spondylolisthesis, diastasis vertebralis, or spondylosis, recent joint replacement surgery or IUD, metal needle or steel plate.

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