RobiComb® Instruction Manual

Facing the facts....

Lice and nits are a pain to deal with. There is no quick solution – it takes time and patience to completely eradicate an infestation. But the good news is the RobiComb® is here to help.

The RobiComb® is an electronic lice zapping comb. It emits a very low energy electric charge that will kill or disable adult lice on contact so they can't reproduce or lay eggs. It is also designed to zap mature lice eggs.

The RobiComb® should replace all chemical treatments; it is effective and safe when used as directed. Most lice have become resistant to the harsh chemicals used in most common treatments but they can never become resistant to the RobiComb®.

What the RobiComb® can do:

- **Kill or disable lice on contact.**
- **Zap large nits.**
- **Detect lice.**

PLEASE READ CAREFULLY BEFORE USING THE ROBICOMB®!!!

FOR BEST RESULTS read the instructions and watch our instructional videos by going to www.lice-guard.com or by scanning the QR code below.

How to use video

Directions for Use:

Step 1: Preparing for use

- Place 1 AA battery in the battery compartment .
- Dry hair completely- RobiComb® must be used on DRY HAIR ONLY.
- Detangle knots- Before using the RobiComb®, long hair should be combed through with a regular brush or comb to untangle any knots.

Step 2: Effectively using the RobiComb®

- Divide hair into small sections: Use hair clips to separate combed hair from the uncombed.
- To turn on: Push the switch to the "ON" position. The RobiComb® makes a continuous high pitched humming sound when working.
- Hold in the correct position: Hold unit with the word "RobiComb®" facing up and with the comb at a 45 degree angle to the roots of the hair.
- Comb close to the scalp: Start at the scalp and pull comb through one inch of hair length (where lice live).
- Humming sound indicator: As you comb the hair, the humming sound will stop momentarily every time a louse, an egg is destroyed or anything such as hair-product build-up or dandruff is large enough to break the connection between the teeth. Brush lice or objects out of combing unit teeth with the special enclosed brush and the humming sound will resume. Resume combing.
- Length of treatment may vary depending on type of hair. Switch off and clean combing unit when treatment is completed.

Step 3: Complete the treatment cycle

- Comb infested hair daily for up to two weeks in order to detect and kill any new lice which have hatched since you began treatment. By treating daily and killing lice before they reach the age of reproduction, you put an end to the infestation cycle.

Warnings:

- Consult with a physician before using RobiComb® on an infant or toddler.
- Avoid contact with eyes, ears or mouth as well as other sensitive areas. Touching both sides of the combing unit may result in tingling sensation which some children may find unpleasant.

- Do not use RobiComb® on sores or broken skin.
- RobiComb® is intended for head hair only.
- RobiComb® should be handled with care
- Persons who suffer from epilepsy or seizure disorders, heart disease, or persons with a pace-maker or other neuro-stimulators should consult with a physician before using RobiComb®.

Clean & Store the RobiComb®:

- Always clean before beginning treatment on another child or storing.
- To clean, remove the combing unit (metal teeth) and gently brush clean. The combing unit can also be washed by dipping the small cleaning brush in isopropyl alcohol and brushing the teeth.
- Let the teeth air dry completely.
- Store in a dry place; preferably not in the bathroom.
- Remove battery if not in use for extended period of time.

How to care video

For any issues or questions please call our USA customer service hotline at 1-888-994-8213, Monday through Friday 8 AM to 8 PM EST or go to www.lice-guard.com.

Troubleshooting:

| PROBLEM | SOLUTION |
|---|---|
| If comb shuts off when touches hair... | Dry hair and combing unit completely |
| If no buzzing is heard when RobiComb® is turned on... | Clean RobiComb® (see cleaning guide) and make sure no debris is caught in between the comb's teeth Ensure that none of the comb's teeth are touching, especially the teeth on the end. If any teeth are touching, gently pry them apart. Insert fresh battery with "+" matching "+" on the RobiComb® Dry combing unit completely |
| If buzzing is silent with combing unit removed... | Remove battery and gently clean metal connections and battery (do not wash) Gently pry metal connections out a little to make better contact with battery Dry combing unit completely |


More Tips and Advice for using the RobiComb®:
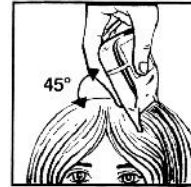
| DO | DON'T |
|---|---|
| Visit our website for videos, tips, and additional information. | Use on wet hair or in humid conditions. |
| Use a brand new battery for best results. | Bring the RobiComb® in contact with eyes, ears, or mouth. |
| Store RobiComb® in a DRY place. | |

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LiceGuard LLC