

ORACURA®

YOUR
ADVANCED
DAILY
ORAL CARE
ROUTINE



PRIMA

OC500

THANK YOU

For Choosing ORACURA® To Protect Your
Beautiful Smile!

Index

Product Specification	3
Configuration	4
Charging The Water Flosser®	4
Steps to Use	5
Steps After Using The Water Flosser®	6
Maintenance	7
Troubleshooting Your ORACURA® Water Flosser®	8
Register Your Warranty	9
Complementary Product Offers	10
How To Dispose?	10
Safety Guidelines	11
ORACURA® Product Range	12

Product Specification



Accessories Included



Operational Manual



450 ml Water Tank



Handle



USB Type-C Cable



2 Standard Nozzles



1 Orthodontic Tip



1 Pocket Tip




1 Tongue Scraper



Protective Pouch

Configuration

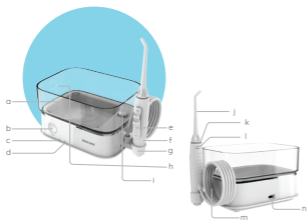

- a] Water Tank
- b] On/Off + Mode Button
- c] Power/Charge Indicator
- d] Mode Indicator
- e] Pressure Control Switch
- f] Water Flosser® Handle
- g] Handle Bracket
- h] Water Volume Mark
- i] Water Tank Rubber Valve
- j] Nozzle
- k] Identification Ring
- l] Nozzle Release Button
- m] Water Pipe
- n] Charging Point

Charging The Water Flosser®


1] The water flosser® is powered by an internal rechargeable battery and charged with a USB Type C port cable.

2] Insert the water flosser® plug securely into the charging point socket. Before inserting the plug, ensure that there are no traces of liquid inside the socket.

3] Make sure that the charging cable is dry and has no damage.

4] The "CH" indicator blinks when the Water Flosser® is charging. The indicator will be steadily on after it is fully charged.

5] If the "CH" indicator blinks during or after water flossing, it means that the water flosser® battery needs to be recharged. Full recharging will take about 4 hours.

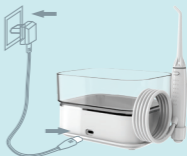

Steps To Use The Water Flosser®

- 1] Attach the water flosser® handle bracket at the bottom of the flosser and rotate clockwise to ensure that it is properly locked.
- 2] Attach the nozzle to the water flosser® handle.
- 3] Fill the water tank with water and place it firmly on the water flosser®.
- 4] Hold the handle and point the nozzle to the sink, then press the Main Button to turn on the water flosser®.
- 5] Press the button again to select the desired mode. If you're using it for the first time, then the clean mode will be selected by default. The massage mode will be activated by pressing the On/Off button again.
- 6] Push the pressure control switch to adjust the water pressure. One can select amongst the 5 pressure settings. If you're using the water flosser® for the first time, we suggest using the lowest pressure setting & gradually increasing the same to suit one's oral care requirements.
- 7] Adjust the nozzle direction by rotating the knob of the nozzle.
- 8] Hold the handle in an upright position with the nozzle tip pointing towards your teeth & gums at a 90° angle.
- 9] Move the water jet slowly along the teeth to ensure a thorough cleanse. The water flosser® will work till the 2-minute cycle is completed, after which it will shut off. You can turn off the flosser on your own as well.
- 10] Remove the nozzle by pressing the nozzle release button.

Note: Avoid using water warmer than 40° C



Scan for Live Demo

You can schedule a live product demonstration by scanning the above QR code.

2 Operation Modes

You can select the desired operation mode by pressing the On/Off button on your water flosser®.

Clean Mode: The water flow is ideal to thoroughly cleanse all the germs & bacteria from the mouth. It is suitable for everyday oral care routine. Blue light is indicated when the clean mode is on.

Massage Mode: It helps in advanced gum stimulation and is used ideally for massaging the gums. It assists in ensuring clean, fresh, and healthy gums and teeth. Green light is indicated when the massage mode is on.



Steps After Using The Water Flosser®

- 1] Dispose of the water in the tank.
- 2] Push the pressure control switch to the 'High' position. Turn on the water flosser® & drain out the remaining water. Push the pressure control switch to the 'Low' position again.
- 3] Turn off the water flosser®.
- 4] Wipe the water flosser® with a dry cloth.
- 5] Remove the nozzle and store it properly.

Disassembly for Travel

- 1] Remove and empty the water tank.
- 2] Disconnect the water flosser® handle bracket from the main unit by rotating it in an anticlockwise direction.
- 3] Point the tip into the sink and extend the water pipe to drain all remaining water.
- 4] To prevent bacteria, make sure there are no water traces inside the water tank.
- 5] Be sure to wipe and dry the water flosser® before preserving it.

To prevent bacteria, make sure there are no water traces inside the water tank. If you're not going to use the flosser for a long time, be sure to wipe and dry the same before preserving it.

Maintenance

Use only water or neutral washing agents to clean the water flosser®. Never use caustic or abrasive agents (For eg. vinegar or limescale removers). Do not use water/ cleaning fluids over 40°C. Keep the water flosser® in a cool place and avoid sunlight.

Main Unit

1] To remove dirt, wipe the unit with a damp cloth (using liquid soap/ water) and then wipe it again with a dry cloth.

2] Do not submerge the main unit in water.



Handle

Wash with water and dry up with a soft cloth. Do not bend, pull or twist the water pipe.


Nozzle

Wash with water and wipe with a soft cloth. Replace the nozzle after every 3-4 months.

Tank

1] Wash with water and wipe with a soft cloth.

2] Keep the tank dry and away from sunlight when not in use for a long period.

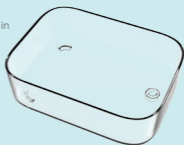

Troubleshooting Your ORACURA® Water Flosser®

Problem	Possible Cause	Solution
It does not work.	You have just purchased the product or it has been left unused for more than 3 months.	Charge the product for at least 4 hours.
It only works for a few minutes even after charging.	The battery has reached the end of its life cycle.	Contact the service center for support.
	Charging time was insufficient.	Charge the product for at least 4 hours.
Water Pressure is very low.	There is dirt stuck in the nozzle.	Replace the nozzle.
	The nozzle is deformed.	
	There is dirt stuck in the filter.	Clean the filter.
Water is not coming out of the product.	The tank is empty.	Fill up the tank.
	The product is tilted too much while in use.	Use the product in the vertical position.
	There must be air caught in the pipe.	Slide the switch to adjust the pressure.
Tank Leaks	The water tank rubber valve was lost	Replace with a new water tank rubber valve.
	The water tank rubber valve isn't installed in place.	Take off the water tank rubber valve and reinstall again.



Please Note: The nozzle tips should be replaced after every 3-4 months.

Register Your Warranty

ORACURA's OC500 Countertop Smart Water Flosser® is covered under 18 + 6 months warranty period.



For the extended 6 months warranty, register your product by scanning the QR code or visiting the link given below.
<https://oracura.in/pages/register>

18-months of limited warranty is applicable for material, manufacturing, or operational defects existing at the time of delivery or for the occurrence of any operational error during regular use. Any mishandling of the product or alteration/repair done by an unauthorized person will void the warranty. Consumer is required to submit the original purchase invoice as proof of purchase date and if requested, the entire product, to support a warranty claim.

If you face any problems while using our product or have any complaints regarding the same, we advise you to visit our return policy page.
<https://oracura.in/pages/return-policy>



Complementary Product Offers

Go Beyond Manual Brushing with Our Sonic Electric Toothbrushes.



If you have purchased an ORACURA® product, then this is your chance to explore various offers complementing your purchased product. Scan the QR code for further details.

Visit: <https://oracura.in/collections/sonic-electric-toothbrush>

DISPOSAL

The rechargeable battery inside your water flosser® is easy to remove for recycling. At the end of the water flosser®'s life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below. During the below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

1] Please discharge the battery power before removing the battery.

2] Remove the water tank, use a sharp hard object to remove the bottom rubber pad (Refer Fig. 1), remove the screws under the rubber pad, and open the shell. (Refer Fig. 2)

3] Open the shell and take out the core components, and use a screwdriver to pry open the battery cover (Refer Fig. 3).

4] Cut off the battery line and pull out the batteries from the inside. (Refer Fig. 4)

If you have recycling questions, please contact your local waste management office.

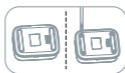


Fig. 1



Fig. 2



Fig. 3



Fig. 4

Safety Guidelines

Read All Instructions Carefully Before Using The Water Flosser®

CAUTION

- Toddlers or infants should not be allowed to use the product. Keep it out of the reach of children.
- Never try to alter, disassemble or repair the product.
- Do not place or store the product where it can fall or be pulled into a tub or sink.
- Do not use the product under wet conditions. (For eg. bathing)
- Be careful while using the product. The strong force of water spray may hurt sensitive organs of your body.

WARNINGS

- This product can be used by children aged 8 years and above. People with reduced physical, sensory or mental capabilities, lack of experience, and knowledge can also use the product if given proper supervision or instruction.
- Use the product only with the charging cable that is provided.
- Make sure to remove the battery before scrapping the product. The battery should be disposed of safely.
- Please refer to the "disposal" section for any concerns regarding the removal of the battery. The product must be disconnected from the supply while removing the battery.
- Do not fill the tank with water warmer than 40°C.
- Never use the adaptor if it has a damaged charging cable or plug.
- Keep the cable away from heated surfaces.
- Do not use the product without a nozzle.
- Do not submerge the product in water. Use the product only for oral cleaning purposes.
- Do not only use mouthwash or any antibacterial rinses in the the water tank.
- The concentrated chemical may damage the product by causing the water tank to crack.
- People who experience tooth or gum pain may not be able to use this product and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums.
- You may experience gum bleeding during the initial stage of using the flosser. However, with regular usage, the product shouldn't hurt anymore.

ORACURA® Product Range

ORACURA® also offers various products such as:

- Water Flossers
- Sonic Electric Toothbrushes
- Combo Packs
- Spare Parts
- Accessories

For the complete product range, Visit: <https://oracura.in/>

Keep Your Product Up & Running, Always!

Accessories & Spare Parts To Buy with OC500



Standard Tip



Orthodontic Tip



Tongue Scraper



Pocket tip



Water Tank
Rubber Valve



Water Tank



Water Flosser®
Handle



USB Type-C
Cable



Protective
Pouch

Service Station:

ORACURA® Solutions Pvt. Ltd.

6th Floor, Unit 1, Plot No 233/2 Vidhani Kalpavruksha, MIDC Industrial Area, T.T.C. Industrial Area, Mumbai - Pune Hwy, Shiravane, Nerul, Navi Mumbai, Maharashtra, 400706

Visit Us: www.oracura.in Email: support@oracura.in



ISO9001 & ISO13485 Certified Manufacturer



Certificated Manufacturer