

ORACURA®

YOUR
ADVANCED
DAILY
ORAL CARE
ROUTINE



SB100

THANK YOU

For choosing ORACURA's SB100 Sonic Lite Electric Battery Operated Toothbrush. Please check if the package is intact before using it.

Index

Package Content	1
Parts of the toothbrush	2
How to attach the toothbrush Head	3
How to Select Modes	4
How to use the toothbrush	5
How to replace the battery	6
FAQs	7
Complementary Product Offers	8
Register Your Warranty	9
After Sales Support	10
Maintenance	11
Important instructions	12
FCC Statement	13
Service station	14

Package Content



Toothbrush Handle



Toothbrush Head



User Manual



Toothbrush Cap



2 AAA Batteries

Parts of the Toothbrush



How to Attach the Toothbrush Head



Set the brushhead into the handle slightly facing towards the right side.


Firmly press the head down on the handle cavity and rotate the brushhead towards the left side to lock it in place.

Getting Started with Brush Modes



- The first time you use the brush, press once for Clean, twice for Sensitive, thrice or Massage mode, and four times to power OFF your toothbrush.
- The second time you use the brush it will remember which mode you used last and it automatically starts on that mode.
- Clean Mode: for daily brushing.
- Sensitive Mode: For sensitive teeth & gums.
- Massage Mode: For gentle brushing & gum massage.

Using the Toothbrush

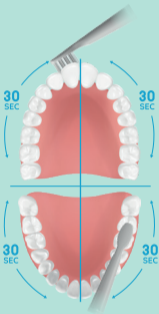

- Wet the bristles and apply toothpaste.
- Place the toothbrush in the mouth before turning it on.
- Turn on and glide the brush slowly along your teeth.
- Spend about 30 seconds per quadrant; the toothbrush will automatically turn off after 2 minutes of cleaning.
- Built-in smart timer pauses every 30 seconds to help you get an even clean.
- After use, rinse your toothbrush and let it air dry.
- Repeat twice daily or as recommended by your dentist.

Replacing the Battery



- Remove the battery cover by rotating it as shown (Rotate the battery cover towards the left)



- Insert 2 AAA batteries into the battery compartment correctly. Make sure the battery polarity matches the label shown)
- Attach the battery cover by rotating it as shown. (Attach and rotate towards the right to lock)

FAQs

1. What should I do if the toothbrush I just received cannot be switched on?

- It may be because of insufficient battery power. Kindly replace the batteries of the toothbrush before using it.

2. Why does the vibration decrease after a period of time?

- More power is consumed with the increase in the use time of the toothbrush and this results in a decrease in the vibration. Simply change the batteries and it will be restored.



Scan for Live Demo

You can schedule a live product demonstration by scanning the above QR code.

Complementary Product Offers

Go Beyond Just Brushing with Our Water Flossers.



If you have purchased an ORACURA® product, then this is your chance to explore various products complementing your purchased product. Scan the QR code for further details.

Visit: <https://oracura.in/collections/all-collections>

Register Your Warranty

WARRANTY
1 YEAR
WARRANTY

SB100 Sonic Lite Electric Battery Operated Toothbrush is covered under 1 year warranty period.



For warranty, register your product by scanning the QR code or visiting the link given below.

<https://oracura.in/pages/register>

1 Year of limited warranty is applicable for material, manufacturing, or operational defects existing at the time of delivery or for the occurrence of any operational error during regular use. Any mishandling of the product or alteration/repair done by an unauthorized person will void the warranty.

Consumer is required to submit the original purchase invoice as proof of purchase date and if requested, the entire product, to support a warranty claim.

If you face any problems while using our product or have any complaints regarding the same, we advise you to visit our return policy page.

<https://oracura.in/pages/return-policy>

After-sales Support

Thank you for choosing ORACURA to upgrade your oral care routine. We provide a one-year free warranty for the product and unopened toothbrush head. In case of any fault within 12 months from the date of purchase, you will be provided with a free repair after a thorough verification by ORACURA after-sales Service Center as a quality problem.

Any damage or defect which is not attributable to the product quality is not covered by the warranty, including but not limited to:

- Failure or damage caused by acts of God like earthquakes, lightning strikes, abnormal voltages, and other environmental factors.
- Dismantling, modifying, or replacing with non-original parts, or
- If the toothbrush is not operated according to the manual or damage is caused due to improper use, such as using the wrong current and voltage, plugging the product into an unqualified outlet, or other man-made damage.

If you need repair within the scope of the warranty, please fill up the warranty card. If we are unable to perform the obligation of free warranty service due to inaccurate information, the consumer shall bear the cost of the repair.

For details of service and warranty, please visit the official website. or send an email to support@oracura.in

Maintenance

- Wash the brush head and bristles after each use and allow them to dry naturally.
- Do not pull the bristles, otherwise, they will be easily deformed.
- Do not immerse the handle of the toothbrush in water.
- Do not clean the toothbrush using water higher than 40°C or an acidic cleaning solution.
- Clean the bottom cover of the toothbrush with a soft cloth.
- It is recommended to replace the toothbrush head every three months.

Model Title	SB100 Sonic Lite Electric Battery Operated Toothbrush
Color Variants	Blue, Peach, Grey, Green
Package contents	1 Toothbrush Handle, 1 Brush head, 1 Brush Head Cap, User Manual, 2 AAA Batteries
Modes	Clean, Soft & Massage
Specifications	36,000 Strokes/Minute IPX7 Waterproof 2 Minutes Auto-Timer 30 Seconds Interval Pause Reminder Dupont Bristles

Important Instructions

- Before discarding the toothbrush, the battery must be removed
- The battery must be removed after switching off the toothbrush.
- Disassembly should be done with care, under the guidance of professionals only
- The battery should be disposed of safely.

Steps to Remove the Battery:

- ➊ Remove the bottom cover
- ➋ Remove the movement support assembly
From the body shell.

FCC Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and radiate radio frequency energy and, if not installed and used in accordance with the instructions, it may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment
And receiver
- Connect the equipment into an outlet on a circuit
different from that to which the receiver is connected
- Consult the dealer or an experienced radio/ TV technician
for help

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Service Station:

ORACURA® Solutions Pvt. Ltd.
6th Floor, Unit 1, Plot No 233/2 Vidhani Kalpavruksha,
MIDC Industrial Area, T.T.C. Industrial Area, Mumbai -
Pune Hwy, Shiravane, Nerul, Navi Mumbai,
Maharashtra, 400706

Visit Us: www.oracura.in

Email: support@oracura.in



ISO9001 & ISO13485 Certificated Manufacturer



Certificated Manufacturer