

ORACURA®

YOUR
ADVANCED
DAILY
ORAL CARE
ROUTINE



SB200

THANK YOU

Thank you for purchasing the ORACURA SB200 Sonic Lite Electric Rechargeable Toothbrush is a professional high-quality device used for oral hygiene. The ORACURA SB200 Sonic Lite Electric Rechargeable Toothbrush has 3 operation modes that allow you to adjust the toothbrush to your individual needs. The toothbrush is made of a high-quality material, as a result of which brushing your teeth is even more enjoyable. In order to safely use the sonic technology offered by ORACURA, carefully read the instruction manual before using the toothbrush!

Index

Package Content	1
Parts of the toothbrush	2
How to attach the toothbrush Head	3
How to Select Modes	4
How to use the toothbrush	5
Charging the toothbrush	6
FAQs	7
Complementary Product Offers	8
Register Your Warranty	9
After Sales Support	10
Maintenance	11
Safety Rules	12
FCC Statement	13
Service station	14

Package Content



Toothbrush Handle



Toothbrush Head



User Manual



Toothbrush Cap



USB Charging Cable

Parts of the Toothbrush



How to Attach the Toothbrush Head



Set the brushhead into the handle slightly facing towards the right side.

Firmly press the head down on the handle cavity and rotate the brushhead towards the left side to lock it in place.

Getting Started with Brush Modes




- The first time you use the brush, press once for Sensitive, twice for Clean, thrice for massage mode, and four times to power OFF your toothbrush.
- The second time you use the brush it will remember which mode you used last and it automatically starts on that mode.
- Sensitive Mode: For sensitive teeth & gums.
- Clean Mode: for daily brushing.
- Massage Mode: For gentle brushing & gum massage.

If you face any problems while using our product or have any complaints regarding the same, we advise you to visit our return policy page.

<https://oracura.in/pages/return-policy>

Using the Toothbrush

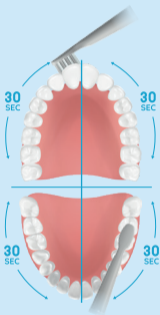

- Brushing your teeth - slowly and gently. Glide the head of the toothbrush on the teeth, above the gum line, in forward and backward motions, so that the longer bristles are able to reach the gaps between the teeth.
- The toothbrush should not be used in the same way as a traditional toothbrush.
- While brushing it is necessary to regularly gently press the toothbrush against the teeth, which will increase the efficiency of brushing.
- The ORACURA's SB200 Sonic Lite Electric Rechargeable Toothbrush has a two-minute timer, triggering an intermittent signal every 30 seconds in order to remind the user to brush the four sections of the mouth.

Charging the toothbrush

Before first use, the toothbrush must be charged for 6 hours.
During the charging on LED will light up.



- **Input Power:** DC5V/1A
- **Charging Time:** max. 6 hours
- **Operating Voltage:** 3.7V
- **Battery Capacity:** 600mAh



USB Charging Cable

FAQs

1. What should I do if the toothbrush I just received cannot be switched on?

- It may be because of insufficient battery power. Kindly charge the toothbrush before using it.

2. Why is the working time of the toothbrush less than described?

- The working time of the toothbrush is subject to usage. 25 days of battery life is endorsed, given the toothbrush is used for a period of no more than 2 minutes, twice a day.

3. Why does the vibration decrease after a period of time?

- More power is consumed with the increase in the use time of the toothbrush and this results in a decrease in the vibration. Simply recharge the battery and it will be restored.



Scan for Live Demo

You can schedule a live product demonstration by scanning the above QR code.

Complementary Product Offers

Go Beyond Just Brushing with Our Water Flossers.



If you have purchased an ORACURA® product, then this is your chance to explore various products complementing your purchased product. Scan the QR code for further details.

Visit: <https://oracura.in/collections/all-collections>

Register Your Warranty

WARRANTY
1 YEAR
WARRANTY

ORACURA's SB200 Sonic Lite Electric Rechargeable Toothbrush.



For warranty, register your product by scanning the QR code or visiting the link given below.

<https://oracura.in/pages/register>

- Damage resulting from random events (electrical fault, fire, inundation, flooding, etc.).
- Mechanical damage caused during the use of the device by the holder of the warranty, and any defects caused by them, Damage resulting from the user's fault.
- Damage resulting from the use of consumables not recommended by the manufacturer in the instruction manual.

The warranty provider ensures the purchaser, that the device to which this warranty applies, has been designed and manufactured on the basis of high-quality standards. ORACURA SB200 Sonic Lite Electric Rechargeable Toothbrush is covered by a 12 months warranty period from the date of purchase.

The warranty is valid from the date of purchase, certified by proof of purchase, which should be attached to the warranty card. The breach of the warranty seals, or any other interference of persons other than authorized service personnel.

After-sales Support

Thank you for choosing ORACURA to upgrade your oral care routine. We provide a one-year free warranty for the product and unopened toothbrush head. In case of any fault within 12 months from the date of purchase, you will be provided with a free repair after a thorough verification by ORACURA after-sales Service Center as a quality problem.

Any damage or defect which is not attributable to the product quality is not covered by the warranty, including but not limited to:

- Failure or damage caused by acts of God like earthquakes, lightning strikes, abnormal voltages, and other environmental factors.
- Dismantling, modifying, or replacing with non-original parts, or
- If the toothbrush is not operated according to the manual or damage is caused due to improper use, such as using the wrong current and voltage, plugging the product into an unqualified outlet, or other man-made damage.

If you need repair within the scope of the warranty, please fill up the warranty card. If we are unable to perform the obligation of free warranty service due to inaccurate information, the consumer shall bear the cost of the repair.

For details of service and warranty, please visit the official website. or send an email to support@oracura.in

Maintenance

- Wash the brush head and bristles after each use and allow them to dry naturally.
- Do not pull the bristles, otherwise, they will be easily deformed.
- Do not immerse the handle of the toothbrush in water.
- Do not clean the toothbrush using water higher than 40°C or an acidic cleaning solution.
- Clean the bottom cover of the toothbrush with a soft cloth.
- It is recommended to replace the toothbrush head every three months.

Model Title	SB200 Sonic Lite Electric Rechargeable Toothbrush
Color Variants	Blue, Peach, Grey, Green
Package contents	1 Toothbrush handle, 1 Brush head, 1 Brush Head Cap, User Manual, USB Charging Cable
Modes	Clean, Sensitive & Massage
Specifications	36,000 Strokes/Minute IPX7 Waterproof 2 Minutes Auto-Timer 30 Seconds Interval Pause Reminder DuPont Bristles

Safety Rules

- Before connecting the device, make sure that the voltage indicated on the device is compatible with the voltage in the electrical network.
- Do not leave the charger cable connected to a power source. Only use the toothbrush in accordance with the intended purpose described in this instruction manual. Do not use any accessories that are not recommended by the manufacturer.
- Do not attempt to open, repair, or modify the product without authorization.
- If the toothbrush is dropped, the head of the toothbrush should be replaced before the next use even if there are no visible signs of damage.
- Do not use the charging base and the cable outdoors or near sources of heat.
- If the device does not work, do not attempt to disassemble it or to replace any part of the toothbrush as this could cause the risk of electric shock. Contact your technical support hotline.
- Keep the charger in a dry and safe place, when it is not used.

FCC Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and radiate radio frequency energy and, if not installed and used in accordance with the instructions, it may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment
And receiver
- Connect the equipment into an outlet on a circuit
different from that to which the receiver is connected
- Consult the dealer or an experienced radio/ TV technician
for help

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Service Station:

ORACURA® Solutions Pvt. Ltd.
6th Floor, Unit 1, Plot No 233/2 Vidhani Kalpavruksha,
MIDC Industrial Area, T.T.C. Industrial Area, Mumbai -
Pune Hwy, Shiravane, Nerul, Navi Mumbai,
Maharashtra, 400706

Visit Us: www.oracura.in

Email: support@oracura.in



ISO9001 & ISO13485 Certificated Manufacturer



Certificated Manufacturer