

ORACURA®

YOUR
ADVANCED
DAILY
ORAL CARE
ROUTINE



SB300

THANK YOU

for choosing ORACURA's SB300 Sonic Smart Electric Rechargeable Toothbrush. Please check if the package is intact before using it.

Index

Package Content

Attaching The Toothbrush Head

Parts Description

How to Select the Mode

Select the Intensity

Power Off

Charging

Service Station

FCC Statement

Important Instructions

Maintenance

After-sales Support

Register Your Warranty

Complementary Product Offers

FAQs

Package Content



Toothbrush Handle



Toothbrush Head



User Manual




USB Type - C Cable

How to Attach the Toothbrush Head



Set the brush head into the handle and ensure both face the same side. Firmly press the head down on the handle cavity and it will be attached.

Parts of the Toothbrush

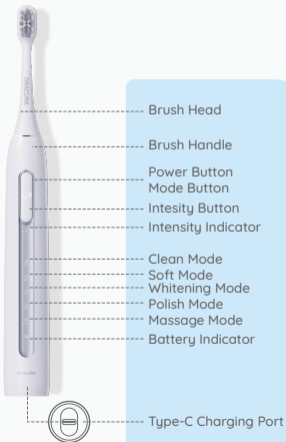

How to Select the Mode



Press the Power button to switch the toothbrush on. Once it's switched on, press the Power button again to change mode, the mode will be selected if there is no further operation in 5 seconds.

How to Select the Intensity

Power on the toothbrush and press the intensity button to choose from three different intensity levels. The first indicator is for low intensity, the second one for medium intensity, and the third indicator is for the highest intensity.

How to Switch Off the Toothbrush



- Press and hold the power button for about 1.5 seconds while brushing.
- After all the modes are switched in a sequence, press the On/Off button once to power off.
- Automatically power off while charging.



Scan for Live Demo

You can schedule a live product demonstration by scanning the above QR code.

Charging



- If the indicator keeps flashing after switching on the toothbrush, it needs to be recharged; if the indicator does not blink, it is fully charged.
- Ensure to wipe the bottom and keep it dry before charging. Please keep the bottom cover closed after charging.



FAQs

1. What should I do if the toothbrush I just received cannot be switched on?

- It may be because of insufficient battery power. Kindly charge the toothbrush before using it.

2. Why is the working time of the toothbrush less than described?

- The working time of the toothbrush is subject to usage. 90 days of battery life is endorsed given the toothbrush is used for a period of no more than 2 minutes, twice a day.

3. Why does the vibration decrease after a period of time?

- More power is consumed with the increase in the use time of the toothbrush and this results in a decrease in the vibration. Simply recharge the battery and it will be repair.

Complementary Product Offers

Go Beyond Just Brushing with Our Water Flossers.



If you have purchased an ORACURA® product, then this is your chance to explore various products complementing your purchased product. Scan the QR code for further details.

Visit: <https://oracura.in/collections/all-collections>

Register Your Warranty

ORACURA's SB300 Sonic Smart Electric Rechargeable Toothbrush is covered under 1 year warranty period.



For warranty, register your product by scanning the QR code or visiting the link given below.

<https://oracura.in/pages/register>

1 year of limited warranty is applicable for material, manufacturing, or operational defects existing at the time of delivery or for the occurrence of any operational error during regular use. Any mishandling of the product or alteration/repair done by an unauthorized person will void the warranty.

Consumer is required to submit the original purchase invoice as proof of purchase date and if requested, the entire product, to support a warranty claim.

If you face any problems while using our product or have any complaints regarding the same, we advise you to visit our return policy page.

<https://oracura.in/pages/return-policy>



After-sales Support

Thank you for choosing ORACURA to upgrade your oral care routine. We provide a one-year free warranty for the product and **unopened toothbrush head**. In case of any fault within 12 months from the date of purchase, you will be provided with a free repair after a thorough verification by ORACURA after-sales Service Center as a quality problem.

Any damage or defect which is not attributable to the product quality is not covered by the warranty, including but not limited to:

- Failure or damage caused by acts of God like earthquakes, lightning strikes, abnormal voltages, and other environmental factors.
- Dismantling, modifying, or replacing with non-original parts, or
- If the toothbrush is not operated according to the manual or damage is caused due to improper use, such as using the wrong current and voltage, plugging the product into an unqualified outlet, or other man-made damage.

If you need repair within the scope of the warranty, please fill up the warranty card. If we are unable to perform the obligation of free warranty service due to inaccurate information, the consumer shall bear the cost of the repair

For details of service and warranty, please visit the official website. or send an email to support@oracura.in.

Maintenance

- Wash the brush head and bristles after each use and allow them to dry naturally.
- Do not pull the bristles, otherwise, they will be easily deformed.
- Do not immerse the handle of the toothbrush in water.
- Do not clean the toothbrush using water higher than 40°C or an acidic cleaning solution.
- Clean the bottom cover of the toothbrush with a soft cloth.
- Wipe the plug with a dry cloth to avoid short circuits.
- It is recommended to replace the toothbrush head every three months.
- When you use it for the first time, make sure the battery is exhausted before charging.
- If you have not used it for more than 6 months, charge it for 3-5 hours before using it again.
- During use, if the charging indicator flashes continuously for a few seconds and then stops, it indicates that the battery is low.
- The indicator flashes slowly when it is being charged, and stays on after it is fully charged. If it has not been charged for a long time, it is normal for the indicator light to take 6-10 minutes to flash again.

Important Instructions

- Before discarding the toothbrush, the battery must be removed
- The battery must be removed after switching off the toothbrush.
- Disassembly should be done with care, under the guidance of professionals only
- The battery should be disposed of safely.

Steps to Remove the Battery:

- Remove the bottom cover
- Remove the movement support assembly from the body shell
- Heat the battery pins with electric ferrochrome to remove the battery

FCC Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and radiate radio frequency energy and, if not installed and used in accordance with the instructions, it may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and
(2) this device must accept any interference received, including interference that may cause undesired operation.

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment and receiver
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected
- Consult the dealer or an experienced radio/ TV technician for help

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Service Station:

ORACURA® Solutions Pvt. Ltd.
6th Floor, Unit 1, Plot No 233/2 Vidhani Kalpavruksha,
MIDC Industrial Area, T.T.C. Industrial Area, Mumbai -
Pune Hwy, Shiravane, Nerul, Navi Mumbai,
Maharashtra, 400706

Visit Us: www.oracura.in

Email: support@oracura.in



ISO9001 & ISO13485 Certificated Manufacturer



Certificated Manufacturer