




QuickPore™ Daily Moisturizing and Lifting Essential NanoPen Lux Device Kit

BD-DM-K707L




Step 1
Charge NanoPen Lux before use

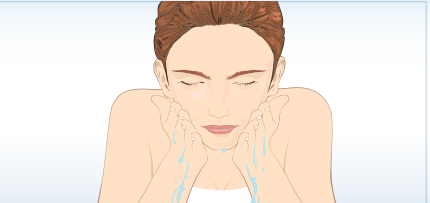


Step 2
Begin with freshly washed skin




Step 3
Install Daily NanoChip to NanoPen Lux

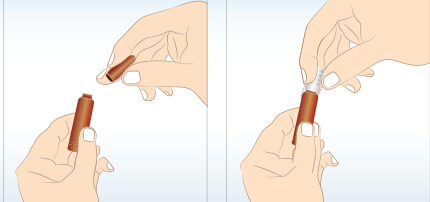

Step 4
Open serum, Place rubber head over the bottle



Step 5
Apply Deminuage Blanq Serum

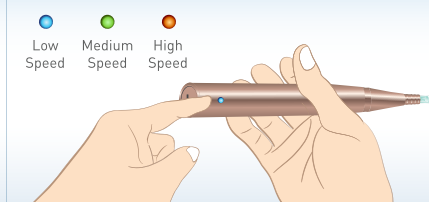

Step 6
Turn on the NanoPen Lux by pressing and holding down the indicator light

Low Speed Medium Speed High Speed



Step 7
Use NanoPen Lux with Daily NanoChip

Gliding Stamping



Step 8
Remove and Throw Daily NanoChip

TRASH CAN