




QuickPore™ NanoPen Prime Acne Scar Removal and Acne Spot Treatment Kit

BD-DM-K709




Step 1
Charge NanoPen Lux before use

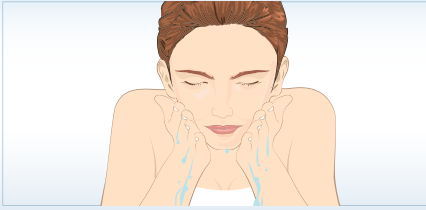


Step 2
Begin with freshly washed skin




Step 3
Install Booster NanoChip to NanoPen Lux

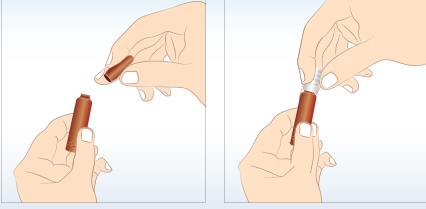

Step 4
Open serum, Place rubber head over the bottle




Step 5
Apply Deminuage Blanq Serum

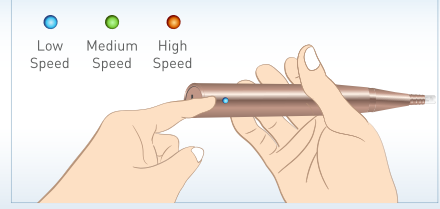

Step 6
Turn on the NanoPen Lux by pressing and holding down the indicator light



Step 7
Use NanoPen Lux with Booster NanoChip

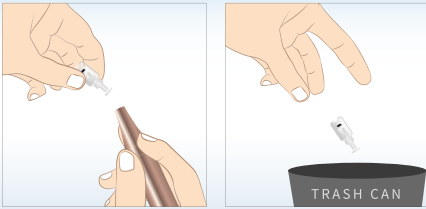

Step 8
Remove and Throw Booster NanoChip