

Deminuage

Skincare



Boost
your
inner
glow




NanoPen Lux Kit
User Guide



1. Be sure to charge your NanoPen before use. The micro USB rechargeable port is located at the bottom of the pen. The indicator light will remain red while charging and will turn off once the battery is completely full. For safety purposes, remove from the charger before use.



2. Begin with freshly washed skin.

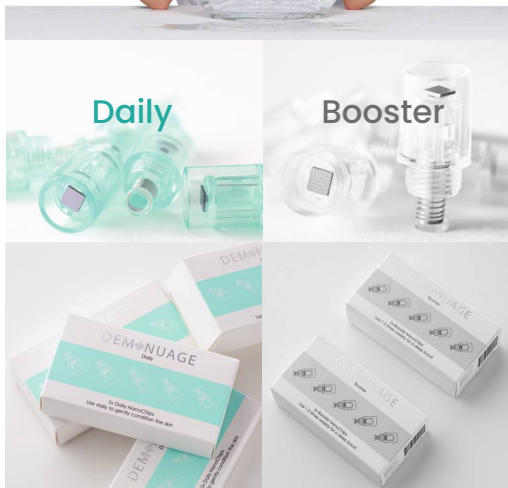

3. Choose NanoChip and apply serum

Daily NanoChip: Use daily with Deminuage Blanq Brightening Serum or your own serum and apply it evenly to your face or the area you'd like to focus on.

Booster NanoChip: Use 1-2 times weekly (in place of Daily NanoChip) with included Deminuage Blanq Brightening Serum.



Note:
How to apply Deminuage Blanq Brightening Serum


a. Open serum by placing hands on each end of the bottle and twisting along the bottleneck.

b. Place rubber head over the bottle, ensuring there's no gap between the dipper and bottle.

c. Apply serum by holding bottle vertically to face and gently tapping the rubber head against skin

4. Use NanoPen with NanoChip

Remove a new NanoChip cartridge from its sealed bag and screw into the NanoPen with the NanoChip facing outward. Turn on the NanoPen by pressing and holding down the indicator light until NanoPen begins to vibrate. To switch speeds, quickly press down the indicator light (do not hold).

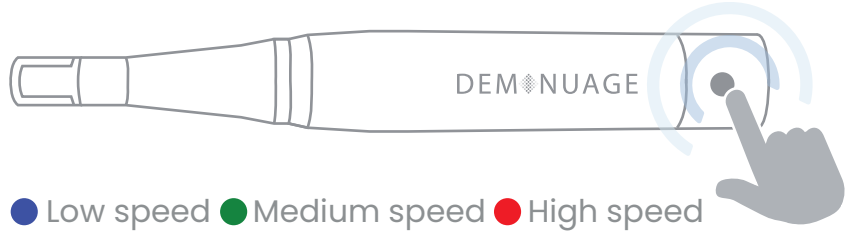

Around the eyes and nose: Use **stamping** motion with **low speed** (blue).



Stamping: Hold the pen perpendicular to the eye and nose area. Lightly press the pen down against the skin for 2 to 3 seconds at a time and release (in a stamping motion).

Over the face: Use **gliding** motion with **medium speed** (green) or **high speed** (red).



Gliding: Hold the NanoPen perpendicular to the face and gently begin sliding the NanoChip onto the skin, moving from inside to outside and bottom to top.

*Please note that high speed is for experienced users only. It doesn't improve performance but shortens treatment time.

5. Continue until the desired area has been treated.

6. Turn off by once again pressing and holding down the light indicator until the NanoPen stops vibrating.

7. Unscrew and dispose of the used NanoChip cartridge. Do not reuse. Then clean the pen using a damp cloth. Do not use harsh chemicals or alcohol.

