

Deminuage

Skincare



Boost
your
inner
glow

NanoPen Prime Kit
User Guide



1. For first use, open the battery compartment by unscrewing the middle portion of the NanoPen. Install one AAA Battery in the direction indicated above.



2. Begin with freshly washed skin.


3. Choose NanoChip and apply serum

Daily NanoChip:

Use daily with included Deminuage Blanq Brightening Serum or your own serum and apply it evenly to your face or the area you'd like to focus on.

Booster NanoChip:

Use 1-2 times weekly (in place of Daily NanoChip) with included Deminuage Blanq Brightening Serum.

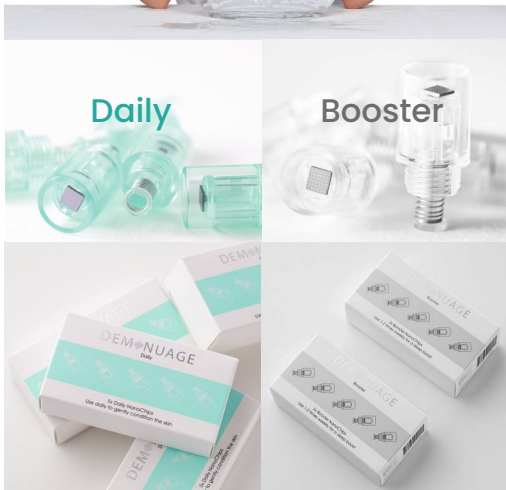

Note:

How to apply Deminuage Blanq Brightening Serum

a. Open serum by placing hands on each end of the bottle and twisting along the bottleneck.

b. Place rubber head over the bottle, ensuring there's no gap between the dipper and bottle.

c. Apply serum by holding bottle vertically to face and gently tapping the rubber head against skin.

4. Use NanoPen with NanoChip

Remove a new NanoChip cartridge from its sealed bag and screw into the NanoPen with the NanoChip facing outward. Turn on the NanoPen by pressing the power button located at the bottom of the NanoPen.



Around the eyes and nose: Use stamping motion.



Stamping: Hold the pen perpendicular to the eye and nose area. Lightly press the pen down against the skin for 2 to 3 seconds at a time and release (in a stamping motion).

Over the face: Use gliding motion.



Gliding: Hold the NanoPen perpendicular to the face and gently begin sliding the NanoChip onto the skin, moving from inside to outside and bottom to top.

5. Continue until the desired area has been treated.

6. Turn off by once again pressing the power button.

7. Unscrew and dispose of the used NanoChip cartridge. Do not reuse. Then clean the pen using a damp cloth. Do not use harsh chemicals or alcohol.