



MONIQUE

INITIAL ASSEMBLY

Initial Assembly

Parts List & Tools

Your Monique comes partially assembled. You will need an M6 Allen to finish assembly. Your ski will come with a frame, seat, footrest, handle, shock, gas struts as well as any additional spare parts you ordered.

For Steps 1-6 please also see the Initial Assembly Video on the product manuals page of our website.

1) Shift the Lever to the “Up” position and pull the A-Arms down.


2) Install the Gas Struts.

The left and right Gas Struts have been labelled with corresponding red and green stickers.

- Remove the stickers and the lock pins from the Ball Sockets.
- Thread the Gas Struts into the Eyelets on the Frame.
- Snap the Ball Sockets onto the Ball Studs.
- Secure the Lock Pins onto the Ball Sockets.

3) Remove the rear Shock Bolt and Shock Spacers.

You will need an M6 Allen.

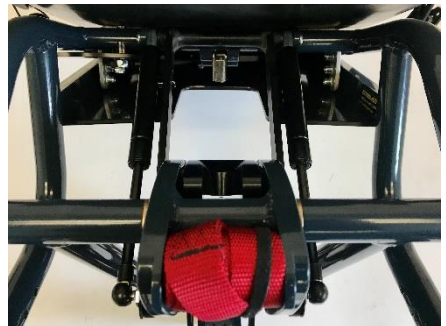

4) Remove the front Shock Bolt.

Shift the Lever to the “Load / Ski” position and push the A-Arms up to expose the front Shock Bolt in the Rocker Link. You will need an M6 Allen.



5) Install the front eyelet of the Shock.

Feed the Shock into the Rocker Link from the back side of the Load Link. The rebound knob must face down. Fasten with an M6 Allen.



6) Install the rear eyelet of the Shock to the Load Link.

Pull the A-Arms down to align the rear eyelet. Install the Shock Spacers between the Shock and the Load Link with the wide side of the flanges facing out.



7) Compress the Monique into the lowered skiing position.

Compress with the Shifter Lever in the “Load / Ski” position.



8) Set up the Seat.

See the image to the right for proper setup.

- The waist and thigh use pads secured with a ratcheting ladders.
- The seat angle adjustment has two extension buckles at the front of the thighs secured to two ratcheting ladders from the seat back.
- On the seat back is a Velcro limiter strap. Loop it under the webbing guide and back over onto itself.

**9) Remove and discard the wooden dowels from the foot and handle tubes.**

The pin does not go through the dowel. Remove the dowel first and then remove the pin.

**10) Install the Footrest.**

Be careful to slide each side in straight and at the same time. Slide / tap the footrest into the tubes. If necessary, a rubber mallet can be used to overcome initial friction. If necessary, tighten the ankle angle adjustment clamps.

**11) Install the handle.**

**ENABLING
TECHNOLOGIES**

enablingtech.com

Enabling Technologies // 303.578.9345
2200 S Jason St Denver, CO 80223
info@enablingtech.com