

eTherm[®]

Infrared Ear & Forehead Thermometer

Thermometer
Quick User Guide




Attention


- 1 The optimum calibrated measuring position of the thermometer is placement of the probe in the center of the forehead and ear canal. Please do not measure on other body parts.
- 2 When changing environments (cold to hot) for the thermometer please wait 30 minutes for unit to reach an equilibrium steady state.
- 3 Baby crying and bathing can cause body temperature fluctuations- for accurate measurements please do not measure immediately after crying or bathing.
- 4 If the forehead is covered with bruises, ointments or perspiration, local skin temperature can be abnormal. Please keep the forehead and ear canal exposed and clean for at least 30 minutes before measuring.
- 5 Please do not measure in a strong electro-magnetic interference environment such as near an active mobile phone.

Please read this Quick User Guide before using the device and read the Instruction Manual for the full and comprehensive instructions of use.


Forehead

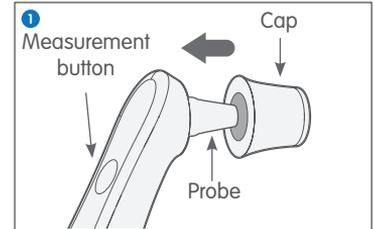

Put the thermometer cap on the unit.

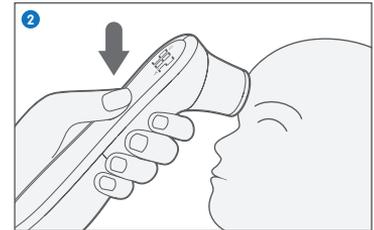


Aim the thermometer probe at the center of the forehead and lightly touch the skin. **Press and release** the "measurement button"

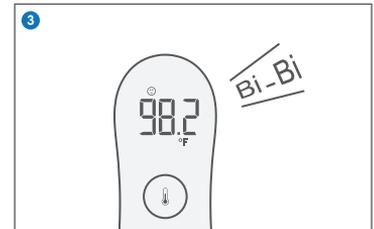


When hearing the beeper sounds, the measurement process is completed and you can check the measurement result.

(The thermometer will automatically shutdown after 8 seconds.)


Ear

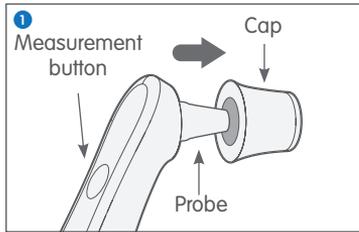

Remove the thermometer cap from the unit.

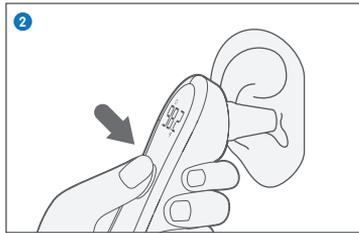

Carefully and gently position the thermometer probe into the ear canal and **press** the "measurement button"

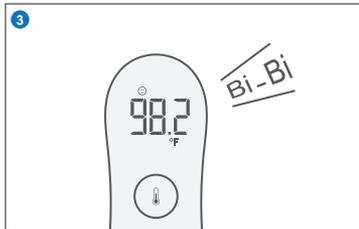


When hearing the beeper sounds, the measurement is completed and you can check the measurement result.


(The thermometer will automatically shutdown after 8 seconds.)

Battery Replacement

- 1 Press the battery door and slide it down to open the battery door.

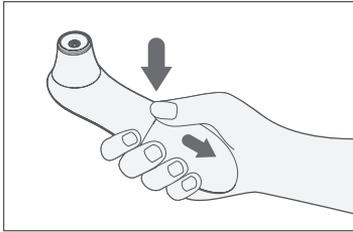

- 2 Take out the old batteries and replace them with new ones.
- 3 Insert the new batteries; pay attention and follow the polarity symbols on the plastic casing to avoid incorrect installing of the batteries.
- 4 Replace back the battery door.

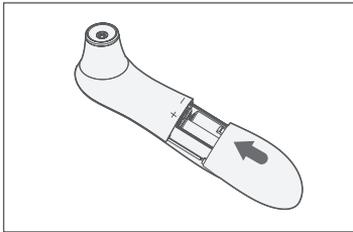

- Please follow local laws for disposing the used batteries.
- Please do not throw the batteries directly to the garbage and dispose of them in accordance with local laws for used batteries.
- Please take out the batteries from the unit if the device is not being used for long periods of time. (say for 3 months).
- Please do not put the batteries into an open fire to prevent an explosion.

Temperature Reading Guide

The temperature of the human body has a normal range but may vary between people. Temperature of individuals may also vary from time to time. The following table is for reference purposes only and indicates the average temperature ranges for the human body.

Ampit	94.5-99.1 °F	34.7-37.3 °C
Oral	95.9-99.5 °F	35.5-37.5 °C
Anus	97.9-100.4 °F	36.6-38.0 °C
Ear	95.9-100.0 °F	35.5-37.8 °C

Self-diagnosis and self-treatment based on the measurement results could be dangerous; please contact and consult a medical professional for advise and required recommended follow-up action.

On-line product support at:
www.elepho.com/help-support/

E-mail support via: customercare@elepho.com

Telephone support via Elepho Consumer Relations at: 1-833-435-3746 (1-833-4-ELEPHO)

General information regarding Elepho full product line via: www.elepho.com

Manufactured for:
Elepho Inc, 466 Central Ave (2nd floor), Cedarhurst, New York 1 516. USA.

Website: www.elepho.com
Item #: FDIR-V4
Product/model name: eTherm

This product design is protected and registered under the EUIPO patent and trademark office