Rhino[®]

180 kg / 100 g
400 lb / 0.2 lb

Capacidad Máxima / División Mínima
Maximum Capacity / Minimum Division
Capacidade Máxima / Divisão Mínima

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- **Conexión bluetooth 4.0**

Bluetooth 4.0 connection / Com conexão bluetooth 4.0

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Download the app scanning the QR code

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US You have purchased a RHINO product, which guarantees quality and durability, since the BABAIN-180NG scale was manufactured according to American quality standards. Keep this instruction manual, it will help you understand the operation of the RHINO BABAIN-180NG will be useful for future reference.

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You have purchased a RHINO product, which guarantees quality and durability, since the BABAIN-180NG scale was manufactured according to American quality standards. Keep this instruction manual, it will help you understand the operation of the RHINO BABAIN-180NG will be useful for future reference.

TECHNICAL SPECIFICATIONS

Max Capacity: 180 kg / 400 lb

Division: 100 g / 0.2 lb

Power: 2 batteries AAA (1,5 V included)

Temperature: 32 °F to 284 °F

QUICK USE GUIDE

1. Remove the scale from its packaging, make sure it's in perfect conditions and that includes all its accessories (manual and 2 AAA batteries).
2. Place the scale on stable and leveled surface, be sure to place batteries correctly any movement no matter how small can affect accuracy.
3. Place your feet in the center of the scale.

PRECAUTIONS

1. The display will show "Lo" when it is necessary to replace the batteries.

2. The display will show "Err" when the scale is overloaded.
3. Remove excess weight to avoid damage.
4. Place the scale in a cool, dry place.
5. For domestic use only.
6. DO NOT USE for commercial purposes.

WARNINGS

-  Do not open the scale yourself or the warranty will be void.
-  To clean the scale use a damp cloth or sponge, never a direct stream of water.
-  Protect the scale from insects, do not spray, only solid products against insects.
-  Do not drop the scale suddenly.
-  Do not direct air flow from fans or heaters directly to the scale during use.
-  Do not expose your scale to dripping or splashing liquids.

If you need technical assistance, go to one of the authorized Rhino service centers.

WARNING

DO NOT use with medical implants as pacemakers.

DO NOT stand on the edge of the scale or jump on it.

DO NOT overload the scale (MAX 180 kg / 400 lb).

DO NOT drop the scale or objects on it, as this could damage the sensors.

NOT recommended for women pregnant.

Any measurement obtained with this device is for reference only and should not be considered as a medical opinion.

Always place the scale on a hard, dry and flat surface before of measurement

Make sure your feet are dry before stepping on the scale.


Remove batteries if scale is not used for an extended period of time.

The user under 6 years old, just use weight function without BMI (body mass index) and body fat, data.

The user over 10 years old, just use only weight and BMI function, no body fat, data etc.

BATTERY PLACEMENT


1. Open the lid underneath the scale as shown in the diagram.
2. Install the batteries taking into account the polarity indicated in the compartment.
3. Replace the battery compartment cover.

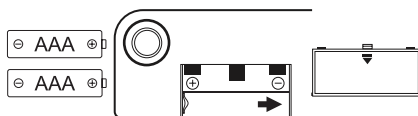OPERATION

1. To ensure the accuracy of the scale, place it on a hard, flat surface.

Note. Avoid placing it on a slippery surface.

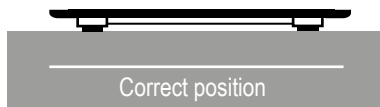
2. Gently step on the scale, it will turn on automatically. Stay on the scale without

move until the weight appears in the display, when the numbers on the display flash 3 times and then remain fixed, it will indicate its weight.

WEIGHT UNIT CHANGE

The unit change button is located underneath the scale **[UNIT]**, as shown in the diagram. With the button **[UNIT]** you can change the unit of measurement, kilograms (**kg**) to pounds (**lb**) and vice versa.

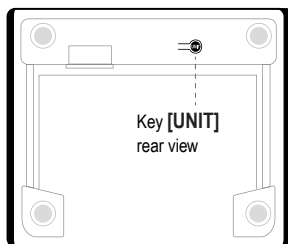INDICATION SYMBOLS

Lo	Low battery; Change batteries.
Err	Overload.
C	An error has occurred during measurement; Get off the scale and repeat the process.



DOWNLOAD THE APP

To use FITDAYS app you must to download the app to your smartphone. You can do it from Google Play or App Store.

NOTE: The app is compatible only with IOS 8.0 or higher Android 5.0 or higher operating systems and only supports one smartphone device at a time

1. After downloading the app, check that your bluetooth is turned on.
2. Enter the app from your phone or tablet and register your email.

3. Read and accept the privacy policies.
4. Add personal data like sex, age, height etc (Enter incorrect data may result in inaccurate measurements).

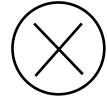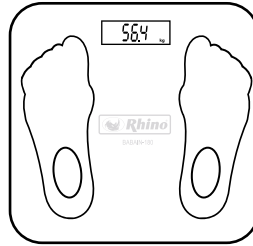
HOW TO USE THE APP

Kick the scale to activate the display. Wait until the display shows 0.0 (kg / lb).

With your app open on the home screen, place your bare feet on top of the metal bars.

NOTE: You must place your bare feet on the metal bars.

If you step on the scale with socks or shoes, only the weight and BMI (body mass index) will be displayed on the main page.

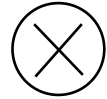The numbers in the screen will be blinking while the scale is analyzing data.

Keep your feet on the scale till the weight showed in the app be stable then the scale will show the obtained information.

you will find the result such as BMI, visceral fat, body water skeletal muscle, muscle mass, body age, protein, subcutaneous fat, muscle frequency, fat loss and BMR. Also you can see the weekly, monthly and year result.

USER


Add or change users by clicking on my profile or click on new user.

The scale and the app recognize the variation of results, and gives you the option to change the user if the data is different from the selected user.

NOTE. Each account can support 24 users.

BABY WEIGHING MODE

1. Click options, then select **Baby Mode** in the app.
2. Get on the scale so the app is determined its weight before carrying the baby.

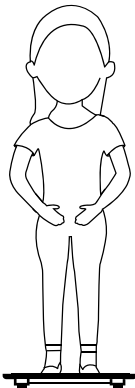3. After I register your weight, hold the baby and wait stabilize the scale to get the baby's weight on the application.

TEAM BODY PROBLEMS AND SOLUTIONS

PROBLEM	SOLUTIONS
Bluetooth not connected	<p>The scale already has bluetooth activated once turned on.</p> <p>Check the IOS version of your smartphone must be IOS higher than 8.0 or higher / Android 5.0 or higher.</p> <p>Restart your mobile device.</p> <p>Remove the batteries from the scale and put them back.</p> <p>Verify that in your phone Bluetooth connection is turned on.</p> <p>Verify that location is activated on phone for a correct connection to the App.</p> <p>NOTE: The scale only supports one smartphone device at a time.</p>
Body fat data is not measured	<p>User must be barefoot to measure body fat.</p> <p>Make sure the scale and your feet are dry.</p>
The scale is producing inconsistent weights	<p>Check if the scale is on a hard flat surface.</p> <p>Check each foot of the sensor to make sure nothing is stuck to the bottom of the sensor.</p>
The app shows the wrong ratio of muscle	<p>The body parameter setting is incorrect, check to make sure I correctly enter the sex, age and height data.</p>
The scale does not turn on.	<p>Remove the batteries and put them back. be sure to insert the batteries into the correct address.</p> <p>Carry out the battery change.</p>

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WARRANTY

Rhino Maquinaria warrants to customers who purchase any Rhino products from an authorized seller that Rhino products will be free of defects in workmanship for the duration of 1 year from the date of purchase.

During the 1 year warranty period, Rhino Maquinaria will replace the product that has a manufacturing defect as long as it has been determined by Rhino Maquinaria.

This warranty does not apply to defects caused by misuse or alterations of the product, nor does it apply if the damage was caused by repairs or service provided or attempted by anyone other than Rhino Maquinaria.

To make this warranty effective, you must first contact our Customer Service by dialing (619) 393-1079 or by mail to info@rhinoglobal.net, once you have contacted our Customer Service they will request from you a warranty claim by email.

Rhino Maquinaria will review the warranty claim. If the warranty claim is found to be valid, please return the product with a copy of the invoice and order ID to the of the invoice and order ID to the indicate address.

Upon receiving the product, Rhino Maquinaria will send a replacement product to customer's address.

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Rhino Maquinaria S.A. de C.V.