


65mm

112mm

正面

GLO24K 6-IN-1 Beauty Therapy Wand for the Face, Eyes & Neck with Dual LED, Thermal, Vibration, and Micro-Current Technologies  
Item #: BD5

**USER MANUAL**

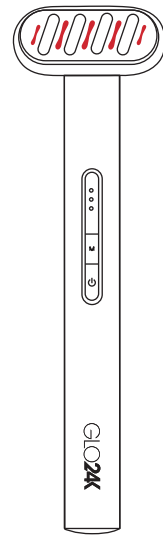

Please read this manual carefully before use and keep it properly for future reference

—1—

Thank you for purchasing this product!

**Product description**

**Main functions:**

**Red Light Therapy**  
Red LED Light Therapy with wavelengths of 600-650 nm, which can help with stimulating blood circulation and can generate anti-aging benefits.

**Blue Light Therapy**  
Blue Light Therapy with wavelengths of 450-500 nm, which can help with multiple skin disorders and assist in detoxifying and balancing the skin.

**Micro-Current Technology**  
Uses weak current to muscles and skin cells to help stimulate the skin, increase movement and flexibility, promote blood circulation, increase skin elasticity, and can help in achieving the visible effect of skin tightening.

**Thermal Therapy**  
The massage head automatically and gently heats up. This can help in accelerating blood circulation and enhancing skin vitality.

—2—

**90°Rotating Massage Head**  
Innovative rotatable and adjustable massage head with up to 90° (T shape) tilt design. Convenient for sliding and massaging on the forehead, around the eyes, cheeks, jaws, neck, and various parts of the face.


**Adjustable Vibration at 3 Gears**  
Built-in high-frequency adjustable vibration and massage generator in 3 gears: low/medium/high to stimulate the skin and help in making it look more firm and radiant.

**Product features**

Micro Current   3 Gears   Rotating/Adjusting Massage Head   Dual mode Light Switching   High Frequency Vibration

—3—

**An illustrative**

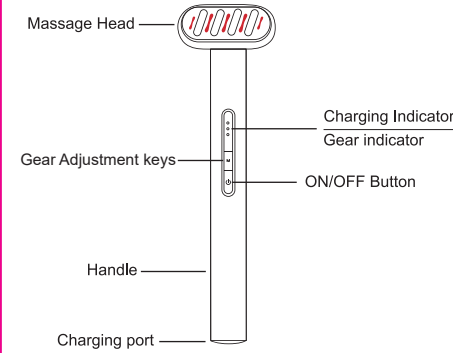

**How to use**

1. Apply a Face or Eye Serum (preferably GLO24K) until full absorption. (This step is not mandatory but recommended for best results).
2. Turn ON: A long press on the bottom ON/OFF button will turn on the wand with Light, Micro-Current, Thermal, and Vibration functions all set to turn ON.
3. Gently slide the massage head over areas of concern (Face, Neck, and around the Eyes).
4. You can rotate and adjust the massage head counter-clock wise up to 90°, so that the massage head is in "T" shape.

—4—

反面

5. To change from Red Light to Blue Light apply a short press on the ON/OFF button while the product is ON. You can change between lights at any time during the application.
6. Click the top Gear Adjustment Key to change intensity of both vibration and micro-current that change together (low/medium/high) indicated by 1/2/3 indicating lights respectively.
7. Enjoy 5-10 minutes of therapeutic relaxing massage.
8. Shut Down: long press on the ON/OFF button to turn off.
9. After use, wipe the massage head with a clean cloth. (Do not rinse with water as it is not water resistant)
10. NOTE: after turning on, if you do not make any key operation for 5 minutes, the product will automatically shut-off. Turn it back ON to resume operation.

**Charging Instructions**

1. Plug the charging USB cable into an output USB charger (charging source) and plug the charging cable end into the charging port at the bottom of this product for charging. The product must be OFF for proper charging and should NOT be turned ON while charging. The charger is designed to charge this product while NOT being all the way into the charging port.
2. When the indicator light flashes (slow flash), it is in a charging state.
3. When the indicator light is steady and stops flashing, it means that the power is full.
4. If the indicator light blinks several times in a row it signals low power. Please charge the product in time, and between applications otherwise, it will stop working.
6. Do not operate this product while being charged. To turn ON disconnect from the charger.

**NOTE: The charger does NOT go all the way into the charging port to charge this product.**

—5—

**Warning**

1. This product is not waterproof or water resistant, do not rinse or wash with water or other liquid. Keep this product clean and dry.
2. Make sure to charge the product between applications and charge it to full at least once in 2 months to keep its battery life.
3. While using, the massage head can be turned and adjusted only with a counter-clock wise rotation, and the maximum rotation range is of 90°. The reverse rotation applies after use (clock-wise).

**Matters needing attention**

1. If you have skin disease, issues, or a very sensitive skin please consult a dermatologist before use.
2. For the first time, please try it on your arm skin to confirm no sensitivities.
3. In the rare case you have any discomfort while or after using this product, stop using immediately and contact the manufacturer.
4. For hygiene reasons, please do not share this product with others.
5. Place the device out of the reach of children to prevent accidents.
6. Do not use corrosive chemicals to wipe the product;
7. Please do not disassemble this product.

—6—

**Product Specifications**

Product material	metal + PC plastic
Battery capacity	350mAh
Input voltage	5V
Charging time	2H
Product size	158*59.5*21.8mm
Net weight	40g
Power	≤3W

**Warranty Policy**

GLO24K offers a standard limited warranty on its beauty tools, massagers, and wands (check website for details)

—7—

**GLO24K**  
UNLEASH THE POWER OF GOLD



TO MAXIMIZE RESULTS SCAN HERE  
For more information about the product please visit  
**WWW.GLO24K.COM**

Dist. by Y&E Ent., LLC  
Atlanta, GA 30338 USA  
71-75 Shelton Street, London, UK  
Hauptstrasse 22 Neugersdorf 02727, Germany

—8—