

# 5

## MAINTENANCE

1. Wipe off with a dry cloth after each use. Do not wash off, product is not water resistant.
2. Turn Off the device before wiping it.
3. Always unplug after use and do not disassemble the product.

## STORAGE

1. Store this product in a cool and dry place avoiding humidity and direct sun light or store it in the box.
2. Do not submerge the product in water or liquid
3. Keep out of reach of children.

## TROUBLESHOOTING


Trouble	Action for trouble
The light is not very bright	Charge the device
Cannot charge	Change the plug/usb
The light flashes	Charge the device
Others	Contact seller

# 6

## SPECIFICATIONS

Product Name	7 Color LED Beauty Mask
Product Size	195*165*125mm
Product Weight	76g
Materials	ABS
Output Voltage	5V / 1A
Accessories	Mask, Manual, Cable Line
Power Supply	Rechargeable Battery
Charging Time	2 Hours

## PRODUCT APPEARANCE

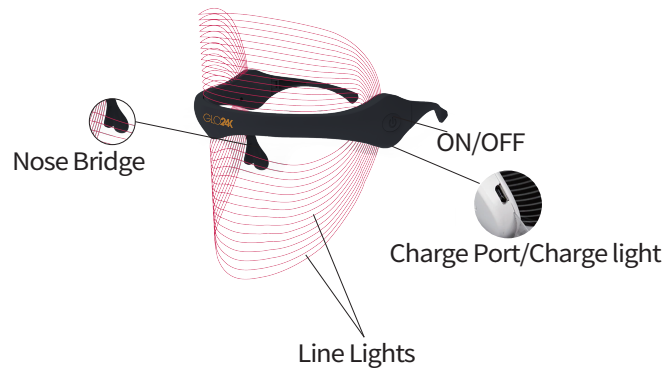

**GLO24K**  
UNLEASH THE POWER OF GOLD

## 7 Color LED Beauty Mask



# 1

## User Manual

Please Read This Manual Carefully Before Using

# 4

**GLO24K**  
UNLEASH THE POWER OF GOLD

## 7 Color LED Beauty Mask



TO MAXIMIZE RESULTS SCAN HERE  
For more information about the device please visit  
[WWW.GLO24K.COM](http://WWW.GLO24K.COM)  
Dist. by Y&E Ent., LLC  
Atlanta, GA 30338 USA  
71-75 Shelton Street, London, UK

## PRODUCT BENEFITS

- ※ This product uses 7 types of light sources including RED, BLUE, GREEN, YELLOW, PURPLE, CLEAR BLUE and LASER colors.
- ※ Assists in reducing wrinkles, fine-lines, age-spots. Helps in balancing the skin. Improves skin's flexibility and elasticity. Purifies Skin and fights acne. Firms skin for a radiant complexion. Accelerates blood circulation.

## FUNCTIONS

Red Light (650-730nm)	Reduces the appearance of fine lines and wrinkles, anti-aging
Blue Light (525-550nm)	Fights Acne/ Anti Inflammatory
Green Light (430-450nm)	Reduces and Prevents spots, soothing and purifying
Yellow Light (565-570nm)	Improves skin conditions, skin redness, and pigmentation
Purple Light (380-440nm)	Battles Acne/ Scars/ Marks, Minimizes Pores
Clear Blue (160-195nm)	Balances skin tone, Accelerates blood circulation
Laser Light (soft light)	Penetrates the skin deep, accelerates active tissue metabolism

# 3

## DIRECTIONS

1. Cleanse your face
  2. Apply your favorite skincare products (preferably GLO24K)
  3. Press ON/OFF button for 3 seconds to Turn On the device
  4. Click ON/OFF button to select desired LED light
  5. Place Mask on face for 10 minutes
  6. To Turn Off press the ON/OFF button for 3 seconds
- ※ Use as needed (auto timer is set for 10 minutes per use)
  - ※ Do NOT use the device while charging

## NOTES

1. This product is only recommended for use by an adult.
2. Mask is designed to distribute lights evenly on the face with no side effects.
3. The linear area is only the light guide and not the real source of light, which is the subtlety of the product.
4. The circuit inside the device should not be pressed or folded and it should be placed flat to protect the product after daily use.
5. In the case of sensitivity, please use product with eyes closed.
6. Using the Mask after applying skincare (preferably GLO24K) will improve results.
7. Use as needed. Each application for up to 10 minutes.
8. In the rare case of any discomfort stop using the device and contact the seller.
9. While charging, red light next to charging port indicates that the mask is being charged. When the red light is turned off, mask is fully charged.