



To Maximize Your Results Scan Here
For more information about the device
please visit WWW.GLO24K.COM



Skin Rejuvenation Beauty Device **User Manual**

To Maximize Your Results Go to Back Side

LED Infrared Photon Therapy

1. When your skin is in contact with the infrared device head, infrared light can penetrate into the subcutaneous tissue, then the tissue temperature will be rising, thereby it can promote blood circulation, revive enzymes, strengthen blood and cell tissue metabolism, and restore the cells for a youthful vitality.

2. After the infrared light is being absorbed by the skin, the skin water molecules will be resonated and activated, which enhances the binding strength between molecules, thereby activating the proteins in the skin cells produce resonance effects, which can transmit the infrared heating energy to the deeper tissue, and then, the skin temperature will be rising, so that the generated heat will be emitted from inside to outside.

3. GLO24K Skin Rejuvenation Beauty Device is an anti-aging facial care product, which provides effective relief for the skin layers. It is equipped with two different wavelengths of LED and far infrared, with wavelengths of 660nm (visible) and 850nm (invisible). The light waves can penetrate deep into the skin tissue cells, it can promote blood circulation, dilute facial spots and enhance skin elasticity.

- 1 40°C+/-3°C heating treatment head
- 2 Lightens facial spots and marks
- 3 Accelerates blood circulation
- 4 Minimizes fine-lines and wrinkles
- 5 Improves skin tone and balance
- 6 Repairs collagen
- 7 Restores skin elasticity

Remember:

1. Please fully charge the device before use.
2. Please recharge the device if not used for a long time.
3. Full charging time is about 2-3 hours.
4. Auto shut-off is set at 15 minutes for safety reasons.



Product Features




LED Infrared Skin Rejuvenation

Deeply heats the skin, dilates the micro-vessels, promotes blood circulation, repairs collagen, improves skin tone, reduces fine lines and wrinkles, tightens pores.

Stimulates Skin Rejuvenation
Balances Skin Texture

Repairs Collagen

Minimizes fine-lines and wrinkles

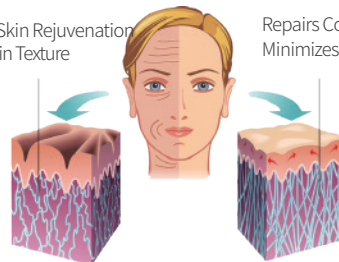

Heating, Vibrating, and Magnet Massage

40°C+/-3°C constant temperature, heating massage, deep cleansing, firming skin, magnet massage.
The metal head is embedded with two magnets, they can improve the body's immunity and regulate body functions.

40°C+/-3°C
heating treatment head

Heating massage

Healthy magnet

Skin Rejuvenation
Beauty Device



Your Device



You can charge this device directly on bottom of device or use the base.

Usage



1. Clean your face, remove make-up, and apply the device only on dry, clean skin.



2. Press A button, Turn on the device, the infrared will enter the standby mode (blue light), clean skin.



3. Press the B Button to start up the Heating, Vibration, and LED functions. Red light will show. Once in touch with your skin, these functions are ON and the device head will warm up gradually.



4. Gently massage your face by gliding the device over areas of concern in up/down and circular motions. A recommended use is 10-15 minutes per application 1-2 times a week.



5. Press the B Button to switch to standby mode.



6. To Turn Off the Device press the A button shortly (Not a long press). For best results, apply a moisturizing cream after each application of the device.

Tips: The above methods of care can be adjusted according to your needs.

Safety Tips

To avoid LED light directly into your eyes, we design the device light as a touch-sensitive sensor switch. Holding the device body, metal massage head touching the skin, the light therapy function will start. When you take it away from the skin, the light therapy function will stop working.

As the device generates heat remember to massage it over the skin of your face and neck and avoid keeping it on one skin area for too long. Move the device around eyes, please do not touch eyes directly and please close eyes to avoid discomfort.

Discontinue use in cases of discomfort or any skin irritation and contact supplier.

This device is non-waterproof, please do not expose it to corrosive substances and water.

Children should not use this product, please keep it away from children.

Product Parameters

Model No: BD1	Material: ABS + Stainless steel
Size: 200*46.5*29.5mm	Net Weight: 240g
Battery Capacity: 1100mAh	Working Voltage: 7.4V
Working Current: <1.5A	Heating Temperature: 40°C+/-3°C
Using Time: 10-15mins/time	Auto shut-off is set at 15 minutes for safety reasons

Charging: Flashing indicator lights when charging, blue light when fully charged.

- * Note: This device is set for an Auto shut-off at 15 minutes for safety reasons.
- * Charging time may take about 2-3 hours when battery is low.
- * The product cannot be turned on when the battery is low. Please make sure to charge device between each use.

Product Accessories: Infrared Beauty Device *1, Charging base *1
Charger* 1, User Manual *1

Repair & Maintenance

To clean your device make sure it is turned off and unplugged.

Clean the device with a soft, dry cloth between uses.

Put the device in a safe place after use. Do not place heavy objects on it.

Store the device in a dry place to avoid moisture.

Recharge device between uses.

If trouble-shooting the device is needed follow the steps below. If the problem is not resolved contact the product service center. Do not disassemble the device.

Breakdown Maintenance	Breakdown Reason	Breakdown Measure
No light therapy or device does not work on the skin	Make sure the device is ON. Make sure the device is touching the skin	Re-adjust gesture to touch skin
Device does not turn ON or does not work	Battery ran out	Please try charging, charging time is about 2-3 hours
Device lights are acting out	Device needs a restart	Turn off the device. Wait a few minutes. Recharge. Turn back ON
Device automatically shuts-off after 15 minutes	Auto shut-off is set at 15 minutes for safety reasons	Wait a few minutes and turn the device back on

WWW.GLO24K.COM GLO24K offers a limited warranty on its Beauty Devices
Dist. by Y&E Ent., LLC (check website for details)

Atlanta, GA 30338 USA

71-75 Shelton Street, London, UK

