



Triple Action Eye Care Therapy Wand (User Manual)



Precautions

Thank you for purchasing this product. Please read this instruction manual carefully before using this product and keep it for future reference.


1. This product should not be used by small children and is for external use only. If you have any medical condition, please consult with your physician before using this product.
2. Stop using this device if you feel any discomfort. Do not use this product while wearing contact lenses. If you feel abnormalities in your eyes, eyelids or skin during use or right after use and it persists, consult your medical specialist.
3. Avoid direct contact with the eyes.
4. Do not use while charging.
5. For first time users or people with sensitive skin: Apply for a short duration to ensure there are no side-effects. In the rare case that you feel uncomfortable or the skin becomes irritated please stop using the product and contact the manufacturer.

6. Do not repair, disassemble or modify this product without authorization.

Function Introduction

This Eye Care Therapy Wand is based on triple action: LED, Vibration, and Temperature Controlled Light Heating Technologies. It was developed specifically for the delicate skin around the eyes. It can assist in reducing the visible signs of the aging process around the eyes for a more rejuvenated, revived appearance. For best results apply some Eye Cream before using the wand.

Product Structure

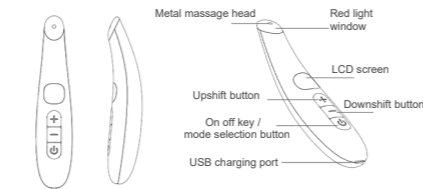

Product Specifications

Product name	Intelligent eye beauty device	Product model no	1203
Net weight	47.5g	Power	0.7W
Rated voltage	DC 3.7V	Material	Zinc alloy+ABS
Temperature	up to 40°C	Working voltage	3.7-4.2v
Vibration frequency	10000RPM		
Charging time	1.5 hours		
Unit size	139×29×28.5mm		

Three Working Modes

Mode 1: Vibration + Red Light + Temperature Controlled Light Heating

Triple effect care. This mode can assist in diminishing the visual signs of the aging process around the eyes such as wrinkles and fine-lines.

Mode 2: Red Light + Temperature Controlled Light Heating

This mode can help in stimulating cell activity and promoting skin elasticity. It can also assist in rebalancing the skin around the eyes for a more radiant appearance.


Mode 3: Vibration + Temperature Controlled Light Heating

This mode of high-frequency massage and gentle temperature-controlled heating can assist in stimulating the muscles in the eye area for a more invigorated appearance. It can also make the absorption of any Eye Cream applied more efficient.

Using Method



1. Apply your favorite Eye Cream around the eyes and massage in circular motions until it is fully absorbed. (Recommended)

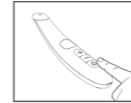

2. Long press (switch button) to turn on product and short press (switch button) to switch the mode, and press "+" / "-" to adjust the temperature as needed.



3. Holding the upper part of the product body, use the massage head to massage as shown.



4. After using each time, long pressing (on/off button) to turn off, please wipe the massage head with a dry towel or paper towel.

Charging Method

1. Charge this Beauty Wand between applications.
2. During the charging process, the display screen battery symbol flashes and the battery symbol stops flashing after charging is completed.
3. Please charge occasionally to avoid aging after no using for a long time. This wand is not water-proof or water-resistant.
4. When charging is complete, please unplug the charger.

Maintenance

1. If applying Eye Cream before application, spread the cream evenly and let it absorb before using the product.
2. Please wipe the product clean with a dry towel after each application.
3. Please don't pour water on the product or wash the body with water. This product is not water-proof or water-resistant.
4. Please don't use corrosive to clean the product.
5. Please store it in a cool dry place.
6. Please don't use the product if it is damaged.
7. Please don't give a strong impact on the body.

Troubleshooting

Frequent Questions	Possible Causes	Treatment
The product can't be turned on	The turning on mode is wrong	Short press the on/off button
	Power is exhausted	Please try charging, if you still can't turn on after charging, please contact the seller
The product does not work after being turned on.	The battery is exhausted	Please try charging, if it is still not available after charging, please contact the seller.
The product's vibration is not strong after touching the skin or does not work.	The resistance is big after touching the skin.	Using together with skin care products is recommended to achieve more significant effect.

Warranty Policy

GLO24K warrants all of its beauty tools. Please check website for details.



TO MAXIMIZE RESULTS SCAN HERE
For more information about the device please visit
WWW.GLO24K.COM
Dist. by Y&E Ent., LLC
Atlanta, GA 30338 USA
71-75 Shelton Street, London, UK