GOZONE

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MASSAGE GUN / PISTOLET DE MASSAGE

⚠ Please read all the instructions in the manual carefully before assembling or using this product, and keep instructions for future reference.

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GZE-2243

Thank you for purchasing our product. To ensure proper operation and use, please read this manual carefully prior to use.

Please keep this manual properly for future reference.

Massage gun is suitable for physical massage before/after exercise to relieve muscle fatigue, soreness, and stiffness.

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IMPORTANT INFORMATION


⚠ Notice:

To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment damage, please read all instructions and warnings in this manual before using this product.

Please consult your physician before using this product.

⚠ Attention:

1. For adults use only. Stop the use of this item if injury is sustained.
2. Do not use this item on non-dry or unclean skin. Always wear a shirt/clothing when using.
3. Do not use this product on the neck, head or any hard areas of the body.
Do not use this product on bones. Stop the use of this item if you experience pain or discomfort.
4. Use this product only on the soft tissue of your body.
5. Do not place any objects in the vents of this product.
6. Do not submerge the product in water or allow water to enter the vents of the product
7. Keep away from sources of fire.
8. Do not drop this product intentionally.
9. Do not use an adapter of any kind other than the adapter included in the attachment.
10. Do not disassemble this product by yourself.
11. Do not allow fingers, hair, or other body parts to enter the joint between the massage node and machine while the machine is running.

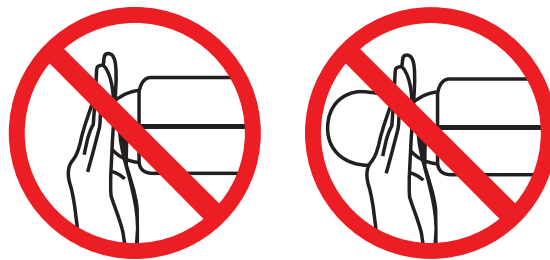Do not touch this part while the machine is running

⚠ Warnings and Cautions:

1. Keep this product away from water or heated surfaces. Keep the vents free from dust.
2. Unplug the charger while charging is completed.
3. The product cannot be used while being charged.
4. Thoroughly inspect products according to the operating instructions in this manual before every use.

⚠ Important:

Do not use this product in any of the following cases without consent of your physician: Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, epilepsy or migraine, disc herniation, Scoliosis, vertebral detachment, or spinal joint stiffness, recent joint replacement surgery or IUD, having metal pins/screws or steel plate, and/or any problems related to your general health. The contraindications mentioned above do not mean that you are not allowed to use this product, we recommend you consult your physician before use.

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

PRODUCT INTRODUCTION

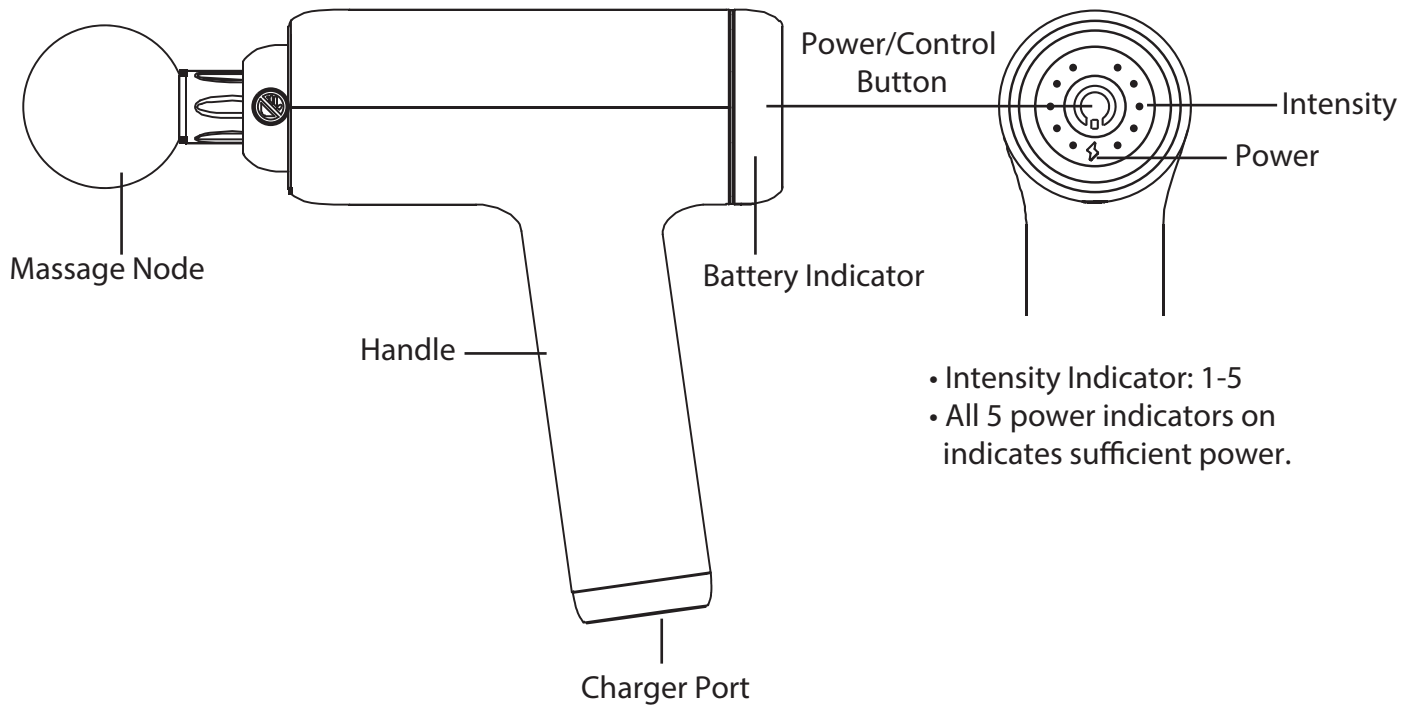- Intensity Indicator: 1-5
- All 5 power indicators on indicates sufficient power.


1. Whisper quiet high-speed brushless balancing motor.
2. EVA soft massage node to provide a more comfortable massage.
3. Rechargeable lithium-ion battery.
4. Delicate handle for a more comfortable grip.
5. Wireless, lightweight and portable.
6. 5-speed adjustable intensity.
7. User-friendly design: sweat-proof and moisture-proof.
8. Integrated design combining the battery with a machine for safe use.

4 REPLACEABLE MESSAGE NODES

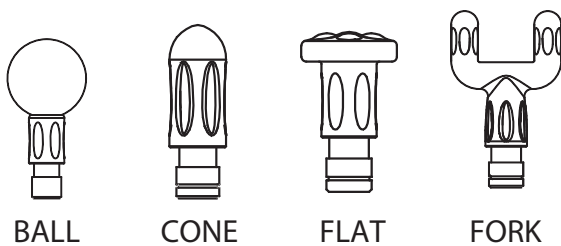
* Only swap massager nodes when the massager is off.

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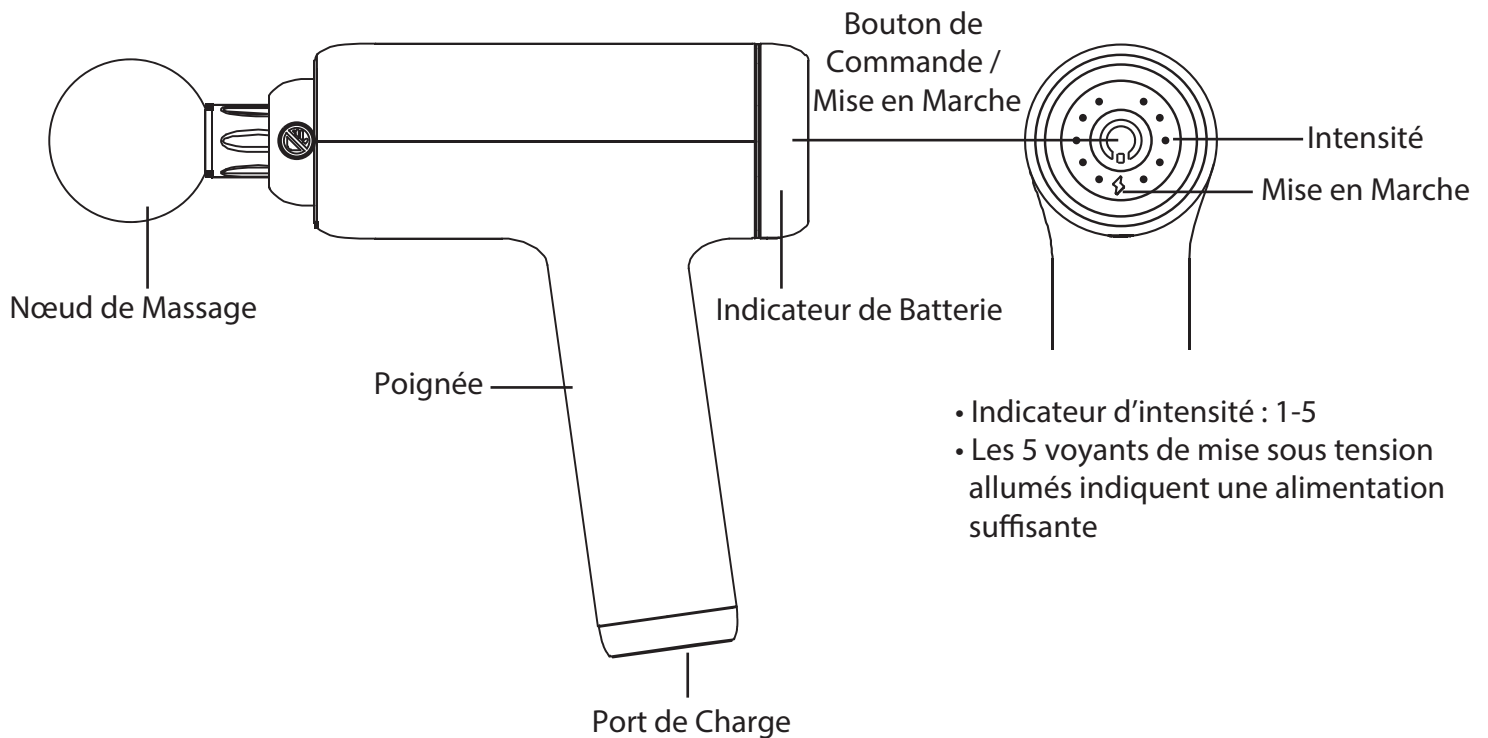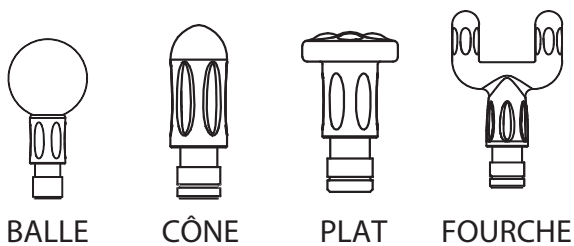CHARGE

1. Charge the battery for at least 6 hours before first use.
 2. Plug the adapter into the massager charging port and then plug into an outlet to charge.
 3. The power indicator will flicker regularly while charging and keep steady when fully charged.
 4. The battery can be charged at any time and under any power level.
- *Do not use the massager while charging.

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OPERATION INSTRUCTIONS

Press the control button for 3 seconds to turn the massager on.
Press it again to cycle through the 5 intensity levels.

When finished using, press the control button for 3 seconds to turn the massager off.

MAINTENANCE, CLEANING, STORAGE, TRANSPORTATION

1. Spot clean only. Unplug and wipe with a damp cloth.
2. Air-dry only. Never place in a dryer.
3. Clean the massager and pack well if not using for a long time and charge every half a year.
4. This product is safe to be transported by automobiles, trains, ships, and aircraft.

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


USER GUIDE

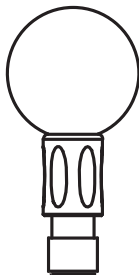
BALL
Mainly used for large muscle groups

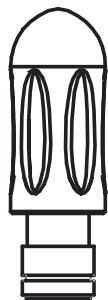
CONE
Mainly used for trigger point and deep muscle groups

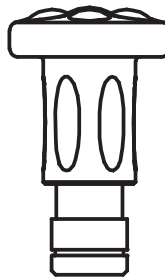FLAT
Suitable for all body parts

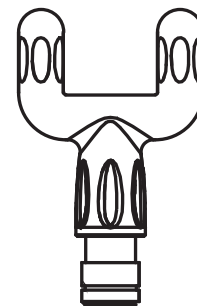
FORK
Use around spine, waist, and large muscle groups

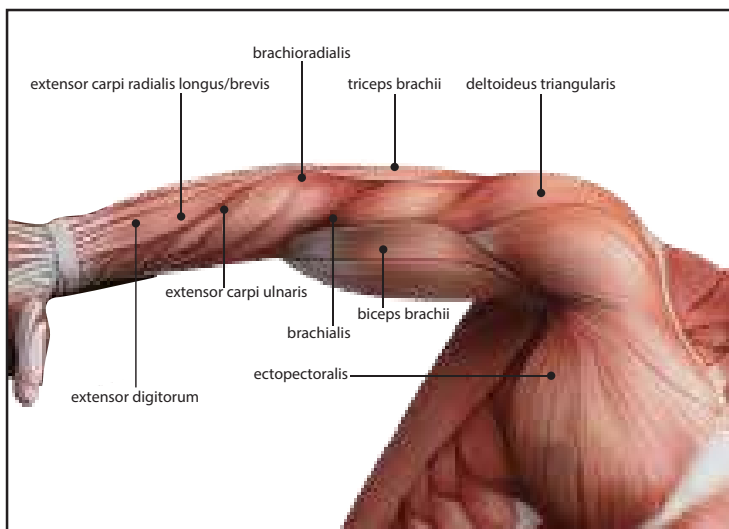ARM

brachioradialis	10-15s x 3 times
flexor carpi radialis	10-15s x 3 times
palmaris longus	12-15s x 3 times
extensor carpi ulnaris	10-15s x 3 times
extensor digitorum	5-10s x 2 times
extensor carpi radialis longus	5-10s x 2 times
extensor carpi radialis brevis	5-10s x 2 times
biceps brachii	12-20s x 3 times
brachialis	12-20s x 3 times
deltoideus triangularis	12-20s x 3 times
ectopectoralis	12-20s x 3 times

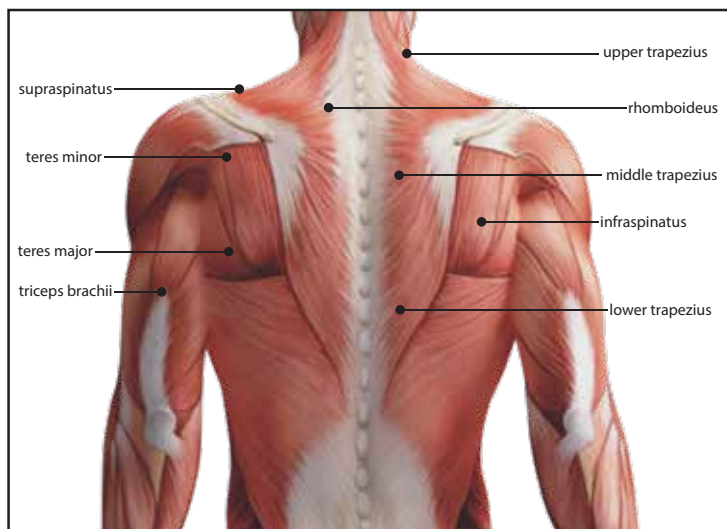
SHOULDERS

triceps brachii	15-20s x 3 times
rhomboides	20-30s x 3 times
upper trapezius	15-20s x 3 times
middle trapezius	15-20s x 3 times
lower trapezius	15-20s x 3 times
supraspinatus	10-15s x 3 times
levator scapulae	10-12s x 2 times
biceps brachii	10-12s x 2 times
teres major	10-12s x 2 times
teres minor	10-12s x 2 times
infraspinatus	10-12s x 2 times

*Suggestion: Massage from the middle to two sides, top to bottom.

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
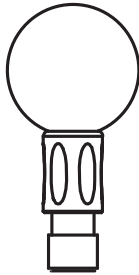
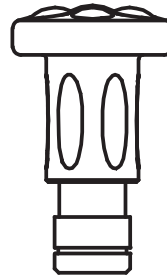
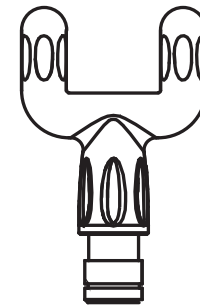
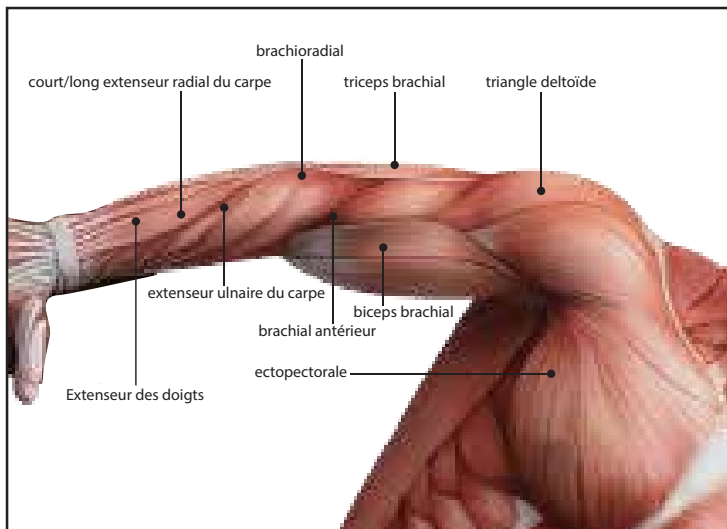
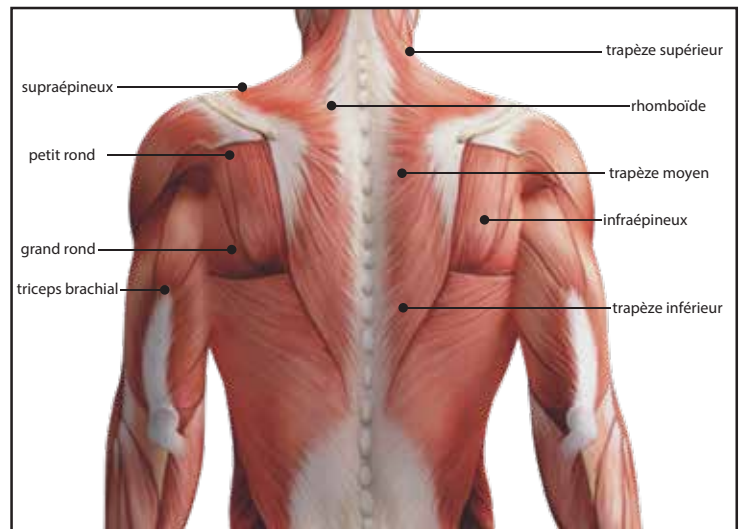
USER GUIDE

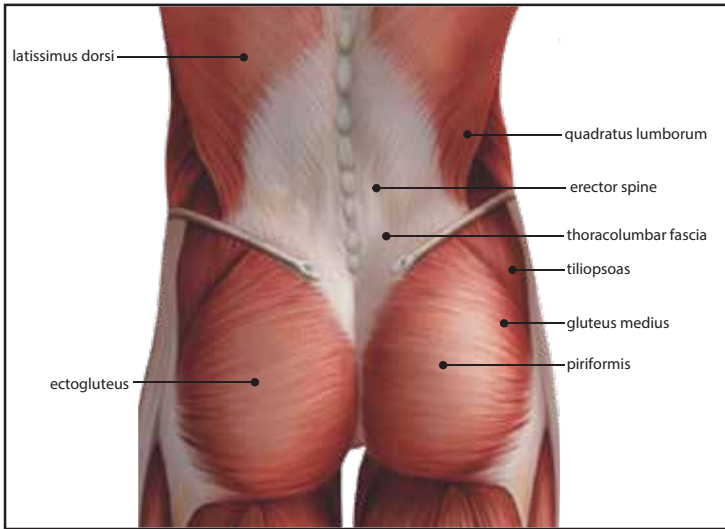WAIST

erector spine	30-60s x 4 times
latissimus dorsi	30-60s x 4 times
thoracolumbar fascia	40-60s x 3 times
psoas major	30-60s x 3 times
obliquus externus abdominis	10-15s x 2 times
erector spine	15-20s x 2 times
quadratus lumborum	20-30s x 2 times
iliopsoas	20-30s x 3 times
ectogluteus	30-60s x 4 times
gluteus medius	30-60s x 4 times
piriformis	60-90s x 1 times


*Suggestion: Massage in the direction of middle to two sides, top to bottom.

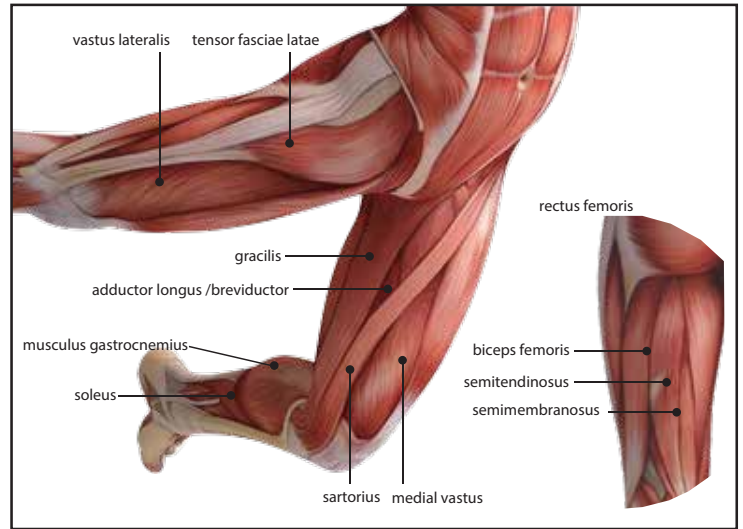THIGH

biceps femoris	20-30s x 3 times
semitendinosus	20-30s x 3 times
semimembranosus	20-30s x 3 times
gracilis	10-15s x 2 times
breviducto	10-15s x 2 times
adductor longus	10-15s x 2 times
rectus femoris	30-45s x 2 times
vastus lateralis	30-45s x 2 times
medial vastus	30-45s x 2 times
sartorius	30-45s x 1 times
tensor fasciae latae	15-20s x 2 times

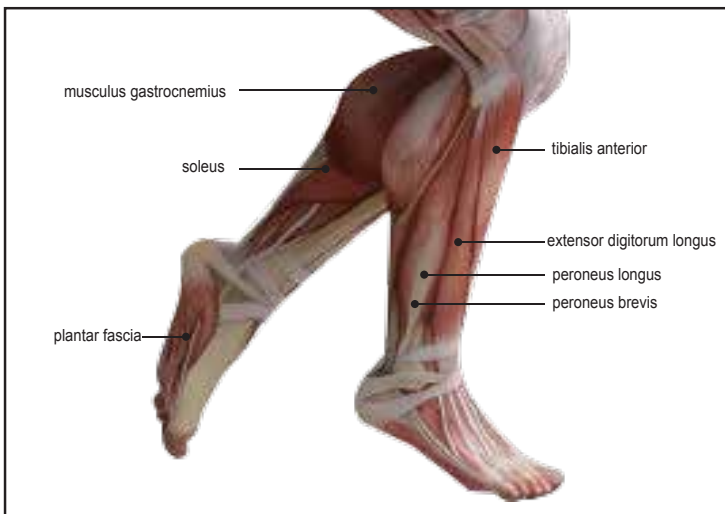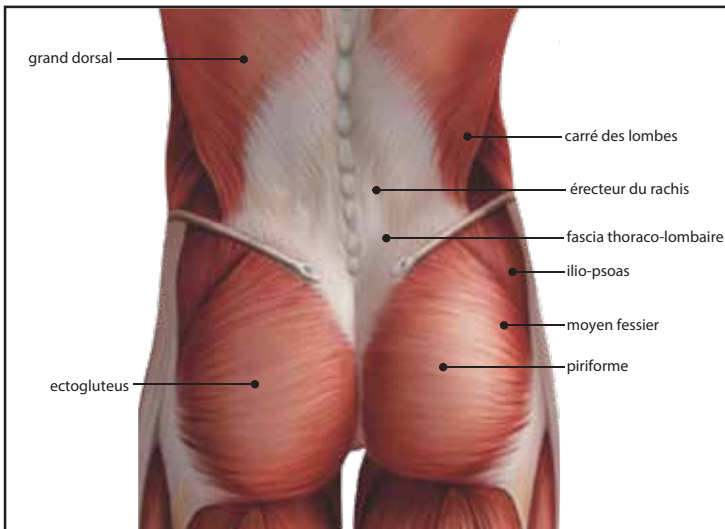
*Suggestion: Massage from top to bottom.

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