


# L' Masque

## LED light therapy mask

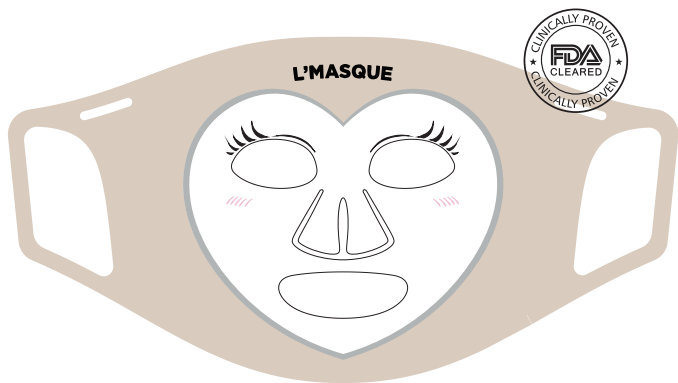

Instruction Manual

# PRODUCT INTRODUCTION

The L'Masque light therapy mask uses LED light to stimulate collagen production under the skin. This helps to speed up the skin's natural processes, reduce pigmentation and redness, and improve the elasticity of the skin. With regular use, the mask can help to make the skin appear brighter, smoother, and an.

A single charge can be used 4-6 times, and just 10-15 minutes three times a week can promote skin metabolism and enhance product absorption.

L'Masque uses two wavelengths of red light: 633nm and 830nm or near infrared light. These wavelengths have been shown in clinical studies to stimulate collagen and elastic production cells and to improve blood flow and tissue oxygenation through a process called photobiomodulation. These effects help skin absorb topical skincare products more efficiently and helps rebuild the damage caused by the aging process.

L'Masque also uses blue light. Blue light therapy on the face has shown several benefits in clinical studies. It is effective in treating acne by killing acne-causing bacteria and regulating oil production. Blue light also has anti-inflammatory effects, promotes wound healing, and contributes to overall skin rejuvenation.

Please consult your dermatologist if under treatment for any concerns.


# CONTENTS

|                                   |    |
|-----------------------------------|----|
| 1.Product parameters .....        | 01 |
| 2.Equipment parts .....           | 02 |
| 3.Product description .....       | 03 |
| 3.1 Equipment list .....          | 03 |
| 3.2 Product features .....        | 03 |
| 3.3 Control marks .....           | 03 |
| 4. Safety .....                   | 04 |
| 4.1 Fair use .....                | 04 |
| 4.2 Safety Instructions .....     | 04 |
| 5. Charging .....                 | 06 |
| 6. Operation .....                | 07 |
| 6.1 Operation selection .....     | 07 |
| 6.2 The function of light.....    | 08 |
| 6.3 Mode selection .....          | 09 |
| 6.4 Sensitive testing .....       | 10 |
| 6.5 Usage steps .....             | 11 |
| 7. Maintenance .....              | 13 |
| 7.1 Maintenance Tips .....        | 13 |
| 7.2 Cleaning steps .....          | 13 |
| 8. Check the product.....         | 14 |
| 9. Discard the product .....      | 14 |
| 10. Warranty .....                | 16 |
| 11.Safety signs and meaning ..... | 18 |
| 12.Distributor .....              | 20 |


# Product parameters

|                      |                |
|----------------------|----------------|
| Brand                | Boujee Beauty  |
| Product name         | L'Masque       |
| Model No.            | FM-02          |
| Rated voltage        | 3.7V           |
| Material             | Silicone       |
| Battery capacity     | 1800mAh        |
| Number of LED lights | 100 pcs        |
| Weight               | 213.3g         |
| Size (L×W×H)         | 383×198×33.5mm |


# Equipment partsv

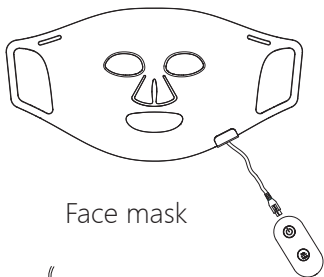


Face mask

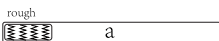


a

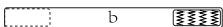

b

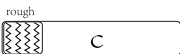


c

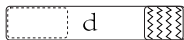


d

4 Head straps

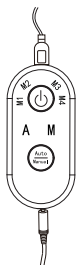


ON/OFF

AUTO/MANUAL

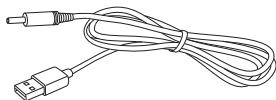

USB cable

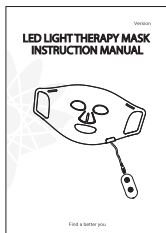

User manual

# Product description

## 1、Equipment list

- ①Silicone mask    ②Remote control    ③ 4 straps    ④ USB cable    ⑤Instruction manual

## 2、Product features:

- ①Top industrial aesthetics and simple fashion design, better fit on the face, and the LED light makes the skin absorb faster.
- ②There are medical-specific silicone eye pads on the eyes to protect the eyes from light stimulation.
- ③You can choose the corresponding working mode and adjust the time according to your skin care needs.
- ④Small and lightweight, easy to carry, you can take care of your skin anytime, anywhere.

## 3、Control marks

- ①The controller turns on/off to control the working state and mode of the mask and control the power supply of the mask.
- ②The controller contains a rechargeable lithium-ion polymer battery.

| Symbol | Function       |
|--------|----------------|
| M1     | Red light      |
| M2     | Blue light     |
| M3     | Infrared Light |
| M4     | Mixed light    |
| A      | Auto mode      |
| M      | Manual mode    |

# Safety

## 1、 Fair use

This product is specially used to improve the skin's excessive oil secretion, pigmentation, acne, large pores, dull skin and other problems, so that the skin can be improved from the inside out, and can effectively restore luster and elasticity. In addition, this product can also enhance metabolism, delay skin aging, inhibit oil secretion, effectively sterilize, treat acne, anti-inflammatory, reduce swelling, and relieve pain. Please note that the mask and USB charging cable are not waterproof, keep them away from water or other liquids.

This product can only be used in the above scenarios and cannot be used for other purposes, otherwise it may cause damage to the product or the person. The manufacturer and seller are not responsible for damages caused by improper or incorrect use.

## 2、 Safety Instructions

Warning! Risk of electric shock

The supplied USB charging cable can only be used to connect the mask and not for other purposes. Connect the USB charging cable to an easy-to-remove plug so you can quickly disconnect it from the power source if something goes wrong.

- If the mask is damaged, do not continue to use it, and contact customer service as soon as possible to replace or repair it to avoid risks. In the case of self-repair,

disassembly, improper connection or improper operation by the user, the merchant is exempt from liability and warranty services.


- During the charging process, the mask should not be used.

During charging, do not touch the mask, USB charging cable and USB adapter with wet hands.

- If you feel any discomfort or any adverse reaction during use, please stop using it immediately.



# Charging

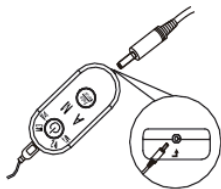

Insert the supplied USB charging cable into the charging hole of the remote control, and insert the other end into the adapter to ensure that the voltage marked on the product matches the power supply voltage. When charging, one indicator light on the remote control will light up in white, and all four indicator lights will light up to indicate that the charging is complete.

## Tips:

- When charging, use only the original USB charging cable provided with the product, and do not connect the mask during charging.
- Using an adapter lower than 5V will cause the charging to fail. Using an adapter higher than 5V may damage the remote control.
- When the remote control voltage is lower than 3.2V, the indicator light on the remote control will flash, and it will automatically shut down after 10 flashes.
- The mask has overvoltage protection (the protection range is higher than 6V and lower than 20V). While charging, the machine cannot be used.

# Operation

## 1、 Operation selection

- ① Press the power button for 1 second, turn on the mask, the function red light is bright, red light mode: wavelength 620nm.
- ② Press the switch again, the function blue light is bright, blue light mode: wavelength 460nm.
- ③ Press the power button for the third time, the function is near infrared light 830nm. This light is not visible to the human eye but the mask will illuminate with 3 red lights so you know it is on.
- ④ Press the power button for the fourth time, the functions mixed light is bright.
- ⑤ Press the power button for the fifth time, the mask will be turned off.

## 2、 The function of light

Red light: Promote and improve blood circulation, anti-inflammatory, reduce swelling, relieve pain, lighten pigmentation, suitable for all skin types.



Blue light: Inhibit oil secretion, balance oil, effectively sterilize, treat acne, suitable for oily skin and acne-prone skin.



NIR light: Improve dullness, lighten spots, accelerate the excretion of pigmentation, and enhance the skin's radiance.



Mixed light: gathers a variety of light, sterilization, anti-inflammatory and anti-aging, suitable for all skin types.



### 3、 Mode selection

The machine is divided into manual mode and automatic mode

Manual mode: The mask is turned on by default in manual mode (default red light). If there is no button operation for 3 seconds, the mask will automatically time and automatically shut down after 10 minutes.

Automatic mode: After the mask is turned on, press the mode button once to switch to automatic mode without any operation, the mask will automatically time after 3 seconds, red light 5 minutes, blue light 3 minutes, yellow light 4 minutes, mixed light 3 minutes, arrival time automatically shuts down after that.

Reminder: When using the mask in manual mode, it should be controlled within 15 minutes. Excessive use may cause damage to the skin.

## 4、 Sensitive testing








Before your first treatment, you must have a sensitization test. Follow the step-by-step guide below:

1. Place the machine on your inner arm.
2. Turn on the controller and let the machine run a full 10-minute treatment program.
3. If your skin feels hot or uncomfortable, stop using the device immediately.
4. Wait at least 6 hours after treatment. If the skin on your inner arm is red and/or itchy, then the 5、

Usage steps

Recommendations:

For best results, use 3-4 times a week for at least 4 weeks. Maintenance use should be performed as required, according to lifestyle and skin condition.

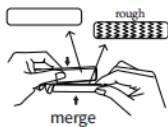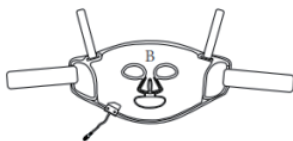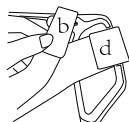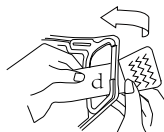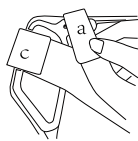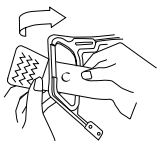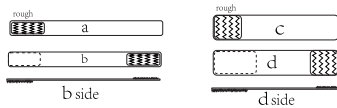

1. Assemble the mask as instructed above.
2. Clean and dry your face and hands thoroughly.
3. Adjust the straps on the mask to make sure the mask fits comfortably.
4. Sit or lie down with the mask over your face.
5. Use the controller to select light color, brightness, blinking or steady light, and treatment time.
6. Use your favorite skin care serum or cream after treatment as absorption will be enhanced after photo-therapy.
7. Continue your skin care routine. 15 minutes per treatment, 3-4 treatments per week, reducing to 1-2 treatments per week once results are available.

### **Tips:**

If the mask is used in manual mode, it should be controlled within 15 minutes. May cause skin damage if used for too long. If you are using this product for the first time, please refer to the following maintenance plan and adjust it according to your skin condition.

Week 1: Use every 3 days, 10 minutes each time.

Week 2: Use every 2 days for 10 minutes each time.

Week 3: Use every 2 days, 10-15 minutes each time.

Week 4: Use once a day for 10-15 minutes each time.

# Maintenance

## 1、 Maintenance Tips

- Clean the surface of the LED light therapy mask after each treatment.
- Store the LED light therapy mask in a dry room temperature condition.
- Avoid storage in dusty environments, direct sunlight or humid environments.
- The LED light therapy mask contains no serviceable parts.

## 2、 Cleaning steps

- (1) When cleaning and arranging the product, please turn off the power.
- (2) Wipe the mask surface twice with 70%-80% ethanol solution for 3 minutes each time, and then wipe with a damp or dry cloth.
- (3) When cleaning this product, do not use flammable materials such as benzene, thinner, butyl acetate, etc. to wipe it.
- (4) When cleaning the product, do not spray water directly on the product. Wipe the silicone surface twice with 70%-80% ethanol solution for 3 minutes each time, then wipe with a damp or dry cloth.
- (5) Do not use any cleaning products, cleaning with other products may damage the mask.



## Check the product

After receiving the product, if you use scissors or knives to open the package, be careful not to scratch the product. In addition, it is necessary to check whether the product accessories are complete and whether the product can be used normally. If there is any omission or the product cannot be used normally, you should contact customer service for replacement in time.

## Discard the product

1. When the appliance is discarded, the battery needs to be taken out from the appliance; take-out steps: use a tool to pry open the remote control and take out the battery.
2. The appliance must be switched off when the battery is removed.
3. After the battery is taken out, it should not be discarded at will, and garbage should be sorted.

# Warranty

Thank you for purchasing the LED Light therapy mask. You can enjoy a one-year warranty from the date of purchase. During the warranty period, the product will be repaired free of charge due to technical defects. Failures caused by improper use, self-repair, improper transportation, etc. are not covered by the warranty.

This warranty does not cover any shipping or handling charges which may incur. This warranty does not include any damages caused by misuse, accident, neglect, abuse, use with a voltage other than what's marked on the product, or for commercial purposes.

If your Boujee Beauty appliance is found to be defective within the time period of this warranty, please contact the Boujee Beauty warranty department at [warranty@boujeebeauty.com](mailto:warranty@boujeebeauty.com). If the product is discontinued or no longer available, a replacement may be made with a similar product of equal or greater value. Any attempt to repair or service this product on your own will void the warranty.

## **MISSING, DAMAGED, OR INCORRECT ITEMS:**

It is the responsibility of the customer to notify Boujee Beauty of any damaged, incorrect, or missing items within 48 hours of delivery for immediate correction.

# Warnings

LED light therapy mask is a household wearable LED light therapy device. It can only be used according to the usage method described in the manual, and cannot be used for other purposes, otherwise it may cause damage to the product or the human body. Manufacturers and distributors are not responsible for damage caused by misuse or incorrect use.

1. It is recommended to close your eyes when using, if you have skin redness, itching, allergies, etc. during use, please stop using it immediately.
  2. Pregnant women, infants, people with epilepsy, skin damage, thyroid patients or thyroid medication users, and people with light allergies should not use this product.
  3. This product is not suitable for people with broken skin, trauma or post-operative wounds.
  4. Do not use if you suffer from Lupus erythematosus, photosensitive eczema or Albinism.
  5. Do not use if are taking any medication that can cause photosensitivity or if you suffer from any photosensitive disorder, you may cause a severe skin reaction. Be aware that other substances can also cause photosensitivity such as St John's Wort, deodorants, antibacterial soaps, mothballs, petroleum products, etc.
- ④The mask cannot be used during charging.

## FAQ

Q: Is the LED light therapy mask waterproof?

A: The mask is not waterproof, it can only be wiped with water or ethanol, and the mask cannot be rinsed.

Q: How long can the LED light therapy mask last on a single charge?

A: In the automatic mode, it can be used 4 times for 15 minutes each time; in the manual mode, if the mixed light is selected, it can be used 3 times for 15 minutes each time. If red light, blue light or yellow light are selected, it can be used 4 times for 15 minutes each time.

Q: The skin will be a little dry after using the LED light therapy mask, how can I improve it?

A: Because the wavelength of light stimulates the dermis layer of the skin and stimulates cell metabolism, the skin's absorption capacity is enhanced, and more moisture and nutrients are needed. Therefore, it is recommended to replenish the skin with a large amount of moisture and essence after using the LED light therapy mask to promote skin absorption and give the skin enough water and nutrients.

## Safety signs and meaning

Q: Will the temperature of the LED light therapy mask increase during use?

A: When operating with red light, blue light and yellow light, the temperature of the LED light therapy mask will not rise during use. When operating with mixed light, the temperature of the mask will rise during use, and the maximum temperature is  $38\pm 1$  °C.

Q: Is the LED light therapy mask effective immediately after use?









A: The effect of using once or twice may not be obvious, because skin metabolism needs a cycle, and a skin metabolism cycle is 28 days, so after a month of consistent use, the effect will be more significant.

Q: How long is the best time to use the LED light therapy mask?

A: We recommend that you use the LED light therapy mask for no more than 15 minutes at a time, because it may damage the skin for too long, so you need to control the time of use.

Q: Can the LED light therapy mask be folded?

A: No, it will cause internal circuit breakage.

| Safety sign   | meaning  |
|---|--|
|    | Electrical appliances for indoor use.  |
|    | Don't throw away at will.  |
|    | Conforms to certain EC Directives on health, safety, and environmental protection standards.   |
|    | This symbol on the device label means: Read the user manual before you start using the device. |
|    | This symbol means keep away fire   |
|    | This symbol means place upwards  |
|   | This symbol means keep dry   |
|  | TYPE BF APPLIED PARTS  |

|                         |   |
|-------------------------|---|
| <p><b>RoHS</b></p>      | <p>This symbol means: Directive on Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment. ROHS stands for Restriction of Hazardous Substances</p>   |
| <p><b>REACH</b></p>     | <p>This symbol means: "Registration, Evaluation, Authorization and Restriction of Chemicals, chemical registration, evaluation, authorization and restriction", which is the EU's regulations for preventive management of all chemicals entering its market.</p> |
| <p><b>UK<br/>CA</b></p> | <p>This symbol means: "United Kingdom Conformity Assessed", electrical and electronic products entering the UK market need to pass standard tests under UK regulations.</p>   |
| <p><b>IP22</b></p>      | <p>IP classification</p>  |

**Distributor:**

Boujee Beauty San Bernardino, CA 92407.