

ageLOC.

LumiSpa™



USER'S MANUAL

SAFETY INSTRUCTIONS

ageLOC® LumiSpa™ was designed for use by competent adults. It can be used by children aged 8 years old and above and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning its use in a safe way and understand the hazards involved. ageLOC LumiSpa is not a toy and children should not play with it. Cleaning and user maintenance should not be performed by children.

WARNING

The detachable power supply must be placed where it cannot fall into water. For the purposes of recharging the battery, only use the power supply unit provided with this device.

PRECAUTIONS

USAGE

- Do not use harsh chemicals or abrasives on your ageLOC LumiSpa device.
- Keep out of reach of children.
- Do not use device if treatment head is damaged
- Only use device as directed and do not apply excessive product to the treatment head.
- Only use device for standard treatment times and do not use for an excessive period.
- Do not use for an excessive period on one area of the skin.
- Only use device as directed and do not apply to compromised skin.
- Do not use on raised moles or compromised skin.
- Do not conduct treatment more frequently than recommended.
- Do not expose charging base to water.
- Do not use device if it is damaged.
- Do not use power cord in a wet environment.
- Do not use power cord if it becomes wet.
- Do not handle charging base with wet hands.
- Regularly clean handgrip.
- Do not share treatment heads.

TEMPERATURE

Never place your device near or over a radiator, heat register, or any other heat source. Avoid leaving your device in an overly hot or cold environment for prolonged periods of time. The storage temperature range for the device is between -10° C and 50° C (14° F and 122° F).

ELECTROMAGNETIC COMPATIBILITY

Conforms to the requirements of the 2004/108/CE directive on electromagnetic compatibility.

For United States Users: The ageLOC® LumiSpa™ device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) the device may not cause harmful interference, and (2) the device must accept any interference received, including interference that may cause undesired operation.


SERVICING

Do not attempt to repair your device yourself. This will void any warranty. There are no user-serviceable parts inside. Please see warranty section for service.

TABLE OF CONTENTS

System Components	5
Treatment Heads and Treatment Cleanser Options	6
Charging Your ageLOC® LumiSpa™ Device	7
Attaching and Removing Your Treatment Head	8
Using Your ageLOC LumiSpa Device	9
Cleaning and Caring for Your ageLOC LumiSpa Device	10
Tips	10
Replacement and Warranty Information	11

SYSTEM COMPONENTS

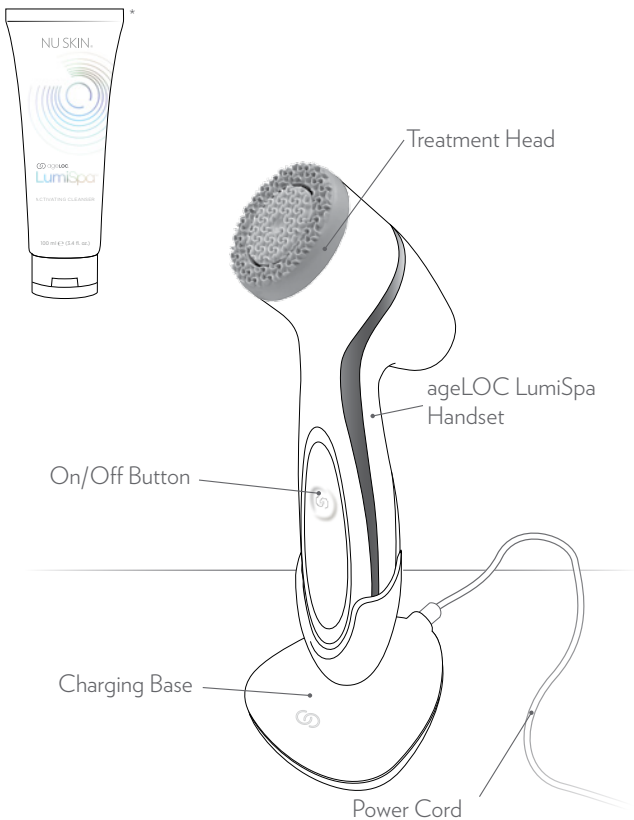

NOTE: Illustrations do not reflect actual size.

*There are three ageLOC LumiSpa treatment cleanser options available.

TREATMENT HEADS AND TREATMENT CLEANSER OPTIONS

ageLOC® LumiSpa™ offers a selection of treatment heads and treatment cleansers, so you can find the best combination for your skin.

Treatment Heads

Each ageLOC LumiSpa treatment head features a gentle, nonabrasive silicone face that is embedded with silver particles. There are two treatment head options available: Normal and Firm.



Normal



Firm

Treatment Cleanser Options

The ageLOC LumiSpa treatment cleansers are designed to be compatible with the ageLOC LumiSpa device. There are three ageLOC LumiSpa treatment cleanser options, which are based on skin type and personal preference.


ageLOC LumiSpa Treatment Cleanser—Dry

ageLOC LumiSpa Treatment Cleanser—Normal/Combo

ageLOC LumiSpa Treatment Cleanser—Oily


CHARGING YOUR AGELOC® LUMISPA™ DEVICE

1

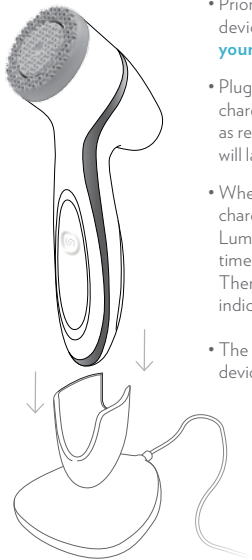


- Prior to placing your ageLOC LumiSpa device in the charging base, **make sure your device is completely dry.**
- Plug the charging base into an outlet and charge your device overnight. When used as recommended, the fully charged device will last at least a week.
- When you place your device in the charging base, the lights on the ageLOC LumiSpa handset will radiate for a short time and the On/Off button will remain lit. Then, the On/Off button will pulse slowly, indicating the device is charging
- The On/Off button will remain lit when the device is fully charged.

WARNING: Rechargeable batteries, like those used in ageLOC LumiSpa, need to carry a minimal charge in order to function properly. If the battery is allowed to fully discharge and then is not charged for an extended period of time, it may adversely impact the battery's ability to be charged. For optimal battery life, we recommend charging your device whenever low battery is indicated and to avoid storing the device unused for months at a time.

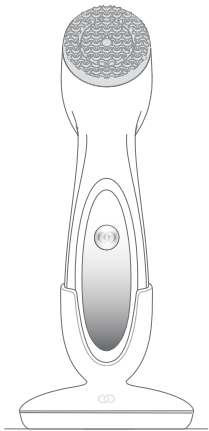

ATTACHING AND REMOVING YOUR TREATMENT HEAD


1

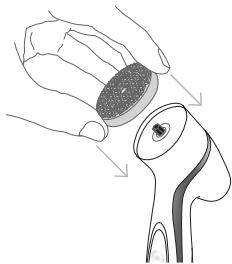

Attaching your ageLOC® LumiSpa™ treatment head

- Grip the sides of the treatment head.
- Align the hole on the back of the treatment surface with the black rotating axle on the LumiSpa handset.
- Gently press the head onto the axle until it clicks.

2

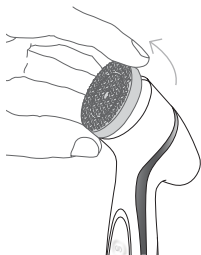

Removing your ageLOC LumiSpa treatment head

- Grip the sides of the treatment head.
- Gently lift the top of the treatment head and pull until it releases.

USING YOUR AGELOC® LUMISPA™


ageLOC LumiSpa is a twice daily regimen. Use it morning and night for best results.



STEP 1

- Moisten face with water.
- Apply ageLOC LumiSpa treatment cleanser to all areas of the face.

STEP 2

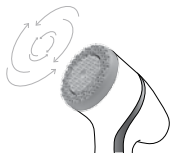

- Wet treatment head by placing it under running water. Then turn on device by pressing the On/Off button.
- Each of the four lights around the center button will flash once and the treatment head will immediately begin to counter-oscillate, first slowly, then up to normal speed.
- Slowly glide the treatment head back and forth over one of the areas of your face using broad strokes. Do not scrub or move across the face quickly.
- The device will do three brief stops every 30 seconds to prompt you to move from one area of the face to the next.



NOTE: If at any time you would like to pause your treatment, press the On/Off button once. To turn off the device, press and hold the On/Off button down for two seconds.

Warning #1: If you press the device too firmly onto your skin, the counter-oscillating motion will stop and the device will vibrate once slowly. Lift the device slightly to resume treatment.

Warning #2: If you scrub your face too aggressively with the device, the counter-oscillating motion will stop and the device will vibrate quickly several times. Glide the device across your face using slow, broad strokes to resume treatment.

NOTE: ageLOC LumiSpa is safe to use in the shower or wet environment. However, the charging base should not be exposed to water.

NOTE: ageLOC LumiSpa is suitable for makeup removal. However, do not use the device around the eye area.



STEP 3

- When the treatment is done, the device will stop.
- The device will turn off automatically.
- Rinse your face with water to remove residual treatment cleanser.


CLEANING AND CARING FOR YOUR AGELOC® LUMISPA™ DEVICE



AFTER EVERY TREATMENT:

- Remove the treatment head from the device. Rinse it with water while rubbing to remove residual treatment cleanser. Dry thoroughly.


- Rinse device under water.

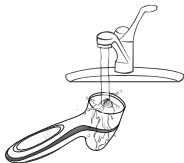

- Wipe device dry with washcloth.

- Only reattach the treatment head to the device after both components are thoroughly dry.

- Set device in charging base when it is completely dry.

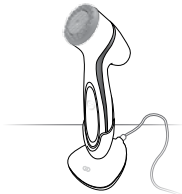

TIPS

- To ensure optimal performance, replace your treatment head every three months.

- The lights on the front of the handset indicate battery power. An amber light filling only the bottom third of the handset signifies a low battery.

DEVICE DISPOSAL

You must dispose of ageLOC LumiSpa properly according to local laws and regulations. Because ageLOC LumiSpa contains electronic components and a NiMH (Nickel-Metal Hydride) battery pack, it must be disposed of separately from household waste. When ageLOC LumiSpa reaches its end of life, contact local authorities to learn about disposal and recycling options.



Dispose of batteries properly for your area.

REPLACEMENT AND WARRANTY INFORMATION

Limited Two-Year Warranty: Nu Skin guarantees your device to be free from defects in materials and workmanship for a period of two years from the original date of consumer purchase. This warranty does not cover damage to the product as a result of misuse or accident. If the product becomes defective within the two-year warranty period, please call your local Nu Skin Support Services to arrange for a replacement.

PATENTS

Numerous US and international patents issued and pending

©2017 NSE PRODUCTS, INC.
75 WEST CENTER STREET, PROVO, UT 84601 NUSKIN.COM
1-800-487-1000



NU SKIN.®