## TABLE OF CONTENTS

ENGLISH	1
FRANÇAIS	15
ESPAÑOL	29
中文简体	43

# ENGLISH

## **CONGRATULATIONS ON PURCHASING YOUR TRIPOLLAR® STOP EYE2 DEVICE!**

This user manual describes the TriPollar STOP EYE2 device and how to use it.

Supplementary information to be found at:  
[www.tripollar.com](http://www.tripollar.com)

## TABLE OF CONTENTS

ABOUT TRIPOLLAR RF TECHNOLOGY	6
TRIPOLLAR STOP EYE2 KIT INCLUDES	7
OVERVIEW OF THE TRIPOLLAR STOP EYE2 DEVICE	7
TREATMENT PROCEDURE	9
TROUBLESHOOTING	11
TECHNICAL SPECIFICATIONS	11
WARRANTY & CUSTOMER SUPPORT	12
APPENDIX	14

## IMPORTANT NOTICE

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice. Manufactured by Pollogen Ltd. Patented technology.

The names Pollogen, TriPollar are registered trademarks of Pollogen Ltd.

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

## SAFETY

**Important:** Be sure to read the entire user manual prior to use of the product. Always adhere to all warnings and cautions. Be sure to follow all safety precautions below when using the TriPollar STOP EYE2 product. To learn more about TriPollar RF technology, visit our website at: [www.tripollar.com](http://www.tripollar.com).

### ⓘ **WARNINGS AND CAUTIONS**

TriPollar STOP EYE2 is an electrical device that must be used properly. To reduce the risk of injury from electric shock, adhere to all warnings and cautions below:

- Do not touch the device should it fall into water.
- Do not place or store the device near any source of water, such as a sink or bathtub.

- Do not submerge the device under water.
- When plugged in, do not leave the device unattended.
- Before use, ensure that the system and its USB cord do not show any visible damage. Do not operate if the electric cord is damaged or if the system is disassembled.
- Only use the device with the electric cord provided.
- Keep the device out of the reach of children.
- Be sure that your hands are completely dry when operating the device.
- Turn the device OFF after use or before cleaning and maintenance.
- Use TriPollar Preparation Balm during treatment according to its instructions for use.

#### **SAFETY INSTRUCTIONS**

- Only use this device for its intended purpose, as described in this manual.
- Do not use the device over upper-eyelids, closed eyes, thyroid area, lips, ears, breasts or genital areas.
- The sensation while using the device should be warm and pleasant. If you experience discomfort, discontinue use immediately.

Failure to follow the instructions in this manual may cause damage or burns to your skin.

#### **DO NOT USE THE DEVICE IF ANY OF THE FOLLOWING APPLY TO YOU:**

- Are under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Have metal implants in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which present no problem.
- Currently have cancer or a history of cancer, especially skin cancer or premalignant moles.
- Have severe concurrent conditions, such as cardiac disorders.
- Are pregnant or nursing.
- Have an impaired immune system due to immunosuppressive diseases such as AIDS, or are using immunosuppressive medications.

- Have any disease that is stimulated by heat, such as recurrent Herpes Simplex (within the treatment area).
- Have any active skin condition within the treatment area, such as sores, psoriasis, eczema or a rash.
- Have a history of keloid scarring, abnormal wound healing or fragile skin.
- Have had any surgical procedure, laser resurfacing or deep chemical peeling within the treatment area within the past three months or are still undergoing a process to complete healing.
- Have permanent fillers within the treatment area. If you have temporary dermal fillers within the treatment area, consult your physician before using the device.
- Have had Botox administered in the treatment area within the past three months.
- Have suspicious lesions within the treatment area.
- Currently use medications, herbal preparations, vitamins or food supplements that may cause fragile or dry skin.

## **POSSIBLE SIDE EFFECTS**

The following side effects are rare. However, should they occur, they are usually temporary. If any of the following conditions appear and last more than 24 hours, consult your physician:

- Pain
- Damage to natural skin texture (crusting, blistering, burns, fragile skin, bruising)
- Mild headache

For other aesthetic procedures or should you have any uncertainty regarding potential side effects, consult your physician before using the device.

## ABOUT TRIPOLLAR MULTI-RF TECHNOLOGY

TriPollar MULTI-RF technology delivers safe and effective skin renewal and tightening. TriPollar home-use devices bring professional level aesthetic results from the professional clinic to the comfort of your home. TriPollar treatments produces real, noticeable immediate and long-lasting results. TriPollar home-use devices focuses four low power beams deep into the skin to stimulate dermal activity, increasing and accelerating the natural production of collagen, (the protein that is responsible for our skins' supple texture). The increase in the skin's collagen results in tighter and suppler skin.

**The TriPollar STOP EYE2 Intended Use is as follows:**

The TriPollar STOP EYE2 system is intended for non-invasive skin rejuvenation at home for delicate facial areas, specifically indicated for skin tightening in the Periorbital, Perioral and Nasolabial folds.

## TRIPOLLAR STOP EYE2 KIT INCLUDES

1. TRIPOLLAR STOP EYE2 DEVICE (FIG.1)
2. TRIPOLLAR STOP EYE2 USB CHARGING CORD (FIG.2)
3. TriPollar PREPARATION BALM (FIG.3)  
After the container has been opened, the gel can be used for up to 24 months.
4. TRIPOLLAR STOP EYE2 USER MANUAL AND WARRANTY POLICY (FIG.4)
5. TRIPOLLAR STOP EYE2 QUICK START GUIDE (FIG.4)

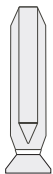
FIG.1

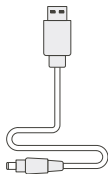
FIG.2


FIG.3

FIG.4

## OVERVIEW OF THE TRIPOLLAR STOP EYE2 DEVICE

- |                   |                                    |
|-------------------|------------------------------------|
| 1. DEVICE BODY    | 6. ON/OFF BUTTON AND INDICATOR LED |
| 2. TREATMENT HEAD | 7. GREEN INDICATOR LED             |
| 3. ELECTRODE      | 8. BLUE INDICATOR LED              |
| 4. CONTROL PANEL  | 9. PROTECTIVE COVER                |
| 5. USB CONNECTION |                                    |

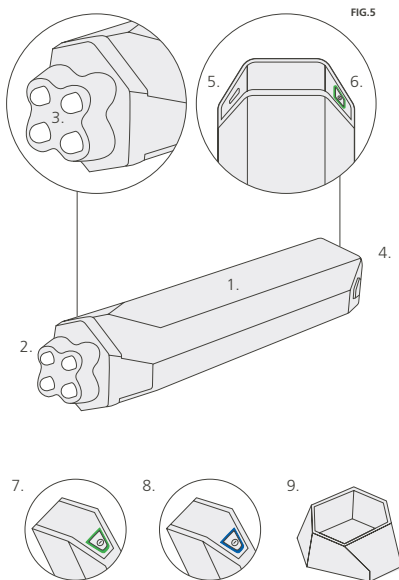**THE INDICATOR LEDS ON THE TRIPOLLAR STOP EYE2 DEVICE WORK AS FOLLOWS:**

**GREEN INDICATOR LED:**

- Indicates that the device is ON and ready for treatment.
- Slow blinking - low battery, approximately 5 minutes remaining.
- Fast blinking - battery soon to be depleted, approximately 2 minutes remaining.

**BLUE INDICATOR LED:**

- Blinking- battery charging.
- Steady- battery is fully charged.

## TREATMENT PROCEDURE

### COURSE OF TREATMENT

1. Treatment duration is approximately 4 minutes per area.
2. For optimal results, use your TriPollar STOP EYE2 device 2-3 times a week for a period of 6-8 weeks.
3. Make sure you have at least one day interval between treatments.
4. In order to maintain long term results after you have finished the full course of treatments, it is recommended that you use your TriPollar STOP EYE2 device 2-4 times a month or as needed.

### TREATMENT PROCEDURE

1. Charge the device so that its battery is fully charged and its LED indicator is steady blue.
2. Cleanse and dry the treatment area.
3. TriPollar Preparation Balm must be used throughout treatment. Spread an even layer over the entire treatment area. In case you feel that the balm has been absorbed during the treatment, spread an additional layer of TriPollar Preparation Balm.
4. Remove the protective cover from the device. Turn TriPollar STOP EYE2 ON by short pressing the Power Button (FIG.6). Make sure the green indicator LED appears and the Micro Massaging effect is felt.

FIG.6

FIG.7

5. Place the device on your skin ensuring that all 4 electrodes are touching the skin. For best results, all four electrodes must be in full contact with the skin during treatment.

6. Start moving the device in slow motions over the treatment area (FIG.7). You will feel pleasant warming and Micro Massaging.

**IMPORTANT:** be sure to keep the device constantly moving on your skin. Do not leave the device on one place without movement, as skin overheating may occur.

7. Treat each area for 4 minutes.

#### **POST TREATMENT**

1. Turn the device OFF by pressing the power button for two seconds.
2. Gently remove the TriPollar Preparation Balm from the skin.
3. Clean the TriPollar STOP EYE2 device with a soft, damp cloth or paper towel and cover it with its protective cover. Do not submerge the TriPollar STOP EYE2 under water!

**IMPORTANT:** It is recommended to perform a patch test 24 hours prior to initial use. Apply a small amount of TriPollar Preparation Balm on the inner arm and treat the area with TriPollar STOP EYE2 device. Ensure that the skin does not have a negative reaction to the gel or treatment.

#### **INDICATIONS OF SUCCESSFUL TREATMENT**

As a natural result of the treatment, the skin may show some redness (erythema) and minor skin swelling (edema) in the treatment area. Both should disappear after a short time.

## TROUBLESHOOTING

In the event that your TriPollar STOP EYE2 device is not working, prior to contacting the service center, take the following precautions:

1. Make sure you have removed the device's protective cover before use.
2. Make sure the device's battery is charged.
3. Make sure power button was pressed and the green indicator LED appears.
4. Make sure you have applied the TriPollar Preparation Balm according to the instructions. TriPollar STOP EYE2 will not work properly if TriPollar Preparation Balm is not applied.

## TECHNICAL SPECIFICATIONS

### Technical Specifications

Weight:	60gr
Physical Dimensions:	H=131mm L=26mm W=26mm
TriPollar RF Alternating Frequencies:	0.9 MHz / 1.1 MHz

### Operating conditions

Temperature:	+5°C to +30°C
Humidity:	20RH to 80RH
Rating:	DC 5V 2A

### Charging

When battery is fully charged (after approximately 3 hours), the LED indicator will emit a steady blue.

### Transport and Storage

Temperature:	0°C to 40°C
Humidity:	10RH to 90RH

## WARRANTY AND CUSTOMER SUPPORT

### WARRANTY

#### **Congratulations on your purchase of this TriPollar STOP EYE2 device.**

This product is powered by TriPollar RF technology, the most advanced RF technology in the market, and has been subject to strict quality controls performed by the manufacturer. Pollogen Ltd. guarantees proper functioning of your product for the specified period outlined in this warranty. If you use the device properly, and follow all instructions outlined in this user manual, you can enjoy the product to your ultimate satisfaction.

This policy warrants the support that Pollogen Ltd. provides for the device in the event of failures attributable to manufacturing defects. In such case, you shall be entitled during the warranty period to the repair and replacement of parts, free of charge (regarding both service and parts) to the extent necessary to be able to make use of the product for its intended use. This warranty is valid for the product marketed by Pollogen Ltd.

#### **The warranty is not valid when:**

- There are amendments in the warranty policy, sales receipt or invoice or the absence of any of these original documents.
- There is damage or deterioration of the device due to improper use or if the device has been modified or repaired by an unauthorized third party.
- There is damage of the product due to defects incurred during transportation, such as bumps, bangs or improper handling.
- There is damage due to humidity or any kind of liquids, foreign bodies inside the device or the use of unapproved, unadvised Accessories with the device contradictory to the user manual.

- The device is not used in accordance with the respective user manual or with its supplied accessories.
- There is partial or total damage or loss of the device, due to natural disasters (such as earthquakes, floods, fire, lightning and so on), force majeure or accident.
- There is damage to the device after being subjected to high dust environments, humidity or application of Excessive voltage from the power supply.
- The product's serial number does not match the one in the policy or it has been amended or erased.
- There is cosmetic damage, including but not limited to scratches and dents.

#### **CUSTOMER SUPPORT**

For service, please refer to the store where you purchased the TriPollar STOP EYE2 or visit our website for a store locator: [www.tripollar.com](http://www.tripollar.com).

When reporting to local service center, please have the device serial number ready. The serial number can be found on the top panel of the device.

Customer satisfaction is our top priority. To help us provide you with the best support possible, please send us your comments and suggestions. Contact details can be found on the back cover of this user manual or via our website at: [www.tripollar.com](http://www.tripollar.com).

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This table is prepared in accordance with the provisions of SJ/T 11364

- Indicates that said hazardous substance contained in all of the homogeneous materials for this part is below the limit requirement of GB/T 26572
- x Indicates said that hazardous substance contained in at least one of the homogeneous materials used for this part may be above the limit requirement of GB/T 26572

The restricted materials are used in accordance with manufacturing needs. It has no impact on device safety, and no harmful to users. The declaration is to remind end users that the device need to be disposed reasonably and environment friendly.

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