ENGLISH

Congratulations on purchasing your TriPollar Stop VX2 device!

You have taken the first step toward flawless beautiful skin!
This user manual describes the TriPollar Stop VX2 product and how to use it.

Supplementary information to be found at:
www.tripollar.com | cn.pollogen.com

TABLE OF CONTENTS

INTRODUCING TRIPOLLAR STOP VX2	6
• OVERVIEW	
• THE TRIPOLLAR STOP VX2 KIT	
USING TRIPOLLAR STOP VX2	8
• THE TRIPOLLAR STOP VX2 DEVICE	
• BEFORE YOU BEGIN	
• TREATING YOUR SKIN	
• RECOMMENDED TREATMENT GUIDELINES	
TROUBLESHOOTING	15
CUSTOMER SUPPORT	15
TECHNICAL SPECIFICATIONS	17
APPENDIX	18

IMPORTANT NOTICE

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice. Manufactured by Pollogen Ltd.

Patented technology.

The names Pollogen, TriPollar and Stop are registered trademarks of Pollogen Ltd.

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

SAFETY

Important: Be sure to read the entire user manual prior to use of the product. Always adhere to all warnings and cautions. Be sure to follow all safety precautions below when using the TriPollar Stop VX2 product. To learn more about TriPollar technology, visit our website at: www.tripollar.com | cn.pollogen.com.

WARNINGS AND CAUTIONS

TriPollar Stop VX2 is an electrical device that must be used properly. To reduce the risk of injury from electric shock, adhere to all warnings and cautions below:

- Do not touch the device should it fall into water. Unplug the device immediately with dry hands.
- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not submerge the device under water.
- When plugged in, do not leave the device unattended.
- Before use, ensure that the system and its power supply do not show any visible damage. Do not operate if the electric cord is damaged or if the system is disassembled.
- Only use the device with the electric cord and power supply provided.
- Keep the device out of the reach of children.
- Be sure that your hands are completely dry when operating the device.
- Turn the device OFF and then unplug the power supply from the electrical outlet after use or before cleaning and maintenance.
- Use the device for personal use only.

- Use only TriPollar Preparation Gel during treatment according to its instructions for use.
- Do not use TriPollar Preparation Gel after its expiration date.

SAFETY INSTRUCTIONS

- Only use this device for its intended purpose, as described in this manual.
- Do not use the device over eyelids, closed eyes, thyroid area, lips, ears, breasts or genital areas.
- The sensation while using the device should be warm and pleasant. If you experience discomfort, discontinue use immediately & consult service.
- The use of this product may cause discomfort or adverse effects on people with allergic constitution, people with skin diseases and allergic dermatitis, people with sensitive skin, etc. If adverse effects occur during use, you should stop using this appliance immediately.
- In case the device or the power supply generates excessive heat disconnect power supply from the electrical outlet & consult service.

Failure to follow the instructions in this manual may cause damage or burns to your skin.

DO NOT USE THE DEVICE IF ANY OF THE FOLLOWING APPLY TO YOU:

- Are under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Have metal implants in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which present no problem.
- Currently have cancer or a history of cancer, especially skin cancer or premalignant moles.
- Have severe concurrent conditions, such as cardiac disorders.
- Are pregnant or nursing.
- Have an impaired immune system due to immunosuppressive diseases such as AIDS, or are using immunosuppressive medications.
- Have any disease that is stimulated by heat, such as, but not limited to, recurrent Herpes Simplex (within the treatment area).
- Have any active skin condition within the treatment area, such as, but not limited to, sores, psoriasis, eczema or a rash.
- Have a history of keloid scarring, abnormal wound healing or fragile skin.
- Have had any surgical procedure, laser resurfacing or deep chemical peeling within the treatment area within the past three months or are still undergoing a process to complete healing.
- Have permanent fillers within the treatment area.
If you have temporary dermal fillers within the treatment area or threads, consult your physician.
- Have had Botox administered in the treatment area within the past three months.
- Have suspicious lesions within the treatment area.
- Currently use medications, herbal preparations, vitamins or food supplements that may cause fragile or dry skin.

POSSIBLE SIDE EFFECTS

The following side effects are rare. However, should they occur, they are usually temporary. If any of the following conditions appear and last more than 24 hours, consult your physician:

- Pain
- Damage to natural skin texture (crusting, blistering, burns, fragile skin, bruising)
- Mild headache

For other aesthetic procedures or should you have any uncertainty regarding potential side effects, consult your physician before using the device.

INTRODUCING TRIPOLLAR STOP VX2

OVERVIEW

TriPollar Stop VX2 presents the **latest in RF technology at the comfort of your home**. TriPollar Stop VX2 device is intended for increasing and accelerating collagen production in the skin to achieve skin tightening and rejuvenating effect.

TECHNOLOGY

The revolutionary TriPollar Stop VX2 combines two effective technologies - the groundbreaking TriPollar Multi-RF technology (RF) and Dynamic Muscle Activation (DMA) that work simultaneously to deliver non-invasive collagen remodeling and muscle activity, as well as skin tightening and treatment of wrinkles.

HOW DOES IT WORK?

TriPollar Stop VX2 affects homogeneously all layers of the dermis and the superficial muscle layer, for faster, safer and more effective treatment. It delivers aesthetic results that are both immediate and long-lasting.

TriPollar Stop VX2 reaches the effective temperature quickly with its proprietary Turbo function. Thanks to this unique function shorter treatment would still have the best results.

TriPollar Stop VX2 users enjoy the combined benefits of a unique non-conductive gel that ensures the delivery of focused RF energy into the right layers of the dermis using the newly designed electrodes. The temperature sensor ensures an effective and safe treatment while stimulating dermal activity that accelerates the natural production of collagen, the protein responsible for the skin's smooth and supple texture.

Simultaneously, DMA targets the SMAS layer (superficial muscular aponeurotic system), toning and lifting it by stimulating gentle muscle contractions inducing a "face lift" effect.

Using TriPollar Stop VX2 regularly according to its recommended treatment can effectively postpone the signs of aging.

THE TRIPOLLAR STOP VX2 KIT

The TriPollar Stop VX2 kit contains the following components:

1. TRIPOLLAR STOP VX2 DEVICE (FIG.1)
2. TRIPOLLAR STOP VX2 POWER SUPPLY AND ELECTRIC CORD (FIG.2)
3. TRIPOLLAR PREPARATION GEL (FIG.3), MAY BE PROVIDED SEPARATELY IN SOME REGIONS.
4. After the container has been opened, the gel can be used for up to 12 months.
Important: It is recommended to perform a patch test 24 hours prior to initial use. See page 10 for details.
5. TRIPOLLAR STOP VX2 USER MANUAL AND INTERNATIONAL WARRANTY (FIG.4)

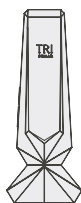
FIG.1

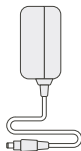
FIG.2

FIG.3

FIG.4

USING TRIPOLLAR STOP VX2

THE TRIPOLLAR STOP VX2 DEVICE

TriPollar Stop VX2 focuses RF (Radio Frequency) energy deep into the dermis. The TriPollar Stop VX2 DMA technology targets the SMAS layer, as shown below.


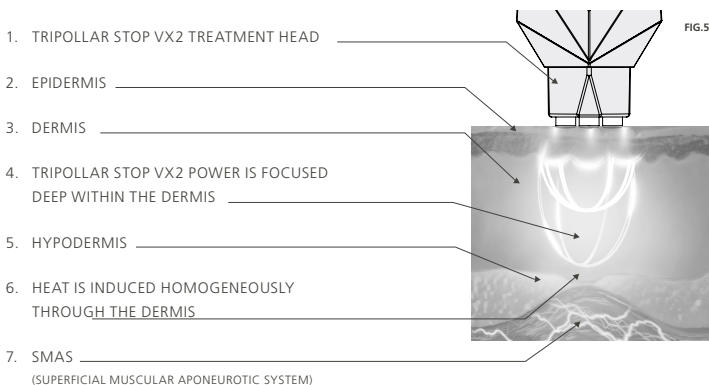
FIG. 6 SHOWS THE COMPONENTS COMPRISING THE TRIPOLLAR STOP VX2 DEVICE

FIG. 6

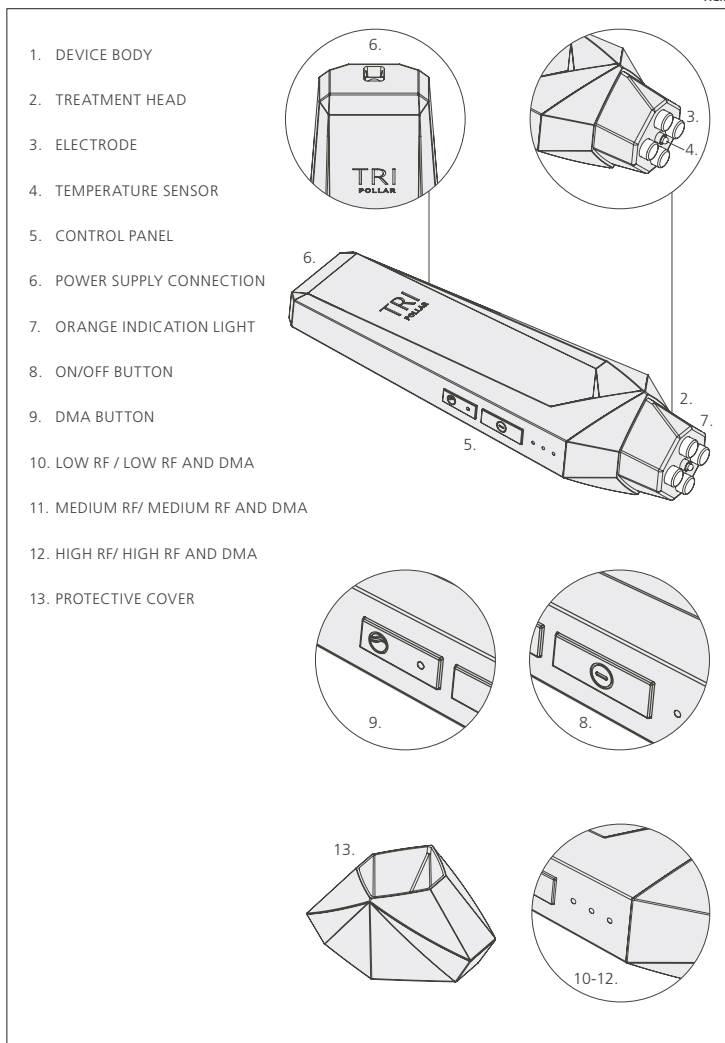BEFORE YOU BEGIN

Before using the TriPollar Stop VX2 device, do the following:

1. Open the box containing the device, and remove the device and its power supply.
2. Connect the power supply to the device (FIG.7)
3. Plug the power supply into an electrical outlet (FIG.8)
4. Remove the TriPollar Preparation Gel from its box.

Important: Before using the device for the first time, it is recommended to perform a patch test 24 hours prior to initial use. To do so, apply a small amount of TriPollar Preparation Gel to the inner arm and treat the area with the TriPollar Stop VX2 device. Ensure that the skin does not have a negative reaction to the gel or treatment.

5. Remove the protective cover from the device.

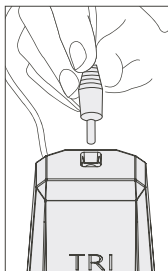
FIG.7

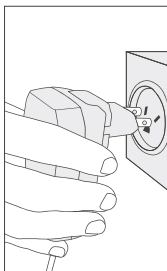
FIG.8

TREATING YOUR SKIN

To treat your skin using the TriPollar Stop VX2 device:

1. Cleanse and dry the treatment area.
2. Apply a very thin, even layer of TriPollar Preparation Gel over the entire treatment areas. This gel must be used during the entire treatment session. Your skin may feel a warm sensation when the gel is applied.
3. Turn the device ON by pressing the Power button (FIG.9). **Important:** All Leds should turn on briefly, followed by a steady single white led – this insures all Leds are operational. **Note:** A long press of the Power button turns the device OFF.
4. Select the RF power setting (Low/Medium/High) that is most comfortable for you. When the white indicator light shows a single white light illuminated (FIG.10), the device is operating at the Low power level. You can press the Power button again to change the power level (Fig.11, Fig.12) the white LEDs will light up according to the level chosen.

Important: During your first treatment session, it's recommended to use the Low power setting and then adjust the power level according to your skin's tolerance for subsequent treatments.

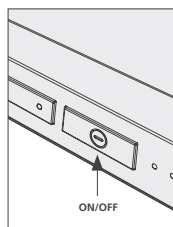
FIG.9

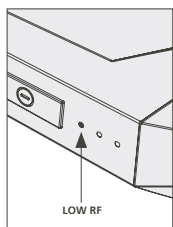
FIG.10

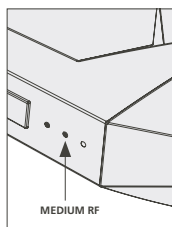
FIG.11

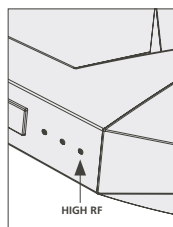
FIG.12

Treatment:

5. Divide the face to 3 treatment areas (FIG.13):
 - Right cheek
 - Left cheek
 - Forehead
6. Place the device on your skin, make sure to apply a slight pressure in order to assure a full contact between the electrodes and the skin.
7. Once the device is placed on the skin, the Turbo RF will start. Turbo RF will be initiated once at the beginning of the treatment or if a new level or mode is chosen.


8. Start moving the device in a very slow circular motion on a small part of your treatment area (approximately twice the size of the treatment head), slowly heating your skin to the optimal temperature (FIG.14).

Important: Be sure to constantly move the device slowly over the skin to ensure proper heating of the treated area. Do not keep the device on one place without movement.

9. Once the orange indicator light illuminates followed by a 'beep' sound, you have reached the optimal temperature. Continue the same slow circular motion for a few moments and then gradually enlarge your movements in the treatment area. Proceed with this technique covering the entire treatment area, making sure the orange indicator light is illuminated as much as possible during the treatment. As you move from one treatment area to another (FIG.13), the orange indicator light may turn-off if the temperature of the skin is not warm enough.
10. Continue working on the entire treatment area in a slow circular motion **for 3-4 minutes**. Make sure to maintain a uniform heating while covering the entire treatment area. Treatment can be augmented by "lifting" movements from the jaw line upwards toward your cheeks.
11. Once you finished the 3-4 minutes of treatment, you can proceed to the next treatment area repeating steps 6-10.
12. For jaw line treatment switch to RF + DMA mode by pressing on the DMA button . Make sure that the Blue indicator light turns ON. The DMA level will start to operate according to the last RF level chosen. You can press the Power button again to change the RF+DMA power level, the white LEDs will light up according to the level chosen. It's recommended to start at the Low RF+DMA.

Important: Use the RF + DMA mode ONLY on the jaw line area (FIG.15).

Do not use on other facial areas.

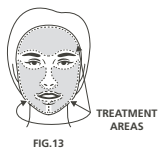13. Start moving the device using a very slow linear motions along the jaw line.
14. To ensure optimal results, use the device for 3-4 minutes on each side. After the entire treatment area has been treated, the treatment session is complete. As a natural result of the treatment, the skin may show some redness (erythema) and minor skin swelling (edema) in the treatment area. Both should disappear after a short time.
15. Press the Power button for two seconds to turn the device OFF.
16. Unplug the power supply from the electrical outlet.
17. Gently remove the TriPollar Preparation Gel from your skin and wash your face thoroughly. You can apply a moisturizer if desired.
18. Clean the electrodes of the device with a damp cloth or paper towel and cover the device tip with its protective cover.

Warning: Do not submerge the TriPollar Stop VX2 device under water.

SAFETY MECHANISM:

The device's automatic temperature-optimization mechanism protects the skin from overheating. When the skin has reached the maximal temperature, the energy source will shut down automatically.

Important:

The orange light is ON when the skin temperature in the treatment area is sufficiently heated and is OFF when it is not.

RECOMMENDED TREATMENT GUIDELINES

The following treatment guidelines are recommended for using the TriPollar Stop VX2 device:

- Treating the entire face takes up to 15 minutes. For shorter treatment times focus on selected areas of concern. Treatment times differ for each person.
- Facial DMA treatment takes approximately 10 minutes.
- For optimal results, use the TriPollar Stop VX2 device 2-3 times a week for a period of 6-8 weeks.
- Make sure that you allow at least one day in between treatments on the same treatment area.
- To maintain long-term results after finishing the full course of treatments, we recommend that you use your TriPollar Stop VX2 device 2-4 times a month or as needed.
- If the orange indicator light and sound indicators still do not appear within 5 minutes, or if you do not feel any heating, try moving the device as slowly as possible without stopping the motion. Slow, circular, constant movements are more effective than rapid movements. The slower the movement, the faster the treatment area heats internally.
- If the orange indicator light and sound indicators do not appear after 3–5 minutes, select a higher power level. If the indicators still do not appear after increasing to the High power setting, try moving the device more slowly, without stopping the motion.
- If the orange indicator light and sound indicators turn ON immediately or if you feel excessive heat, select a lower power level.

TROUBLESHOOTING

If your TriPollar Stop VX2 device is not working, take the following steps before contacting the service center:

- Ensure that you have **removed** the device's protective cover.
- Make sure that the electric cord is securely connected to the device and the Power supply adapter is firmly connected.
- Press the Power button and verify that the white indicator light illuminates.
- Ensure that you applied TriPollar Preparation Gel according to instructions described in this manual. The TriPollar Stop VX2 device does not work properly when the preparation gel is not applied.
- If orange, white or/and blue indicator lights do not illuminate when power button is pressed, contact the customer support.
- In case headlight LEDs blinking white or orange, discontinue use immediately & consult service.

CUSTOMER SUPPORT

For service, refer to the store where you purchased the TriPollar Stop VX2 device or visit our website at: www.tripollar.com | cn.polloegen.com. When contacting a local service center, have the device serial number ready. The serial number is located on the top panel of the device. Customer satisfaction is our top priority. To help us provide you with the best support possible, please send us your comments and suggestions. Contact information is provided on the back cover of this manual, or you can contact us via our website (see above).

WARRANTY TERMS

This product has been subjected to strict quality control performed by the manufacturer. For this reason, the distributor in your country guarantees proper function of your product for 24 months from date of actual purchase. If you use the device properly, follow all instructions outlined in the user manual, you can enjoy the product to your ultimate satisfaction.

This policy warrants the support that the distributor provides for the device, should it happen that failures attributable to manufacturing defects be found. In such case, you shall be entitled to make full use of the benefits of this warranty. This implies the reason and replacement of parts, free of charge (regarding both service and parts) at any Authorized Service Center. Within its validity term, you are entitled to make use of the warranty if necessary to have a fully enjoyable product, provided the above-indicated conditions are fully met. This warranty is valid for the device marketed by the distributor in your country.

The warranty is not valid when:

1. Amendments in the warranty policy, sales receipt or invoice, or the absence of any of these original documents.
2. Damage or deterioration of the device due to improper use or if the device has been modified or repaired by an unauthorized third party.
3. Damage of the product due to defects incurred during transportation, such as bumps, bangs or improper handling.
4. Damage due to humidity or any kind of liquids, foreign bodies inside the device or the use of unapproved, unadvised accessories with the device contradictory to the user manual.
5. When the device is not used in accordance with the respective user manual or with its supplied accessories.
6. Partial or total damage or loss of the device, due to natural disasters (such as earthquakes, floods, fire, lightning, etc.), force majeure or accident.
7. Damage to the device after being subjected to high dust environments, humidity or application of excessive voltage from the power inlet.
8. If the product's serial number has been amended or erased.
9. To cosmetic damage, including but not limited to scratches and dents.

TECHNICAL SPECIFICATIONS

TABLE 1

Technical Specifications

Weight:	94gm
Physical Dimensions:	H=154mm L=40mm W=60mm
TriPollar RF Frequency:	0.85, 1.0, 1.15 \pm 10% MHZ

TABLE 2

Operating conditions

Temperature:	+5C to +30C
Humidity:	20RH to 80RH
Rating:	100~240VAC, 0.4A, 50/60Hz VDC 8V 1.5A

TABLE 3

Transport and Storage

Temperature:	0°C to 40°C
Humidity:	10RH to 90RH

TABLE 4

Electrical Safety

Power supply is classified as Class II double insulation

APPENDIX

Restriction of Hazardous Substance Material Declaration

IN ACCORDANCE WITH THE PROVISIONS OF SJ/T 11364

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- x Indicates that the named hazardous substance contained in at least one of the homogeneous materials used for this component may be above the limit requirement of GB/T 26572. The restricted materials are used in accordance with manufacturing needs. Use of such materials has no impact on the device's safety and no harmful effect on device users. **This declaration is included to remind device users that the device must be disposed of reasonably and in accordance with the environmental policy of your locale.**
- Indicates that the named hazardous substance contained in all of the homogeneous materials for this component is below the limit requirement of GB/T 26572.

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