INTELLIGENT
BIOLOGY™

TRI

POLLAR®

USER MANUAL

ENGLISH

FRANÇAIS

ESPAÑOL

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DESIRE

FACIAL RENEWAL AND
REJUVENATION DEVICE

ENGLISH

Congratulations on purchasing your TriPollar® Desire device!

You have taken the first step toward flawless beautiful skin!
This user manual describes the TriPollar® Desire product and how to use it.

Supplementary information to be found at:
www.tripollar.com

TABLE OF CONTENTS

INTRODUCING TRIPOLLAR® DESIRE	5
• OVERVIEW	
• THE TRIPOLLAR® DESIRE KIT	
USING TRIPOLLAR® DESIRE	7
• THE TRIPOLLAR® DESIRE DEVICE	
• BEFORE YOU BEGIN	
• TREATING YOUR SKIN	
• RECOMMENDED TREATMENT GUIDELINES	
TROUBLESHOOTING	14
CUSTOMER SUPPORT	14
TECHNICAL SPECIFICATIONS	15
APPENDIX	16

IMPORTANT NOTICE

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice. Manufactured by Pollogen Ltd. Patented technology.

The names Pollogen and TriPollar are registered trademarks of Pollogen Ltd.

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

SAFETY

Important: Be sure to read the entire user manual prior to use of the product. Always adhere to all warnings and cautions. Be sure to follow all safety precautions below when using the TriPollar® Desire product.

To learn more about TriPollar technology, visit our English website at: www.tripollar.com

WARNINGS AND CAUTIONS

TriPollar® Desire is an electrical device that must be used properly. To reduce the risk of injury from electric shock, adhere to all warnings and cautions below:

- Do not touch the device should it fall into water. Unplug the device immediately with dry hands.
- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not submerge the device under water.
- When plugged in, do not leave the device unattended.
- Before use, ensure that the system and its power supply do not show any visible damage. Do not operate if the electric cord is damaged or if the system is disassembled.
- Only use the device with the electric cord and power supply provided.
- Keep the device out of the reach of children.
- Be sure that your hands are completely dry when operating the device.
- Turn the device OFF and then unplug the power supply from the electrical outlet after use or before cleaning and maintenance.
- Use TriPollar® Preparation Gel during treatment according to its instructions for use.

SAFETY INSTRUCTIONS

- Only use this device for its intended purpose, as described in this manual.
- Do not use the device over eyelids, closed eyes, thyroid area, lips, ears, breasts or genital areas.
- The sensation while using the device should be warm and pleasant. If you experience discomfort, discontinue use immediately.

Failure to follow the instructions in this manual may cause damage or burns to your skin.

DO NOT USE THE DEVICE IF ANY OF THE FOLLOWING APPLY TO YOU:

- Under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Have metal implants in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which present no problem.
- Currently have cancer or a history of cancer, especially skin cancer or premalignant moles.
- Have severe concurrent conditions, such as cardiac disorders.
- Pregnant or nursing.
- Have an autoimmune disease or impaired immune system due to immunosuppressive diseases such as AIDS, or are using immunosuppressive medications.
- Have any disease that is stimulated by heat, such as recurrent Herpes Simplex (within the treated area).
- Have any active skin condition within the treated area, such as sores, psoriasis, eczema or a rash.
- Have a history of keloid scarring, abnormal wound healing or fragile skin.
- Have had any surgical procedure, laser resurfacing or deep chemical peeling within the treated area within the past three months or are still undergoing a process to complete healing.
- Have permanent fillers within the treated area.
If you have temporary dermal fillers within the treated area, consult your physician before using the device.
- Have had Botox administered in the treated area within the past three months.
- Have suspicious lesions within the treated area.
- Currently use medications, herbal preparations, vitamins or food supplements that may cause fragile or dry skin.

POSSIBLE SIDE EFFECTS

The following side effects are rare. However, should they occur, they are usually temporary. If any of the following conditions appear and last more than 24 hours after treatment, consult your physician:

- Pain
- Damage to natural skin texture (crusting, blistering, burns, fragile skin, bruising)
- Mild headache

For other aesthetic procedures or if you have any uncertainties regarding potential side effects, consult your physician before using the device.

INTRODUCING TRIPOLLAR® DESIRE

This chapter introduces the TriPollar® Desire product and describes its components.

OVERVIEW

TriPollar® Desire – The Next Generation of Turbo RF Home Treatments

TriPollar® Desire presents the latest in RF (Radio Frequency) technology at the comfort of your home. Based on the groundbreaking TriPollar® Multi-RF technology, TriPollar® Desire affects all layers of the dermis homogeneously for faster, safer, more focused and effective treatment. It delivers professional-level aesthetic results that are both immediately visible and long-lasting. TriPollar Desire reaches the effective temperature quickly with its proprietary Turbo function. Thanks to this unique function shorter treatment would still have the best results. TriPollar® Desire users enjoy the combined benefits of a unique gel that ensures the delivery of focused RF energy into the right layers of the dermis using the newly designed electrodes. A proprietary and innovative temperature sensor ensures an effective and safe treatment while stimulating dermal activity that accelerates the natural production of collagen, the protein responsible for the skin's smooth and supple texture. Using TriPollar® Desire regularly according to its treatment protocol can effectively postpone the signs of aging.

THE TRIPOLLAR® DESIRE KIT

The TriPollar® DESIRE kit contains the following components:

1. TRIPOLLAR® DESIRE KIT (FIG.1)
2. TRIPOLLAR® DESIRE DEVICE (FIG.2)
3. TRIPOLLAR® DESIRE POWER SUPPLY AND POWER CORD (FIG.3)
4. TRIPOLLAR® DESIRE ELECTRICAL OUTLET ADAPTERS (FIG.4)
5. TRIPOLLAR® PREPARATION GEL (FIG.5)

After the container has been opened, the gel can be used for up to 12 months after first use.

Important: It is recommended to perform a patch test 24 hours prior to initial use. See page 9 for details.

6. TRIPOLLAR® DESIRE USER MANUAL AND WARRENTY POLICY (FIG.6)

FIG.1

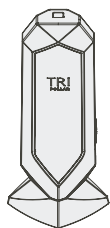
FIG.2

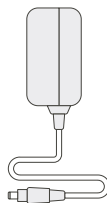
FIG.3

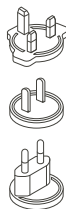
FIG.4

FIG.5

FIG.6

USING TRIPOLLAR® DESIRE

This chapter describes the TriPollar® Desire device and how to use it.

THE TRIPOLLAR® DESIRE DEVICE

TriPollar® Desire is a facial renewal and rejuvenation device. The device is used to increase and accelerate collagen production in your skin for rejuvenation effect.

TriPollar® Desire focuses heat deep within the dermis of your skin, as shown below.


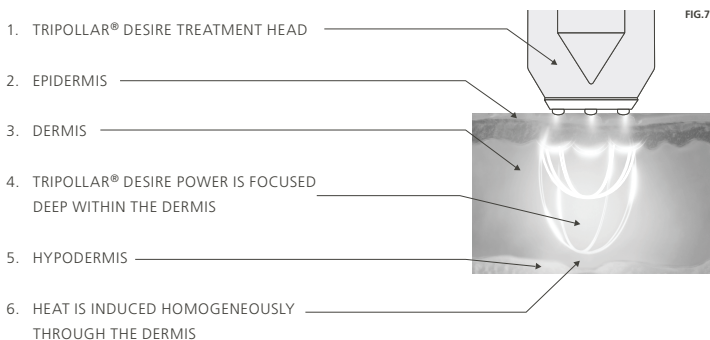
FIG.8 SHOWS THE COMPONENTS COMPRISING THE TRIPOLLAR® DESIRE DEVICE.

1. DEVICE BODY
2. TREATMENT HEAD
3. ELECTRODE
4. TEMPERATURE SENSOR
5. CONTROL PANEL
6. POWER SUPPLY CONNECTION
7. ORANGE INDICATOR LIGHT
8. ON/OFF BUTTON
9. LOW RF
10. MEDIUM RF
11. HIGH RF
12. PROTECTIVE COVER

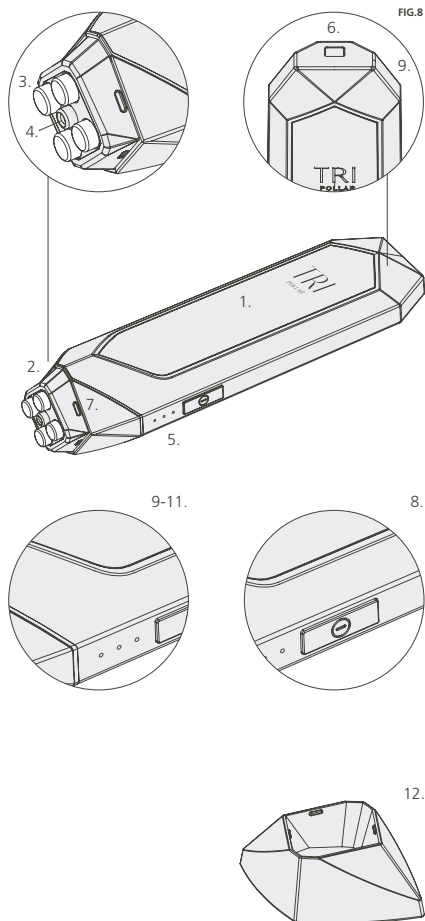BEFORE YOU BEGIN

Before using the TriPollar® Desire device, do the following.

1. Open the box containing the device, and remove the device and its power supply.
2. Insert the adapter into the slot of the power supply (Fig.9)
3. Slightly push and turn the adapter in the slot of the PS, until it securely "clicks into place" and cannot be removed easily (Fig. 10)
4. Connect the power supply to the device (Fig.11)
5. Plug the power supply into an electrical outlet (Fig.12)
6. Remove the TriPollar® Preparation Gel from its box.
Important: Before using the device for the first time, it is recommended to perform a patch test 24 hours prior to initial use. To do so, apply a small amount of TriPollar® Preparation Gel to the inner arm and treat the area with the TriPollar® Desire device. Ensure that the skin does not have a negative reaction to the gel or treatment.
7. Remove the protective cover from the device.

You are now ready to treat your skin using TriPollar® Desire, as described on page 10.

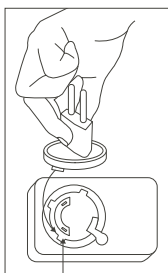
FIG.9

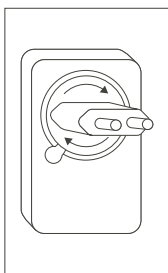
FIG.10

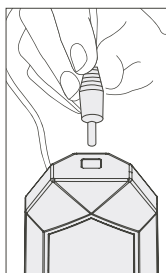
FIG.11

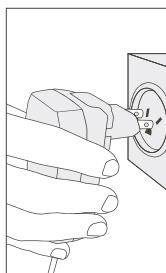
FIG.12

TREATING YOUR SKIN

To treat your skin using the TriPollar® Desire device:

1. Cleanse and dry the treated area.
2. Apply a very thin, even layer of TriPollar® Preparation Gel over the entire treatment area. This gel must be used during the entire treatment session. Your skin may feel a warm sensation when the gel is applied.
3. Turn the device ON by pressing the Power button (Fig.13). Make sure that the white led turns ON.
Tip: A long press of the Power button turns the device OFF.
4. Select the power level (Low/Medium/High) that is most comfortable for you. When the white led shows a single light illuminated (Fig.14), the device is operating at the Low power level. You can press the Power button again to change to the power level (Fig.15, Fig.16) the white LEDs will light up according to the level chosen.

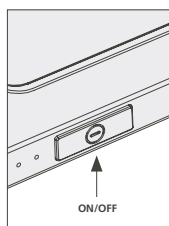
FIG.13

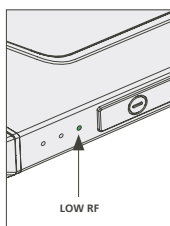
FIG.14

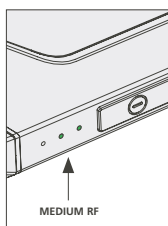
FIG.15

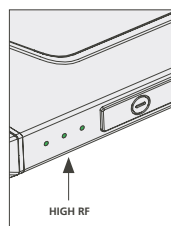
FIG.16


Important: During your first treatment session, we recommend that you use the Low power setting and then adjust the power level according to your skin's tolerance for subsequent treatments.

Treatment:

5. Divide the face into 3 optional treatment areas (Fig.17):
 - Right cheek
 - Left cheek
 - Forehead
6. Place the device on your skin, make sure to apply a slight pressure in order to assure a full contact between the electrodes and the skin and then start moving the device slowly in a circular motion.
7. Once the device is placed on the skin, the Turbo RF will start. Turbo RF will be initiated once at the beginning of the treatment or if a new level is chosen (Low RF, Medium RF, High RF).

Important: Be sure to constantly move the device slowly over the skin to ensure proper heating of the treated area. Do not leave the device in one place without movement.

8. Once the orange indicator light illuminates followed by a 'beep' sound, you have reached the optimal temperature. Continue the same slow circular motion for a few moments and then gradually enlarge your movements to cover a wider treatment area.
Proceed with this technique covering the entire treatment area, making sure the orange indicator light is illuminated as much as possible during the treatment.
9. Continue working on the entire treatment area in a slow circular motion **for 3-4 minutes**. Make sure to maintain a uniform heating while covering the entire treatment area. Treatment can be augmented by "lifting" movements from the jaw line upwards toward your cheeks (Fig.19).
10. Once you finished the treatment area, you can proceed to the next treatment area repeating steps 6-9.

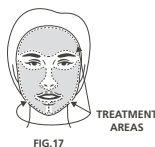11. As a natural result of the treatment, the skin may show some redness (erythema) and minor skin swelling (edema) in the treatment area. Both should disappear after a short time.
12. Press the Power button for two seconds to turn the device OFF.
13. Unplug the power supply from the electrical outlet.
14. Gently remove the TriPollar® Preparation Gel from your skin and wash your face. You can apply a moisturizer if desired.
15. Clean the electrodes of the device with a damp cloth or paper towel and cover the device tip with its protective cover.

Warning: Do not submerge the TriPollar® Desire device under water.

SAFETY MECHANISM:

The device's automatic temperature-optimization mechanism protects the skin from overheating. When the skin has reached the maximal temperature, the energy source will shut down automatically. The indication is - power level LEDs turn off. The device automatically turns the energy source back ON, when the skin temperature drops below the maximal temperature.

Important:

The orange light turns ON when the temperature of the skin in the treatment area is sufficiently heated.

DESCRIPTION	LED STATUS
Device is on/Low RF level	1 white lights
Device is on/Mid RF level	2 white lights
Device is on/High RF	3 white lights
Optimal temperature reached	Constant orange light
Overheat protection	White leds turn off + constant orange light

RECOMMENDED TREATMENT GUIDELINES

The following treatment guidelines are recommended for using the TriPollar® Desire device:

- Treating the entire face takes 12-15 minutes. For shorter treatment times focus on selected areas of concern. Treatment times differ for each person.
- For optimal results, use the TriPollar® Desire device 2-3 times a week for a period of 6-8 weeks.
- Make sure that you allow at least one day in between treatments on the same treatment area.
- To maintain long-term results after completing the full course of treatments, we recommend that you use your TriPollar® Desire device 2-4 times a month or as needed.
- If the orange light/ sound indicators do not appear within five minutes, or if you do not feel any heating, try moving the device as slowly as possible without stopping the motion. Slow, circular, constant movements are more effective than rapid movements. The slower the movement, the faster the treatment area heats internally.
- If the orange indicator light turn ON immediately or if you feel excessive heat, select a lower power level.
- If the orange indicator light and sound indicator do not appear after 3–5 minutes, select a higher power level. If the indicators still do not appear after increasing to the High power setting, try moving the device slower, without stopping the motion.

TROUBLESHOOTING

This chapter describes how to troubleshoot problems that may arise when using the TriPollar® Desire device.

If your TriPollar® Desire device is not working, take the following steps before contacting the service center:

- Ensure that you have removed the device's protective cover.
- Make sure that the electric cord is securely connected to the device and the Power supply adapter is firmly connected.
- Press the Power button and verify that the white LED indicator light illuminates.
- Ensure that you applied TriPollar® Preparation Gel according to instructions described in this manual. The Tripollar® Desire device does not work properly when the preparation gel is not applied.

CUSTOMER SUPPORT

This chapter provides customer support contact information for the Tripollar® Desire device.

For service, refer to the store where you purchased the TriPollar® Desire device or visit our website at: **www.tripollar.com**. When contacting a local service center, have the device serial number ready. The serial number is located on the top panel of the device. Customer satisfaction is our top priority. To help us provide you with the best support possible, please send us your comments and suggestions. Contact information is provided on the back cover of this manual, or you can contact us via our website (see above).

TECHNICAL SPECIFICATIONS

This chapter presents the specifications for the TriPollar® Desire system.

TABLE 1

Technical Specifications

Weight:	85gm
Physical Dimensions:	H=146mm L=40mm W=58mm
TriPollar RF Frequency:	0.82/1/1.115 MHz

TABLE 2

Operating conditions

Temperature:	+5°C to +30°C
Humidity:	20RH to 80RH
Rating:	100~240V AC, 1.5A, 50/60Hz DC 8V 1.5A

TABLE 3

Transport and Storage

Temperature:	-10°C to 60°C
Humidity:	MAX 95RH

TABLE 4

Electrical Safety

Power supply is classified as Class II double insulation

APPENDIX

This appendix provides the hazardous materials declaration for the TriPollar® Desire device.

Restriction of Hazardous Substance Material Declaration

TABLE 5 IS IN ACCORDANCE WITH THE PROVISIONS OF SJ/T 11364

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- x Indicates that the named hazardous substance contained in at least one of the homogeneous materials used for this component may be above the limit requirement of GB/T 26572. The restricted materials are used in accordance with manufacturing needs. Use of such materials has no impact on the device's safety and no harmful effect on device users. **This declaration is included to remind device users that the device must be disposed of reasonably and in accordance with the environmental policy of your locale.**
- Indicates that the named hazardous substance contained in all of the homogeneous materials for this component is below the limit requirement of GB/T 26572.

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