ENVIG

EDGE

Facial device  
for luminous  
and smooth skin

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**ENGLISH**

**Congratulations on purchasing your Edge device!**

You have taken the first step toward flawless beautiful skin!  
This user manual describes the Edge product and how to use it.

Additional information can be found at [ENVIG.COM](http://ENVIG.COM)

## TABLE OF CONTENTS

<b>SAFETY &amp; CAUTIONS</b>	<b>4</b>
<b>POSSIBLE SIDE EFFECTS</b>	<b>5</b>
<b>INTRODUCTION EDGE DEVICE</b>	<b>6</b>
• OVERVIEW	
• THE EDGE KIT	
• GETTING TO KNOW YOUR DEVICE	
<b>USING THE EDGE DEVICE</b>	<b>9</b>
• CHARGING INSTRUCTIONS	
• TREATING YOUR SKIN	
• RECOMMENDED TREATMENT GUIDELINES	
<b>TROUBLESHOOTING</b>	<b>12</b>
<b>WARRANTY AND CUSTOMER SUPPORT</b>	<b>12</b>
<b>TECHNICAL SPECIFICATIONS</b>	<b>14</b>
<b>APPENDIX</b>	<b>15</b>

## IMPORTANT NOTICE

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice. Manufactured by Pollogen Ltd. Patented technology.

**The names Pollogen and ENVIG are registered trademarks of Pollogen Ltd.**

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

## SAFETY & CAUTIONS

### READ ALL INSTRUCTIONS BEFORE USE

- Operate the device only according to its intended use.
- Use the device only for personal use.
- Keep out of the reach of children. Children, or youth under the age of 18, should not use or play with the device. The device has not been tested on children.
- DO NOT use the device if you have severe allergies, sensation disorders, infected skin, or a heart condition.
- DO NOT use if you have a pacemaker or other electrical medical device without consulting with your physician.
- DO NOT use the device if you are pregnant, lactating, or had given birth less than six months ago.
- Be careful using the device around facial jewelry or metal implants in the treatment area (Not including dental fillings, implants, and crowns).
- Do not immerse the device in water.
- If the product appears damaged, disassembled, broken or has been submerged in water, discontinue use.

### IMPORTANT FOR OPTIMUM SAFETY:

- If you have medical concerns consult your dermatologist before using the device. Also consult your dermatologist, if you have recently had a cosmetic procedure, acid or chemical peels, extreme sunburn, or thread lift.
- Apply particular care when using around or the under-eye area, and do not bring the device into contact with the eyelids or eyes themselves. DO NOT USE THE DEVICE on unintended regions of the body. The device is designed ONLY for facial use.
- DO NOT use the device on active acne, psoriasis, eczema, rosacea, open sores; skin with a rash, sunburned or distressed skin, and on skin surfaces with cuts and abrasions or any healing is present.

- DO NOT use the device on Tattoo or permanent makeup in the treated area.
- DO NOT use the device on moles, cancerous areas, skin lesions or on skin which has recently been treated with a laser or chemicals.
- DO NOT use the device if you use medications, herbal preparations, food supplements, or vitamins that might cause fragile skin or impaired skin healing.
- You may feel a tingling sensation while using the device and after the treatment.
- If you experience extreme discomfort, discontinue use immediately and consult a physician. DO NOT use the device if you experience a worsening of the skin (skin feels uncomfortable, irritated, burning, etc.).
- DO NOT use the device while charging.
- If the device or charging cable is damaged or appears damaged, do not use it. This product contains no serviceable parts.
- DO NOT use, place, or store the device or charging components where it can fall or be pulled into a tub or sink or in water.
- CAUTION: Before charging, make sure the charging socket and cable are free from water.
- For any questions or concerns, contact us at ENVIG.COM.
- Make sure the skin is completely dry before use. Do not treat on wet or sweaty skin.
- Do not use the device right after cosmetic peeling/exfoliation treatment.

## POSSIBLE SIDE EFFECTS

The following side effects are rare. However, should they occur, they are usually temporary. If any of the following, conditions appear and last more than 24 hours after treatment, consult your physician:

- Damage to natural skin texture
- Discomfort or pain

For other aesthetic procedures or if you have any uncertainties regarding potential side effects, consult your physician before using the device.

## INTRODUCING EDGE DEVICE

### OVERVIEW

The EDGE device is intended for facial exfoliation to improve the appearance of skin and give the skin a luminous and smoother look and feel.

Treatment time is short, allowing best results with minimal effort.

The technology has been proven to be safe and effective.

The technology is based on Pollogen's fractional resurfacing technology which was originally developed for professional clinics. Pollogen has adapted the technology for effective and safe use at the convenience of home setting, promoting healthy skin appearance.

## THE EDGE KIT

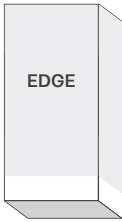
FIG.1

The Edge kit contains the following components:

1. EDGE DEVICE (FIG.2)
2. USB CABLE (FIG.3)
3. USER MANUAL AND warranty POLICY (FIG.4)

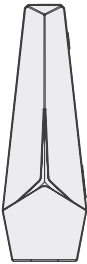
FIG.2

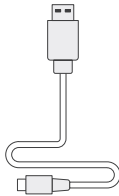
FIG.3


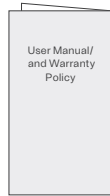
FIG.4

GETTING TO KNOW YOUR DEVICE

- 1. DEVICE BODY
- 2. TREATMENT HEAD
- 3. HEADLIGHT
- 4. ON/OFF & POWER BUTTON
- 5. TRIGGER BUTTON
- 6. CHARGING PORT + CHARGING INDICATOR LED
- 7. LOW POWER LEVEL
- 8. MEDIUM POWER LEVEL
- 9. HIGH POWER LEVEL
- 10. PROTECTIVE COVER

FIG.5
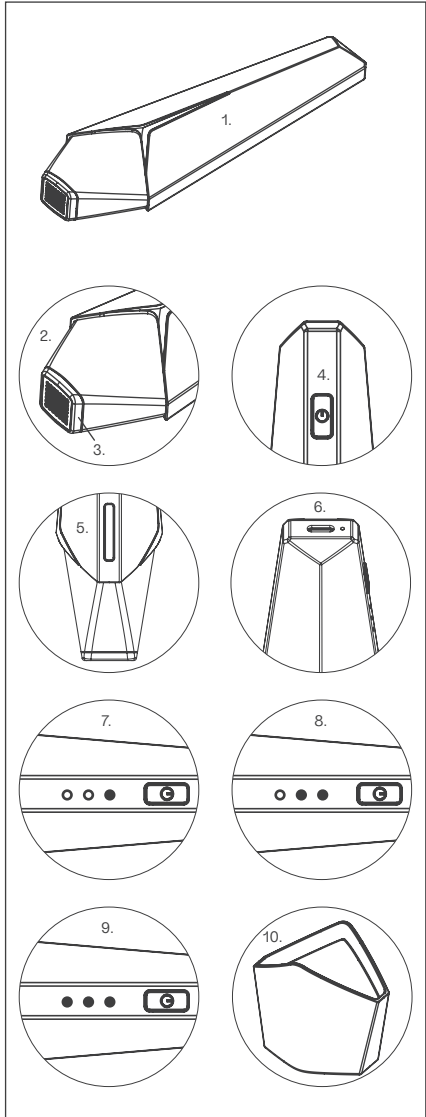
## USING THE EDGE DEVICE

### CHARGING INSTRUCTIONS

We recommend charging the device before each use.

1. Open the box containing the device and remove the device and its USB cable.
2. Connect your Device to a USB outlet to begin charging (Fig.6).
3. Once charging is completed disconnect the USB cord.

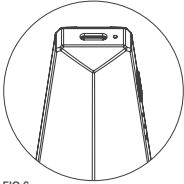
FIG.6

### LED INDICATORS

The LED next to the USB port indicates the current state of the battery.

LED COLOR	INDICATOR STATE
Blinking green	Battery Charging
Steady green	Battery Full
Blinking Red (slow)	Battery Low
Blinking Red (fast)	Battery Depleted

#### Notes:

- It should take approximately 1-2 hours to fully charge the device.
- A fully charged device will last approximately 30 minutes of continuous use.
- The device will not work, while charging.

## TREATING YOUR SKIN

To treat your skin using the Edge device:

1. Remove all jewelry
2. Cleanse and properly dry the treatment areas. Ensure the skin is free of topical products (such as makeup, moisturizers etc.).
3. Divide your face into 4 treatment areas (Fig.7)
  - Right cheek
  - Left cheek
  - Forehead
  - Perioral (around the mouth)

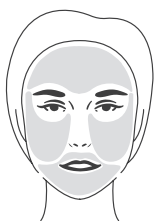
FIG.7

4. Turn the device ON by pressing the Power Button (Fig.8). Both the headlight and the low power level led (located next to the Power Button) will turn white.

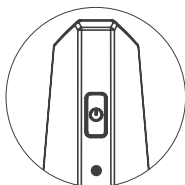
FIG.8

5. Press the Power button to change the power level (Fig.9), the white LEDs will light up according to the level chosen.
6. Select the power level (Low/Medium/High) that is most comfortable for you. It is recommended to start at the lowest level and increase per personal feeling.

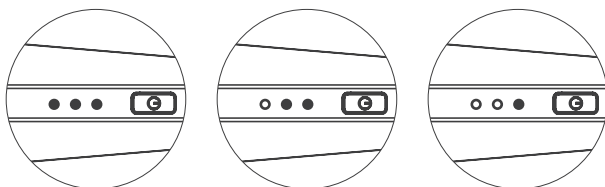
FIG.9

7. Place the Edge's treatment head flat and in full contact against your skin and then press the Trigger Button to emit energy. During the energy pulse, the headlight will briefly turn blue, accompanied by a vibration sensation. At the end of the pulse, the device will beep.
8. After the beep you can move to the adjacent location.
9. Proceed with this technique covering the entire treatment area, use methodological pattern to ensure the entire face is covered (Fig.10). There is no need to treat the same area more than once.

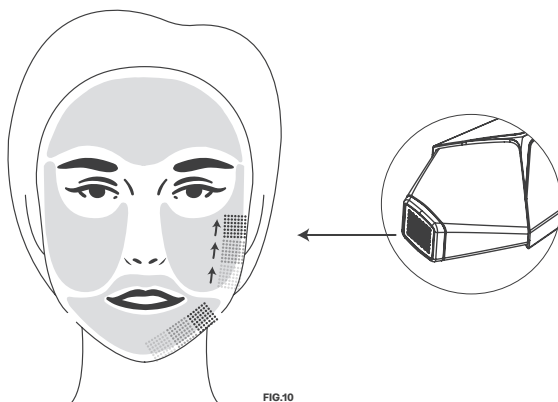
FIG.10

**Note:** you may use different power settings on different parts of your face depending on the sensitivity of the treatment area. Such as perioral or periorbital (around the mouth and the eyes).

10. The treatment is complete once you covered the entire face.
11. Press the Power Button for two seconds to turn the device OFF.
12. Clean the device treatment head with alcohol wipes and cover the device with its protective cover.

#### NOTES:

- You can easily change the power setting by pressing on the Power Button throughout the treatment.
- The skin may show some redness (erythema) and minor skin swelling (edema) in the treatment area. Both should disappear after a short time.
- You may feel tingling sensation while and after using the device.

## RECOMMENDED TREATMENT GUIDELINES

The following treatment guidelines are recommended for using the Edge device:

- Treating the entire face takes 3-4 minutes. For shorter treatment times, focus on selected areas of concern. Treatment times differ for each person.
- For optimal results, use the Edge device 2-3 times a week.
- After treatment, maintain your regular skincare routine.
- Make sure that you allow at least one day in between treatments on the same treatment area.
- Response to Edge treatment may differ between individuals depending on environmental and biological factors, heredity and aging.
- Please avoid direct sun exposure right after treatment.

## TROUBLESHOOTING

This chapter describes how to troubleshoot problems that may arise when using the Edge device.

**If your Edge device is not working, take the following steps:**

- Ensure that you have removed the device's protective cover.
- Ensure that the Power button is not stuck in and verify both headlight and low level led will turn white.
- Make sure the treatment level LED are working.
- Make sure that the device is fully charged (the led color is steady green)
- Ensure that the Trigger button is not stuck in and verify the Blue LED is illuminating followed by a soft vibration and "Beep" sound.

**If one of the above is not working properly contact the service center**

## WARRANTY AND CUSTOMER SUPPORT

### WARRANTY

#### CONGRATULATIONS ON YOUR PURCHASE OF THIS EDGE DEVICE.

This product has been subjected to strict quality control performed by the manufacturer. For this reason, the distributor in your country guarantees proper function of your product for 24 months from date of actual purchase. If you use the device properly, follow all instructions outlined in the user manual, you can enjoy the product to your ultimate satisfaction. This policy warrants the support that Pollogen Ltd. provides for the device in the event of failures attributable to manufacturing defects. In such case, you shall be entitled during the warranty period to the repair and replacement of parts, free of charge (regarding both service and parts) to the extent necessary to be able to make use of the product for its intended use. This warranty is valid for the product marketed by Pollogen Ltd.

**The warranty is not valid when:**

1. There are amendments in the warranty policy, sales receipt or invoice or the absence of any of these original documents.
2. There is damage or deterioration of the device due to improper use or if the device has been modified or repaired by an unauthorized third party.
3. There is damage of the product due to defects incurred during transportation, such as bumps, bangs or improper handling.
4. There is damage due to humidity or any kind of liquids, foreign bodies inside the device or the use of unapproved, unadvised accessories with the device contradictory to the user manual.
5. The device is not used in accordance with the respective user manual or with its supplied accessories.
6. There is partial or total damage or loss of the device, due to natural disasters (such as earthquakes, floods, fire, lightening and so on), force-majeure or accident.
7. There is damage to the device after being subjected to high dust environments, humidity or application of excessive voltage from the power supply.
8. The product's serial number does not match the one in the policy or it has been amended or erased.
9. There is cosmetic damage, including but not limited to scratches and dents.

**Customer Support**

This chapter provides customer support contact information for the Edge device. For service, refer to the store where you purchased the Edge device or visit our website at **ENVIG.COM**

When contacting a local service center, have the device serial number ready. The serial number is located on the top panel of the device. Customer satisfaction is our top priority.

To help us provide you with the best support possible, please send us your comments and suggestions. Contact information is provided on the back cover of this manual or via our website.

## TECHNICAL SPECIFICATIONS

This chapter presents the specifications for the Edge device

**TABLE 1**

**Technical Specifications**

Weight:	120gm
Physical Dimensions:	H=38mm L=157mm W=52mm

**TABLE 2**

**Operating conditions**

Temperature:	+5°C to +30°C
Humidity:	30RH to 80RH
Rating:	5V, max 2A

**TABLE 3**

**Transport and Storage**

Temperature:	-10°C to 60°C
Humidity:	30RH to 80RH

## APPENDIX

This appendix provides the hazardous materials declaration for the Edge device.

### Restriction of Hazardous Substance Material Declaration

TABLE 4 IS IN ACCORDANCE WITH THE PROVISIONS OF SJ/T 11364

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- x Indicates that the named hazardous substance contained in at least one of the homogeneous materials used for this component may be above the limit requirement of GB/T 26572. The restricted materials are used in accordance with manufacturing needs. Use of such materials has no impact on the device's safety and no harmful effect on device users. **This declaration is included to remind device users that the device must be disposed of reasonably and in accordance with the environmental policy of your locale.**
- Indicates that the named hazardous substance contained in all of the homogeneous materials for this component is below the limit requirement of GB/T 26572 & RoHS 3 (EU 2015/863).

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© POLLOGEN LTD., ALL RIGHTS RESERVED.
THE NAMES POLLOGEN AND ENVIG ARE REGISTERED TRADEMARKS OF POLLOGEN LTD.
**POLLOGEN LTD.**
**OBELIS UK LTD.**
**OBELIS S.A.**
MANUFACTURED BY POLLOGEN LTD.