ENGLISH

Congratulations on purchasing the TriPollar Pose VX home-use device.

You are on your way towards a smoother, firmer-looking body shape!
This manual describes how to use the TriPollar Pose VX device.

Supplementary information to be found at: www.tripollar.com

TABLE OF CONTENTS

INTRODUCING TRIPOLLAR POSE VX	5
• OVERVIEW	5
• TECHNOLOGY	5
• HOW IT WORKS	5
• TRIPOLLAR POSE VX KIT	6
USING TRIPOLLAR POSE VX	7
• GET TO KNOW THE DEVICE	7
• BEFORE YOU BEGIN	9
• START TREATING YOUR SKIN	10
• RECOMMENDED TREATMENT GUIDELINES	13
TROUBLESHOOTING	14
CUSTOMER SUPPORT	14
TECHNICAL SPECIFICATIONS	15
APPENDIX	16

IMPORTANT NOTICE: TRADEMARKS, COPYRIGHTS & PATENTS

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice.

Manufactured by Pollogen Ltd.

Patented technology.

The names Pollogen, TriPollar and Pose VX are registered trademarks of Pollogen Ltd.

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

SAFETY

Important: Be sure to read the entire user manual before using the device.

Always adhere to all the warnings, cautions and safety precautions below when using the TriPollar Pose VX device.

To learn more about TriPollar technology visit www.tripollar.com

ⓘ WARNINGS AND CAUTIONS

TriPollar Pose VX is an electrical device that must be used properly.

To reduce the risk of injury, adhere to all warnings and cautions below:

- Do not place or store the device near water, such as a sink or bathtub.
- Do not submerge the device in water.
- Do not touch the device if it is dropped in water. Unplug the device immediately with dry hands.
- When plugged in, do not leave the device unattended.
- Before use ensure that the system and its power supply do not show any visible damage. Do not operate if power cord is damaged or if system is disassembled.
- Use the device only with electric cord and power supply provided.
- Keep device out of the reach of children.
- Ensure that your hands are completely dry when operating the device.
- Turn device OFF and then unplug power cord from electrical outlet after use or before cleaning.
- Use the preparation gel during treatment according to its instructions.
- Do not use TriPollar® Preparation Gel after its expiration date.

Important: Do not use TriPollar Pose VX on any area of your face, neck, thyroid area, breasts, groin, genitals, or armpit as this device is not designed to treat these areas.

SAFETY INSTRUCTIONS

Use this device only for its intended purpose, as described in this manual. Failure to follow instructions may cause damage or burns to your skin.

- The use of this product may cause discomfort or adverse effects on people with allergic constitution, people with skin diseases and allergic dermatitis, people with sensitive skin, etc. If adverse effects occur during use, you should stop using this appliance immediately.

DO NOT USE THE DEVICE IF ANY OF THE FOLLOWING APPLIES TO YOU:

- Under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Have metal implants in treated area.
- Have cancer or a history of cancer, especially skin cancer or premalignant moles.
- Have severe concurrent conditions, such as cardiac disorders.
- Pregnancy or nursing or on abdomen during the menstruation.
- Have an autoimmune disease or impaired immune system due to immunosuppressive diseases such as AIDS, or are using immunosuppressive medications.
- Have a disease that is stimulated by heat, such as recurrent herpes simplex within treated area.
- Have any active skin condition within the treated area, such as sores, psoriasis, eczema or a rash.
- Have a history of keloid scarring, abnormal wound healing or fragile skin.
- Any surgical, invasive, ablative procedure in the treatment area before complete healing.
- Have suspicious lesions within the treated area.
- Currently use medications, herbal preparations, vitamins or food supplements that may cause fragile or dry skin.

POSSIBLE SIDE EFFECTS

The following side effects are rare, and usually temporary. Consult your physician if any of the following conditions appear and last more than 24 hours after treatment:

- Pain
- Damage to natural skin texture (crusting, blistering, burns, fragile skin, bruising)
- Mild headache

If you are considering other aesthetic procedures in the treated area, or have any uncertainties regarding potential side effects, consult your physician before using the device.

INTRODUCING TRIPOLLAR POSE VX

OVERVIEW

The revolutionary TriPollar Pose VX skin rejuvenation device, is clinically-proven. The device is intended for body skin tightening, improvement of cellulite appearance and circumference reduction.

With this safe, non-invasive and easy-to-use device, women and men alike can enjoy professional aesthetic treatments at home, with long-lasting results.

TECHNOLOGY

TriPollar Pose VX combines two patented technologies: TriPollar Radio Frequency (RF) and Dynamic Muscle Activation (DMA), which work synergistically and non-invasively-to tighten the skin and improve its appearance.

HOW DOES IT WORK?

During treatment TriPollar Pose VX - delivers focused RF, or combined RF & DMA energies into the skin, utilizing 4 electrodes. The low power high frequency RF energy transmitted to the dermis causes the collagen fibers to contract, which creates an immediate skin tightening effect noticeable from the first treatment. In the long term, this increases the production of new collagen - the protein that gives the skin its firmness.

Additionally, the heat increase the permeability of fat cells membrane and thus allow liquid fat to secrete from the cells. This process results: improves the appearance of cellulite and body circumference.

Simultaneously, the action of the RF and DMA technologies creates a bi-directional effect on the skin – the DMA contributes to lymphatic drainage, an overall tightening effect and muscle toning.

The result is a healthier, and more contoured look.

TRIPOLLAR POSE VX KIT

The kit contains the following components:

1. TRIPOLLAR POSE VX DEVICE (FIG.1)
2. TRIPOLLAR POSE VX POWER SUPPLY AND POWER CORD (FIG. 2)
3. PREPARATION GEL (FIG.3) CAN BE USED UP TO 12 MONTHS AFTER FIRST USE.
Important: Perform a patch test 24 hours before first use (See BEFORE YOU BEGIN page 9).
4. TRIPOLLAR POSE VX USER MANUAL AND WARRANTY POLICY (FIG. 4)

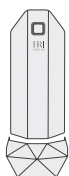
FIG.1

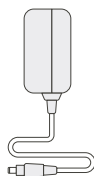
FIG.2

FIG.3

FIG.4

USING TRIPOLLAR POSE VX

GET TO KNOW YOUR DEVICE

The TriPollar Pose VX device focuses RF (Radio Frequency) energy deep into the dermis and hypodermis, and utilizes DMA technology to target the muscle layer. (FIG. 5).








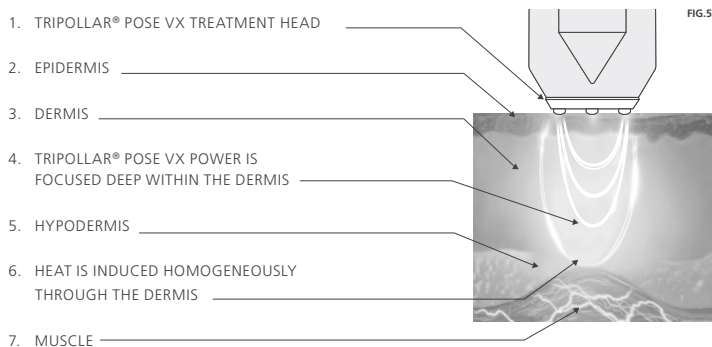
FIG. 6 SHOWS THE COMPONENTS OF THE TRIPOLLAR POSE VX DEVICE:

1. DEVICE BODY
2. TREATMENT HEAD
3. ELECTRODES
4. TEMPERATURE SENSOR
5. POWER SUPPLY CONNECTION
6. ORANGE INDICATOR LIGHT
7. ON/OFF BUTTON
8. LOW RF
9. MEDIUM RF
10. HIGH RF
11. LOW RF AND LOW DMA
12. MEDIUM RF AND MEDIUM DMA
13. HIGH RF AND HIGH DMA
14. PROTECTIVE COVER

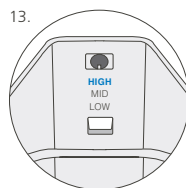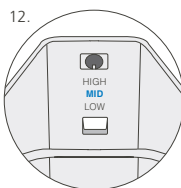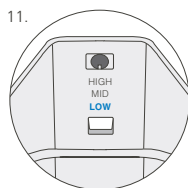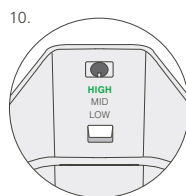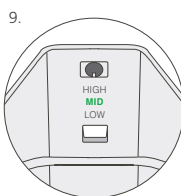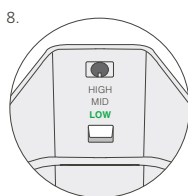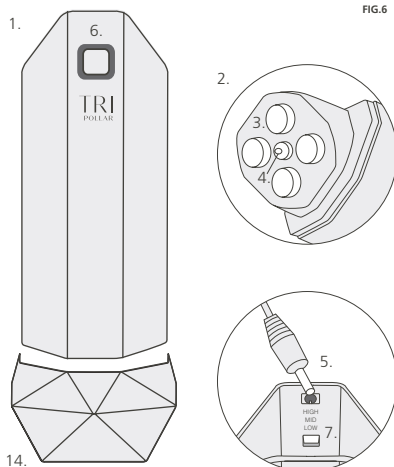BEFORE YOU BEGIN

Before using the TriPollar Pose VX device, perform the following:

1. Open the box containing the device, remove the device and power supply.
2. Connect the power supply to the device (FIG.7).
3. Plug the power supply into an electrical outlet (FIG.8).
4. Remove the preparation gel from its box.
5. Remove the protective cover from the device.

Important: Before using the device for the first time perform a patch test 24 hours prior to use: Apply a small amount of preparation gel to the inner arm, and treat the area with TriPollar Pose VX device. Ensure that the skin does not have a negative reaction to the gel or to the treatment.

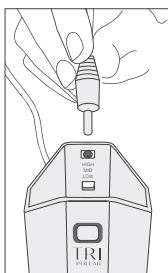
FIG.7

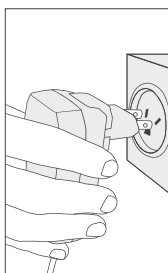
FIG.8

START TREATING YOUR SKIN

TriPollar Pose VX is intended for treating the following areas: thighs, buttocks, abdomen, waist, hips, and arms (front and back).

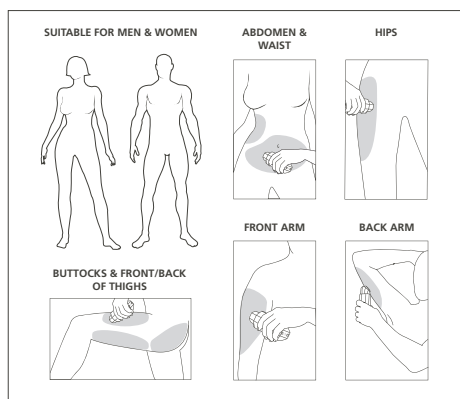
FIG.9

TREATMENT PROCEDURE

Treat your skin by following the steps below:

1. Choose a treatment area: Thigh, buttock, abdomen, waist, hip, or arm. Do not use this device on your face, neck, ears, breasts, genital, chest or armpit areas.
2. Clean and dry the treated area.
3. Apply a thin, even layer of preparation gel over the entire treated area. You may feel a warm sensation.
4. Turn on the device by pressing once on the POWER button (FIG. 10). The Green indicator light turns on to the LOW RF power level.


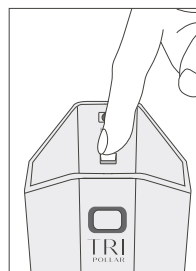
FIG.10

Tips:

- Every additional press on the power button switches the power level from LOW RF, to MID RF, to HIGH RF, then to LOW RF+DMA, MID RF+DMA and HIGH RF+DMA (FIG.10).
- A long press for 2-3 seconds will turn the device OFF.
- Green indicator light represent RF mode, blue indicator light represent RF + DMA mode.

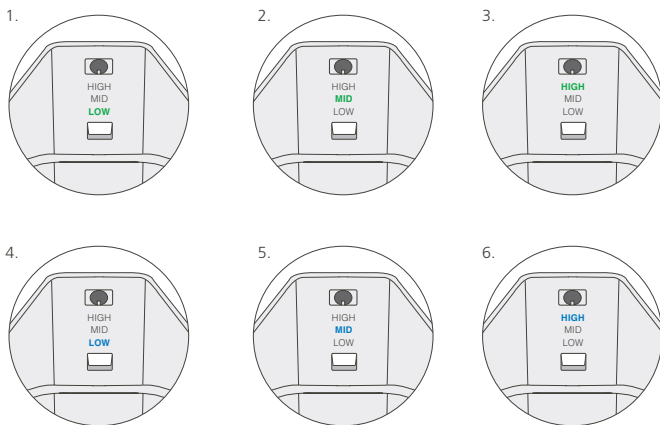- | | |
|------------|-----------------------|
| 1. Low RF | 4. Low RF + Low DMA |
| 2. Mid RF | 5. Mid RF + Mid DMA |
| 3. High RF | 6. High RF + High DMA |

5. Select the power level most comfortable to you. You can choose RF only or RF+DMA.

Important: Start your first treatment on the LOW level, and adjust the power level according to your comfort. It is recommended to use higher-level settings in the areas of the thighs and buttocks, and lower level settings for the arms.

6. Divide the abdomen / buttocks / thighs into treatment areas according to table on FIG.12. Treat the front, back arm or waist as single treatment area.

7. Place the device on treated area and start moving it slowly and constantly on your skin in small circular or figure eight motions (FIG.11).

Slow movements will heat the area internally, whereas fast motions will not achieve the effect to the same extent. As a result, the treated area will start to heat up. If you choose the RF+DMA mode, you will feel muscle contraction accompanied by internal heating of the treated area.

If you feel excessive, intolerable heat, reduce the power level. The orange indicator light will turn on followed by a 'beep' sound when the skin reaches to optimal temperature.

Important: During treatment, all 4 electrodes must remain in contact with the skin at all times. Do not stop the device movement at any time while it is touching your skin, as it may cause the skin to overheat, resulting in unwanted side effects.

For optimal treatment, it is recommended to use the LOW power setting first and to follow the guidelines below:

- If the orange indicator light does not turn on within 5 minutes, or if you do not feel any heating, try moving the device slowly possible without stopping the motion. Slow circular, constant movements are more effective than rapid movements. The slower the movement, the faster the treated area heats internally.

- If the orange indicator light still does not turn on after 3-5 minutes, select a higher power level. If the indicators still do not illuminate after increasing a higher power setting, try moving the device slower, without stopping the motion.
- If the orange indicator light turns ON immediately or if you feel excessive heat, select a lower power level.

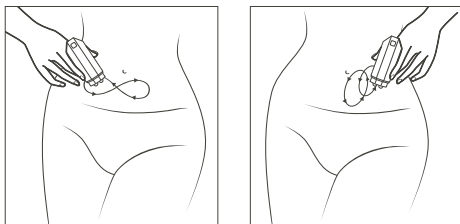
FIG.11

8. Once the orange indicator light illuminates followed by a 'beep' sound, you have reached the optimal temperature. Continue the same slow circular motion for a few moments and then gradually enlarge your movements in the treated area. Proceed with this technique covering the entire treated area, making sure the orange indicator light is illuminated as much as possible during the treatment. As you move from one treated area to another, the orange indicator light may turn-off if the temperature of the skin is not warm enough.

The device's automatic temperature sensor protects the skin from overheating. After the skin has reached the maximal treatment temperature, the GREEN or the BLUE indicator light turns off, indicating that the energy source has been automatically shut down. The device will turn on automatically when the skin temperature drops below the maximal temperature, or when you move the device to an adjacent area where the skin temperature is lower.

9. Continue working on the entire treated area in a slow circular motion for **5 minutes**. Make sure to maintain constant heating while covering the entire treated area.
10. Once you have finished the treatment time according to the body area, you can proceed to the next treatment area repeating steps 5-9.
Important: After treatment, some redness (erythema) may be visible, and slight skin swelling (edema) may occur in the treated area. Both should disappear after a short time.
11. Turn OFF the device by holding down the Power button for two seconds.
12. Unplug the power supply from electrical outlet.
13. Remove the preparation gel gently from your skin and rinse the area. You can apply moisturizer.
14. Clean the 4 electrodes with a damp cloth or paper towel and cover the device's tip with its protective cover.
Warning: Do not submerge the device in water.

RECOMMENDED TREATMENT GUIDELINES

Typically, each treated area takes 5 minutes (recommended treatment area is 12x12 cm). However, treatment times can vary from one person to another, according to body size, and individual treatment pace, etc. In addition, you may choose to treat specific areas of interest using the same working technique.

For optimal results use the TriPollar Pose VX device 2-3 times a week for 6-8 weeks. Allow at least 24 hours between treatments in the same treated area. For long-term results, after completing a course of treatment, it is recommended to repeat it 2-4 times a month, or as needed.

For best results divide the treated area according to the following table:

TriPollar Pose VX Treated area	Recommended treatment
Abdomen	4-6 areas
Thigh	2-3 areas
Arm (single back or front)	1 area
Buttock	2 areas
Hip	1-2 areas
Waist	1 area

FIG.12

TROUBLESHOOTING

If your TriPollar Pose VX device is not working, take the following steps before contacting the service center or store where you purchased the device:

- Ensure that you have removed the device's protective cover.
- Ensure the power cord is securely connected to the device, and the power supply adapter is connected to the electrical outlet.
- Press the Power button and verify that the green indicator light illuminates.
- Ensure that you applied the preparation gel according to the instructions in this manual. The device will not work properly without proper application of the preparation gel.

CUSTOMER SUPPORT

Customer satisfaction is our top priority. For service, refer to the store where you purchased the TriPollar Pose VX device or visit www.tripollar.com. When contacting your local service center or store where you purchased the device, have the serial number ready. The serial number is located on the top panel of the device. To help us provide you with the best support possible, please send us your comments and suggestions. Contact information is provided on the back cover of this manual, or you can contact us via our website (see above).

TECHNICAL SPECIFICATIONS

This chapter presents the specifications for TriPollar Pose VX device.

TABLE 1

Technical Specifications

Weight:	125gr \pm 10% (without stand)
Physical Dimensions:	H=146mm L=50mm W=40mm
TriPollar RF Frequency:	1MHz

TABLE 2

Operating conditions

Temperature:	+5°C to +30°C
Humidity:	20RH to 80RH
[AB] input:	AC 100Vac to 240Vac 0.6 A 50/60Hz
[AB] output:	DC 12V 2.0A

TABLE 3

Transport and Storage

Temperature:	0°C to 40°C
Humidity:	10RH to 90RH

TABLE 4

Electrical Safety

Power supply is classified as
Class II double insulation

APPENDIX: HAZARDOUS MATERIALS DECLARATION

This appendix provides the hazardous materials declaration for TriPollar Pose VX device.

Restriction of Hazardous Substances (RoHS) Declaration

TABLE 5 IS IN ACCORDANCE WITH THE PROVISIONS OF SJ/T 11364

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LEGEND

- X Indicates that the hazardous substance in that column is contained in at least one of the homogeneous materials used for this component, and may be above the limit requirement of GB/T 26572. The restricted materials are used in accordance with manufacturing needs. Use of such materials has no impact on the device's safety and no harmful effect on device users. **This declaration is included to remind users that the device must be disposed of reasonably and in accordance with the environmental policy of your country.**
- Indicates that the named hazardous substance contained in all of the homogeneous materials for this component is below the limit requirement of GB/T 26572.

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