# STOP V<sup>X</sup> GOLD<sup>2</sup>

Facial Treatment Gel:

 **IK CTEC**

**ENGURES 5, JURMALA,**

**LV-2016, LATVIA**

Tel. +371(0)2 9993414

E-mail. office@ctec.lv

Facial Treatment Device:

 **OBELIS S.A.**

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**B-1030 BRUSSELS, BELGIUM**

Tel. +32 (0) 2 732 5954

Fax. +32 (0) 2 732 6003

E-mail. mail@obelis.net

 **POLLOGEN LTD.**

**6 KAUFMAN ST.**

**GIBOR HOUSE P.O.BOX 50320**

**TEL AVIV, ISRAEL 6801298**

TEL: +972(0)3 5104110

FAX: +972(0)3 5174112

EMAIL: INFO@POLLOGEN.COM

WEBSITE: WWW.TRIPOLLAR.COM

© POLLOGEN LTD., ALL RIGHTS RESERVED.

DESIGN AND TECHNICAL SPECIFICATION ARE SUBJECT TO CHANGE WITHOUT PRIOR NOTICE.

MANUFACTURED BY POLLOGEN LTD. PATENTED ISRAELI TECHNOLOGY, MADE IN CHINA.

THE NAME STOP IS TRADEMARK OF POLLOGEN LTD.

THE NAMES POLLOGEN AND TRIPOLLAR ARE REGISTERED TRADEMARKS OF POLLOGEN LTD.

FOR MORE INFORMATION AND FOR USER-MANUAL FILE, PLEASE VISIT: WWW.TRIPOLLAR.COM

PN 14700501CH REV 1

INTELLIGENT  
BIOLOGY™

# STOP V<sup>X</sup> GOLD<sup>2</sup>

FACIAL RESHAPING AND REJUVENATION

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## ENGLISH

**Congratulations on purchasing your TriPollar® Stop VX Gold<sup>2</sup> device!**

You have taken the first step toward flawless beautiful skin!  
This user manual describes the TriPollar® Stop VX Gold<sup>2</sup> product and how to use it.

Supplementary information to be found at:  
[www.tripollar.com](http://www.tripollar.com)

## TABLE OF CONTENTS

|  |           |
|--|-----------|
| <b>INTRODUCING TRIPOLLAR® STOP VX GOLD<sup>2</sup></b> | <b>8</b>  |
| • OVERVIEW   |           |
| • THE TRIPOLLAR® STOP VX GOLD <sup>2</sup> KIT         |           |
| <b>USING TRIPOLLAR® STOP VX GOLD<sup>2</sup></b>       | <b>11</b> |
| • THE TRIPOLLAR® STOP VX GOLD <sup>2</sup> DEVICE      |           |
| • BEFORE YOU BEGIN                                     |           |
| • TREATING YOUR SKIN                                   |           |
| • RECOMMENDED TREATMENT GUIDELINES                     |           |
| <b>TROUBLESHOOTING</b>                                 | <b>20</b> |
| <b>CUSTOMER SUPPORT</b>                                | <b>20</b> |
| <b>WARRANTY TERMS</b>                                  | <b>21</b> |
| <b>TECHNICAL SPECIFICATIONS</b>                        | <b>23</b> |
| <b>APPENDIX</b>  | <b>24</b> |

## IMPORTANT NOTICE

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice.

Manufactured by Pollogen Ltd.

Patented technology.

The names Pollogen, TriPollar and Stop are registered trademarks of Pollogen Ltd.

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

## SAFETY

**Important:** Be sure to read the entire user manual prior to use of the product. Always adhere to all warnings and cautions. Be sure to follow all safety precautions below when using the TriPollar® Stop VX Gold<sup>2</sup> product. To learn more about TriPollar technology, visit our English website at: [www.tripollar.com](http://www.tripollar.com).

### ① WARNINGS AND CAUTIONS

TriPollar® Stop VX Gold<sup>2</sup> is an electrical device that must be used properly. To reduce the risk of injury from electric shock, adhere to all warnings and cautions below:

- Do not touch the device should it fall into water. Unplug the device immediately with dry hands.
- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not submerge the device under water.
- When plugged in, do not leave the device unattended.
- Before use, ensure that the system and its power supply do not show any visible damage. Do not operate if the electric cord is damaged or if the system is disassembled.
- Only use the device with the electric cord and power supply provided.
- Keep the device out of the reach of children.
- Be sure that your hands are completely dry when operating the device.
- Turn the device OFF and then unplug the power supply from the electrical outlet after use or before cleaning and maintenance.

- Use TriPollar® Preparation Gel during treatment according to its instructions for use.

### **SAFETY INSTRUCTIONS**

- Only use this device for its intended purpose, as described in this manual.
- Do not use the device over eyelids, closed eyes, thyroid area, lips, ears, breasts or genital areas.
- The sensation while using the device should be warm and pleasant. If you experience discomfort, discontinue use immediately.
- The use of this product may cause discomfort or adverse effects on people with allergic constitution, people with skin diseases and allergic dermatitis, people with sensitive skin, etc. If adverse effects occur during use, you should stop using this appliance immediately.

Failure to follow the instructions in this manual may cause damage or burns to your skin.

### **DO NOT USE THE DEVICE IF ANY OF THE FOLLOWING APPLY TO YOU:**

- Are under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Have metal implants in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which present no problem.
- Currently have cancer or a history of cancer, especially skin cancer or premalignant moles.
- Have severe concurrent conditions, such as cardiac disorders.

- Are pregnant or nursing.
- Have an impaired immune system due to immunosuppressive diseases such as AIDS, or are using immunosuppressive medications.
- Have any disease that is stimulated by heat, such as, but not limited to, recurrent Herpes Simplex (within the treatment area).
- Have any active skin condition within the treatment area, such as, but not limited to, sores, psoriasis, eczema or a rash.
- Have a history of keloid scarring, abnormal wound healing or fragile skin.
- Have had any surgical procedure, laser resurfacing or deep chemical peeling within the treatment area within the past three months or are still undergoing a process to complete healing.
- Have permanent fillers within the treatment area.  
If you have temporary dermal fillers within the treatment area, consult your physician before using the device.
- Have had Botox administered in the treatment area within the past three months.
- Have suspicious lesions within the treatment area.
- Currently use medications, herbal preparations, vitamins or food supplements that may cause fragile or dry skin.

## **POSSIBLE SIDE EFFECTS**

The following side effects are rare. However, should they occur, they are usually temporary. If any of the following conditions appear and last more than 24 hours, consult your physician:

- Pain
- Damage to natural skin texture (crusting, blistering, burns, fragile skin, bruising)
- Mild headache

For other aesthetic procedures or should you have any uncertainty regarding potential side effects, consult your physician before using the device.

## INTRODUCING TRIPOLLAR® STOP VX GOLD<sup>2</sup>

### OVERVIEW

TriPollar® Stop VX Gold<sup>2</sup> presents **the latest in RF technology at the comfort of your home**. The revolutionary TriPollar® Stop VX Gold<sup>2</sup> combines two effective technologies – the groundbreaking TriPollar® Multi-RF technology (RF) and Dynamic Muscle Activation (DMA) that work simultaneously to deliver non-invasive collagen remodeling and muscle activity, as well as skin tightening and treatment of wrinkles.

TriPollar® Stop VX Gold<sup>2</sup> affects homogeneously all layers of the dermis and the superficial muscle layer, for faster, safer and more effective treatment. It delivers professional-level aesthetic results that are both immediate and long-lasting.

TriPollar® Stop VX Gold<sup>2</sup> reaches the effective temperature quickly with its proprietary Turbo function. Thanks to this unique function shorter treatment would still have the best results.



TriPollar® Stop VX Gold<sup>2</sup> users enjoy the combined benefits of a unique non-conductive gel that ensures the delivery of focused RF energy into the right layers of the dermis using the newly designed electrodes. The temperature sensor ensures an effective and safe treatment while stimulating dermal activity that accelerates the natural production of collagen, the protein responsible for the skin's smooth and supple texture.

Simultaneously, DMA targets the SMAS layer, toning and lifting it by stimulating gentle muscle contractions inducing a "face lift" effect. The TriPollar® Stop VX Gold<sup>2</sup> system is intended for non-invasive skin renewal at home, specifically indicated for skin tightening and treatment of wrinkles.

Using TriPollar® Stop VX Gold<sup>2</sup> regularly according to its recommended treatment can effectively postpone the signs of aging.

2. TRIPOLLAR® STOP VX GOLD<sup>2</sup> DEVICE (FIG.2)
3. TRIPOLLAR® STOP VX GOLD<sup>2</sup> POWER SUPPLY AND ELECTRIC CORD (FIG.3)
4. TRIPOLLAR® PREPARATION GEL (FIG.4)  
After the container has been opened, the gel can be used for up to 12 months.  
**Important:** It is recommended to perform a patch test 24 hours prior to initial use.  
See page 14 for details.
5. TRIPOLLAR® STOP VX GOLD<sup>2</sup> USER MANUAL AND INTERNATIONAL WARRANTY (FIG.5)
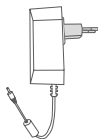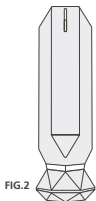
## USING TRIPOLLAR® STOP VX GOLD<sup>2</sup>

### THE TRIPOLLAR® STOP VX GOLD<sup>2</sup> DEVICE

TriPollar® Stop VX Gold<sup>2</sup> is a facial skin renewal and tightening device. The device is used to increase and accelerate collagen production in the skin for a rejuvenating effect.

TriPollar® Stop VX Gold<sup>2</sup> focuses heat deep within the dermis. The TriPollar® Stop VX Gold<sup>2</sup>'s DMA technology targets the SMAS layer, as shown below.


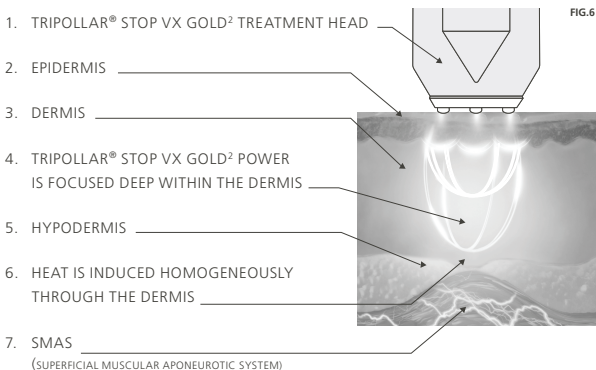
FIG.7 Shows the components comprising the TriPollar® Stop VX Gold<sup>2</sup> device.

1. DEVICE BODY
2. TREATMENT HEAD
3. ELECTRODE
4. TEMPERATURE SENSOR
5. CONTROL PANEL
6. POWER SUPPLY CONNECTION
7. GREEN/BLUE INDICATION
8. ORANGE INDICATION LIGHT
9. ON/OFF BUTTON
10. LOW RF
11. HIGH RF
12. LOW RF AND DMA
13. HIGH RF AND DMA
14. PROTECTIVE COVER

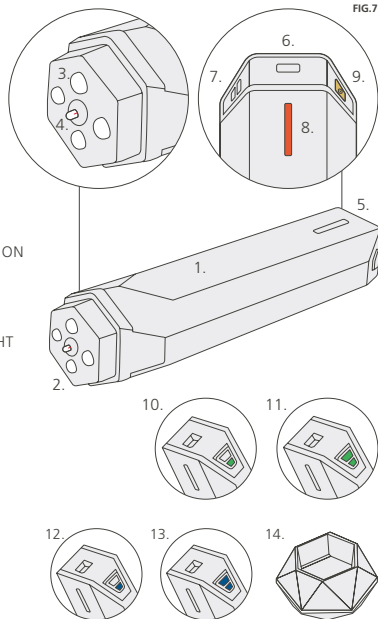The indication lights on the TriPollar® Stop VX Gold<sup>2</sup> device works as follows:

**GREEN LIGHTS INDICATE RF MODE**

- A single green light means that the device is ON and is operating in low-level RF mode.
- Two green lights mean that the device is operating in high-level RF mode.

**BLUE LIGHTS INDICATE RF 2 DMA MODE**

- A single blue light means that the device is operating in low-level RF and DMA mode.
- Two blue lights mean that the device is operating in high-level RF and DMA mode.

## BEFORE YOU BEGIN

Before using the TriPollar® Stop VX Gold<sup>2</sup> device, do the following.

1. Open the box containing the device, and remove the device and its power supply.
2. Connect the power supply to the device (FIG.8)
3. Plug the power supply into an electrical outlet (FIG.9)
4. Remove the TriPollar® Preparation Gel from its box.

**Important:** Before using the device for the first time, it is recommended to perform a patch test 24 hours prior to initial use. To do so, apply a small amount of TriPollar® Preparation Gel to the inner arm and treat the area with the TriPollar® Stop VX Gold<sup>2</sup> device. Ensure that the skin does not have a negative reaction to the gel or treatment.

5. Remove the protective cover from the device.

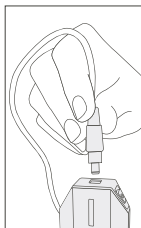
FIG.8

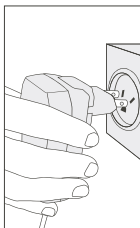
FIG.9

## TREATING YOUR SKIN

To treat your skin using the TriPollar® Stop VX Gold<sup>2</sup> device:

1. Cleanse and dry the treatment area.
2. Apply a very thin, even layer of TriPollar® Preparation Gel over the entire treatment area (face and/or neck). This gel must be used during the entire treatment session. Your skin may feel a warm sensation when the gel is applied.
3. Turn the device ON by pressing the Power button (FIG.10). Make sure that the green indicator light turns ON.

**Note:** A long press of the Power button turns the device OFF.

4. Select the RF power setting (Low/High) that is most comfortable for you. When the green indicator light shows a single green light illuminated (FIG.11), the device is operating at the Low power level. You can press the Power button again to change to the High power level (FIG.12).

**Important:** During your first treatment session, it's recommended to use the Low power setting and then adjust the power level according to your skin's tolerance for subsequent treatments.

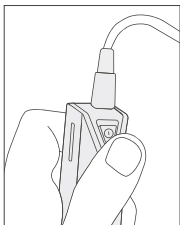
FIG.10

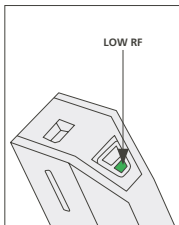
FIG.11

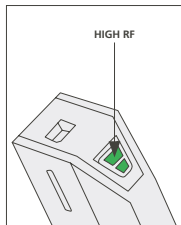
FIG.12

Treatment:

5. Divide the face to 3 treatment areas (FIG.13):
  - Right cheek
  - Left cheek
  - Forehead
6. Place the device on your skin, make sure to apply a slight pressure in order to assure a **full** contact between the electrodes and the skin.
7. Once the device is placed on the skin, the Turbo RF will start. Turbo RF will be initiated once at the beginning of the treatment or if a new level or mode is chosen.
8. Start moving the device in a very slow circular motion on a small part of your treatment area (approximately twice the size of the treatment head), slowly heating your skin to the optimal temperature (FIG.14).

**Important:** Be sure to constantly move the device slowly over the skin

to ensure proper heating of the treated area. Do not keep the device on one place without movement.

9. Once the orange indicator light illuminates followed by a 'beep' sound, you have reached the optimal temperature. Continue the same slow circular motion for a few moments and then gradually enlarge your movements in the treatment area.

Proceed with this technique covering the entire treatment area, making sure the orange indicator light is illuminated as much as possible during the treatment.

As you move from one treatment area to another (FIG.13), the orange indicator light may turn-off if the temperature of the skin is not warm enough.

10. Continue working on the entire treatment area in a slow circular motion **for 3-4 minutes**. Make sure to maintain a uniform heating while covering the entire treatment area. Treatment can be augmented by "lifting" movements from the jaw line upwards toward your cheeks.
11. Once you finished the 3-4 minutes of treatment, you can proceed to the next treatment area repeating steps 6-10.
12. Switch to RF 2 DMA mode. It's recommended to start at the Low power level. **Important:** Use the RF 2 DMA mode **ONLY** on the jaw line area (FIG.15). Do not use on other facial areas.

FIG.13

TREATMENT  
AREAS

FIG.14

FIG.15

13. Start moving the device using a very slow linear motions along the jaw line.
14. To ensure optimal results, use the device for 3-4 minutes on each side. After the entire treatment area has been treated, the treatment session is complete. As a natural result of the treatment, the skin may show some redness (erythema) and minor skin swelling (edema) in the treatment area. Both should disappear after a short time.
15. Press the Power button for two seconds to turn the device OFF.
16. Unplug the power supply from the electrical outlet.
17. Gently remove the TriPollar® Preparation Gel from your skin and wash your face. You can apply a moisturizer if desired.
18. Clean the electrodes of the device with a damp cloth or paper towel and cover the device tip with its protective cover.

**Warning: Do not submerge the TriPollar® Stop VX Gold<sup>2</sup> device under water.**

**Safety Mechanism:**

The device's automatic temperature-optimization mechanism protects the skin from overheating. When the skin has reached the maximal temperature, the energy source will shut down automatically. The indication is - green LED off. The device automatically turns the energy source back ON, when the skin temperature drops below the maximal temperature or when you move the device to an adjacent treatment area, where the skin temperature is lower. The indication is - green LED on.

**Important:**

The orange light is ON when the skin temperature in the treatment area is sufficiently heated and is OFF when it is not.

## RECOMMENDED TREATMENT GUIDELINES

The following treatment guidelines are recommended for using the TriPollar® Stop VX Gold<sup>2</sup> device:

- Treating the entire face takes 12-15 minutes. For shorter treatment times focus on selected areas of concern. Treatment times differ for each person.
- Facial DMA treatment takes approximately 10 minutes.
- For optimal results, use the TriPollar® Stop VX Gold<sup>2</sup> device 2-3 times a week for a period of 6-8 weeks.
- Make sure that you allow at least one day in between treatments on the same treatment area.
- To maintain long-term results after finishing the full course of treatments, we recommend that you use your TriPollar® Stop VX Gold<sup>2</sup> device 2-4 times a month or as needed.
- If the orange indicator light and sound indicators still do not appear within five minutes, or if you do not feel any heating, try moving the device as slowly as possible without stopping the motion. Slow, circular, constant movements are more effective than rapid movements. The slower the movement, the faster the treatment area heats internally.
- If the orange indicator light and sound indicators do not appear after 3–5 minutes, select a higher power level. If the indicators still do not appear after increasing to the High power setting, try moving the device more slowly, without stopping the motion.
- If the orange indicator light and sound indicators turn ON immediately or if you feel excessive heat, select a lower power level.

## TROUBLESHOOTING

**If your TriPollar® Stop VX Gold<sup>2</sup> device is not working, take the following steps before contacting the service center:**

- Ensure that you have **removed** the device's protective cover.
- Make sure that the electric cord is securely connected to the device and the Power supply adapter is firmly connected.
- Press the Power button and verify that the green indicator light illuminates.
- Ensure that you applied TriPollar® Preparation Gel according to instructions described in this manual. The Tripollar® Stop VX Gold<sup>2</sup> device does not work properly when the preparation gel is not applied.

## CUSTOMER SUPPORT

For service, refer to the store where you purchased the TriPollar® Stop VX Gold<sup>2</sup> device or visit our website at: **[www.tripollar.com](http://www.tripollar.com)**. When contacting a local service center, have the device serial number ready. The serial number is located on the top panel of the device. Customer satisfaction is our top priority. To help us provide you with the best support possible, please send us your comments and suggestions. Contact information is provided on the back cover of this manual, or you can contact us via our website (see above).

## WARRANTY TERMS

This product has been subjected to strict quality control performed by the manufacturer. For this reason, the distributor in your country guarantees proper function of your product for 24 months from date of actual purchase. If you use the device properly, follow all instructions outlined in the user manual, you can enjoy the product to your ultimate satisfaction.

This policy warrants the support that the distributor provides for the device, should it happen that failures attributable to manufacturing defects be found. In such case, you shall be entitled to make full use of the benefits of this warranty. This implies the reason and replacement of parts, free of charge (regarding both service and parts) at any Authorized Service Center. Within its validity term, you are entitled to make use of the warranty if necessary to have a fully enjoyable product, provided the above-indicated conditions are fully met. This warranty is valid for the device marketed by the distributor in your country.

The warranty is not valid when:

1. Amendments in the warranty policy, sales receipt or invoice, or the absence of any of these original documents.
2. Damage or deterioration of the device due to improper use or if the device has been modified or repaired by an unauthorized third party.
3. Damage of the product due to defects incurred during transportation, such as bumps, bangs or improper handling.

4. Damage due to humidity or any kind of liquids, foreign bodies inside the device or the use of unapproved, unadvised accessories with the device contradictory to the user manual.
5. When the device is not used in accordance with the respective user manual or with its supplied accessories.
6. Partial or total damage or loss of the device, due to natural disasters (such as earthquakes, floods, fire, lightening, etc.), force majeure or accident.
7. Damage to the device after being subjected to high dust environments, humidity or application of excessive voltage from the power inlet.
8. If the product's serial number has been amended or erased.
9. To cosmetic damage, including but not limited to scratches and dents.

## TECHNICAL SPECIFICATIONS

### Technical Specifications

|                         |                       |
|-------------------------|-----------------------|
| Weight:                 | 85gm                  |
| Physical Dimensions:    | H=134mm L=51mm W=32mm |
| TriPollar RF Frequency: | 1.0MHz-1.25MHz        |

### Operating conditions

|              |  |
|--------------|--|
| Temperature: | 5°C to 30°C                                    |
| Humidity:    | 20RH to 80RH                                   |
| Rating:      | 100~240V AC, 0.6A (max), 50/60Hz<br>DC 8V 1.5A |

### Transport and Storage

|              |              |
|--------------|--------------|
| Temperature: | 0°C to 40°C  |
| Humidity:    | 10RH to 90RH |

### Electrical Safety

Power supply is classified as Class II double insulation

## APPENDIX

### Restriction of Hazardous Substance Material Declaration

IN ACCORDANCE WITH THE PROVISIONS OF SJ/T 11364

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- x Indicates that the named hazardous substance contained in at least one of the homogeneous materials used for this component may be above the limit requirement of GB/T 26572. The restricted materials are used in accordance with manufacturing needs. Use of such materials has no impact on the device's safety and no harmful effect on device users. **This declaration is included to remind device users that the device must be disposed of reasonably and in accordance with the environmental policy of your locale.**
- Indicates that the named hazardous substance contained in all of the homogeneous materials for this component is below the limit requirement of GB/T 26572.

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