## BS103 Bathroom Scale

PLEASE READ ALL THE INSTRUCTIONS BEFORE USE

[duronik.com](http://duronik.com)

[duronik](https://www.youtube.com/channel/UCduronik)

[@duronikmedia](https://www.instagram.com/duronikmedia)

[Duronik](https://www.facebook.com/Duronik)

[@duronikmedia](https://www.twitter.com/duronikmedia)

## **Contents**

<b>English Manual</b>	<b>4</b>
<b>Manuel d'utilisation en Français</b>	<b>8</b>
<b>Gebrauchsanleitung auf Deutsch</b>	<b>12</b>
<b>Manual de Instrucciones en Español</b>	<b>16</b>
<b>Manuale D'istruzioni In Italiano</b>	<b>20</b>
<b>Instrukcja obsługi. Polski</b>	<b>24</b>
<b>Gebruikershandleiding Nederlands</b>	<b>28</b>


**English Manual**

## Body Scale

### 1. Features & Specifications

- Weight Capacity: 180kg / 397lb / 28st: 5lb
- Minimum weight: 5kg
- Graduation: 0.1kg / 0.2lb / 1lb
- Weight units: kg / lb / st:lb
- Function key: UNIT
- Display: LCD
- Batteries: 2x AAA
- Step on or press UNIT to switch on
- Auto-off function to save battery power

### 2. LCD display

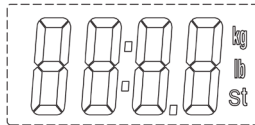
### 3. Changing the Unit

Press the “UNIT” button on the back of the scale to select desired unit in weighing mode.


### 4. Operating Instructions

4-1. Before First Use

- Before first use, pull the plastic strip out of the battery compartment on the back of the scale. If needed, check that the batteries are the correct position and +/- terminals are positioned correctly.
- Wait until 0.0kg is shown on the screen.

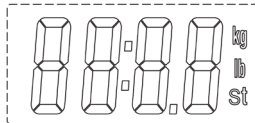
c) Press UNIT to select the desired unit.

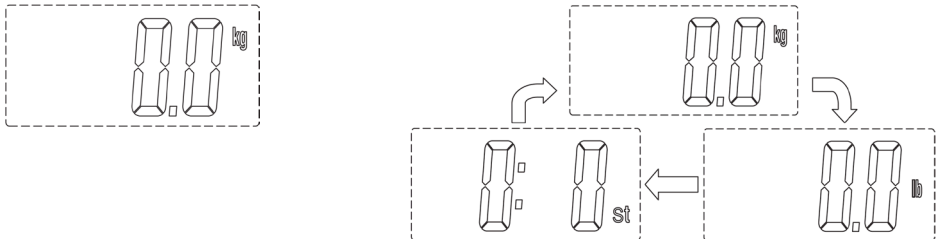
**Note:** After selecting the unit, place the scale on flat hard floor, then you can start to stand on the platform to measure according to the measuring procedure below. If there is nothing placed on the scale whilst on in weighing mode, the scale will turn off automatically after about 10s.

**4-2 Measuring**

- a) Place the scale on flat hard ground.
- b) Stand on the platform to switch on the scale.

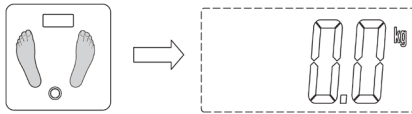- c) Step off of the scale and wait for 0.0 to show on the display.

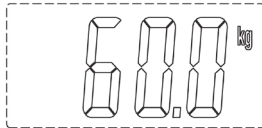
- d) Step back on the scale and stand still. The weight data will blink first, then lock in place when stable. The scale will turn off automatically after about 10 seconds.


**Note:** If you want to measure once again, step off of the scale, then step on the platform again.

**5. Error Messages**

5-1. Low battery indication. When “Lo” is shown, it means the battery is low. The scale will turn off after a few seconds. You will then need to replace the batteries.

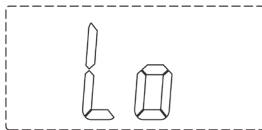5-2 Overload indication. When “Err” is shown, remove the weight from the scale to protect the sensors.

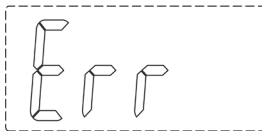**6. Safety Instructions**

- 6-1 Ensure the scale is on a flat, even and firm ground surface. Avoid carpet or soft surfaces.
- 6-2 The surface of the scale will be slippery if wet. Please make sure both the glass and your feet are dry before use.
- 6-3 This scale is not waterproof. Do not expose it to moisture.
- 6-4 Clean the scale with a slightly damp cloth, then wipe again with a dry cloth. Do not use any chemical/abrasive cleaning agents to clean the scale.
- 6-5 If the scale is not being used for a long period of time, remove the batteries before storing away.

# Warranty

1 Year Guarantee from Shine-Mart Ltd, Trademark owners of Duronic Products.

NOTE: THESE STATEMENTS DO NOT AFFECT YOUR STATUTORY RIGHTS AS A CONSUMER

This product has been manufactured under the strictest or quality control procedures, and using the highest quality of materials, to ensure excellent performance and reliability. It will give very good and long lasting service, provided it is properly used and maintained. The product is guaranteed for 1 years from the date of original purchase. If any defect arises due to a faulty materials or poor workmanship, the faulty product must be returned to the original place of purchase. Refund or replacement is at the discretion of that company.

Duronic Products are offered with a 1 year guarantee under the following conditions:

1. The product must be returned to the retailer with original proof of purchase.
2. The product must be installed and used in accordance with the instructions contained in this operating manual.
3. It must be used only for domestic purposes.
4. It does not cover wear and tear, malicious damage, misuse, neglect, inexpert repairs or consumable parts.
5. Shine-Mart Ltd has no responsibilities for accidental or consequential loss or damage.
6. Shine-Mart Ltd are not liable to carry out any type of servicing work, under the guarantee
7. Valid in the UK only

This product is fitted with CE:1363 plug and fuse  
Imported by Shine-Mart Ltd, RM3 8SB

## Information on waste Disposal for Consumers of Electrical & Electronic Equipment.

This mark on a product and/or accompanying documents indicates that when it is to be disposed of, it must be treated as Waste Electrical & Electronic Equipment, (WEEE).

Any WEEE marked products must not be mixed with general household waste, but kept separate for the treatment, recovery and recycling of the materials used.

For proper treatment, recovery and recycling; please take all WEEE marked waste to your Local Authority Civic waste site, where it will be accepted free of charge.

If customers dispose of Waste Electrical & Electronic Equipment correctly, they will be helping to save valuable resources and preventing any potential negative effects upon human health and the environment, of any hazardous materials that the waste may contain.

## Follow Us

[duronic.com](http://duronic.com)

[duronic](https://www.youtube.com/duronic)

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5. Shine-Mart Ltd heeft geen verantwoordelijkheid voor accidentele schade of verlies.
6. Shine-Mart Ltd is niet aansprakelijk voor onderhoud onder de garantie.
Geïmporteerd door Shine-Mart Ltd, RM3 8SB