**DURONIC**  
Creative Consumer Products

## BS201 - Bathroom Scales

PLEASE READ ALL THE INSTRUCTIONS BEFORE USE

[duronic.com](http://duronic.com)

[duronic](https://www.youtube.com/duronic)

[@duronicmedia](https://www.instagram.com/duronicmedia)

[Duronic](https://www.facebook.com/Duronic)

[@duronicmedia](https://twitter.com/duronicmedia)

**Contents**

<b>English Manual</b>	<b>5</b>
<b>Manuel d'utilisation en Français</b>	<b>13</b>
<b>Gebrauchsanleitung auf Deutsch</b>	<b>21</b>
<b>Manual de Instrucciones en Español</b>	<b>29</b>
<b>Manuale D'istruzioni In Italiano</b>	<b>37</b>

**English Manual**

## Specification

1. Capacity: 180kg / 397lb / 28st: 5lb
2. Units: kg / lb / st: lb
3. Unit Conversion: Unit conversion button
4. Division: 0.1kg / 0.2lb/1 lb
5. Power: 3V (2 x battery AAA 1.5V)
6. Step-on
7. Auto switch-off function
8. User profile No: 8 persons (P1-P8)
9. Main feature: Current weight compare with last weight
10. Automatic identification function
11. Automatic identification range: 2kg (weight is below 100kg), 3kg (weight is above 100kg)
12. Full display when power is on:


## Install Or Replace The Battery

This appliance works with 2 x 1.5V AAA batteries. To install batteries, remove the cover of the battery housing on the bottom of the scale, install the battery in under correct positive and negative plate and then slide the cover back on. The scale will then be ready for use.


## Unit Conversion

Press the UNIT button on the back of the scale to select unit of measurement. The default unit is Kg. Key: kg=kilogram, lb=pound, st =stone
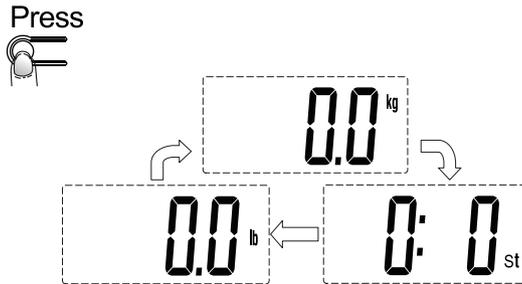
## First Weighing


When using the scale for the first time, place the scale on flat, even and stable ground/surface. Please avoid uneven and carpeted surfaces. Once the user steps on to the scale, the weight will calculate and then be displayed in blue, this will indicate the weight has been recorded.

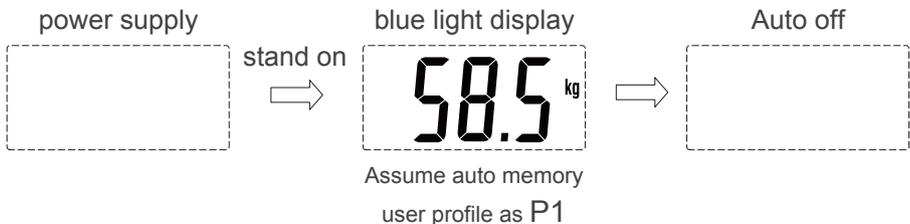**NOTE:** For the first use, an individual's current weight will be displayed and then recorded as a user.

## Weight & Comparison


A) When a user's current weight is less than the last recorded weight, the display will be shown in green to indicate weight loss.

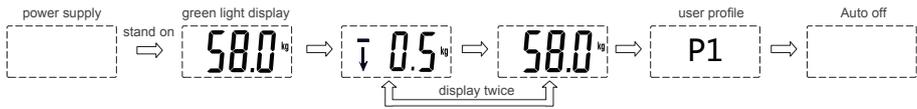B) When a user's current weight is the same as the last recorded weight, the display will be shown in blue.

C) When a user's current weight is greater than their last recorded weight, a red light will display. It will also indicate the exact difference in weight between the last recorded weight and the current weight.
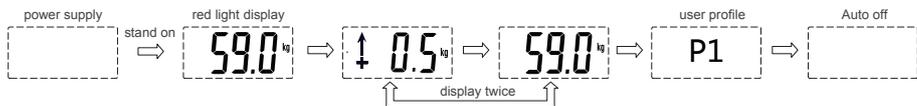
## Automatic Identification Function

The automatic identification range functions as follow: if the weight is below 100 kg, the device will identify any weight difference of up to 2kg. If the weight is above 100kg, the device will identify any weight difference of up to 3kg. Automatic identification allows an individual to be recognized automatically within a 2kg weight range. For example if P1 is set to 50kg and an individual weighing 50-52kg steps on to the scale, the user will automatically be recognized as P1. However, if an individual weighs 52.1kg, a new user profile will be created.

## Delete User Profile

If you want to remove all user weight profiles, press and hold the UNIT key on the bottom, until it displays the symbol 'dEL' on the screen, this means all current user profiles have been deleted.


## Error Indication

A) Overload indication:

If the "Err" symbol appears on the LCD screen, it means that the maximum weight accepted by your scale is exceeded. Please remove the object from the platform to prevent any damage to the scale.

B) Low battery indication

If the "LO" symbol is displayed on the LCD screen, it means that the battery of your scale is low and should be replaced soon. Please dispose of the old batteries in an environmentally friendly way
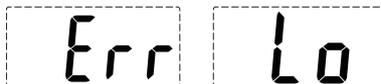
## Reminder & Tips

1. Please ensure before use that the scale is on a flat, even and stable surface.
2. The surface of the scale will be slippery if wet. Please make sure both the glass and your feet are dry before use.
3. This scale is not waterproof, therefore should not be exposed to moisture or submerged under water.
4. Clean the scale with a slightly damp cloth. Do not use any chemical/abrasive cleaning agents to clean the scale.
5. Please remove the batteries if the scale is not to be used for a long period of time.
6. The appliance is intended solely for domestic use within the home.
7. Children should be supervised to ensure that they do not cause any harm to themselves or the appliance.
8. Before you use the scales for the first time, please carefully unpack the appliance and remove all packaging material. Keep the material (plastic bag and cardboard) out of reach of children. Remove the batteries from the protective packaging and insert with the correct polarity.
9. Self-calibration: When using the scales for the first time, or when changing its position or location, the scale should be restarted twice and must be free of any objects and materials to enable accurate measurements. The scales must be placed on a flat, even surface for the calibration to be successful.

# Warranty

1 Year Guarantee from Shine-Mart Ltd, Trademark owners of Duronic Products.

NOTE: THESE STATEMENTS DO NOT AFFECT YOUR STATUTORY RIGHTS AS A CONSUMER

This product has been manufactured under the strictest or quality control procedures, and using the highest quality of materials, to ensure excellent performance and reliability. It will give very good and long lasting service, provided it is properly used and maintained. The product is guaranteed for 1 years from the date of original purchase. If any defect arises due to a faulty materials or poor workmanship, the faulty product must be returned to the original place of purchase. Refund or replacement is at the discretion of that company. Duronic Products are offered with a 1 year guarantee under the following conditions:

1. The product must be returned to the retailer with original proof of purchase.
2. The product must be installed and used in accordance with the instructions contained in this operating manual.
3. It must be used only for domestic purposes.
4. It does not cover wear and tear, malicious damage, misuse, neglect, inexpert repairs or consumable parts.
5. Shine-Mart Ltd has no responsibilities for accidental or consequential loss or damage.
6. Shine-Mart Ltd are not liable to carry out any type of servicing work, under the guarantee
7. Valid in the UK only

This product is fitted with BS:1363 plug and fuse  
Imported by Shine-Mart Ltd, RM3 8SB

#### Information on waste Disposal for Consumers of Electrical & Electronic Equipment.

This mark on a product and/or accompanying documents indicates that when it is to be disposed of, it must be treated as Waste Electrical & Electronic Equipment, (WEEE).

Any WEEE marked products must not be mixed with general household waste, but kept separate for the treatment, recovery and recycling of the materials used.

For proper treatment, recovery and recycling; please take all WEEE marked waste to your Local Authority Civic waste site, where it will be accepted free of charge.

If customers dispose of Waste Electrical & Electronic Equipment correctly, they will be helping to save valuable resources and preventing any potential negative effects upon human health and the environment, of any hazardous materials that the waste may contain.

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1 anno di garanzia da Shine-Mart Ltd, proprietari del marchio commerciale Duronic.
5. Shine-Mart Ltd non è responsabile per Danni accidentali o consequenziali.
6. Shine-Mart Ltd non è tenuta a eseguire qualsiasi tipo di riparazione, ai sensi della garanzia.
Questo prodotto è inclusivo di spina certificata BS importata da Shine-Mart Ltd, RM3 8SB