Bathroom Scales - BS501

PLEASE READ ALL THE INSTRUCTIONS BEFORE USE

duronic.com

[duronic](https://www.youtube.com/channel/UCduronic)

[@duronicmedia](https://www.instagram.com/duronicmedia)

[Duronic](https://www.facebook.com/Duronic)

[@duronicmedia](https://www.twitter.com/duronicmedia)


Contents

English Manual	5
Manuel d'utilisation en français	15
Gebrauchsanleitung auf Deutsch	25
Manual de Instrucciones en Español	35
Manuale D'istruzioni In Italiano	45

English Manual


Specification

1. Maximum Capacity: 183.0kg
2. Overload indication: Err
3. Low battery indication: Lo
4. Zero indication: - - - -
5. High body fat indication: ErrH
6. Low body fat indication: ErrL
7. Number of memory: P0-P9 10 Persons
8. Age range: 10 - 80 years old
9. Height range: 100 - 240cm/3ft 3inch - 7ft 10inch
10. Impedance error: Err2
11. Measuring body fat range: 5% - 50%
12. Measuring water range: 10% - 80%
13. Measuring muscle range: 10% - 80%
14. Measuring bone range: 0.hkg - 8.0kg
15. Automatic identification range: 2.0kg
16. BMI parameter
17. Normal mode and sport mode 

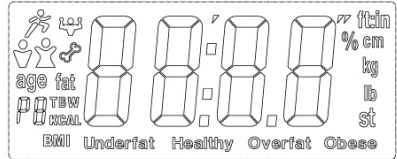
Unit Conversion

1. Default unit is kg.
2. Press UNIT button located on the back of the scale to convert units.

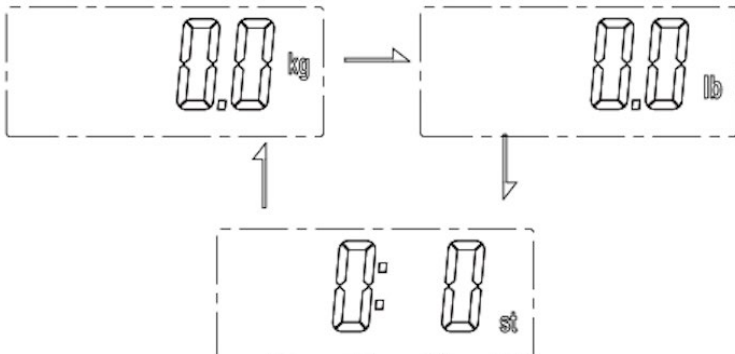
User Information


1. Press SET to enter setting mode.
2. Press ▲ or ▼ to change the user No. Press SET for confirmation.

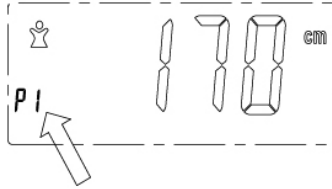
3. Press ▲ or ▼ to choose male or female, then choose normal mode or sport mode (for users who regularly exercise). Press SET to confirm.

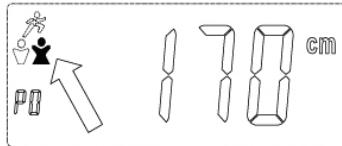4. Press ▲ or ▼ to change the height. Press SET to confirm.

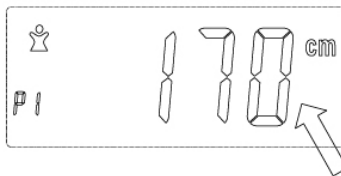
4. Press ▲ or ▼ to change the age. Press SET to confirm.

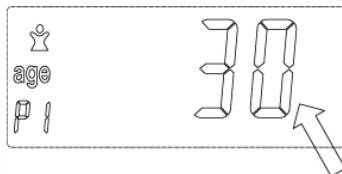
Weighing and Measuring Fat


1. After completing the settings, the scale will indicate “0.0”

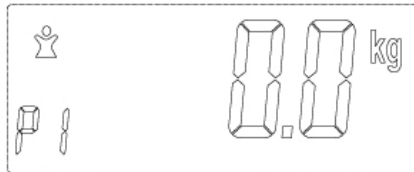
2. Stand with bare feet on the scale and wait for a stable weight to display.


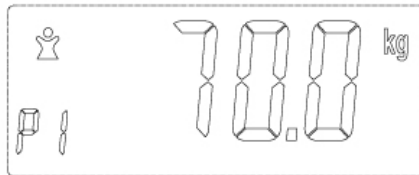
3. Four “O” will appear and disappear one by one from left to right.

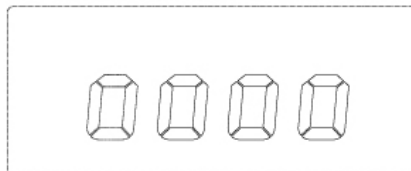4. The measurements will be on display one by one. Body fat, water, BMI, muscle mass, bone mass, calorie and weight.

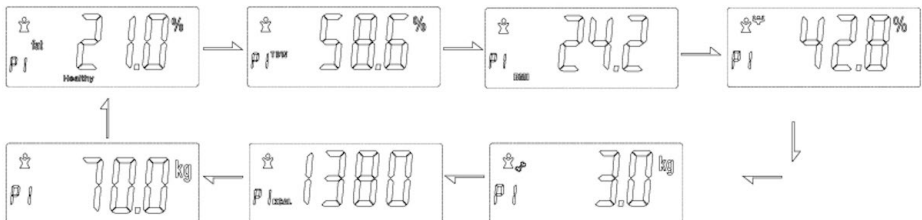5. The scale will display the user information which will let the user confirm the information again. After the display, the scale will turn off automatically.

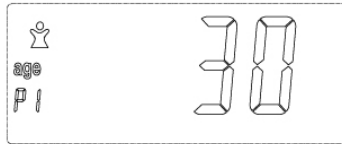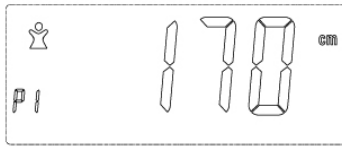
Important (Err2):


To analyse body mass, stand on the scale barefoot. If your feet are cold, if you have poor blood pressure, or if you step onto the scale too quickly, the measurement may be inaccurate. The Err2 message will display to indicate that there has been a problem. Make sure your feet are warm and aligned to the electrodes. Stand still with weight evenly balanced on each foot, and measure again.

Automatic Identification Range


1. If the user has correctly measured the fat weight, the scale will automatically memorize the user's weight and the user can directly stand on the scale to measure the fat weight next time without entering any information again.
2. Stand on the scale and wait for a stable weight to display.

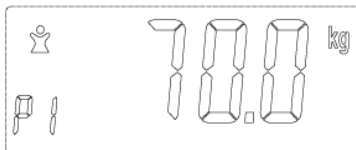
3. Four "O" will appear and disappear from left to right when measuring.

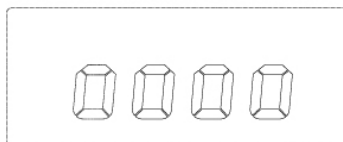4. The LCD display will indicate the body fat, water, muscle, bone, calorie and weight in a sequence.

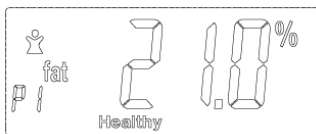5. If the scale identifies the weight is similar to 2 other users, the user should use UP or DOWN to choose the specific user.

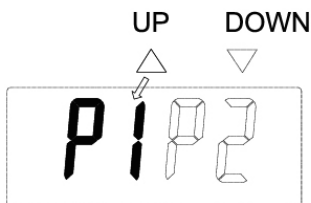Grade Standards

Body fat analysis (Unit %)								
Age	Female				Male			
	Underweight	Healthy	Overweight	Obese	Underweight	Healthy	Overweight	Obese
10-39	<21	21-33	33-39	>39	<8	8-20	20-25	>25
40-59	<23	23-34	34-40	>40	<11	11-22	22-28	>28
60-80	<24	24-36	36-42	>42	<13	13-25	25-30	>30

Body water analysis (Unit %)						
Age	Female			Male		
	Underweight	Healthy	Overweight	Underweight	Healthy	Overweight
10-15	<57	57-67	<67	<58	58-72	<72
16-30	<47	47-57	<57	<53	53-67	<67
31-60	<42	42-52	<52	<47	47-61	<61
60-80	<37	37-47	<42	<47	42-56	<56

BMI parameter (WHO)				
	Underweight	Normal	Overweight	Obese
BMI	<18.5	18.5-25	25-30	>30

Body muscle analysis (UNIT: %)		
Gender	Female	Male
Normal	>34%	>40%

Body bone analysis (UNIT: %)						
Gender	Female			Male		
Weight	<45kg	45-60kg	>60kg	>60kg	60-75kg	>75kg
Bone	1.8kg	2.2kg	2.5kg	2.5kg	2.9kg	3.2kg

Calorie analysis (Unit:KJ)						
	Age	<17	18-29	30-50	>50	
Male	<50kg	1350	1200	1115	1075	
	50-65kg	1485	1320	1227	1183	
	65-80kg	1890	1680	1561	1505	
	80-95kg	2295	2040	1896	1828	
	>95kg	2700	2400	2230	2150	
Female	<50kg	1265	1180	1085	1035	
	50-65kg	1392	1298	1194	1139	
	65-80kg	1771	1652	1519	1449	
	80-95kg	2151	2006	1845	1760	
	>95kg	2530	2360	2170	2070	

Reminder:

1. Ensure the scale is on a flat, even surface. Avoid carpet or soft surfaces.
2. The surface of the scale will be slippery when wet. Please make sure both the glass and your feet are dry before use.
3. This scale is not waterproof. Do not expose to moisture.
4. Clean the scale with a slightly damp cloth. Do not use any chemical cleaning agents to clean the scale.
5. Please remove the batteries, if the scale is not being used for a long period.
6. When using the scales for the first time, or when changing their position or location, they should be restarted twice and should be free of interferences to ensure accurate measurements. This is for calibrating your scale. If you do not perform these measures, the measurements may be incorrect.

Warranty

1 Year Guarantee from Shine-Mart Ltd, Trademark owners of Duronic Products.
NOTE: THESE STATEMENTS DO NOT AFFECT YOUR STATUTORY RIGHTS AS A CONSUMER

This product has been manufactured under the strictest or quality control procedures, and using the highest quality of materials, to ensure excellent performance and reliability. It will give very good and long lasting service, provided it is properly used and maintained. The product is guaranteed for 1 years from the date of original purchase. If any defect arises due to a faulty materials or poor workmanship, the faulty product must be returned to the original place of purchase. Refund or replacement is at the discretion of that company.

Duronic Products are offered with a 1 year guarantee under the following conditions:

1. The product must be returned to the retailer with original proof of purchase.
2. The product must be installed and used in accordance with the instructions contained in this operating manual.
3. It must be used only for domestic purposes.
4. It does not cover wear and tear, malicious damage, misuse, neglect, inexpert repairs or consumable parts.
5. Shine-Mart Ltd has no responsibilities for accidental or consequential loss or damage.
6. Shine-Mart Ltd are not liable to carry out any type of servicing work, under the guarantee
7. Valid in the UK only

Imported by Shine-Mart Ltd, RM3 8SB

Information on waste Disposal for Consumers of
Electrical & Electronic Equipment.

This mark on a product and/or accompanying documents indicates that when it is to be disposed of, it must be treated as Waste Electrical & Electronic Equipment, (WEEE).

Any WEEE marked products must not be mixed with general household waste, but kept separate for the treatment, recovery and recycling of the materials used.

For proper treatment, recovery and recycling; please take all WEEE marked waste to your Local Authority Civic waste site, where it will be accepted free of charge.

If customers dispose of Waste Electrical & Electronic Equipment correctly, they will be helping to save valuable resources and preventing any potential negative effects upon human health and the environment, of any hazardous materials that the waste may contain.

Follow Us

duronic.com

[duronic](https://www.youtube.com/duronic)

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1 anno di garanzia da Shine-Mart Ltd, proprietari del marchio commerciale Duronic.
5. Shine-Mart Ltd non è responsabile per Danni accidentali o consequenziali.
6. Shine-Mart Ltd non è tenuta a eseguire qualsiasi tipo di riparazione, ai sensi della garanzia.
importata da Shine-Mart Ltd, RM3 8SB